

Corsano App User Guide


[Home](#) » [Corsano](#) » Corsano App User Guide 

Contents

- 1 Corsano App
- 2 Product Usage Instructions
- 3 FAQ
- 4 Start Monitoring New Patient
- 5 Documents / Resources
 - 5.1 References

Corsano

Corsano App

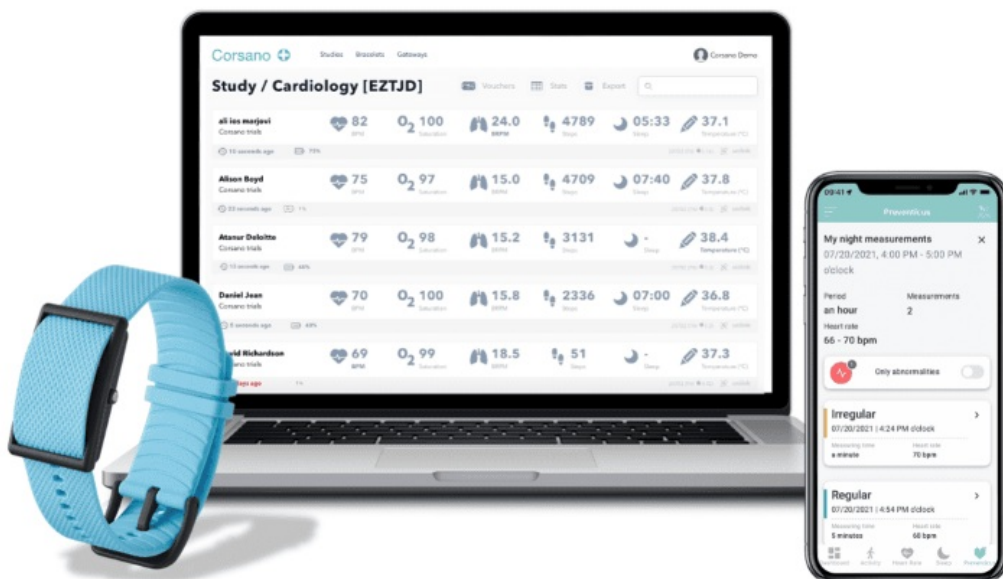

Product Information

This product is a monitoring system that allows users to track and monitor patients using a mobile app and a bracelet. It provides features for registering new patients, entering voucher codes, setting passwords, and initiating the monitoring process.

Specifications

- **Compatibility:** Works with the mobile app
- **Functionality:** Monitoring new patients, entering voucher codes, setting passwords
- **Accessories:** Bracelet for patient identification

Product Usage Instructions

Starting Monitoring for a New Patient

1. Open the app main menu and navigate to Advanced settings.
2. Select “New Patient” from the options.
3. Read the instructions provided and follow them to continue.
4. Start the process as directed.
5. Wait for the process to run.
6. Press “Continue” to register the new patient.

Registering a New Patient with a Voucher Code

1. Choose whether you want to enter the Voucher Code manually or by scanning the QR Code.
2. Enter the Voucher Code in the designated field.
3. Set a password for the new patient.
4. Press “Sign Up” to complete the patient registration process.
5. Your system is now ready to monitor the new patient.

FAQ

Q: Can I monitor multiple patients with the same app?


A: Yes, you can monitor multiple patients using the same app by following the steps outlined in the user manual for registering new patients.

Q: How do I reset a patient’s password?

A: To reset a patient’s password, navigate to the patient’s profile in the app, select the option to reset password, and follow the on-screen instructions.


Start Monitoring New Patient

If you want to start monitoring a new patient with same app and same bracelet.

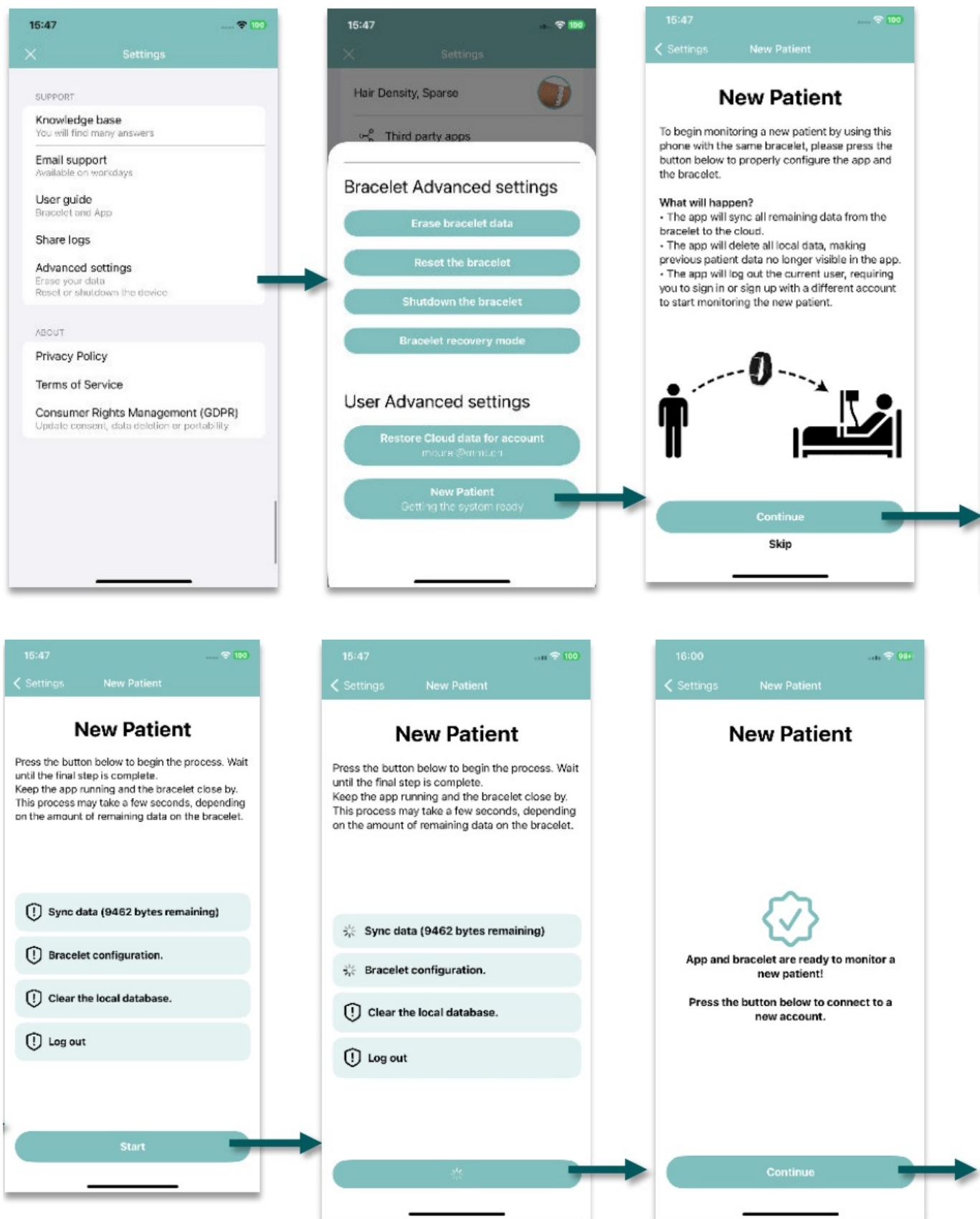


1. Open "Advanced settings" in app main menu
2. Select "New Patient"
3. Read instructions and continue
4. Read instructions and start the process
5. Process is running
6. Press "Continue" to register a new patient


Register a New Patient with a voucher code

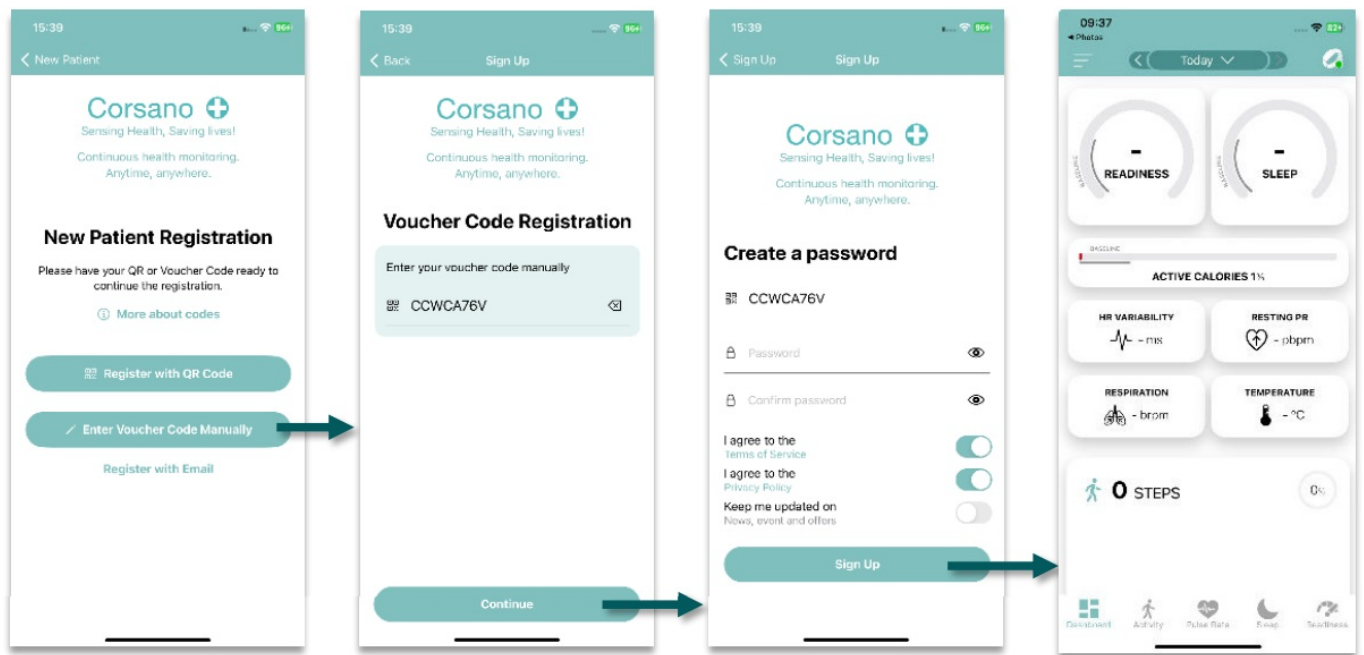

1. Select if you want to enter the Voucher Code manually or by scanning the QR Code
2. Enter the Voucher Code
3. Set a password and press “Sign Up” to finish the Patient registration.
4. Ready to monitor new patient


If you want to start monitoring a new patient with same app and same bracelet.

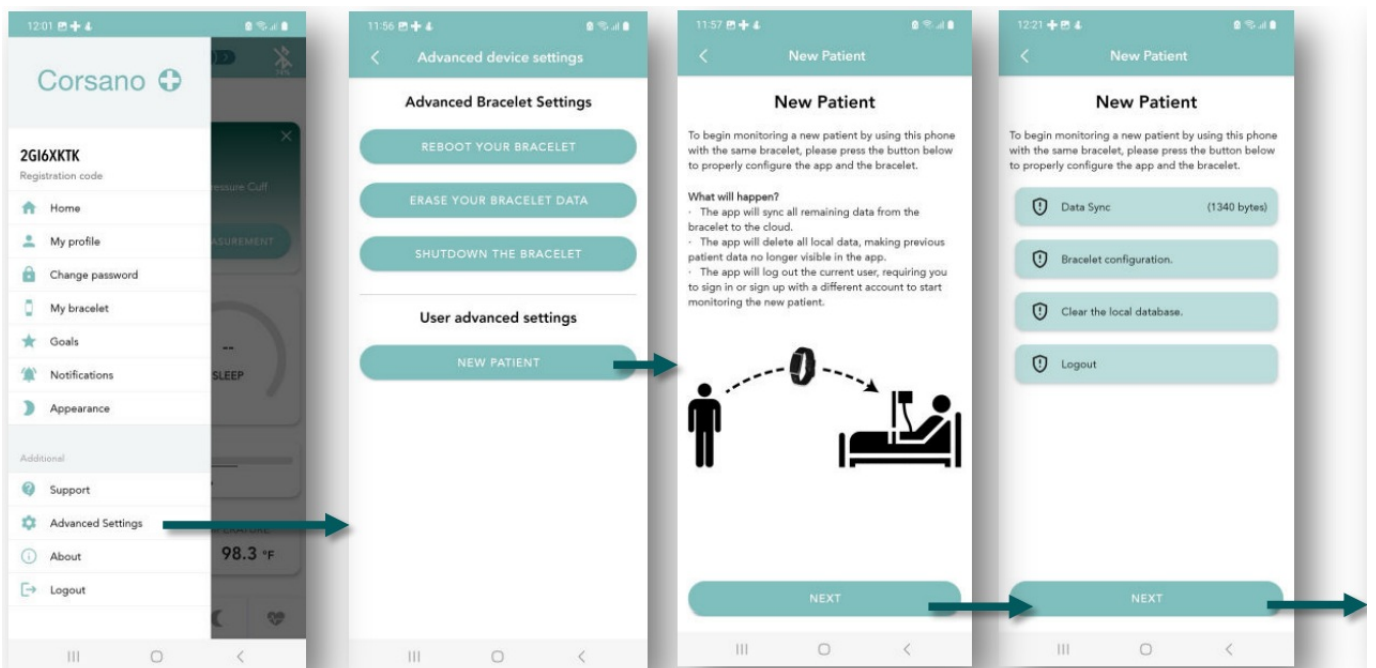

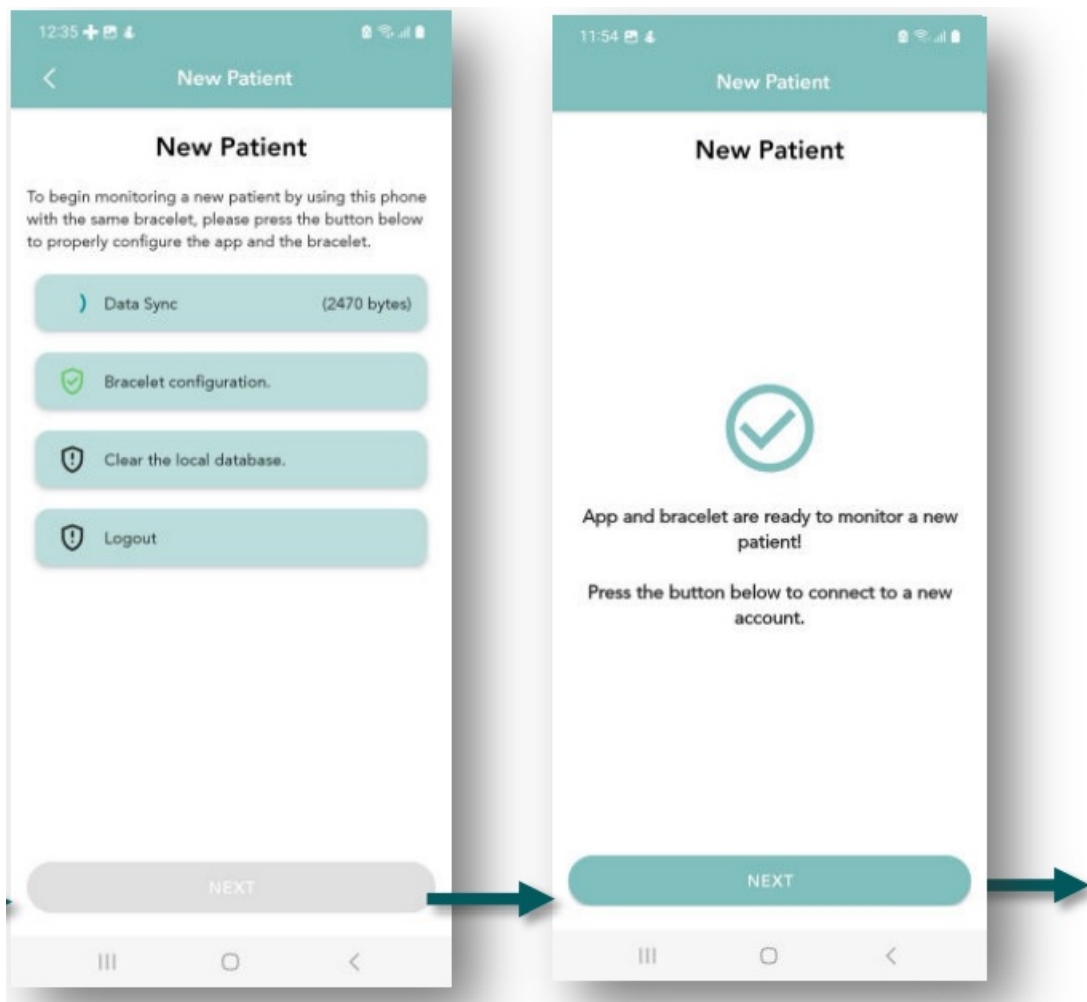


1. Open "Advanced settings" in app main menu
2. Select "New Patient"
3. Read instructions and continue
4. Read instructions and start the process
5. Process is running
6. Press "Continue" to register a new patient

Register a New Patient with a voucher code

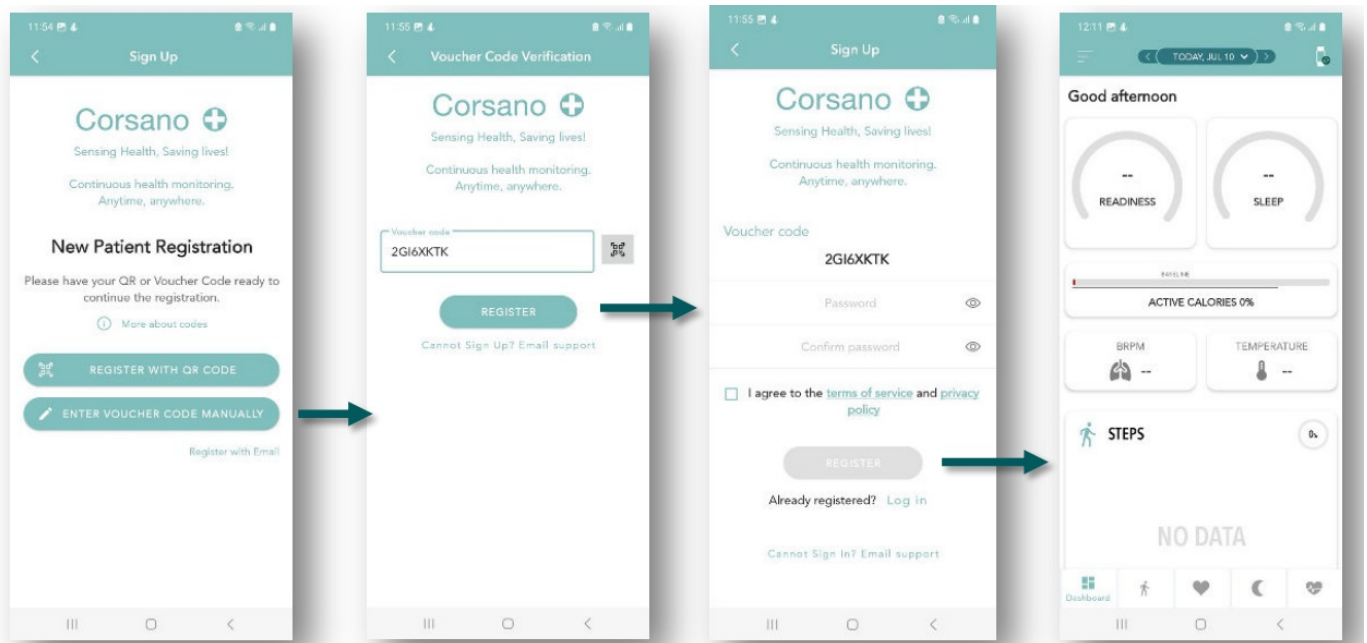

1. Select if you want to enter the Voucher Code manually or by scanning the QR Code
2. Enter the Voucher Code
3. Set a password and press “Sign Up” to finish the Patient registration.
4. Ready to monitor new patient

Documents / Resources

	corsano Corsano App [pdf] User Guide Corsano App, App
---	--

References

- [User Manual](#)

Manuals+, Privacy Policy

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.