



FERROX 260216 Parallel Bars Easy Click Klassic Instruction Manual

[Home](#) » [FERROX](#) » FERROX 260216 Parallel Bars Easy Click Klassic Instruction Manual 

Contents

- [1 FERROX 260216 Parallel Bars Easy Click Klassic](#)
- [2 Product Information](#)
- [3 Product Usage Instructions](#)
- [4 WHATS IN THE BOX](#)
- [5 REQUIRED TOOLS](#)
- [6 INSTALLATION INSTRUCTIONS](#)
- [7 WEIGHT AND CAPACITY](#)
- [8 DECLARATION OF CONFORMITY](#)
- [9 Documents / Resources](#)
 - [9.1 References](#)



FERROX 260216 Parallel Bars Easy Click Klassic



Product Information

Specifications:

- **Art.№:** 3447000 , Art.260216
- **Fabbricante:** FERROX S.r.l.
- **Nome:** PARALLELA EASY CLICK KLASSIC
- **Codice:** 260216.PN (01)08053470432602164U
- **Classe:** I (Regola 1)

Description:

Foldable parallel bars for rehabilitation use. Designed for patients immobilized by simple fractures or those undergoing motor recovery after minor trauma. Suitable for adults who require short periods of walking rehabilitation. Can be used in medical clinics with limited spaces or at home.

Intended Users:

All healthcare professionals, physiotherapists, and related industries. The device is to be used under the instructions given by the reference healthcare professional.

Contraindications:

Follow the safety instructions in the user manual. Not suitable for sports use.

Other Indications:

The device is not considered an essential item.

Manufacturer Information: FERROX S.r.l.

Procedure for the Assessment of Conformity:

Technical documentation as per annexes II and III of the MDR 2017/745.

Safety Work Load:

Please consult the user manual for the recommended safe working load.

Product Usage Instructions

Step 1: Unfolding the Parallel Bars

Start by ensuring that the parallel bars are securely placed on a flat surface. Unlock the folding mechanism by pressing the release buttons located on each side of the bars. Gently unfold the bars until they are fully extended.

Step 2: Adjusting the Height

The parallel bars can be adjusted to different heights to accommodate the user's needs. Locate the height adjustment knobs on each side of the bars. Loosen the knobs and slide the bars up or down to the desired height. Tighten the knobs to secure the bars in place.

Step 3: Using the Parallel Bars

Stand facing the parallel bars with your feet shoulder-width apart. Place your hands on the bars at a comfortable height. Slowly and carefully step forward, using the bars for support. Take small steps and maintain a steady pace. Keep your back straight and engage your core muscles for stability.

Step 4: Storing the Parallel Bars

Once you have finished using the parallel bars, make sure to clean them if necessary. Fold the bars back into their

compact position by pressing the release buttons and gently folding them together. Secure any locking mechanisms in place to prevent accidental unfolding.

Step 5: Maintenance and Care

To ensure the longevity of your parallel bars, regularly inspect them for any signs of damage or wear. Clean them with a mild detergent and warm water, and dry them thoroughly before storing. Avoid exposing the bars to excessive moisture or extreme temperatures.

FAQ:

Can these parallel bars be used for sports activities?

No, these parallel bars are not suitable for sports use. They are designed specifically for rehabilitation purposes.

What is the maximum weight capacity of the parallel bars?

Please consult the user manual for the recommended safe working load and weight capacity of the parallel bars.

Can these parallel bars be used by children?




These parallel bars are intended for use by adults who require short periods of walking rehabilitation. Please consult a healthcare professional for guidance on the suitability of these bars for children.

IT IS IMPERATIVE TO FOLLOW THE INSTRUCTIONS BELOW TO ENSURE CORRECT ASSEMBLY, INSTALLATION AS WELL AS CORRECT USE OF THE PARALLEL BARS 260216 OBLIGATION TO ANCHOR TO THE FLOOR, ACCORDING TO DIN EN 12182.






WHATS IN THE BOX




REQUIRED TOOLS

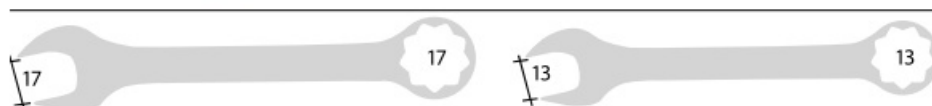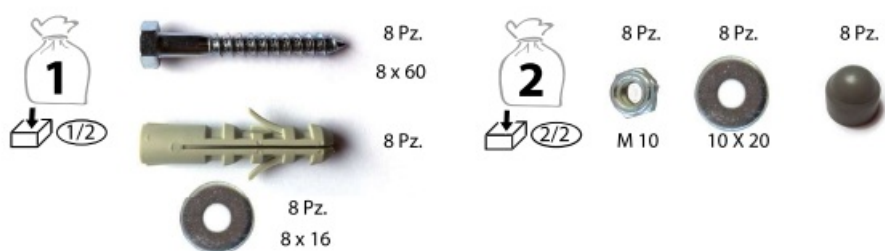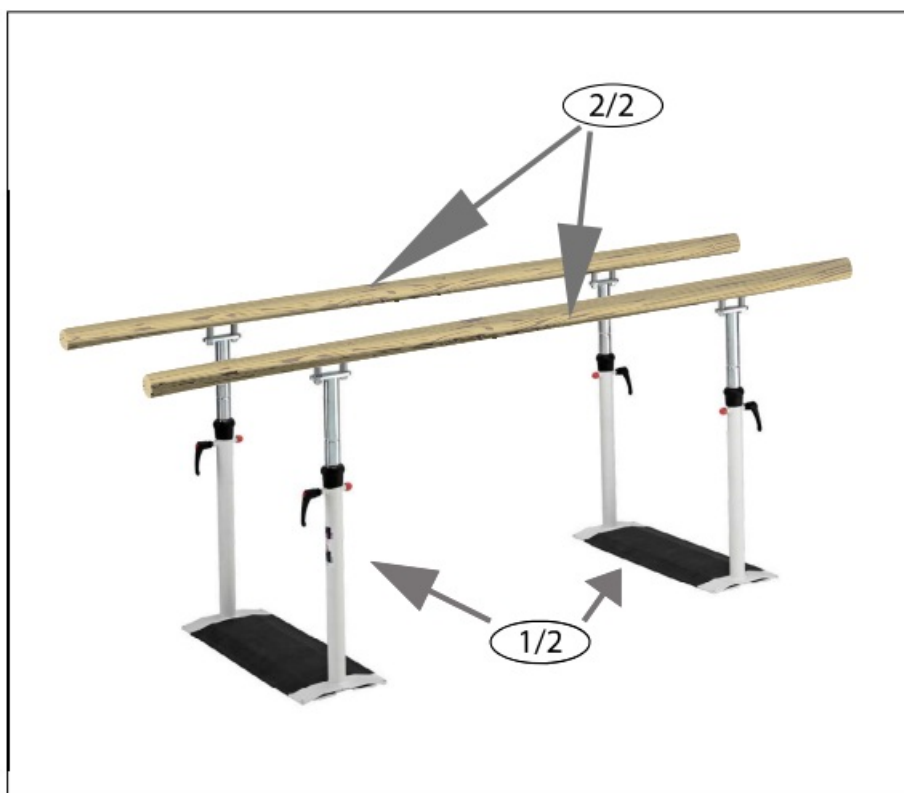




INSTALLATION INSTRUCTIONS

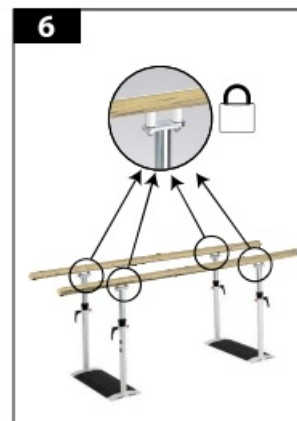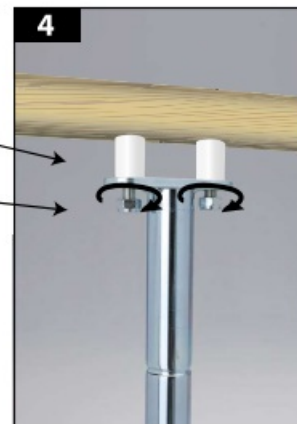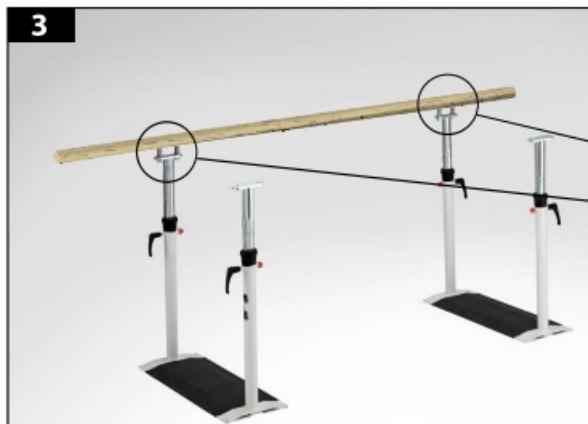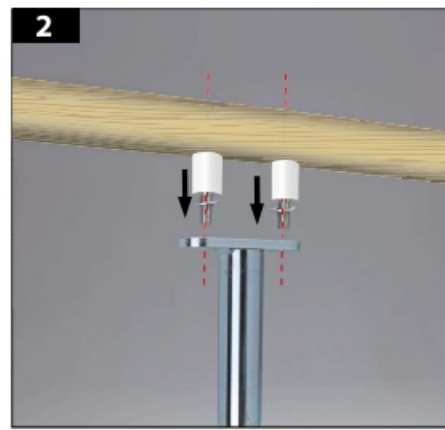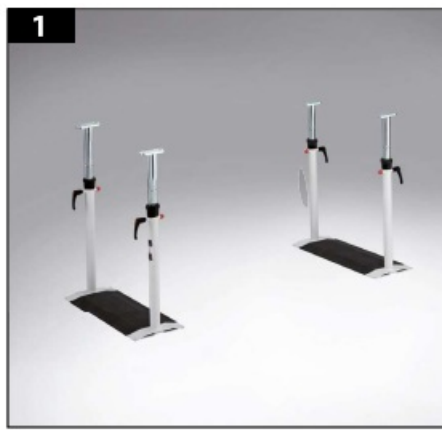


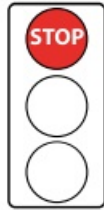


THE NOTCH POSITIONING FOR
HEIGHT ADJUSTMENT IS: **WRONG**.

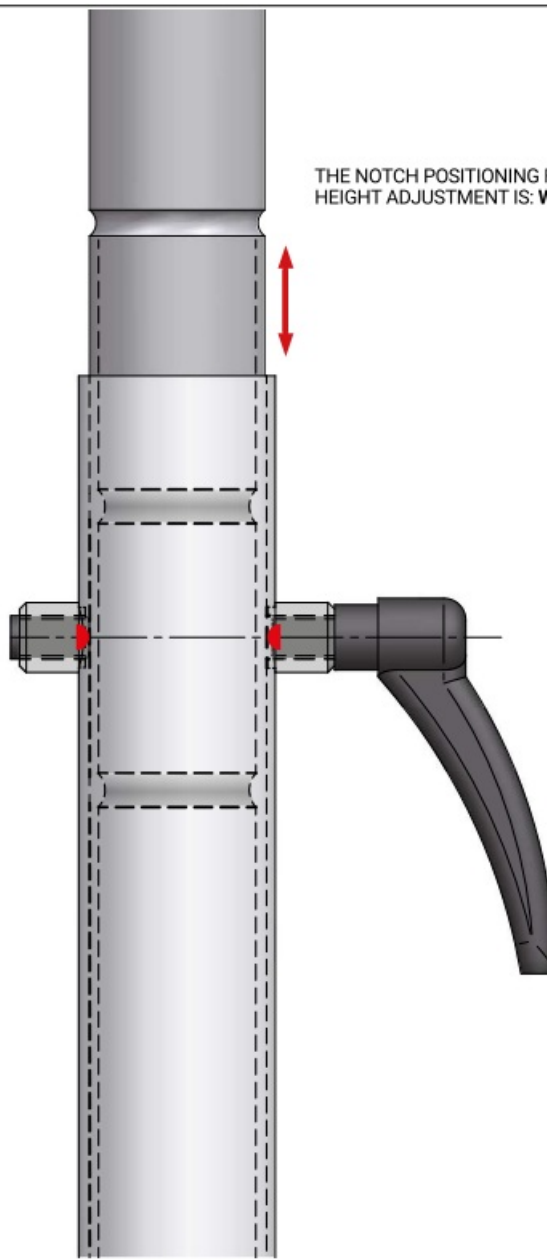

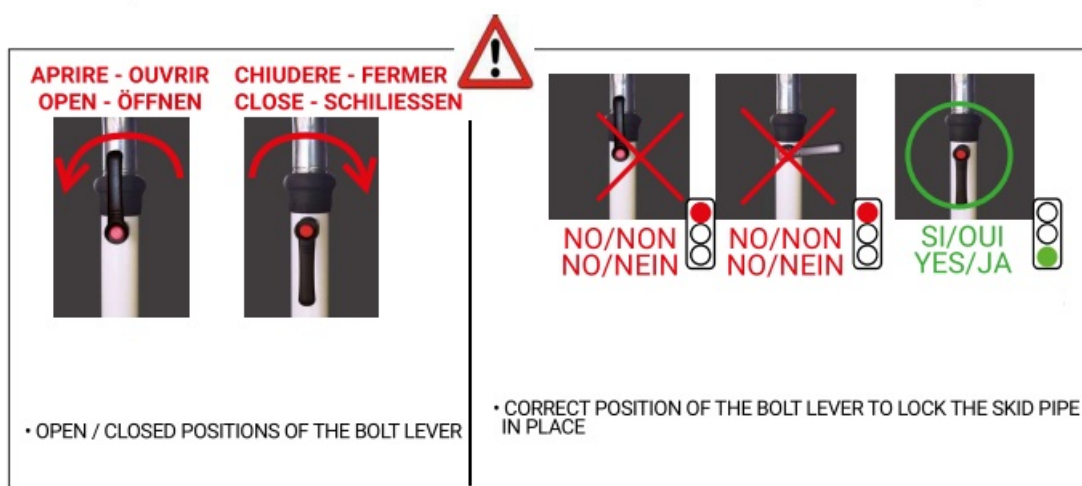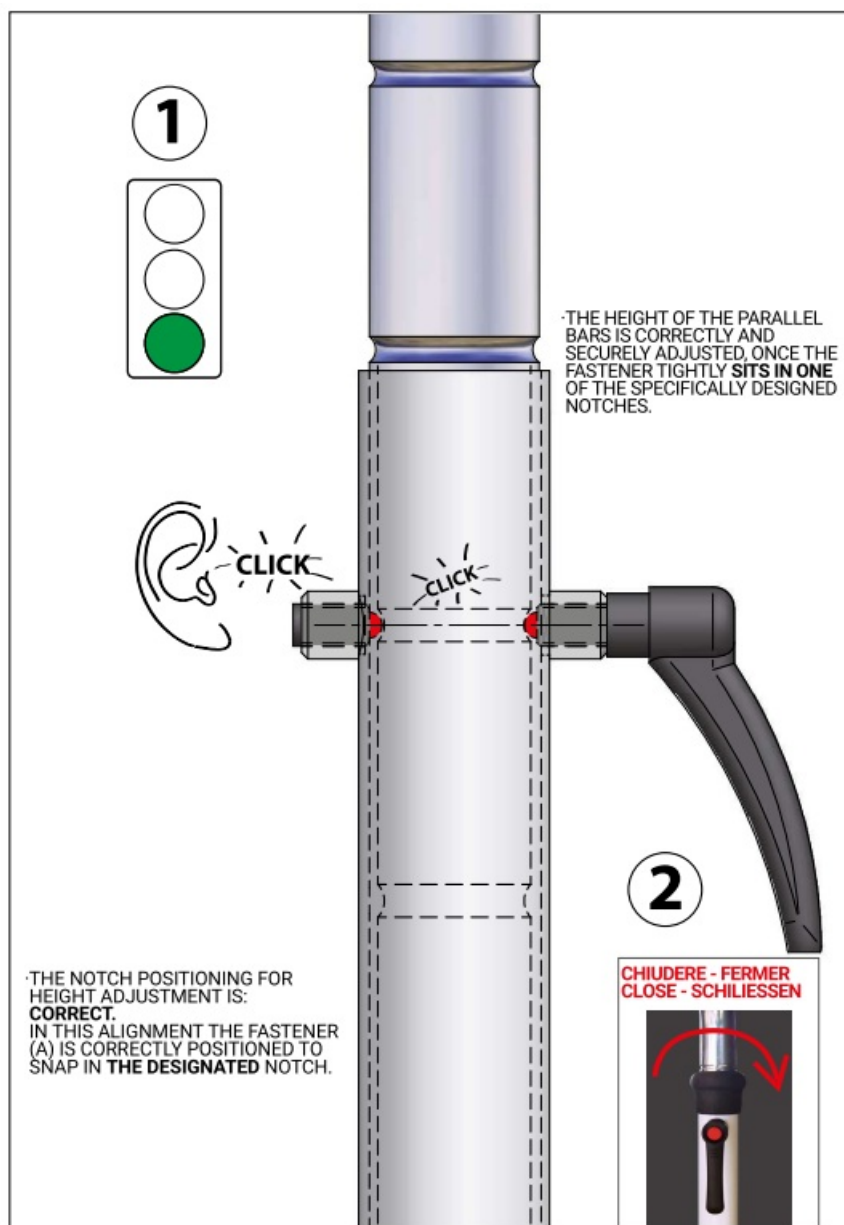

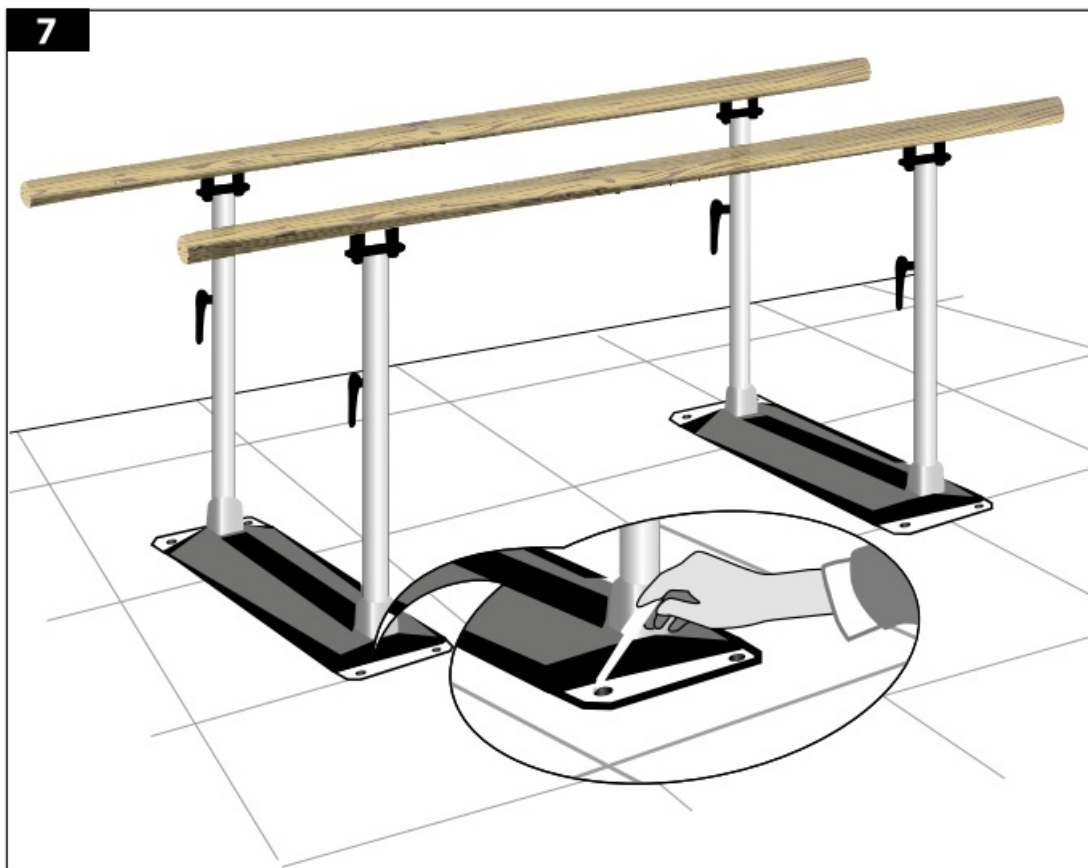



8 A

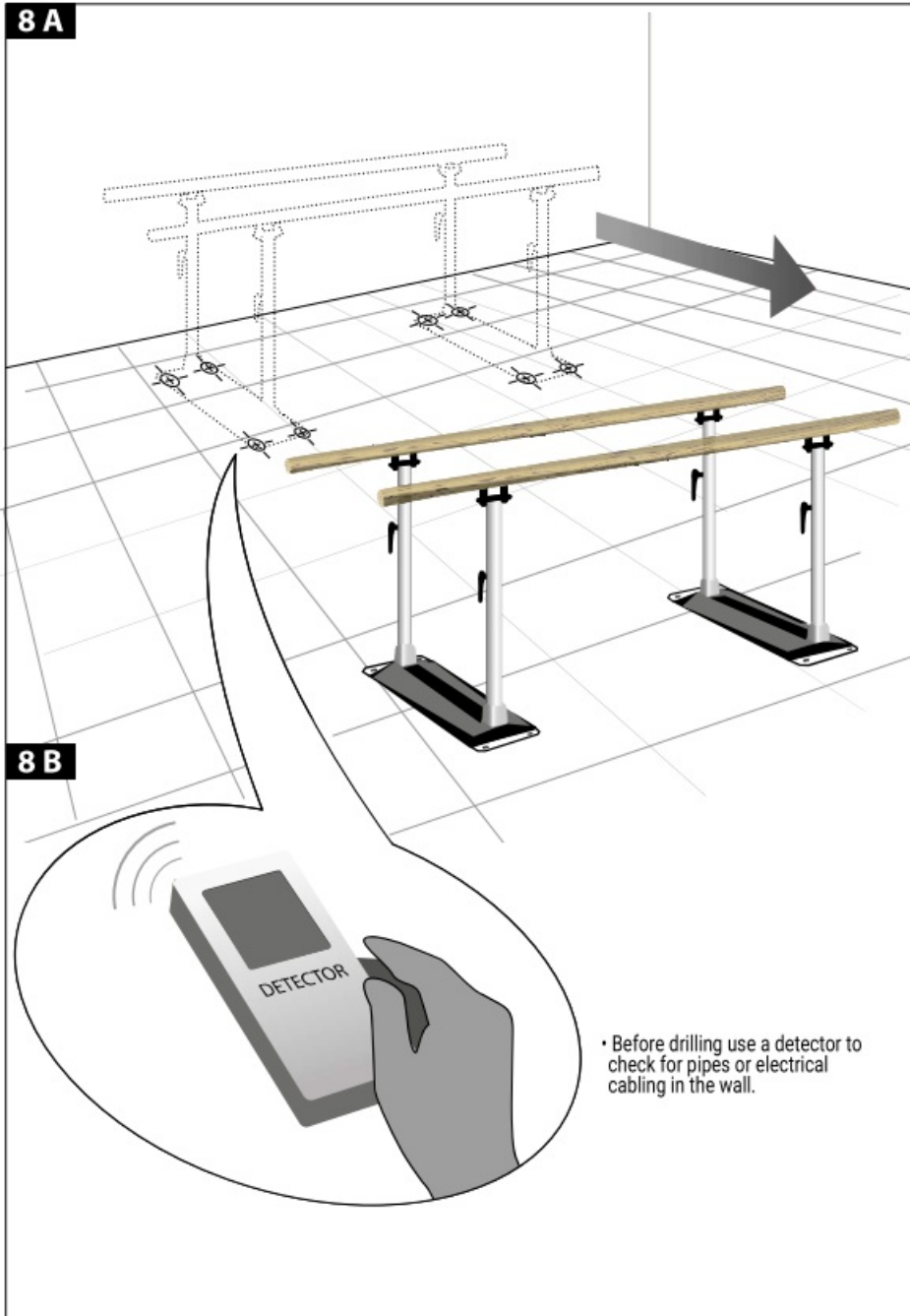


8 B



- Before drilling use a detector to check for pipes or electrical cabling in the wall.

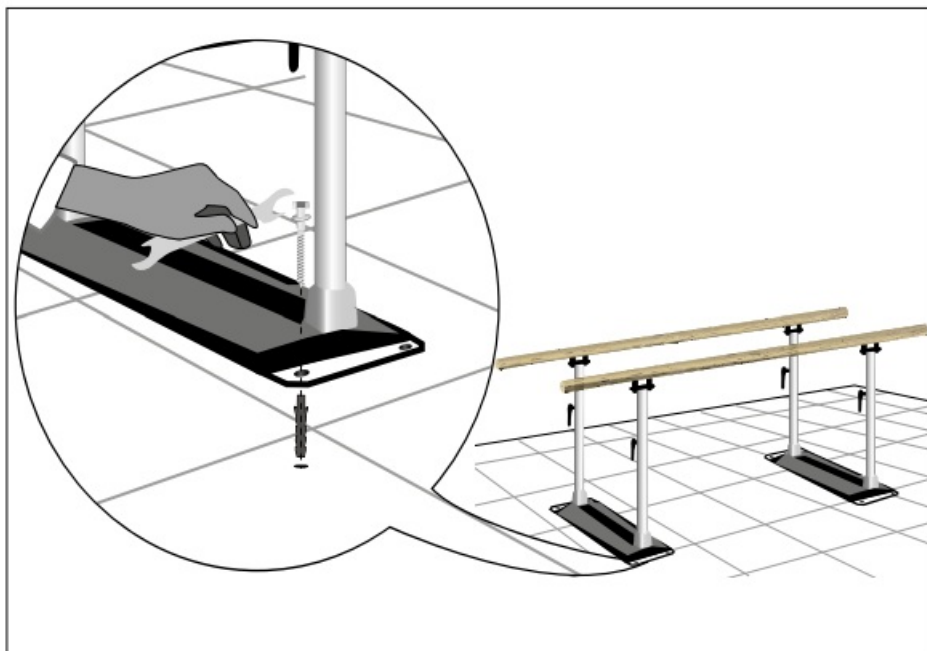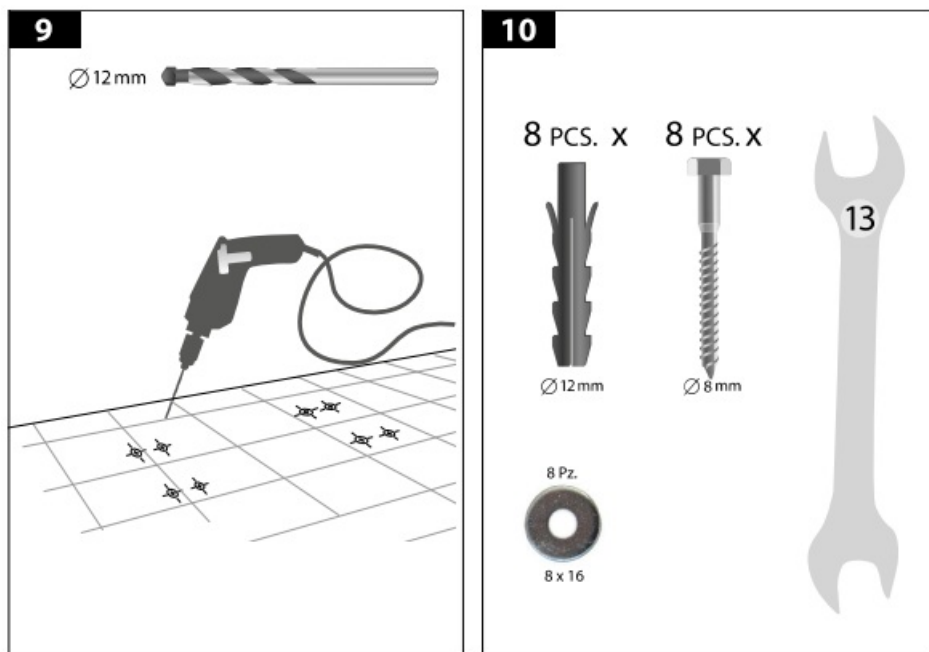

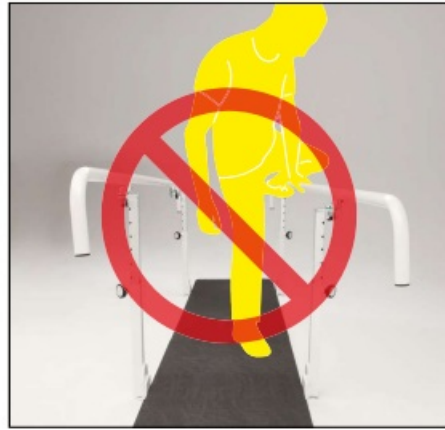



WARNING:



WEIGHT AND CAPACITY



SYMBOLS

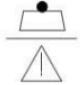SYMBOL.	Description
MD	Medical Device
REF	Catalogue Number
SN	Serial Number
	Manufacturers name
	Date of manufacture
	Safety work load
	Read the manual
CE	CE mark MDR 2017/745

DECLARATION OF CONFORMITY

MDR 2017/745

Manufacturer: FERROX S.r.l. Via del Lavoro, 9 Z. I. – 31013 ITALIA +39 0438 777099 +39 0438 777091 C. F. P.IVA R. I. IT 03613220262 info@ferrox.it

Under its sole responsibility state that the products

Name	Code	
PARALLEL BARS “EASY CLICK” KLASSIC	260216. PN	(01)08053470432602164U

Meet all the essential applicable essential requirements, of Annex I, implemented by the (UE) REGULATION 2017/745 OF THE EUROPEAN PARLAMENT AND COUNCIL on Apr 5 th 2017.

- Classified (Annex VIII, MDR 2017/745): Class I (Rule 1)
- Procedure for the assessment of conformity:
- Date:
- Technical documentation as per annex II and III.
- Responsible FERROX S.r.l.
- Signature:

CLASSIFICATION AND USE FOLDING PARALLEL BARS EASY CLICK:

- 260215,
- 260216,
- 260217,
- 260218,
- 260219

INTENDED USE

Foldable parallel bars for rehabilitation use. For patients immobilized by simple fractures, motor recovery after minor trauma.

TYPE OF PATIENT

Adults who need short periods of walking rehabilitation.

PLACE OF USE

Medical clinics with limited spaces, home.

OPERATORS QUALIFIED TO USE THE PRODUCT

All healthcare professionals, physiotherapists, and related industries. The instructions given to the patient by the reference healthcare professional suffice.

CONTRAINDICATIONS

Follow the safety instructions in the user manual. Not suitable for sports use.

OTHER INDICATIONS

The device is not an essential item.


Identification Lable

 <p>FERROX s.r.l. Via del Lavoro, 9 Z.I. 31013 Codognè (TV) - Italia Tel. +39 0438 777099 Fax +39 0438 777091 www.ferrox.it - info@ferrox.it</p>	MD MODEL	 MADE IN ITALY
	REF YYYY	
	SN XXXXXX	
	 YY/MM/DD	
	 INPUT: / - / - IP: /	
UDI 	 = 140Kg	 (01)0805347043YYYY(11)YYMMDD(21)XXXXXX Rev.02 del 03/20

Lavoro Z.I. 31013 Codogne (TV) Italia Tel. +39 0438 777099 · 778370 Fax +39 0438 777091 P.IVA: IT 03613220262 e-mail: info@ferrox.it www.ferrox.it



Documents / Resources

	<p>FERROX 260216 Parallel Bars Easy Click Klassic [pdf] Instruction Manual 260215, 260216, 260217, 260218, 260219, 260216 Parallel Bars Easy Click Klassic, Parallel Bars Easy Click Klassic, Bars Easy Click Klassic, Easy Click Klassic, Click Klassic, Klassic</p>
-------------------------------------------------------------------------------------	-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

References

- [User Manual](#)