




HANNA instruments HI9810322 Wireless pH Tester for Cheese with Built-In Specialized Electrode Instruction Manual

[Home](#) » [HANNA instruments](#) » HANNA instruments HI9810322 Wireless pH Tester for Cheese with Built-In Specialized Electrode Instruction Manual 



HI9810322

Wireless pH Tester for Cheese with a built-in specialized electrode
Foodcare
INSTRUCTION MANUAL



Dear Customer,

Thank you for choosing a Hanna Instruments product.

Please read this instruction manual carefully before using this instrument.

This manual will provide you with the necessary information for the correct use of this instrument, as well as a precise idea of its versatility.

If you need additional technical information, do not hesitate to e-mail us at tech@hannainst.com or view our contact list at www.hannainst.com.

All rights are reserved. Reproduction in whole or in part is prohibited without the written consent of the copyright owner, Hanna Instruments Inc., Woonsocket, Rhode Island, 02895, USA.

Contents

1	PRELIMINARY EXAMINATION
2	SPECIFICATIONS
3	GENERAL DESCRIPTION & INTENDED USE
4	FUNCTIONAL DESCRIPTION & LCD DISPLAY
5	GENERAL OPERATIONS
6	SETUP
7	BLUETOOTH
8	HANNA LAB APPLICATION
9	CALIBRATION
10	CARE & MAINTENANCE
11	WARNING & ERROR MESSAGES
12	ABBREVIATIONS
13	ACCESSORIES
14	RECOMMENDATIONS FOR USERS
15	WARRANTY
16	CERTIFICATION
17	Documents / Resources
17.1	References
18	Related Posts

PRELIMINARY EXAMINATION

Remove the tester and accessories from the packaging and examine them carefully. For further assistance, please contact your local Hanna Instruments office or email us at tech@hannainst.com. Each HI9810322 is supplied with a starter kit consisting of:

- pH 4.01 buffer solution, 20 mL sachet (2 pcs.)
- pH 7.01 buffer solution, 20 mL sachet (2 pcs.)
- Cleaning solution for cheese deposits, 20 mL sachet (2 pcs.)
- Electrode storage solution, 13 mL dropper bottle (1 pc.)
- 3V Lithium battery – CR2032
- Instrument quality certificate and Instruction manual

Note: Save all packing material until you are sure that the tester works correctly. Any damaged or defective item must be returned in its original packing material with the supplied accessories.

SPECIFICATIONS

Range	pH	0.00 to 12.00 pH
	mV *	pH/mV conversion
	Temperature **	0.0 to 60.0 °C (32.0 to 140.0 °F)
Resolution	pH	0.01 or 0.1 pH
	mV *	0.1 or 1 mV
	Temperature	0.1 °C; 0.1 °F
Accuracy	pH	±0.05 pH
	Temperature	±0.5 °C; ±0.9 °F
Calibration	Up to three points or four points * Automatic buffer recognition with Standard buffers Hanna (pH 1.68 *, 4.01, 7.01, 10.01) or NIST (pH 1.68 *, 4.01, 6.86, 9.18)	
Temperature compensation	Automatic (ATC) or Manual (MTC) *	
Electrode	Body material	Polyvinylidene Fluoride (PVDF)
	Glass	Low Temperature (LT)
	Junction	Open
	Reference cell	Double, Ag/AgCl
	Electrolyte	Viscolene
	Tip / Shape	Conic, Ø 6 x 10 mm (Ø 0.23 x 0.39")
	Outer diameter	13 mm (0.51")
	Length	75 mm (2.95")
Battery type	3V Lithium – CR2032	
Battery life	Approximately 1000 hours (500 hours with Bluetooth enabled)	
Environment	0 to 50 °C (32 to 122 °F)	
Casing	IP65 ingress protection	
Dimensions	51 x 145 x 21 mm (2.0 x 5.7 x 0.8")	
Weight	45 g (1.6 oz.)	

* Available with Hanna Lab App

** Measuring outside the recommended operating temperature range may damage the gel electrolyte and void product warranty.

Note: The tester can display measurements from –2.00 to 16.00 pH. Measurements outside of the pH range will flash. In this case, assess the integrity of the tester and the type of measured sample.

GENERAL DESCRIPTION & INTENDED USE

HI9810322 is a professional wireless pH tester, part of the Hanna Instruments HALO2 family.

- The integrated Bluetooth module allows the tester to be connected to a compatible smart device with the

Hanna Lab App.

- The tester has a compact, waterproof casing, and automatic pH calibration at up to three points, or four points when used with the Hanna Lab App.
- Automatically compensated temperature readings are displayed on a large LCD.
- Accurate and easy to use, the tester is designed to measure and monitor the pH levels in all steps of the cheese manufacturing process, while meeting the food-hygiene and Hazard Analysis Criteria Control Point (HACCP) regulations.

Operating Modes

The HI9810322 can be used as a stand-alone pH tester or connected to the Hanna Lab App.

The Hanna Lab App turns a compatible smart device into a full-featured pH meter. Features include electrode condition, GLP with time-stamp, live readings, mV resolution, manual temperature compensation, stability criteria, calibration reminder, pH (mV) and temperature alarms, tester ID, and data sharing.

Probe Features

The food-grade PVDF body material is easy to clean and disinfect. Resistant to most chemicals (e.g. solvents, sodium hypochlorite), ultraviolet light, and fungal growth, the PVDF body also has high-abrasion resistance and mechanical strength.


The conical tip allows for easy penetration into solids and semisolids.

The double junction design presents a silver-free electrolyte solution interacting with the sample, making the electrodeless susceptible to clogging and guaranteeing a fast response and stable reading.


Built-in temperature sensor at the tip of the pH electrode allows for rapid determination of the sample temperature and a high-accuracy temperature reading.

FUNCTIONAL DESCRIPTION & LCD DISPLAY

Front & Back view

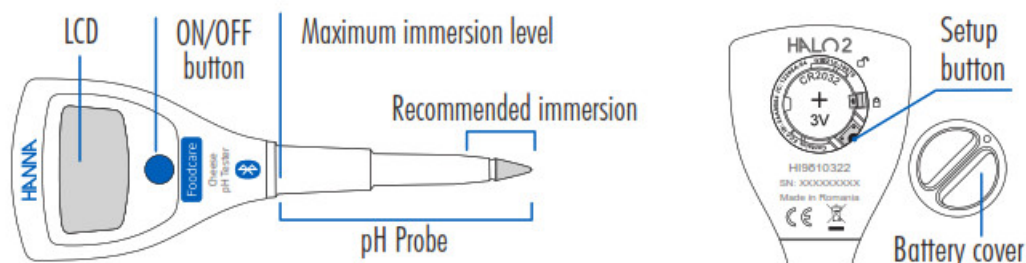

LCD display

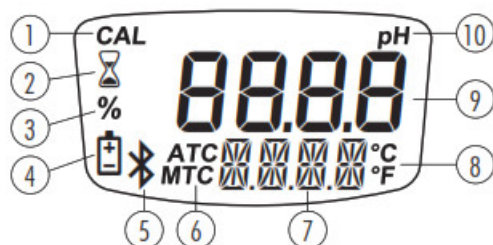



1	Calibration indicator
2	Stability indicator
3	Battery percentage indicator
4	Battery icon
5	Bluetooth icon
6	Temperature compensation
7	Second LCD line
8	Temperature unit
9	First LCD line
10	Measurement unit

GENERAL OPERATIONS


Turning the Tester On & Off

- Turn the tester face down. Rotate the battery cover on the back of the tester counterclockwise and set it aside. Remove the battery insulation film.
- Press the ON/OFF button to turn the tester on or off. The initialization screen displays all LCD segments, followed by the battery percentage. The tester enters measurement mode. Before the tester turns off, “OFF PWR” is displayed briefly.
- Press and hold the ON/OFF button to turn the tester off when connected to Bluetooth.

Battery Replacement

1. Turn off the tester. Turn the tester face down and rotate the battery cover counterclockwise.
2. Set the battery cover aside. Press the metallic pin to push old battery out.
3. Place the new battery with the positive (+) sign facing out.
4. Reset the date and time in setup, or connect to the Hanna Lab App to update it automatically.
5. Align the mark on the cover with the open lock icon () on the case. Rotate the cover clockwise until the mark on the cover aligns with the closed lock icon ().

Note: Only use specified battery type. Dispose of old battery in accordance with local regulations.

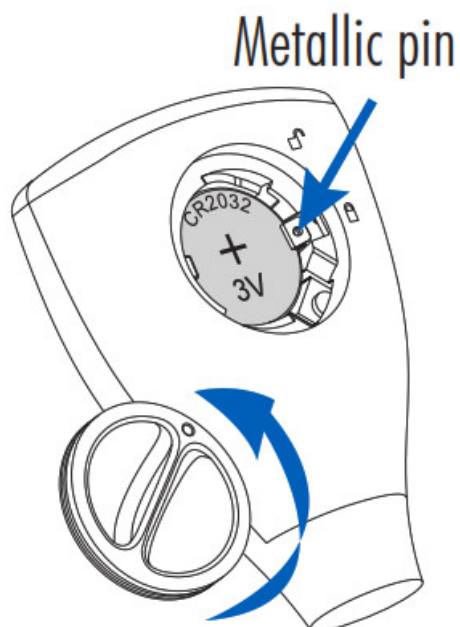

SETUP

The Setup button is located inside the battery compartment. After Setup configuration, replace the cover.

Setup Menu Navigation

- Press the Setup button to enter Setup mode and navigate menu items.
- To exit Setup Mode, press the Setup button after the “SEt TIME” option is displayed.
- Press the ON/OFF button to configure menu item options.

Temperature Unit

Option: °C or °F

Press ON/OFF button to select desired temperature unit.

Auto-Off Interval

Option: 8, 60 min., or “—” (disabled)

Press ON/OFF button to select the desired interval.

To save battery life, after the selected auto-off interval has elapsed, the tester will automatically turn off.

Calibration Points

Option: 2P or 3P

Press ON/OFF button to select between up to two- or up to three-point calibration.

Buffer Set

Option: 7.01 pH (Hanna) or 6.86 pH (NIST)

Press ON/OFF button to select the calibration buffer set (Hanna or NIST).

pH Resolution

Option: 0.01 pH or 0.1 pH


Press the ON/OFF button to select resolution.

Bluetooth Wireless Mode

Option: On, PAIr, or OFF

Press ON/OFF button to select

Bluetooth option at start-up.

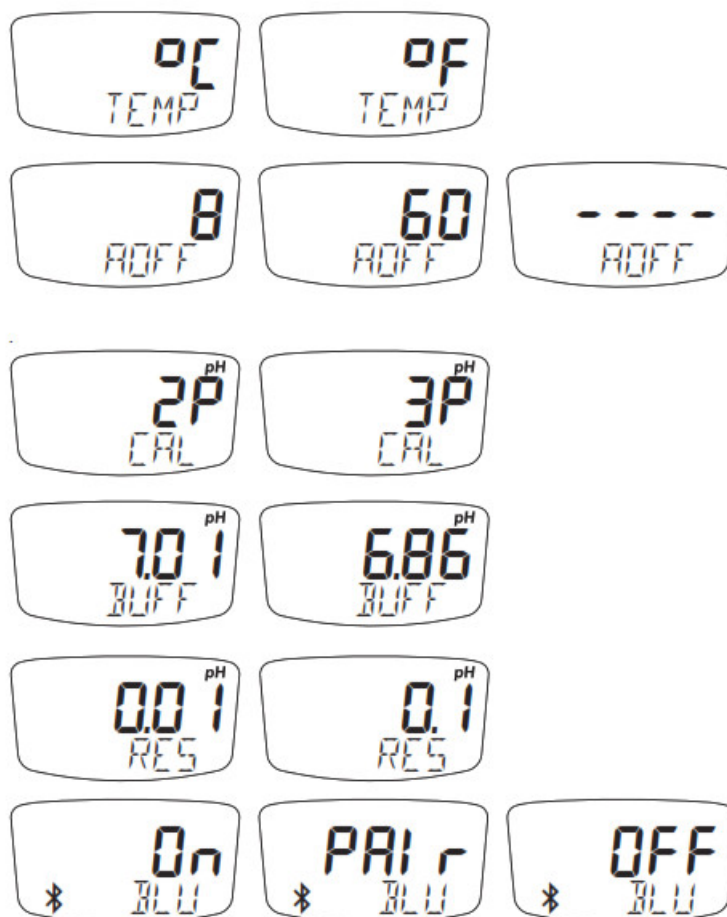

Bluetooth Pairing

Option: dEL PAIr

Press ON/OFF button to delete the saved paired device.

Date & Time

Option: SEt TIME


Press ON/OFF button to set the date and time.

Option: YEAR, MO, DAY, HOUR, and MIN

Use the Setup button to select the option and press the ON/OFF button to change the selected option.



BLUETOOTH

With “PAIr BLU” or “On BLU” selected in Setup, the Bluetooth icon () will blink for up to 45 seconds, indicating the tester is in discoverable mode. Once connected, the icon will stop blinking. If not connected, the icon is not displayed.

- Select “On BLU”, to enable Bluetooth without bonding.
- Select “PAIr BLU”, to enable Bluetooth with bonding. A 6-digit bonding pin is displayed the first time the tester

and smart device are paired. Once the devices are paired, the bonding pin is not required when reconnected.



- Select “OFF BLU”, to disable Bluetooth.
- Select “dEL PAIR”, to delete all paired devices. If PAIr BLU is enabled, a bonding pin will need to be re-entered.



HANNA LAB APPLICATION

- The Hanna Lab App is available on the App Store and on Google Play.
- Consult the Help section of the application for information on calibration, measurement, data logging, and sharing.
- When the tester is in discoverable mode, it will appear in the list of “Available Devices”.
- Within the application, tap “Connect” to pair the tester with the device. All readings are transmitted directly to the application.

Tag a Measurement

Once connected to the application, the ON/OFF button can be pressed to tag the current reading.

- Press the ON/OFF button from measurement mode. The display shows “SEt TAG”, followed by “- TAG”.
- The reading on the application will flash green and the note icon () will be displayed. Tap the note icon () to add an annotation.

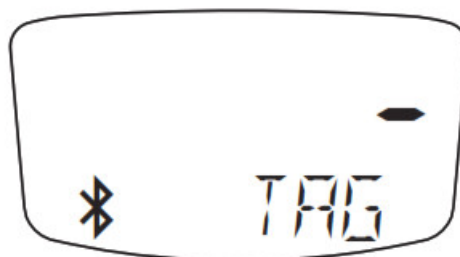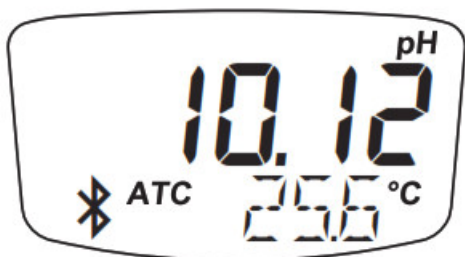

CALIBRATION

Preparation & Guidelines

1. Remove the probe from the plastic storage cap. Save the cap for probe storage.
2. Rinse off any storage solution or salts that may be on the body.
3. Verify there is a solution inside the pH bulb by shaking the probe down to restore continuity as the solution may have moved up the stem during shipping.
4. For best results, use a rinse beaker and a separate calibration beaker for each buffer. Discard rinsing buffers after use.

Procedure

For most applications, it is recommended to start with pH 7.01 buffer (or pH 6.86).

To restore factory defaults, press and hold the ON/OFF button. “CAL CLR” is displayed.


Note: It is recommended to calibrate the electrode with buffers at the temperature it will be used at.

App Store is a trademark of Apple, Inc.



Google Play and the Google Play logo are trademarks of Google LLC.

Stand-Alone Tester (Up to Three-Point Calibration)

1. Rinse the electrode tip with purified water and blot dry. Then rinse with the buffer being used for calibration.
2. Press and hold ON/OFF button until "CAL MODE" is displayed.
3. When "7.01 USE" or "6.86 USE" is displayed with "CAL" tag blinking, place the tip of the electrode in the correct buffer.

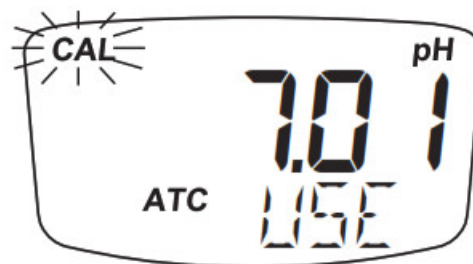

4. When the buffer has been recognized, "REC" is displayed.
"WAIT" is displayed with the stability indicator (#) blinking until the reading is stable.
Wait until the measurement is stored and the stability indicator disappears.

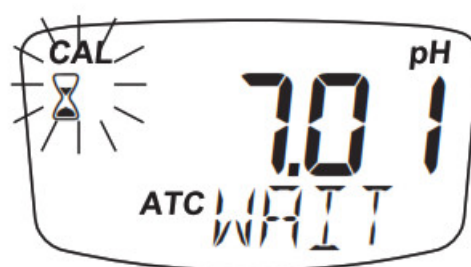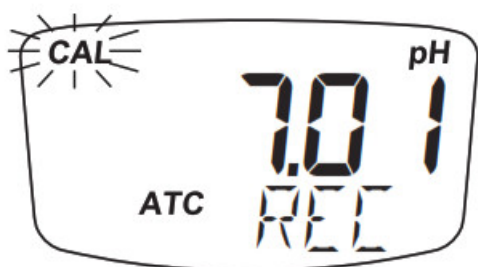

5. To save a one-point calibration and return to measurement mode, press the ON/OFF button.
"CAL SAVE" is displayed briefly.

With 2P option selected

6. Place the tip of the electrode in the second buffer rinse beaker, then in the second calibration buffer. Wait until the measurement is stored and the stability indicator disappears.
7. After the second point has been stored, "CAL SAVE" is displayed briefly and the tester will return to measurement mode automatically.

With 3P option selected

8. Place the tip of the electrode in the second buffer rinse beaker, then in the second calibration buffer. Wait until the measurement is stored and the stability indicator disappears.
9. Place the tip of the electrode in the third buffer rinse beaker, then in the third calibration buffer. Wait until the measurement is stored and the stability indicator disappears.
10. After the third point has been stored, "CAL SAVE" is displayed briefly and the tester will return to measurement mode automatically.

With Hanna Lab App (Up to Four-Point Calibration)

Connect the tester to the Hanna Lab App and follow the calibration procedure.
See the App Help section for the calibration procedure.

CARE & MAINTENANCE

- Never immerse the tester over the maximum immersion level.
- Fresh buffers should be used for each calibration. Once the sachets are opened the buffer value can change over time.
- If the electrode is slow or sluggish, soak it in the cleaning solution for 20 minutes. Rinse with water and hydrate

the electrode in a storage solution for a minimum of 30 minutes before calibrating.

- If measurements are taken successively, rinse the probe thoroughly in distilled or deionized water to eliminate cross-contamination between measurements.

Storage

To ensure a quick response, the glass tip and junction should be kept hydrated.

When not in use, add a few drops of storage solution to the storage cap. If storage solution is not available, pH 4.01pH 7.01 buffer can be used.

Note: Do not store the electrode in distilled or deionized water.

WARNING & ERROR MESSAGES

	"--- WRONG" displayed during calibration. Invalid buffer. Check the buffer value and use a fresh buffer.
	The measured value displayed blinking. The measured value is out of electrode range. Clean the electrode to improve the condition.
	Battery icon (🔋) displayed blinking. The battery is below 10 %, replace the battery.
	"BAtt DEAD" and the tester turns off. Depleted battery, replace the battery.

ABBREVIATIONS

ATC	Automatic Temperature Compensation
GLP	Good Laboratory Practice
MTC	Manual Temperature Compensation
NIST	National Institute of Standards and Technology

ACCESSORIES

Ordering Information	Product Description
pH Buffer Solution	
HI50016-02	pH 1.68 buffer solution, 20 mL sachet (25 pcs.)
HI70004P	pH 4.01 buffer solution, 20 mL sachet (25 pcs.)
HI70006P	pH 6.86 buffer solution, 20 mL sachet (25 pcs.)
HI70007P	pH 7.01 buffer solution, 20 mL sachet (25 pcs.)
HI70009P	pH 9.18 buffer solution, 20 mL sachet (25 pcs.)
HI70010P	pH 10.01 buffer solution, 20 mL sachet (25 pcs.)
HI77400P	pH 4.01 & 7.01 buffer solution, 20 mL sachet (10 pcs., 5 each)
HI770710P	pH 10.01 & 7.01 buffer solution, 20 mL sachet (10 pcs., 5 each)
Electrode Cleaning Solution	
HI700601P	Electrode cleaning solution for general use, 20 mL sachet (25 pcs.)
HI700630P	Meat, grease, and fats acid cleaning solution, 20 mL sachet (25 pcs.)
HI700642P	Cleaning solution for cheese deposits, 20 mL sachet (25 pcs.)
Electrode Storage Solution	
HI70300L	Electrode storage solution, 500 mL
HI70300M	Electrode storage solution, 230 mL
HI70300S	Electrode storage solution, 30 mL dropper
HI9072	Electrode storage solution, 13 mL dropper

RECOMMENDATIONS FOR USERS

Before using this product, make sure it is entirely suitable for your specific application and for the environment in which it is used. Any variation introduced by the user to the supplied equipment may degrade the tester's performance. For your and the tester's safety do not use or store it in hazardous environments.

WARRANTY

HI9810322 is warranted for a period of one year against defects in workmanship and materials when used for its intended purpose and maintained according to instructions. This warranty is limited to repair or replacement free of charge. Damage due to accidents, misuse, tampering, or lack of prescribed maintenance is not covered. If service is required, contact your local Hanna Instruments office. If under warranty, report the model number, date of purchase, serial number, and the nature of

the problem.

If the repair is not covered by the warranty, you will be notified of the charges incurred. If the tester is to be returned to Hanna Instruments, first obtain a Returned Goods Authorization (RGA) number from the Technical Service department and then send it with shipping costs prepaid. When shipping any product, make sure it is shipped back in compliance with shipping regulations, thoroughly cleaned and free of any chemicals, and is properly packaged for complete protection.

CERTIFICATION



All Hanna Instruments conform to the CE European Directives.

Disposal of Electrical and Electronic Equipment. The product should not be treated as household waste.

Instead, hand it over to the appropriate collection point for the recycling of electrical and electronic equipment, which will conserve natural resources.



Disposal of waste batteries. This product contains batteries, do not dispose of them with other household waste. Hand them over to the appropriate collection point for recycling. Ensuring proper product and battery disposal prevents potential negative consequences for the environment and human health. For more information, contact your city, your local household waste disposal service, or the place of purchase.



Regulatory Notices for the stand-alone, Bluetooth, low-energy module.

United States (FCC) FCC ID: 2AA9B04. This device complies with FCC Rules, Part 15 Subpart C "Intentional Radiators" and Subpart B, Chapter §15.105. This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case, users are required to correct the interference at their own expense.

Canada (ISED) IC: 12208A-04. This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

Australia / New Zealand (RCM) BMD-300 complies with the AS/NZS 4268:2017.

Japan (MIC)  R 210-106799

South Korea (KCC)  R-CRM-Rgd-BMD-300

Brazil (ANATEL): Contains ANATEL approved module # 00820-21-05903.

Hanna Instruments reserves the right to modify the design, construction or appearance of its products without advance notice.







World Headquarters

Hanna Instruments Inc.
Highland Industrial Park
584 Park East Drive
Woonsocket, RI 02895 USA

www.hannainst.com

MAN9810322

Printed in ROMANIA

 	<p>HANNA instruments HI9810322 Wireless pH Tester for Cheese with Built-In Specialized Electrode [pdf] Instruction Manual</p> <p>HI9810322, Wireless pH Tester for Cheese with Built-In Specialized Electrode, HI9810322 Wireless pH Tester for Cheese with Built-In Specialized Electrode, Built-In Specialized Electrode, HI9810322 Wireless pH Tester for Cheese, Wireless pH Tester for Cheese, pH Tester for Cheese</p>
 	<p>HANNA instruments HI9810322 Wireless pH Tester [pdf] Instruction Manual</p> <p>HI9810322, Wireless pH Tester, HI9810322 Wireless pH Tester, pH Tester, Tester</p>

References

-  [pH Meters, Photometers, Titrators, Controllers - Hanna Instruments](#)