


5X BUFFER B PH 7.2

running buffer for **proFIRE**

Dynamic Biosensors GmbH & Inc.
BU-P-1000-5 v2.1

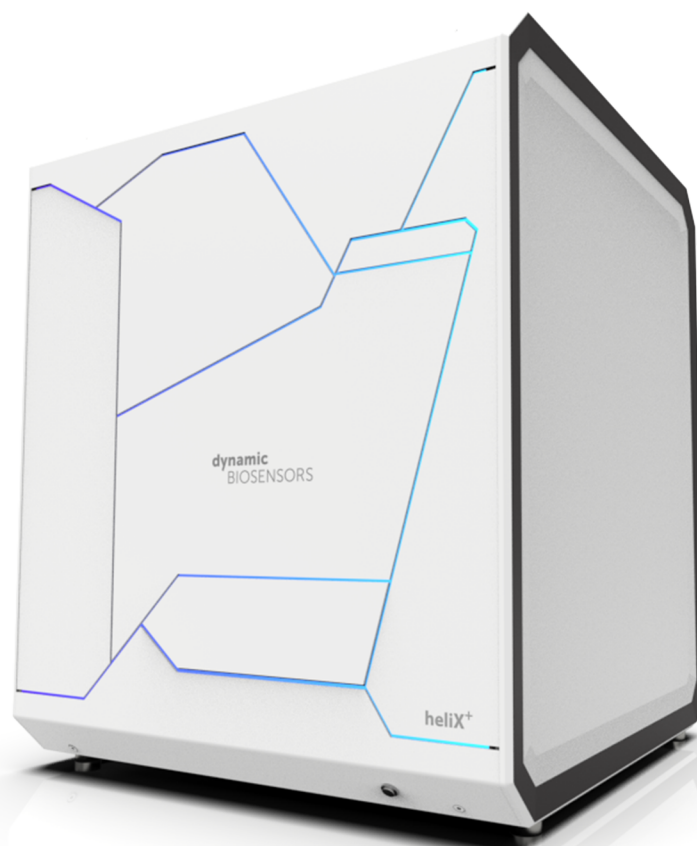

Product Description

Order Number: **BU-P-1000-5**

Table 1. Contents and Storage Information

Material	Composition	Amount	Storage
5x Buffer B pH 7.2	250 mM $\text{Na}_2\text{HPO}_4/\text{NaH}_2\text{PO}_4$, 5 M NaCl; 0.2 μm sterile filtered	50 mL	RT

For research use only.

This product has a limited shelf life, please see expiry date on label.

At temperatures lower than 18°C precipitation may occur.

Preparation

Dilute the complete solution 5x Buffer B pH 7.2 (50 mL) by mixing with 200 mL ultrapure water.

After dilution Buffer B is ready for use (50 mM $\text{Na}_2\text{HPO}_4/\text{NaH}_2\text{PO}_4$, 1 M NaCl).

The diluted buffer should be stored at 2-8°C.

Contact

Dynamic Biosensors GmbH
Perchtinger Str. 8/10
81379 Munich
Germany

Dynamic Biosensors, Inc.
300 Trade Center, Suite 1400
Woburn, MA 01801
USA

Order Information order@dynamic-biosensors.com
Technical Support support@dynamic-biosensors.com

www.dynamic-biosensors.com

Instruments and chips are engineered and manufactured in Germany.
©2024 Dynamic Biosensors GmbH | Dynamic Biosensors, Inc. All rights reserved.