
# User manual

## Molift EvoSling Hygiene

***molift***<sup>®</sup>  
by Etac

BM28899 - 2024-08-13 - Rev. 5.0

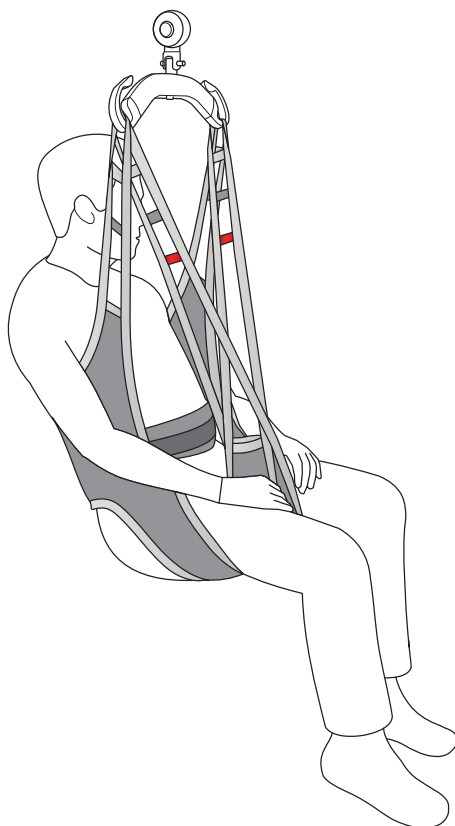[www.etac.com](http://www.etac.com)

|                                                                                     |                           |           |
|-------------------------------------------------------------------------------------|---------------------------|-----------|
|    | <b>Symbols</b>            | <b>4</b>  |
|    | <b>Device overview</b>    | <b>8</b>  |
|    | <b>General</b>            | <b>10</b> |
|                                                                                     | English                   | 10        |
|                                                                                     | Svenska                   | 14        |
|                                                                                     | Norsk                     | 18        |
|                                                                                     | Dansk                     | 22        |
|                                                                                     | Suomi                     | 26        |
|                                                                                     | Deutsch                   | 30        |
|                                                                                     | Nederlands                | 34        |
|                                                                                     | Français                  | 38        |
|                                                                                     | Italiano                  | 42        |
|                                                                                     | Español                   | 46        |
|                                                                                     | <b>Handling</b>           | <b>50</b> |
|    | Sitting application       | 51        |
|    | Sitting removal           | 53        |
|    | Warnings                  | 55        |
|    | Information               | 56        |
|   | <b>Combination list</b>   | <b>58</b> |
|  | <b>Sling measurements</b> | <b>58</b> |
|  | <b>Size guide</b>         | <b>58</b> |

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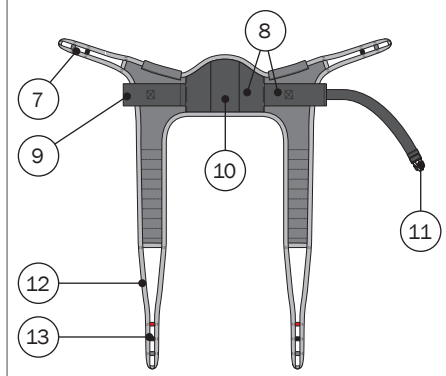
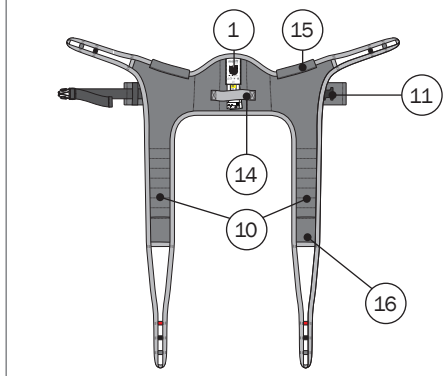
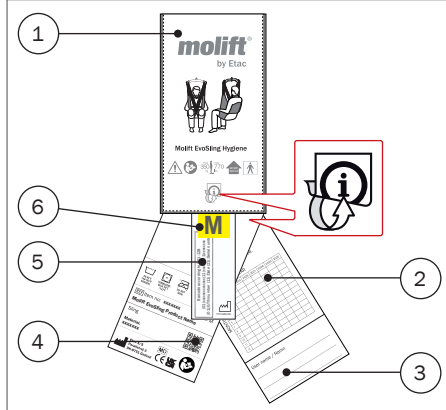
## General

Thank you for choosing an Etac device.

This manual must be read thoroughly before you start using the device as it contains important safety information and gives instructions on how to use the device correctly. Visit [www.etac.com](http://www.etac.com) for all relevant product documentation, for example the user manual, periodic inspection document and recycling instructions. You can select your language via the “International” and “Local websites” link.

You can also consult the QR code on the product label for further information and instructions.

### Copyright

This document is the property of Etac A/S. The content shall not be modified in whole or part without prior approval of Etac A/S.

### Disclaimer

Etac continuously makes improvements to our devices and information. We therefore reserve the right to modify our devices and instructions for use without prior notice. Consult [www.etac.com](http://www.etac.com) to find the latest version of the documentation.

## CE Declaration of Conformity

The device mentioned in this manual and the different models hereof comply with the Medical Device Regulation (EU) 2017/745.

The device is classified as a Class I medical device.

### Product standards

The device has been tested and complies with the requirements applicable for this type of device in the below mentioned standards:

- ISO 10535:2021 Assistive products – Hoists for the transfer of persons – Requirements and test methods.

### Notice

Any serious incident that has occurred in relation to the device should be reported to your distributor/representative and the national competent authority within your region in a timely manner. The distributor/representative will forward the information to the manufacturer.

### Pre-sales information

Please visit [www.etac.com](http://www.etac.com).

## Device description

Etac supplies a wide selection of slings for different types of transfers. The Molift EvoSling series is designed to provide high comfort for users by evenly distributing weight and pressure. Molift EvoSlings are available in different sizes, models and materials. See [www.etac.com](http://www.etac.com) for a complete overview of possible variants.

Molift EvoSling Hygiene can be used for hoisting to and from a sitting position. The sling is primarily designed for a 2-point sling bar, but it also works with a 4-point sling bar.

The wide split-leg support opening facilitates dressing and undressing.

### Benefits

The sling is padded to distribute the weight, which provides a comfortable support for the user through the entire transfer.

The loop connections give multiple options for attachment of the sling.

The sling has a waist belt with non-slip lining inside which helps keep it in place.

Using the sling provides a high level of safety and security for the user.

Application pockets on the leg supports facilitate application and guide handles help with lowering the user into an optimal sitting position.

### Materials

Fabric: Polyester

Padding: Polyethylene

## Introduction

The person referred to as “the user” is the person who is sitting in the product. “The carer” is the person who applies the product.

### Intended purpose

The sling is an assistive device intended for alleviation of or compensation for a functional impairment due to an injury or disability. The device is designed for an individual lacking the ability to stand up and transfer themselves with easy access to toileting to/from a bed, a wheelchair, a commode chair, a toilet or similar, due to reduced mobility or physical strength.

### Intended user

The device can be used for children, adults and elderly people with a need to be lifted.

### Intended environment

For indoor use in acute care, long term care and home care.

### Indications

For users who are in need of support to be lifted for a short time with easy access to toileting; e.g. people with limited ability to move and who require body support during transfers. Not an exhaustive list. The user group for the device is based on individual health and mobility function, and not on a specific diagnosis or age.

### Contraindications

Not suitable for users with hypotonia or reduced core stability. Users with reduced head stability should not use this type of sling.

### Precautions

Make sure the sling fits the user and the sling bar.

### Training requirements

Only carers who have read and understood the user manual are allowed to use the device.

## Safety information

The safety of both user and carer must never be jeopardized – if you suspect that the device has a malfunction you must put it in quarantine marked “out of order” and call your local distributor/representative.

Safety signs, safety measures and warnings symbols are given in the context relevant for the actual handling and/or operating situation – it is mandatory the given instruction is followed to avoid dangerous/hazardous situations.

Signs, symbols and instructions placed on this device must never be covered or removed and must remain present and clearly legible throughout the entire lifetime of the device. Immediately replace or repair illegible or damaged signs, symbols, and instructions.

Contact your distributor/representative for instructions.

### Patient transfer

Transfer of a user will always pose a certain risk. Therefore, a risk assessment must always be carried out before the transfer is performed.

Never leave a user unattended when hoisted.

Never attach the guide handles to the sling bar hooks.

Molift slings shall only be used to hoist persons. Never use the sling to hoist objects of any kind.

### Risk assessment

Do a risk assessment of the sling, the individual user and the intended transfer, and document your analysis. You as the carer are responsible for the safety of the user.

Make a decision on whether one or more carers are required.

Plan the transfer in advance to ensure that it is performed as safe and smooth as possible.

Remember to work ergonomically according to local recommendations.

Before using the sling, mark the date of first use on the periodic inspection label.

## Slings

Always check – before each use:

- that the user has the required ability for the sling.
- that the model, size and fabric are suitable for the user and the transfer situation.
- that the sling is suitable for the hoist and the sling bar.
- that the sling does not have visible damage or frays to seams, fabric, straps and loops.
- that the loop connections are correctly attached to the sling bar hooks before hoisting the user.
- that the wheels on the wheelchair or bed are locked during the transfer. If using a mobile hoist the wheels on the hoist must not be locked.

If the inspection reveals any safety risks, the sling must immediately be taken out of service and be disposed of.

Never use a faulty or damaged sling as it can break and cause personal injury. Destroy and discard damaged and old slings.

## Materials

All used materials which are to be in direct contact with skin are biologically assessed and evaluated according to relevant standards. However, users/carers who are allergic to any of the materials are recommend not to use the device.

## Accessories

User manuals supplied with accessories must be read and understood before use. This includes any limitations in use and combinations. Information about accessories can be found on [www.etac.com](http://www.etac.com).

## Combinations

This device can be combined with Molift sling bars and hoists. See the combination list on page 58.

Statement on sling and hoist compatibility for other manufactures of hoists can be found on [www.etac.com](http://www.etac.com).

When combining the device with other devices, the combination must not change the intended purpose of the devices or modify the devices in such a way that compliance with the applicable requirements may be affected. The person/company responsible for the combination of the devices must ensure that the applicable requirements are fulfilled.

If safe working load (SWL) differs between hoist, sling bar and sling, the lowest load should always be used.

## Troubleshooting

| Problem                                                                             | Action                                                                                                                                                                                                                                                                                                             |
|-------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| There is more than 10 cm between the left and right side padding of the waist belt. | Try a larger size.                                                                                                                                                                                                                                                                                                 |
| The waist belt is too loose after tightening it.                                    | Try a smaller size.                                                                                                                                                                                                                                                                                                |
| The user's head comes too close to the sling bar.                                   | Try another loop connection to extend the distance to the sling bar. Otherwise, try a larger size.                                                                                                                                                                                                                 |
| The straps are not loaded at the same time.                                         | Make sure the straps are attached symmetrically.<br>Make sure the centre of the sling is at the midline of the user.                                                                                                                                                                                               |
| The user's buttocks slide down the opening of the sling.                            | Make sure the sling is placed in waist height and the leg supports are not too close to the knee pits. Otherwise, try a smaller size.<br>Make sure that the user has the physical ability for the sling. Otherwise, try a sling with a more narrow opening.<br>Make sure the user is not wearing slippery clothes. |

## Service, cleaning and maintenance

Maintenance and periodic inspections of the device are wholly the responsibility of the owner of the device. Failure to maintain a device in accordance with the instructions may compromise the safety of the user and/or carer and may affect the expected lifetime of the device.

### Expected lifetime

The device has an expected lifetime of 1 to 5 years under normal use. The lifetime of the device may vary depending on usage frequency, loads and washing/drying frequency and conditions.

Washing and drying at high temperatures wear out the material faster.

### Cleaning

60-85°  
140-185F

Max. 60°  
/140F

### Periodic inspection

Periodic inspection must be carried out at least every 6 months. More regular inspection may be required if the sling is used or washed heavily. See the periodic inspection label on the sling and the periodic inspection document at [www.etac.com](http://www.etac.com).

## Transportation, storage and disposal

### Transportation and storage

When the device is not in use, it can be stored under the following conditions:

### Disposal

The device contains no hazardous substances that require special disposal. Dispose of the device according to the national provisions that apply.

### Recycling

See the recycling instructions on [www.etac.com](http://www.etac.com), and follow the national provisions that apply.

## Technical data

### Size guide

Measurements in the table on page 58 should be used as a reference. The correct sling size depends on the user's weight, function and body shape. Consult the size guide, and measure the user's relevant measurements in a sitting position. Try the sling with the most correct measurements to make sure it fits.

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**EN:** Information

- Slide the sling down behind the user's back until it reaches the user's waist.
- Make sure the straps are fastened symmetrically on left and right side.
- Before performing the transfer - make sure that:
  - Sling is level with the coccyx
  - Leg supports are placed under the thighs without creases
  - Lower straps are crossed and symmetrically applied
- Perform a safety stop when the straps are stretched tense. Check that all straps are correctly and symmetrically attached.
- Lower the user to sitting position using either method A or B.

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**Etac Ltd.**