

Get connected

The Beltone HearMax™
and Beltone HearPlus™ apps



Made for
iPhone | iPad | iPod

Available on
android

Introduction

Thank you for choosing our products. We have created this app to make it easy for you to use and control your hearing aids from your smartphone or mobile device. For information on which smartphones and mobile devices are compatible with the app, visit: beltone.com/compatibility.

Before you can start using the app, you must pair your smartphone or mobile device with your hearing aids. This is easily done – just follow the steps below.


How to pair your hearing aids with your iOS device* (Instructions reflect pairing flow for iOS 13 or later)

1. Turn your hearing aids off.
2. Confirm **Bluetooth** is on by going to the **Settings** menu on your iOS device. The slider turns green when you activate Bluetooth.


*For pairing Beltone Rely 2 with iOS, follow the steps for pairing with Android




3

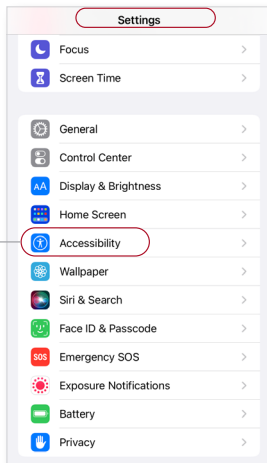

4

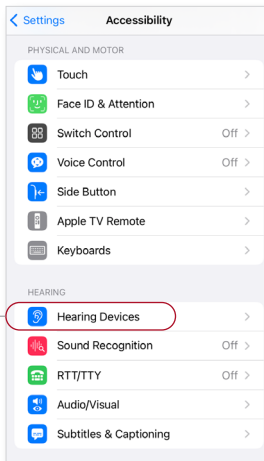

3. Go back to the main **Settings** menu and scroll down to tap on **Accessibility**.
4. Scroll down to **Hearing Devices**.
5. Turn your hearing aids on. Now, your iOS device will be able to detect your hearing aids.

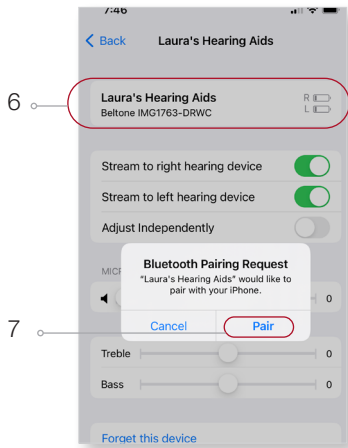

6. Once connected, select your hearing aids under Hearing Devices.
7. Tap on **Pair**.
8. Go to the **App Store**, and in the **Search** field on top enter the name of the app.
9. Download the app, open it and tap **Get Started**. Follow the instructions on the screen to connect to the app.
10. You are now ready to go.

If you need to re-pair your iOS device

1. Start by removing the old pairing by going to the **Settings** menu.
2. Tap **Accessibility**.
3. Select **Hearing Devices** and tap **Forget this device**.
4. Turn your iOS device off and on again and re-pair the hearing aids to the device. If you need more information, contact your hearing care professional.

Hands-free calls for Beltone Imagine users with compatible Apple devices*


If you have Beltone Imagine hearing aids with compatible Apple devices, your hearing aids will become a hands-free headset as default once paired with the Apple devices.

How to take/end calls:

- Short press your hearing aid button (either side) or use your iPhone or iPad.

How to decline/ignore calls:

- Long press your hearing aid button (either side) or use your iPhone or iPad.

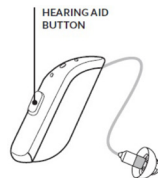

When you receive another call while you are on the phone, simply follow the instructions above.




Note: if the incoming call is declined, it automatically returns to the first/previous call.

How to change the volume of the calls:

- Use the up and down volume buttons on your iPhone or iPad.

*Hands-free calls are compatible with iPhone 11 or later, iPad Pro 12.9-inch (5th generation), iPad Pro 11-inch (3rd generation), iPad Air (4th generation), and iPad mini (6th generation) with software updates iOS 15.3 and iPadOS 15.3 or later.

How to pair your hearing aids with your Android™ device*

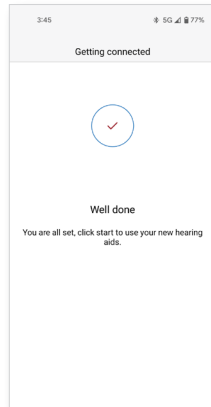
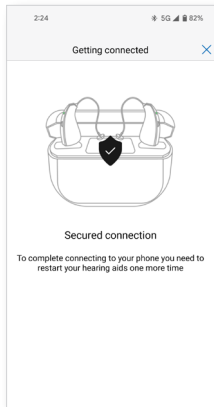
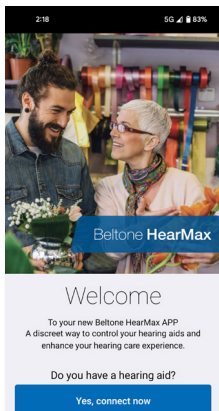


1. Download the Beltone HearMax app, open it and tap **Yes, connect now**.
2. If prompted, provide access to your mobile device location by selecting **Allow**.
3. The app will guide you through the process. Reboot your hearing aids as required by turning them off and on again, following the steps outlined in the app.
4. Pairing is complete.

NOTE: If you have previously paired hearing aids, you need to unpair them in your phone's Bluetooth menu before you pair them again.

*These instructions are also applicable for pairing Beltone Rely 2 with a compatible iOS device

How to disable hands-free calls:

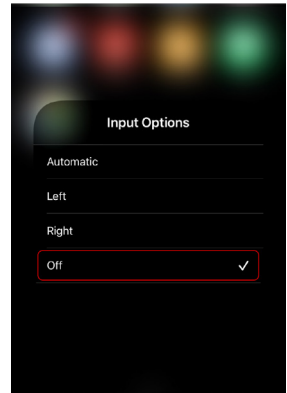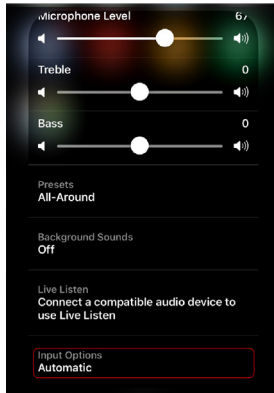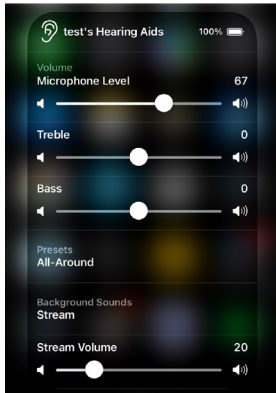

1. Triple click the iPhone or iPad push button to access the iOS hearing aid settings.

2. Scroll down and select "Input Options".

3. Select "Off".

Learn more about hands-free calls at:
<https://www.beltone.com/en-us/support/apps-support/hearmax-app-support>

and





Manufacturer according to FDA:

Beltone

8001 E Bloomington Freeway
Bloomington, MN 55420
USA
1-800-BELTONE
beltone.com

Manufacturer according to Health Canada:

Beltone Canada

2 East Beaver Creek Road, Building 3
Richmond Hill, ON L4B 2N3
Canada
1-800-387-3744
beltone.ca

© 2022 GN Hearing A/S. All rights reserved. Beltone is a trademark of GN Hearing A/S. Apple, the Apple logo, iPhone, iPad, iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries. Android, Google Play and the Google Play logo are trademarks of Google LLC. The Bluetooth word mark and logos are registered trademarks owned by the Bluetooth SIG, Inc.