BLOOD PRESSURE UNIT USER'S MANUAL

REF MDS4001

REF MDS4001LA

REF MDS4001PLUS

REF MDS4001LAT

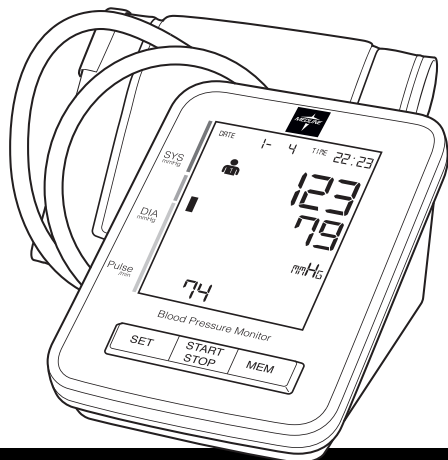
TABLE OF CONTENTS

INTRODUCTION AND INTENDED USE	3
SYMBOL IDENTIFICATION	4
IMPORTANT SAFETY INFORMATION	4
MEASUREMENT ACCURACY	5
HELPFUL INFORMATION ABOUT BLOOD PRESSURE	6
CORRECT METHOD OF MEASUREMENT	10
PARTS IDENTIFICATION	11
PREPERATION BEFORE TAKING A MEASUREMENT	14
SETTING TALKING FEATURE	17
TAKING A MEASUREMENT	18
MEMORY FUNCTION	20
ERROR INDICATORS AND TROUBLESHOOTING	21
TROUBLESHOOTING	22
MEDLINE BLOOD PRESSURE CUFF TABLE	23
CARE AND MAINTENANCE	24
SENSOR CALIBRATION MODE	24
SPECIFICATIONS	25
BLOOD PRESSURE LOG	26
WARRANTY	27

INTRODUCTION AND INTENDED USE

This instruction manual provides important information about the MDS4001 series automatic digital blood pressure monitors. These devices are intended to measure the systolic and diastolic blood pressures and pulse rate of an adult individual by using a non-invasive technique, in which an inflatable cuff is wrapped around the upper arm.

Please read all of these instructions carefully, before using the monitor.

IMPORTANT NOTES

- Only a physician is qualified to interpret changes in your blood pressure. This device is not intended to replace regular medical examinations.
- It is recommended to have your physician review your procedure for using this monitor. Never make adjustments to your medication unless recommended by your physician.
- This Blood Pressure Monitor is intended for use by adults. Children should not use this monitor unless it is under the supervision of an adult.
- Only use the cuff(s) included with this device or offered as approved accessories. Unapproved cuffs can affect measurement accuracy.

SYMBOL IDENTIFICATION

The following symbols may appear in this manual, on the blood pressure monitor, or on its accessories.

	Warning
	Type B: Identifies applied part complying with IEC 60601-1
	Class II
	Refer to instruction manual
	Keep Dry: Protect from moisture.
	Disposal: Do not dispose this product with household waste.

IMPORTANT SAFETY INFORMATION

Warning

- Consult with your physician before using this device on an arm where intravascular access or therapy, or an arteriovenous shunt, is present because of temporary interference to blood flow which could result in injury.
- Consult your physician before using this device if you have had a mastectomy.
- Too frequent measurements can cause injury, due to blood flow interference.
- Do not apply the cuff over a wound, as it can cause further injury.
- Do not kink the cuff hose, as this may cause the cuff pressure to continuously increase, which can interfere with blood flow and result in injury.
- Do not inflate the cuff on the same limb to which other monitoring medical electrical (ME) equipment is applied. This can cause temporary loss of function of the monitoring ME equipment.
- Make sure that operation of this device does not result in prolonged impairment of blood circulation.

MEASUREMENT ACCURACY

The quality of this device has been verified and conforms to the provisions of the EC-Council directive 93/42/EEC as well as the EMC directive 89/336/EEC:

EN 1061-1

Non-invasive blood pressure measuring equipment. General requirements

EN 1060-3

Non-invasive blood pressure measuring equipment. Supplementary requirements for electro-mechanical blood pressure measuring system

EN 60601-1

Safety requirements for medical electrical equipment

EN 60601-1-2

Electromagnetic compatibility and safety for medical electrical equipment

EN 14971

Risk analysis for medical devices in accord with requirement in section 21 of CMDR'S

HELPFUL INFORMATION ABOUT BLOOD PRESSURE

What is Blood Pressure?


Blood pressure is the force that blood exerts on the arteries. This force is constantly changing as the heart beats. When the heart contracts, blood pressure reaches its highest value. This is called systolic blood pressure. When the heart relaxes between beats, the value of blood pressure is lower. This is called the diastolic blood pressure. The unit of measure for blood pressure is millimeter of mercury, abbreviated mmHg. For example, an individual's blood pressure may be measured as 120 mmHg (systolic) and 80 mmHg (diastolic). This would be spoken as "120 over 80" and written as "120/80".

What Affects Blood Pressure?

Blood pressure can be affected by the following factors: eating, drinking alcohol or caffeine, smoking, stress, exercise, physiologic condition, measurement site, body position, temperature, humidity, altitude, and many other factors.

Remember that blood pressure can also vary throughout the day. It is important to take measurements at the same time each day, and under the same conditions.

Typical Daily Blood Pressure Fluctuations

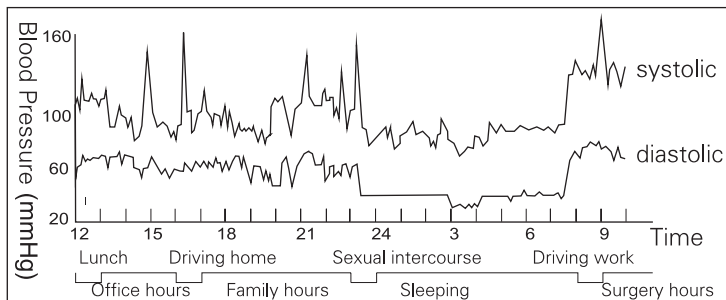Example: 35-yr old male

What is High Blood Pressure?

Hypertension, or high blood pressure, is a condition where an individual's blood pressure remains high over a long period of time. Consult with your physician with any questions or concerns you may have about hypertension.

Blood Pressure Classification

The World Health Organization (WHO) and the International Society of Hypertension (ISH)


Blood Pressure Classification	Systolic (mmHg)	Diastolic (mmHg)
Optimal	Less than 120	Less than 80
Normal	120-129	80-89
High - Normal	130-139	85-89
Hypertension - Stage 1	140-159	90-99
Hypertension - Stage 2	160-179	100-109
Hypertension - Stage 3	180 or Higher	110 or Higher

Monitor Pressure Gauge and Indicator

The pressure gauge allows quick detection of your blood pressure classification. The ■ indicator displays a segment, based on your measurement, corresponding to the WHO and ISH classification.

Note: The pressure gauge indicator is approximate. Please consult your physician to interpret your blood pressure measurements.

Pressure Gauge

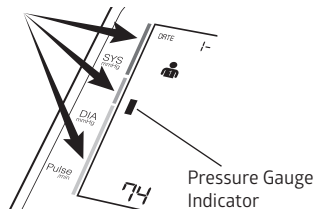Why is it Important to Measure Your Blood Pressure at Home?

Having your blood pressure taken at the doctor's office may cause you to become nervous, thus artificially raising your reading. Having the ability to take your blood pressure at home makes it easy to record a log of your daily readings. This will help you gain a greater understanding of your blood pressure reading and the factors that affect it. Be sure that you share your information with your physician.

Hints for Accurate Measurement

1. Relax and try to remain still for 5 to 10 minutes before a measurement.
2. Remove any clothing on the upper arm so that the cuff can be placed directly on the skin. Constricting garments may cause an inaccurate reading.
3. Perform measurements on the same arm (preferably the LEFT arm).
4. Refrain from eating, smoking and drinking (especially alcoholic beverages) before a measurement. These activities can affect your blood pressure.
5. Remember that blood pressure varies continuously throughout the day. Try to take your blood pressure at the same time each day.
6. Do not be concerned with the results of one measurement. Many measurements, recorded over an extended period of time, will provide a better indication of your blood pressure.
7. Relax for 5 to 10 minutes before remeasuring your blood pressure.
8. Because many factors can affect your blood pressure, daily fluctuations of 25 to 50 mmHg are common.

Please note: Some individuals with hypertension, diabetes, kidney disorders, arteriosclerosis or poor circulation may see a significant difference in blood pressure readings taken from the wrist as compared to readings from the upper arm. It is recommended that you consult with your physician concerning the use of this monitor.

CORRECT METHOD OF MEASUREMENT

To obtain the most accurate blood pressure measurement, please follow these important directions:


1. Sit comfortably in a chair with your back supported. Place both feet flat on the floor and keep your legs uncrossed.
2. Rest your left arm on a table and position your palm upwards.

IMPORTANT: The arm cuff must be at the same level as your heart, or accurate measurement will not be possible. (Your heart is located slightly below your left armpit.)


3. Remain still, avoid moving or talking during the measurement.

Feet Flat on Floor

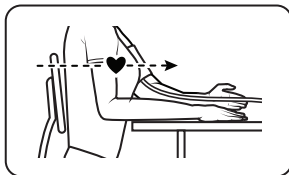Position Arm at Heart Level

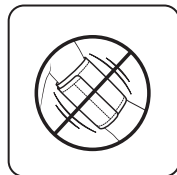
Avoid Moving or Talking

PARTS IDENTIFICATION


Contents

Monitor, arm cuff, four (4) "AA" alkaline batteries, instruction manual, quick start guide


Monitor

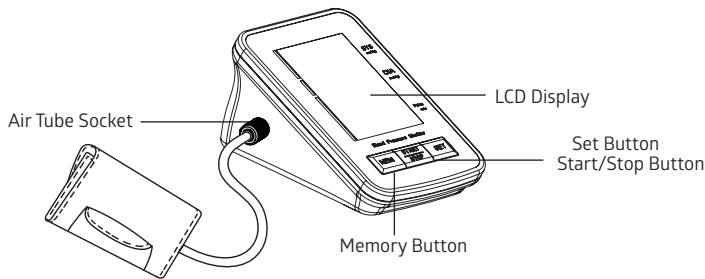Accessories

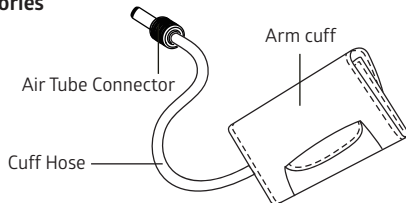AC Adapter

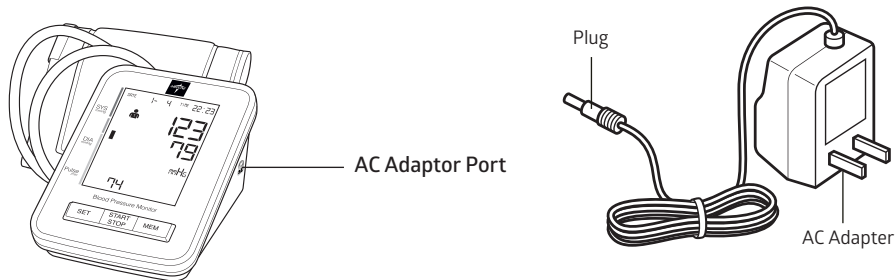The MDS4001 series has a 120V AC adaptor port to allow you to supply power from an outlet in your home. We recommend you only use the AC adapter (item number: M-IADE01) intended for the MDS4001 series items. Please contact Medline at 1-855-441-8737 if you would like to purchase an AC adapter.


Connecting the AC Adapter to the Monitor:

1. Insert the AC adapter into an outlet.
2. Connect the AC adapter plug into the AC adapter port on the **right** side of the device.

Disconnecting the AC Adapter from the Monitor:

1. Turn the device off by pressing the **START/STOP** button.
2. Disconnect the AC adapter plug from the device.

Display Marks

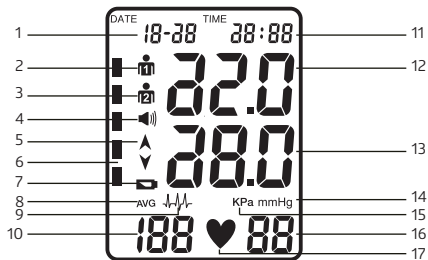- | | |
|---------------------------------|---|
| 1.Date | 10.Pulse Rate |
| 2.User 1 | 11.Time |
| 3.User 2 | 12.Systolic Blood Pressure |
| 4.Volume | 13.Diastolic Blood Pressure |
| 5.Inflating/Deflating | 14.Millimeters of mercury (mmHg) Measurement Unit |
| 6.Pressure Gauge Indicator | 15.Kilopascals (kPa) Measurement Unit |
| 7.Battery Display | 16.Memory Frequency |
| 8.Average (AVG) | 17.Measuring |
| 9.Arrhythmia Detection Function | |

Arrhythmia Detection Function:

The arrhythmia detection symbol  will appear on the display, if pulse irregularities were detected during the measurement. This result may deviate from your normal blood pressure, in which case repeat the measurement. It is important to remain still and avoid talking during the measurement.

Note: We recommend you consult your physician if you see the arrhythmia detection symbol frequently.

PREPERATION BEFORE TAKING A MEASUREMENT



Battery Installation/Replacement

1. Slide the battery cover off in the direction of the arrow.
2. Install or replace the four alkaline "AA" batteries, noting the proper orientation of positive (+) and negative (-) terminals in the battery compartment.
3. Replace the battery cover.

Note:

- Replace all batteries whenever  appears on the display or when the display does not light up.
- Remove the batteries if the device will not be used for an extended period of time.
- Removing the batteries will erase all readings stored in the memory.
- Batteries are hazardous waste. Do not dispose of them with household garbage.

Battery compartment

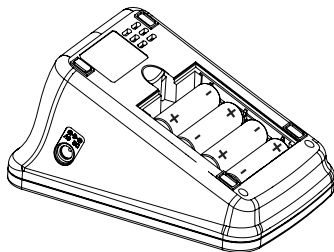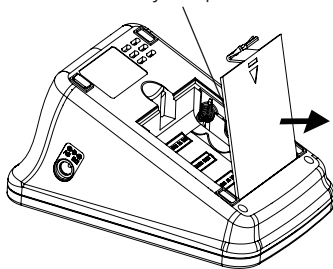Setting Date and Time

1. When the display is off, hold the **SET** button for 5 seconds, or until the year in the upper right corner begins to flash.

2. Press the **MEM** (Memory) button to advance the year, until the proper value is displayed.
3. Once the year is selected, press the **SET** button to switch between month, day, hour, and minute.


4. Press the **MEM** (Memory) button to advance the digit, one value at a time, for the selected date or time setting.

Selecting mmHg or kPa Measurement Units

1. When the display is off, hold **START/STOP** button for 10 seconds, or until the digits on the screen start to flash.
2. Press the **MEM** (Memory) button to select mmHg or kPa.

How to Apply the Arm Cuff

1. Place the cuff directly on your bare upper arm.
2. Adjust the cuff so that the bottom edge lies about 1 inch above the elbow on the inside of the arm.
3. Align the Artery Mark (Φ) over the brachial artery, which runs down the inner arm. Position the cuff hose straight down your arm, in line with the middle finger.
4. Pull the end of the cuff and secure it snugly with the hook and loop closures. You should be able to fit two fingers between your arm and the cuff.
5. Confirm the ▲ indicator at the top, outer edge of the cuff lies within the fit range (↔).
6. Plug the air tube connector into the socket on the left side of the monitor.

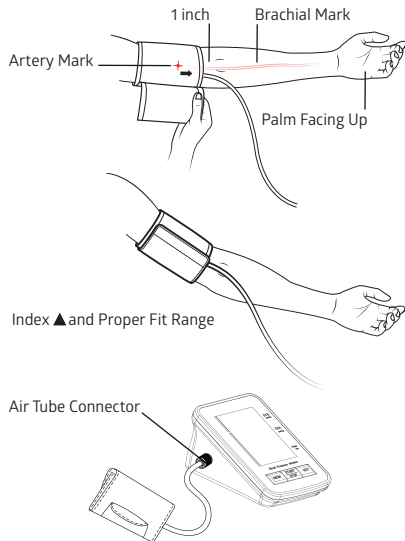IMPORTANT: Make sure the cuff hose remains straight along the arm and unkninked during measurement. Avoid compression and restriction of tubing.

SETTING TALKING FEATURE

The MDS4001LAT Talking Model can be set to announce your measurement results, according to the following method:

NOTE: The MDS4001LAT has a gray faceplate. If your device does not have a gray faceplate, the following instructions are not relevant for your device.

1. When the display is off, hold the **SET** button for about 5 seconds, or until the year in the upper right corner begins to flash.
2. Press the **SET** button again to move through the date and time settings, until **SP** in the lower left corner begins to flash.
3. Press the **MEM** (Memory) button to switch between languages: **SP 01** is English, **SP 02** is Spanish, and **SP 0F** is off.
4. Press the **SET** button again to confirm the language selected and return the device to the initial display off screen.

TAKING A MEASUREMENT

Please read the previous section on proper cuff placement and body position.

1. Press the **SET** button and select [] or [], in the upper left hand corner of the LCD screen, to indicate which user is taking the measurement.
2. Press the **START/STOP** button. The arm cuff will begin to inflate automatically. The cuff will be quite snug for a short while, this is normal.

Note: TO END a measurement for any reason, press the **START/STOP** button again.

3. After the inflation, the cuff will slowly deflate to take the measurement. The  will flash, indicating that the measurement is in process.


Note: If the device determined that the initial inflation pressure was insufficient, it will reinflate at a higher pressure.

4. The device will completely deflate when the measurement is finished. The systolic and diastolic pressure, pulse rate, and pressure gauge indicator will be displayed

Start


Inflating

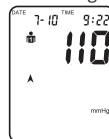Deflating

Complete

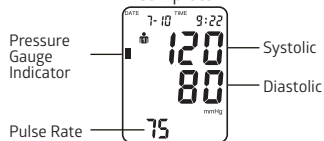Note:

- The device will automatically shut off after 1 minute, or it can be powered off by pressing the **START/STOP** button.
- If during a measurement the power drops to an insufficient level, the device will end the measurement and  (low battery) will be displayed. Please replace the batteries and repeat the measurement.
- If the device cannot detect your pulse, it will end the measurement. Please wait a few minutes, make sure that the cuff is positioned properly, and try again.

MEMORY FUNCTION

Memory Input



When a valid measurement is made, the result (systolic, diastolic, pulse rate, date and time) will be stored in the memory automatically as soon as the device is turned off or when another measurement is initiated. The device can store up to 90 readings for each of 2 users, and also compute the average of the latest 3 readings. When more than 90 readings are stored, the oldest ones will be replaced.

Memory Recall

1. When the display is off, press the **MEM** button to show the average (AVG) of all recorded results.
2. Press the **MEM** button again and 01 will appear in the lower right hand corner, indicating the latest measurement in the memory.
3. Continue pressing the **MEM** button to display additional past readings, sequentially.

Memory Erase

1. When the display is off, press the **MEM** button.
2. Press the **MEM** button again and hold it down for more than 3 seconds to erase all stored results.

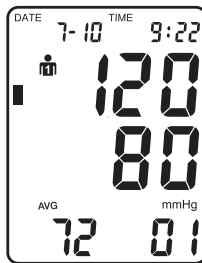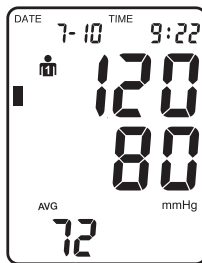Memory
Frequency

ERROR INDICATORS AND TROUBLESHOOTING

ERROR INDICATORS	CAUSE	CORRECTION
E_{r_u}	Incomplete inflation of arm cuff	• Check if the cuff is firmly connected to the device • Make sure the cuff hose remains straight and unkinked during measurement • Reposition the cuff correctly (see pg. 16) and try the measurement again • If the Er message persists, call product support
E_{r_u}	Cuff over-inflated and released automatically	• Take a break, remove all clothing from upper arm, reposition the cuff correctly (see pg. 16) and try the measurement again • Make sure the cuff hose remains straight and unkinked during measurement
E_{r_1}	Could not obtain the pulse rate	• Reposition the cuff correctly (see pg. 16) and try the measurement again • Please refrain from talking or moving during measurement
E_{r_2}	Strong electromagnetic interference (mobile phone or computer)	• Make sure the cuff hose remains straight and unkinked during measurement • Remove sources of interference (mobile phone or computer) • Take a break, reposition the cuff correctly (see pg. 16) and try the measurement again
E_{r_3}	Irregular pulse rate or blood pressure	• Take a break, reposition the cuff and your body position correctly (see pg. 16) and try the measurement again • Please refrain from talking or moving during measurement

TROUBLESHOOTING

PROBLEM	CAUSE	CORRECTION
Nothing appears in the display, after START/STOP button is pressed	Batteries are drained	Replace the batteries with new ones
	Batteries are inserted incorrectly	Reinstall the batteries with positive and negative terminals in correct position (see pg. 14)
Cuff does not inflate	Tube connector is not fully inserted into device	Firmly insert connector into air socket on left side of device
	Tube connector is broken	Call product support
	Cuff hose is compressed or restricted	Make sure hose remains straight and unkinked during measurement
	Cuff is not fastened properly	Reposition the cuff correctly (see pg. 16)
	Battery voltage is too low	Replace all batteries with new ones
Blood pressure readings are too high or too low	Cuff is not fastened properly	Reposition the cuff correctly (see pg. 16)
	Body and arm are not positioned properly	Reposition your body and arm correctly (see pg. 10)
	Movement or talking during measurement	Remain seated, and please refrain from moving or talking during measurement
	Repeated measurement on same arm within short time period	Wait 5 to 10 minutes between measurements
Blood pressure readings are variable	Possible factors: time of day, emotional state, daily activities or exercise, smoking, drinking alcohol or caffeine, eating, certain medications (see pg. 6)	Take measurement at same time each day, and under the same conditions

MEDLINE BLOOD PRESSURE CUFF TABLE

If you need a different-sized cuff, please consult the table below for information about the different sizes available. Replacement cuffs fit all Medline-branded automatic digital blood pressure monitors. Contact Medline Industries, Inc. at 1-800-MEDLINE or visit www.medline.com if you need to purchase a new cuff. Provide the appropriate item number from the table below.

MEDLINE ITEM NUMBER	SIZE	MEASUREMENT RANGE (cm)
MDS9970	Small Adult	17-22 cm
MDS9971	Adult	22-30 cm
MDS9972	Large Adult	30-42 cm
MDS9973	Extra Large Adult	42-48 cm
MDS9974	Universal	22-42 cm

BLOOD PRESSURE MONITOR ITEM NUMBER	INCLUDED ARM CUFF ITEM NUMBER	ARM CUFF SIZE	MEASUREMENT RANGE (cm)
MDS4001	MDS9971	Adult	22-30 cm
MDS4001LA	MDS9972	Large Adult	30-42 cm
MDS4001LAT	MDS9972	Large Adult	30-42 cm
MDS4001PLUS	MDS9971 & MDS9972	Adult & Large Adult	22-30 cm & 30-42 cm

CARE AND MAINTENANCE

- Keep the device away from high temperatures, high humidity and direct sunlight.
- Keep sharp objects away from the cuff and tubing.
- Do not press the **START/STOP** button without the cuff around your arm to avoid excess pressure build up that may damage the cuff.
- Use the device at a sufficient distance from devices with strong electrical fields, such as television, microwave ovens, X-ray equipment, etc.
- Do not subject the device to strong impact or drop it on the floor.
- Do not disassemble or modify the device or the cuff.
- Remove the batteries when the device is to be stored for extended periods of time.

Use only a soft dry cloth to clean the unit. Do not use solvents or other petroleum-based cleaners.

SENSOR CALIBRATION MODE

The blood pressure device should only be calibrated by the manufacturer.

SPECIFICATIONS

Manufacturer Model:	BA-803
Medline Item Numbers:	MDS4001, MDS4001LA, MDS4001PLUS, MDS4001LAT
Memory Function:	Storage and recall 2x90 measurements with average
Display:	LCD
Measuring Range:	Pressure: 0-280 mmHg Pulse: 40-199/minute
Sensor Accuracy:	Pressure: ± 3 mmHg Pulse: $\pm 5\%$
Inflation System:	Electro-pneumatic pump
Arm Size Ranges:	17-42 cm (6.7-16.5 inches)
Batteries:	1.5V alkaline (LR6/AA)X4
Automatic Power-off:	Approx. 1 minute after measurement
Reference Method for Clinical Trials:	Auscultatory measurement
Weight:	Approximately 420g(without batteries)
Operating Condition:	41°F ~ 104°F (5°C ~ 40°C) 10%~90%RH
Storage and Transport Condition:	-4°F ~ 149°F (-20°C ~ 65°C) 10%~95%RH
Standard Cited:	EN 1060-1 EN 60601-1-2 EN 1060-3 EN 14971 EN 60601-1

BLOOD PRESSURE LOG

Date	1/1	1/1	1/1							
Time	7:00	13:30	20:00							
mmHg	240									
	220									
	200									
	180									
	160									
	140									
	128	134								
	120		123							
	100									
	80	90								
	84		76							
	60									
Pulse	70	73	69							
Body Condition										

Warranty

Medline Industries, Inc. warrants this digital blood pressure monitor against any defects in material and workmanship for a period of one (1) year after the date of purchase. The associated parts, specifically the cuff, are warranted for a period of six (6) months after the date of purchase or for ten thousand (10,000) uses, whichever comes first. These uses include, but are not limited to, inflation of the bladder and attachment of the cuff hook and loop closure. Excluding cases of freight damage, tampering, clear abuse, misuse, or accidents, Medline Industries, Inc. will, at its discretion, repair or replace this sphygmomanometer and/or its parts during the warranty period without charge. No representative or person is authorized to assume for us any liability in connection with the sale of the products, of Medline Industries, Inc. This warranty gives you specific legal rights, and you may also have other rights that may vary from state to state. Prior to shipping any product, please contact either Medline Industries, Inc. at 1-800-MEDLINE or an authorized Medline sales representative with any warranty concerns.

www.medline.com ©2021 Medline Industries, Inc.

Manufactured for Medline Industries, Inc., Three Lakes Drive, Northfield, IL 60093 USA.

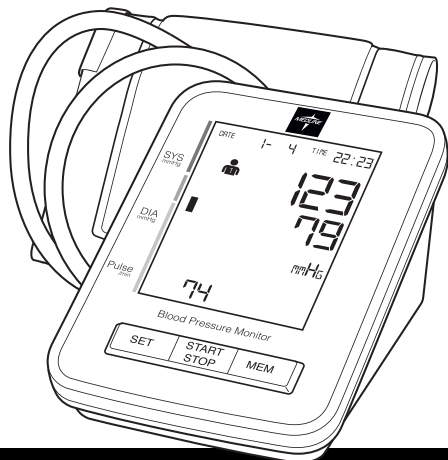
Product of Taiwan, Finished in China. **1-800-MEDLINE** V2 RD21DGA

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Fabricado para: Medline Industries, Inc., Three Lakes Drive,