



Contents [[hide](#)]

- [1 OMNIFIT OCW-H30 Enhancing Brain Training Solution](#)
- [2 Included in the Package](#)
- [3 Description of Parts](#)
- [4 Charging](#)
- [5 Charging Status LED Indicator](#)
- [6 Wearing](#)
- [7 Paring](#)
- [8 LED Indicator](#)
- [9 Specification](#)
- [10 Precautions](#)
- [11 Limited Warranty](#)
- [12 Battery](#)
- [13 FCC Information](#)
- [14 Frequently Asked Questions](#)
- [15 Documents / Resources](#)
 - [15.1 References](#)

OMNIFIT


OMNIFIT OCW-H30 Enhancing Brain Training Solution



Included in the Package



Headset (OCW-H30)

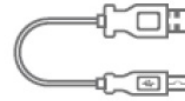

USB C type Cable





User Manual




Portable Pouch

Description of Parts

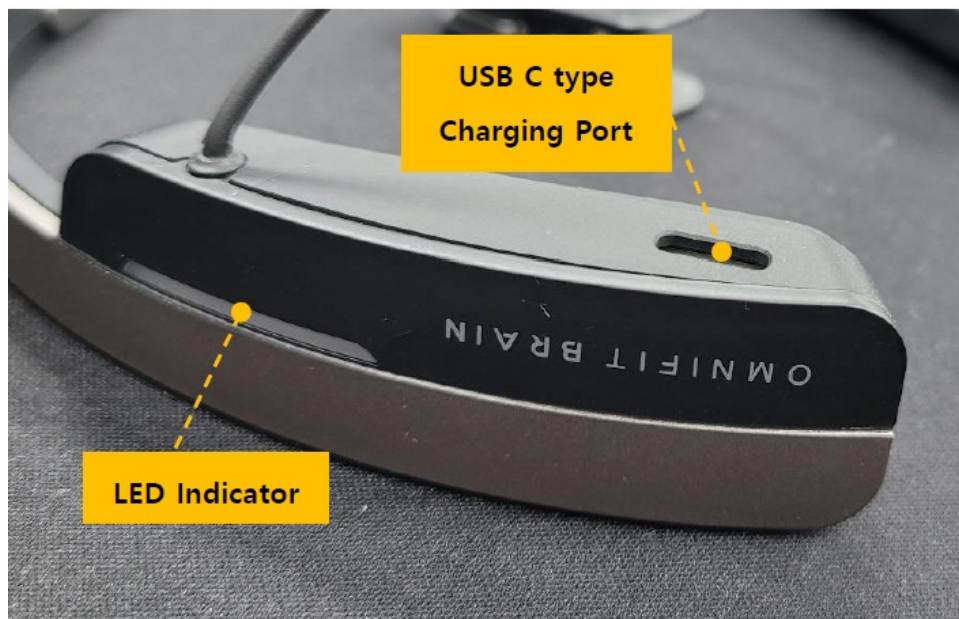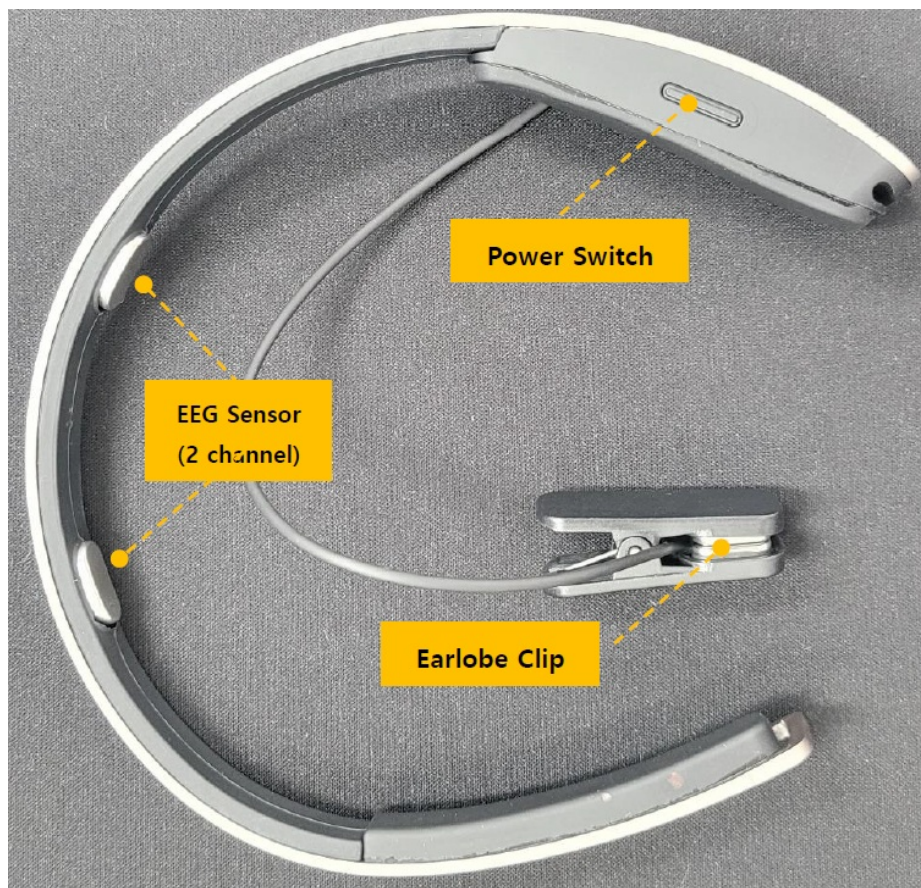

Charging

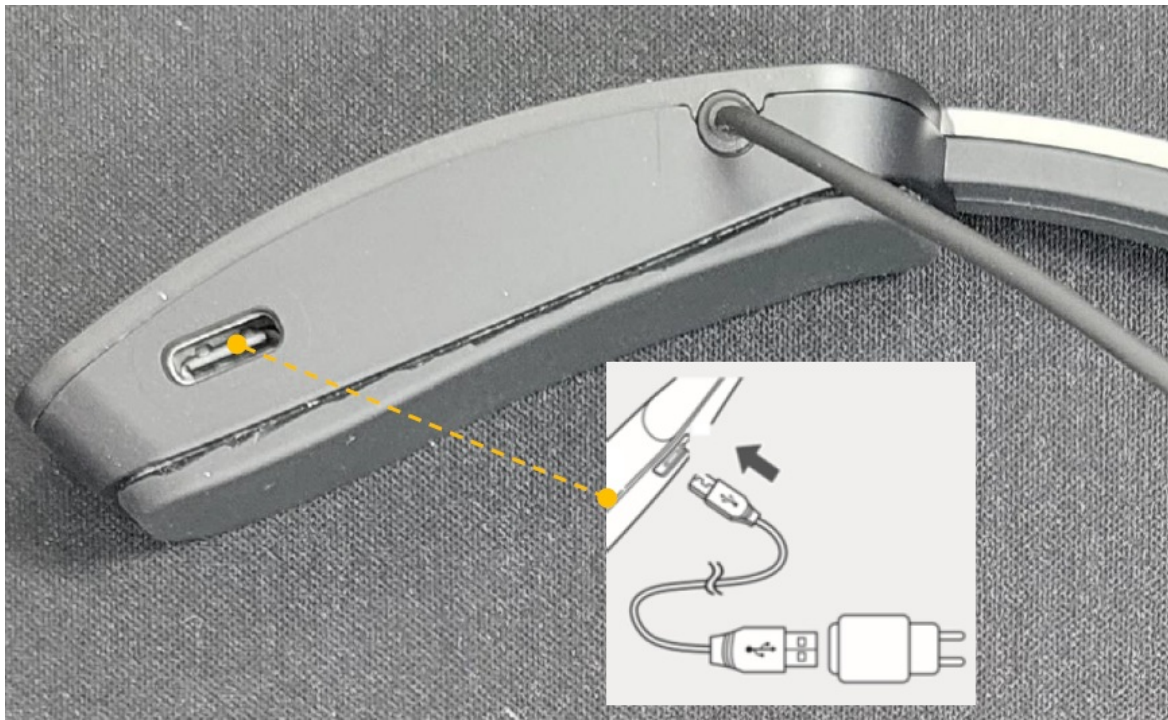

1. For charging, do connect USB cable to right under USB charging port
2. When the charger is connected, the Bluetooth pairing will be disconnected.
3. While charging, the LED indicator will show a solid red light.

Charging Status LED Indicator



Charging(RED)



Full Charged(GREEN)

- OCW-H25 Battery status will appear on OMNIFIT Application
- To prevent power outage while transmitting data to an application, please fully charge the battery before use.
- Please charge the battery below 40°C.

Wearing




1. Wear OCW-H30 where LED indicator must be located on right-side.
2. Contact firmly EEG sensors onto forehead.
3. Put in earlobe clip firmly. earlobe clip have a function of grounding.

Paring

Headset Pairing





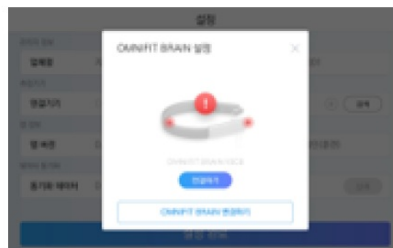


Press a power switch for 3sec for power ON Brain Training App > Select the [Settings] icon on the right side of the screen > Select Device [Search] > (If it is automatically connected, use it as it is) > Select [Connect] > Connection completed.

(Caution: Do not connect directly from the tablet's Bluetooth settings, but must connect through the App.)




Check the 4-digit device number on the side of the headset

LED Indicator

Status	LED Indicator			Action
	Solid Light	Bl i nk	Di mmi ng	
Charging (Red)				Charger connected
Fully Charged (Green)				Charger connected
Power on (White)				Power button pressed for 3sec
Low Battery (Red)				When battery left under 5%
BT disconnected (Blue)				When Power is on and BT not connected
BT connected (Blue)				When BT connection succeeded, the blue light stays on for 3 sec and then turn off.
Measuring (Yellow)				When measuring

Specification

Category	Item	Specification
Main Part	Main Chipset	EFM32BG1
	EEG 	2 Channel
	Memory	N/A
	Bluetooth	BLE 4.2
	Connection	BLE
	LED	RGB LED
	USB	USB C type Connector
	Power on / off	Power Switch
Earlobe Clip	Earlobe	Black housing(PC)Metal ground
Battery	Li-polymer battery	3.7V 135mAh (37x15x5t mm, 7g)
	Play time	Up to 7 Hrs
	Standby time	Up to 220Hrs

	Charging time	< 1.5Hrs
Dimension/Weight	Dimension / Weight	142 X 142 X 22mm (band 10mm) / 43g

Storage condition

- After use, wipe the metal electrode clean and store it in a place that is not inclined, except for a place with dust, liquid, temperature, moisture, and humidity that exceed the usage range.
- Storage temperature: -10 ~ 50 °C
- Storage Humidity: 15 ~ 90 %RH

Precautions

- Do not drop the device.
- Do not modify, repair, or disassemble.
- Do not expose directly to water, alcohol, benzene, etc for cleaning.
- Do not expose directly to flammables.
- Do not place or keep the device near flammables.
- Keep the device away from excessive humidity and dust.
- Only charge the headset with certified charger, otherwise it may be damaged.
- Do not expose the battery to extreme temperatures, never above +60°C (+140°F) or under -20°C (-4°F)
- The product should be connected to power supply not exceeding LPS or PS2 complying with EN IEC 62368-1:2020 + A11:2020

Limited Warranty

- The warranty is valid only if the original receipt issued to the original purchaser by the dealer, specifying the date of purchase and serial number, is presented with the product to be repaired or replaced. OMNIC&S reserved the right to refuse warranty service if this information has been removed or changed after the original purchase of the product from the dealer.
- This Warranty does not cover any failure of the product due to normal wear and tear,

or due to misuse, including but not limited to use in other than the normal and customary manner, in accordance with instructions for use and maintenance of the product. Nor does this warranty cover any failure of the product due to accident, modification or adjustment, acts of God or damage resulting from liquid.

- This warranty does not cover product failures due to improper repair installations, modifications or service performed by a non-OMNIC&S authorized person. Tampering with any of seals on the product will void the warranty.
- OMNIC&S warrants this product to be free from defects in design, material, and workmanship at the time of its original purchase by a consumer, and for a subsequent period of one (1) year.

Battery

- Provided (in instruction manual)
- Replacement of Battery is limited only by service personnel.
- This product contains a battery. During using, storing, or transporting this product, do not place it in extremely hot or cold places, or at high altitude and low atmospheric pressure.
- Do not disposal of battery into fire or a hot oven.
- Do not leave the product including a battery in an extremely high temperature surrounding environment or the leakage of flammable liquid or gas.

OMNIFIT

www.omnifit.co.kr

2014 OMNIC&S Inc. All rights reserved

This product has been manufactured by, and is sold under the responsibility of OMNIC&S Inc., and OMNIC&S Inc. is the warrantor in relation to this product.

OMNIFIT® is a registered trade mark of OMNIC&S Inc.

- Address: #301, 288 Digital-ro, Guro-gu, Seoul, Republic of Korea

FCC Information

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions :

1. This Device may not cause harmful interface, and
2. This device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for CLASS B digital device, pursuant to Part 15 of FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which receiver is connected.
- Consult the dealer or experienced radio/TV technician for help.

WARNING

Changes or modifications not expressly approved by the manufacturer could void the user's authority to operate the equipment.

IMPORTANT NOTE:

FCC RF Radiation Exposure Statement:

This equipment complies with FCC RF radiation exposure limits set forth for an uncontrolled environment.

Frequently Asked Questions

How should I store the device?

After use, wipe the metal electrode clean and store it in a place that is not inclined, free from dust, liquid, and extreme temperatures.

What should I do if the device is damaged?

Do not attempt to modify, repair, or disassemble the device. Contact customer support for assistance.


Can I clean the device with water or alcohol?

Do not expose the device directly to water, alcohol, or other liquids for cleaning. Use a dry cloth for cleaning purposes.

How long is the warranty valid?

The warranty is valid based on the original receipt issued by the dealer. Ensure the receipt is presented for any warranty service.

Documents / Resources

	OMNIFIT OCW-H30 Enhancing Brain Training Solution [pdf] User Manual OCW-H30, OCW-H30 Enhancing Brain Training Solution, Enhancing Brain Training Solution, Brain Training Solution, Training Solution
---	--

References

- [User Manual](#)

OMNIFIT

Brain Training Solution, Enhancing Brain Training Solution, OCW-H30, OCW-H30 Enhancing Brain Training Solution, OMNIFIT, Training Solution

Leave a comment

Your email address will not be published. Required fields are marked *

Comment *

Name

Email

Website

☐ Save my name, email, and website in this browser for the next time I comment.

Post Comment

Search:

e.g. whirlpool wrf535swhz

Search

[Manuals+](#) | [Upload](#) | [Deep Search](#) | [Privacy Policy](#) | [@manuals.plus](#) | [YouTube](#)

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.