



AMD LASERS Monet Laser Curing Light User Guide

[Home](#) » [AMD LASERS](#) » AMD LASERS Monet Laser Curing Light User Guide 

Contents

- [1 AMD LASERS Monet Laser Curing Light](#)
- [2 QUICK START GUIDE](#)
- [3 LIGHT STATUS](#)
 - [3.1 Radiometer Check](#)
- [4 Documents / Resources](#)
 - [4.1 References](#)
- [5 Related Posts](#)



AMD LASERS Monet Laser Curing Light



QUICK START GUIDE

Laser Aperture

Rotating Head

Status Light

Switch

Magnetic removable battery

1. Charge the batteries prior to use

2. Utilize the personal protective equipment (Eyewear, Paddle, and Shield)

3. Connect the charged battery to the Monet device

4. Switch:
 1st Push: Standby mode
 2nd Push: Aiming beam
 3rd Push: Laser on (1 Sec.)
 4th Push: Stop Laser

5. To Emergency Stop the Monet- Remove the battery.

Radiometer Indicator

Battery Ports


Charge indicator

Charging Station

Radiometer

MONET

LIGHT STATUS

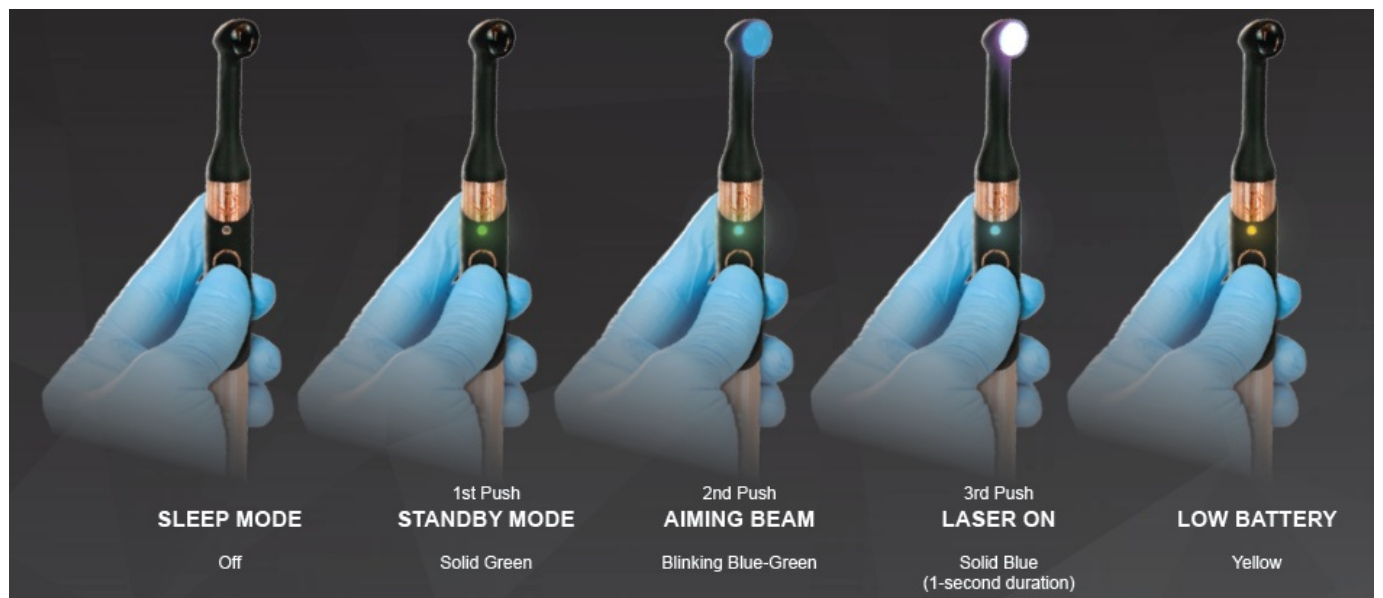

This example of cure time is using an A2 Composite material.

Composite Depth	Curing Cycle	With Attenuator
2mm or less	1 sec	Yes
2 – 3mm	1 sec	No
3 – 5mm	2 sec	No
5mm or more	3 sec	No

Use QR code to see a full matrix of recommended cure times with many brands of composite

Item	Aperture	Usage
	Attenuator	Close to pulp or gingival for adhesive and 2 mm or less composite layer
	2mm	Tacking, interproximal curing
	4mm	Narrow, deep preps, endo caps, interproximal curing, class V
	6mm	Small size, deep preps, interproximal curing

Radiometer Check

1.



Place Monet lens directly on sensor

2. Initiate Monet Laser
3. See results: Low (orange) Borderline
(Yellow) Good (Green)

Note: Radiometer is tuned for the full, unobstructed output of the unit.

Scan & Watch Video

Hosted by Mike DiTolla amdlasers.com/pages/monet-setup-guide

For more information see Owners Manual.

For additional Help contact Customer Service:

1.866.999.2635 | support@amdlasers.com

© 2021 AMD Lasers, Inc. All Rights Reserved. 05019.v06.21

Documents / Resources



[AMD LASERS Monet Laser Curing Light](#) [pdf] User Guide
Monet Laser Curing Light, Curing Light

References

- [🔥 Monet Setup Guide - amdlasers](#)