

[Manuals.plus](#) /

› [Waterpik](#) /

› Waterpik Ultra Water Flosser Instruction Manual

Waterpik WP-100

Waterpik Ultra Water Flosser Instruction Manual

Model: WP-100 | Brand: Waterpik

PRODUCT OVERVIEW

The Waterpik Ultra Water Flosser (Model WP-100) is a premium dental water jet designed for effective oral hygiene. It features an improved design, a large selection of accessories, and is clinically proven to be more effective than traditional flossing for improving gum health. This device is ideal for anyone seeking a thorough clean, especially those with braces, implants, crowns, bridges, or periodontal pockets.



The Waterpik Ultra Water Flosser unit, showcasing its compact design and clear water reservoir.

SETUP GUIDE

Before first use, ensure all packaging materials are removed. Place the unit on a flat, stable surface near an electrical outlet and a sink.

Filling the Reservoir

Lift the hinged lid of the reservoir and fill it with lukewarm water. You may also add mouthwash as desired. The reservoir has markings for ounces and milliliters, holding up to 19 ounces (600 ml) of water.



The removable water reservoir with clear measurement indicators.

Attaching and Changing Tips

The Waterpik Ultra Water Flosser comes with 6 color-coded tips, stored conveniently in a covered compartment on top of the unit. To attach a tip, insert it into the center of the handle until it clicks into place. To remove, press the eject button on the side of the handle and pull the tip out.



The convenient tip storage compartment located within the reservoir lid.



A selection of the 6 included tips: 2 Classic Jet Tips, Plaque Seeker Tip, Orthodontic Tip, Pik Pocket Tip, and Toothbrush Tip.

OPERATING INSTRUCTIONS

Power and Pressure Settings

Plug the unit into a suitable electrical outlet. The Waterpik Ultra Water Flosser features 10 adjustable water pressure settings, ranging from 10 to 90 PSI, delivering 1,400 pulses per minute. Use the dial on the front of the unit to select your desired pressure level. Start with a low setting and gradually increase as you become accustomed to the water pressure.



The ergonomic handle with a convenient pause button for temporary water flow stoppage.

Flossing Technique

Once the water reservoir is filled and a tip is attached, place the tip firmly in your mouth, leaning over the sink. Turn on the unit using the power switch. To begin water flow, release the water stop button on the handle. Direct the water jet at a 90-degree angle to your gumline, moving along the teeth and between them. Keep your lips slightly closed to prevent splashing while allowing water to flow into the sink. The device effectively removes debris and bacteria lodged deep between teeth and below the gum-line, massaging and stimulating the gums for improved oral health.



Proper usage of the Waterpik Ultra Water Flosser for effective cleaning.

CARE AND MAINTENANCE

To ensure the longevity and optimal performance of your Waterpik Ultra Water Flosser, regular cleaning is essential. Rinse the item thoroughly after each usage to avoid calcification of sugar content in the reservoir and hose. Periodically clean the exterior of the unit with a soft, damp cloth.

Tip Care

After each use, rinse the tip under running water. For a more thorough cleaning, tips can be washed with mild soap and water. Replace tips every 3-6 months, or sooner if they show signs of wear.

TROUBLESHOOTING


If you encounter issues with your Waterpik Ultra Water Flosser, consider the following common solutions:

- **No Power:** Ensure the unit is securely plugged into a working electrical outlet.
- **Weak Water Flow:** Check that the reservoir is filled with water and the tip is firmly seated in the handle. Ensure there are no kinks in the hose.
- **Leaking:** Verify that the reservoir is properly seated on the base and that the tip is securely attached.
- **Noise without Water Flow:** This may indicate an empty reservoir or an airlock. Turn off the unit, refill the reservoir, and restart.

For persistent issues, refer to the official user guide or contact Waterpik customer support.

PRODUCT SPECIFICATIONS

Feature	Detail
Model Number	WP-100
Product Dimensions	7 x 9 x 6 inches
Item Weight	2.2 Pounds
Power Source	Corded Electric
Pressure Range	10-90 PSI
Pulses Per Minute	1400
Reservoir Capacity	19 oz (650 mL) / 90+ seconds
Included Tips	6 (2 Classic Jet Tip, Plaque Seeker Tip, Orthodontic Tip, Pik Pocket Tip, Toothbrush Tip)

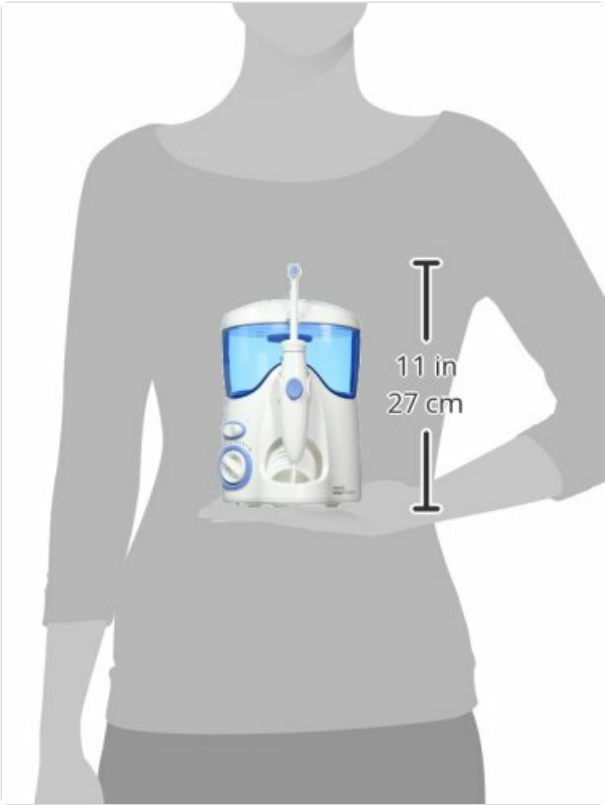

Product dimensions for easy placement and storage.

IMPORTANT SAFETY INFORMATION

Please read all instructions before using this appliance. When using electrical products, especially when children are present, basic safety precautions should always be followed.

- Keep out of reach of children.
- Unplug from outlet when not in use.
- Avoid submerging in water or other liquids.
- Use only as directed to avoid injury.

WARRANTY AND SUPPORT

For detailed warranty information and customer support, please refer to the official user guide provided with your product. You can also download the User Guide (PDF) [here](#).

Please note that the manufacturer’s return policy states: ALL SALES ARE FINAL UNLESS THE ITEM SOLD DOES NOT REFLECT THE DESCRIPTION.

© 2024 Waterpik. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question