

[Manuals.plus](#) /

› [Tria Beauty](#) /

› Tria Beauty Hair Removal Laser User Manual

Tria Beauty Tria Beauty Hair Removal Laser

Tria Beauty Hair Removal Laser User Manual

Brand: Tria Beauty

Model: Tria Beauty Hair Removal Laser

1. PRODUCT OVERVIEW

The Tria Beauty Hair Removal Laser is designed to provide permanent hair reduction results in the comfort of your home. Utilizing advanced laser technology, it targets hair follicles to inhibit future hair growth, offering an alternative to traditional methods like shaving and waxing.





Figure 1: The Tria Beauty Hair Removal Laser device, showcasing its ergonomic design and primary components.

This device is engineered for ease of use, allowing individuals to achieve professional-level results safely and effectively. It is intended for use on various body areas including legs, bikini line, underarms, feet, arms, hands, chest, and stomach.

2. IMPORTANT SAFETY GUIDELINES

Before using the Tria Beauty Hair Removal Laser, please read and understand all safety information. Improper use can lead to adverse effects.

- Must be used as directed.
- Intended for fair to medium skin tones.
- Do not use on moles.
- This product is not intended to diagnose, treat, cure, or prevent any disease.

Hair Removal Laser

The Tria Laser is designed to work on fair to light brown skin tones and light brown to black body hair. This chart is only representative of actual skin tones and hair colors.




Figure 2: Skin Tone and Hair Color Chart. This chart indicates suitable skin tones and hair colors for effective and safe use of the Tria Laser. Not intended for dark skin tones or light hair colors (blonde, white, gray, red).

Refer to the skin tone and hair color chart (Figure 2) to ensure your suitability for this device. The laser targets melanin in the hair follicle, making it most effective on dark hair and less effective or unsafe on very light hair or dark skin tones.

3. SETUP AND FIRST USE

Upon receiving your Tria Beauty Hair Removal Laser, follow these steps for initial setup:

1. **Unboxing:** Carefully remove all components from the packaging. Ensure all items are present as listed in the packing guide.
2. **Charging:** Before first use, fully charge the device. Connect the power adapter to the device and a power outlet. The charging indicator will show the charging status.
3. **Skin Sensor Activation:** The device is equipped with a skin sensor to ensure safe operation. Before each treatment, place the device's treatment window against the skin area you intend to treat. The device will only activate if your skin tone is suitable.

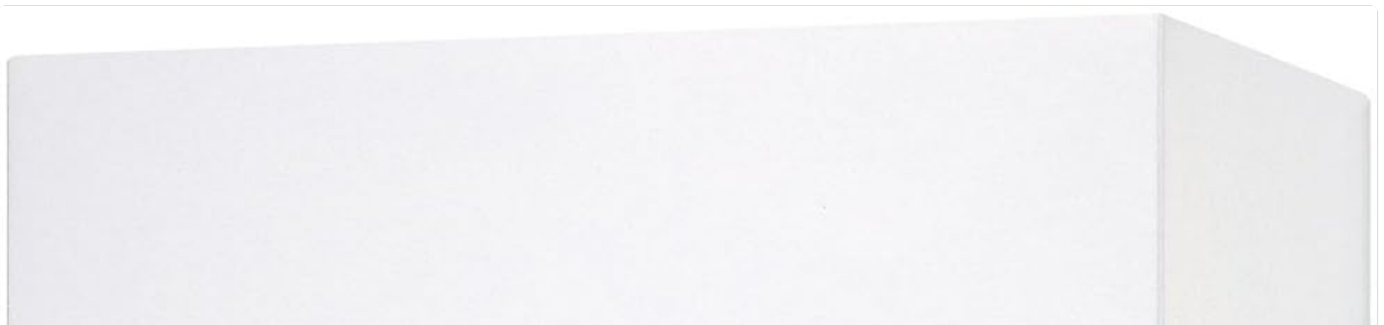

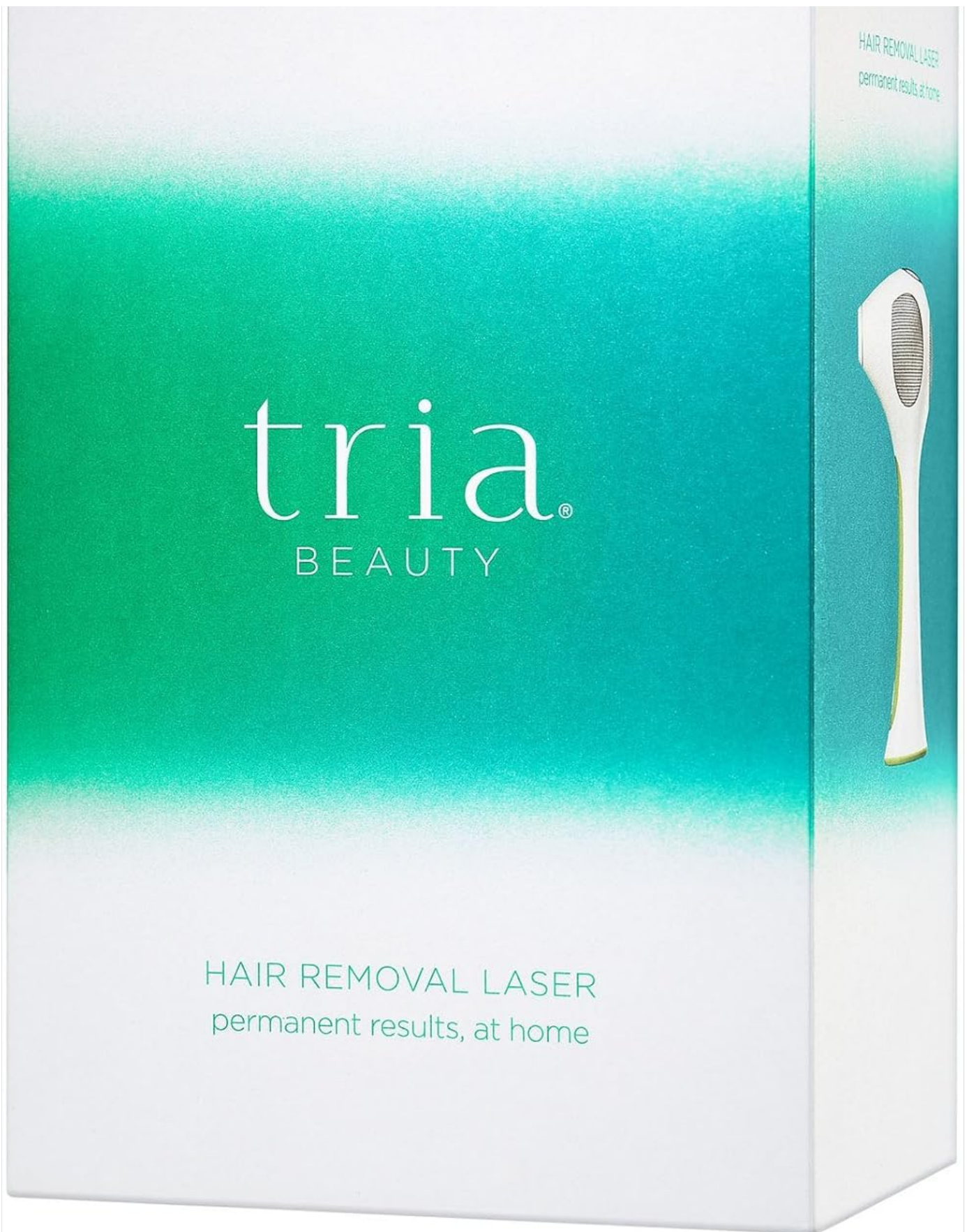


Figure 3: The Tria Beauty Hair Removal Laser packaging, indicating the contents and branding.

It is recommended to perform a patch test on a small area of skin 24 hours prior to full treatment to assess your skin's reaction.

4. OPERATING INSTRUCTIONS

Follow these steps for effective and safe hair removal treatment:

1. **Prepare the Skin:** Shave the area to be treated. The laser works best when hair is at the surface of the skin, not above it. Cleanse the skin to remove any lotions, oils, or makeup.
2. **Select Treatment Level:** Turn on the device and select your desired treatment level. Start with a lower setting and gradually increase as tolerated.
3. **Apply Treatment:** Place the treatment window firmly against your skin. The device will emit a pulse of laser light. Lift the device and move it to an adjacent area, ensuring slight overlap with the previous spot.
4. **Treat All Areas:** Continue this process until the entire desired area has been treated. The device is suitable for use on legs, bikini line, underarms, feet, arms, hands, chest, and stomach.
5. **Post-Treatment Care:** After treatment, you may apply a soothing lotion or gel to the treated area. Avoid sun exposure and use sunscreen on treated areas.





Figure 4: Side view of the Tria Beauty Hair Removal Laser, highlighting the laser aperture and ergonomic grip for precise application.

Consistency is key for optimal results. Follow the recommended treatment schedule provided in the full user manual for best outcomes.

5. MAINTENANCE AND CARE

Proper maintenance ensures the longevity and performance of your Tria Laser:

- **Cleaning:** After each use, gently wipe the treatment window and the device exterior with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case when not in use.
- **Battery Care:** To maintain battery health, avoid completely draining the battery frequently. Charge the device regularly, even if not in active use for extended periods.

6. TROUBLESHOOTING

If you encounter issues with your Tria Beauty Hair Removal Laser, refer to the common problems and solutions below:

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or power issue.	Ensure the device is fully charged. Check power adapter connection.
Laser does not activate.	Skin sensor not detecting suitable skin tone or improper contact.	Ensure the treatment window is in full contact with the skin. Verify your skin tone is suitable using the provided chart.
Device feels warm during use.	Normal operation or extended use.	This is normal. If it becomes uncomfortably hot, pause treatment and allow the device to cool down.
No noticeable hair reduction.	Inconsistent use, incorrect setting, or hair not in active growth phase.	Ensure consistent treatment schedule. Increase treatment level if comfortable. Hair removal requires multiple sessions as hair grows in cycles.

For issues not listed here, please contact customer support.

7. TECHNICAL SPECIFICATIONS

- **Package Dimensions:** 10.5 x 10.4 x 3.4 inches
- **Weight:** 3.7 Pounds
- **Manufacturer:** Tria Beauty
- **ASIN:** B0029OSN4U
- **UPC:** 891753001303, 891753001129, 891753001181

8. WARRANTY AND CUSTOMER SUPPORT

Tria Beauty products are manufactured to high-quality standards. For specific warranty details and terms, please refer to the warranty card included with your product or visit the official Tria Beauty website.

For further assistance, technical support, or to inquire about product registration, please visit the official Tria Beauty store or contact their customer service department:

[Visit the Tria Beauty Store on Amazon](#)

© 2023 Tria Beauty. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question