

Respironics Wisp Nasal Mask Replacement Cushion Small/Medium

Philips Respironics Wisp Nasal Mask Replacement Cushion User Manual

Model: Small/Medium Cushion

1. INTRODUCTION

This manual provides instructions for the proper use, maintenance, and replacement of the Philips Respironics Wisp Nasal Mask Replacement Cushion (Small/Medium). This cushion is designed for use with the Philips Respironics Wisp Nasal Mask system. Please read these instructions thoroughly before use to ensure safe and effective operation.

2. SAFETY INFORMATION

Important: This device is intended for single-patient use in the home and multi-patient use in the hospital/institutional environment. Always consult your healthcare provider for specific instructions regarding your CPAP therapy.

- Do not use the cushion if it appears damaged or discolored.
- Regularly inspect the cushion for wear and tear. Replace as needed to maintain an effective seal and comfort.
- Clean the cushion as instructed to prevent bacterial growth and maintain hygiene.
- Keep out of reach of children.

3. PACKAGE CONTENTS

Your package should contain the following item:

- One (1) Philips Respironics Wisp Nasal Mask Replacement Cushion (Small/Medium)

REF 1094087



S/M, WSP CUSHION, RP


IT: Piccolo/medio, materiale ammortizzante

Wisp, ricambio

12NC



884314600010

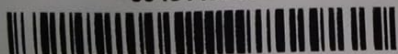

LOT 0302977816 CR 2023-08-22



HCPCS A7032


Made in Costa Rica

UDI



(01)00606959407428(10)0302977816

Respironics Inc.
1001 Murry Ridge Lane
Murrysville, PA 15668 USA

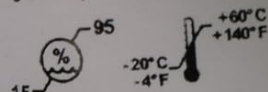

UPC: 606959407428



6 06959 40742 8

For Internal Use Only



PN: 300000949171

PN: 300004449111.E

Image: The Philips Respironics Wisp Nasal Mask Replacement Cushion, Small/Medium. This image displays the clear silicone cushion, which is the primary component of the Wisp nasal mask that makes contact with the user's face.



Image: The Philips Respironics Wisp Nasal Mask Replacement Cushion, still sealed in its clear plastic packaging, with the "PHILIPS RESPIRONICS" logo visible on the packaging. This shows the product as it would appear when first unboxed.

4. SETUP AND INSTALLATION

To replace your Wisp Nasal Mask cushion:

1. **Remove Old Cushion:** Gently pull the existing cushion away from the mask frame. It should detach with minimal effort.
2. **Align New Cushion:** Take the new Small/Medium replacement cushion. Align the cushion with the opening on the Wisp mask frame. Ensure the cushion is oriented correctly, with the narrower part typically pointing towards the top (nose bridge).
3. **Attach Cushion:** Press the edges of the new cushion firmly onto the mask frame until it snaps securely into place all around. Run your finger along the seam to ensure a complete and even seal.
4. **Inspect:** Verify that the cushion is fully seated and there are no gaps between the cushion and the mask frame.



Image: A detailed view of the Philips Respironics Wisp Nasal Mask Replacement Cushion. This image emphasizes the contours and the attachment mechanism of the cushion, illustrating how it connects to the mask frame.

5. OPERATING INSTRUCTIONS

Once the replacement cushion is installed, reassemble your Wisp Nasal Mask with the headgear. Place the mask on your face and adjust the headgear straps to achieve a comfortable and effective seal. The cushion should rest gently around your nose without excessive pressure. Connect the mask to your CPAP machine and turn on the airflow to check for leaks. Adjust the headgear as necessary to minimize leaks while maintaining comfort.

For optimal therapy, ensure the cushion forms a good seal without being overtightened, which can cause discomfort and skin irritation.

6. MAINTENANCE

6.1 Cleaning Instructions

Regular cleaning of your cushion is essential for hygiene and to prolong its lifespan.

1. **Daily Cleaning:** Detach the cushion from the mask frame. Wash the cushion daily in warm, soapy water using a

mild dish soap.

2. **Rinse Thoroughly:** Rinse the cushion thoroughly under running water to remove all soap residue.
3. **Air Dry:** Allow the cushion to air dry completely out of direct sunlight before reattaching it to the mask frame. Do not use a clothes dryer.

Do not use bleach, alcohol, or cleaning solutions containing strong chemicals, as these can damage the silicone material and reduce the cushion's effectiveness.

6.2 Replacement Schedule

The Wisp Nasal Mask Replacement Cushion should be replaced regularly to maintain an effective seal and optimal therapy. General guidelines suggest replacing the cushion every 1 to 3 months, or sooner if:

- The cushion becomes discolored or stiff.
- You notice persistent leaks that cannot be resolved with adjustments.
- The material shows signs of cracking, tearing, or deterioration.
- You experience discomfort or skin irritation from the cushion.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Air leaks around the cushion	• Cushion not properly seated • Headgear too loose or too tight • Cushion worn out or damaged • Incorrect cushion size	• Re-seat the cushion firmly onto the mask frame. • Adjust headgear straps for a snug, but comfortable, fit. • Inspect cushion for damage; replace if necessary. • Ensure you are using the correct size (Small/Medium).
Discomfort or skin irritation	• Headgear too tight • Cushion not clean • Cushion worn out	• Loosen headgear straps slightly. • Clean the cushion thoroughly. • Replace the cushion if it's old or deteriorated.

8. SPECIFICATIONS

- **Product Name:** Wisp Nasal Mask Replacement Cushion
- **Size:** Small/Medium
- **Manufacturer:** Respironics Inc.
- **ASIN:** [B013CXALMI](#)
- **Reference Number (REF):** [1094087](#) (information from packaging)
- **UPC:** [0606959407428](#) (information from packaging)
- **Lot Number:** [0302977816](#) (information from packaging)
- **Material:** Silicone (typical for CPAP cushions)
- **Package Dimensions:** 7.52 x 6.22 x 1.97 inches
- **Weight:** 0.63 ounces
- **Date First Available:** March 13, 2015

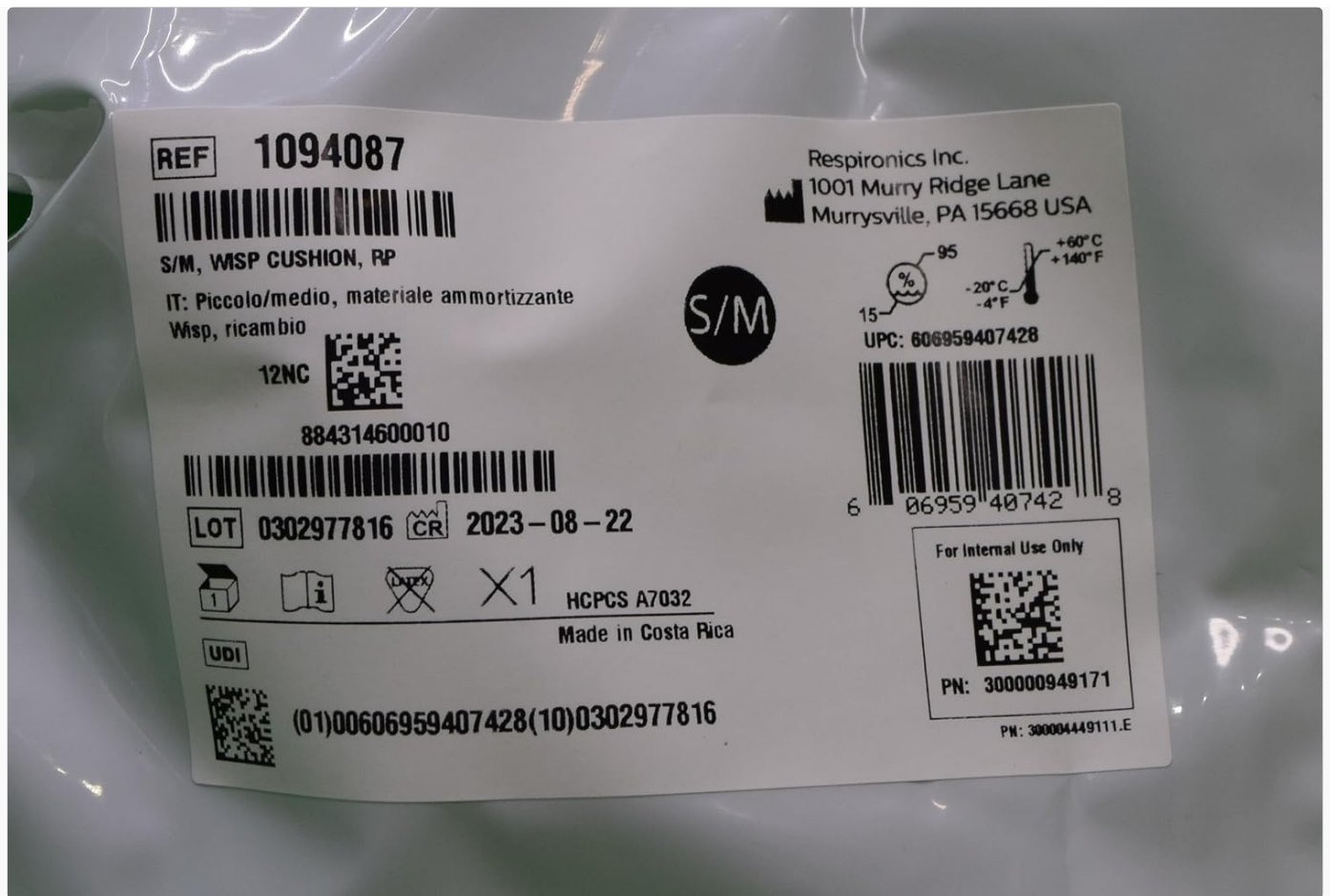


Image: A Philips Respironics Wisp Nasal Mask Replacement Cushion, used here to visually represent the product from which specifications and identifiers like REF, UPC, and Lot numbers are typically found on its packaging.

9. WARRANTY AND SUPPORT

For specific warranty information, please refer to the documentation provided with your original Philips Respironics Wisp Nasal Mask system or contact Philips Respironics directly. Replacement cushions typically fall under the general warranty terms for accessories.

For technical support or further assistance, please contact your medical equipment provider or Philips Respironics customer service. Contact information can usually be found on the manufacturer's official website or on the original product packaging.

© 2026 Philips Respironics. All rights reserved. Information subject to change without notice.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question