Saved from manuals.plus.
Open the live page for links to other documents and user manuals.

Accu-Answer UA01

Accu-Answer UA01 Uric Acid Test Strip User Manual

For use with Accu-Answer Multi-Monitoring Systems

1. Product Overview

The Accu-Answer UA01 Uric Acid Test Strips are designed for quantitative measurement of uric acid in fresh capillary whole blood. These strips are intended for self-testing and professional use to monitor uric acid levels at home or in a clinical setting. Each pack contains 25 test strips.



Image 1.1: Two boxes and two vials of Accu-Answer Uric Acid Test Strips. Each vial contains multiple strips, and the packaging indicates compatibility with Accu-Answer Multi-Monitoring Meters.

Key Features:

- Provides a rapid and effective method for measuring uric acid.
- Designed for ease of use in home monitoring.
- CE certified and approved for use in the EU.

2. Important Safety Information

Caution: These test strips are only compatible with Accu-Answer Multi-Monitoring devices. Do not attempt to use them with other brands or models of monitoring systems, as this may lead to inaccurate results or device malfunction.
Always read the user manual for your Accu-Answer Multi-Monitoring device before performing a test. Keep test strips out of reach of children.

3. Package Contents

- 25 Accu-Answer Uric Acid Test Strips (UA01)
- Instruction leaflet (within the product packaging)

Note: The Accu-Answer Multi-Monitoring device and lancing device are sold separately.

4. Setup and Preparation

Before performing a test, ensure you have your Accu-Answer Multi-Monitoring device, a lancing device, lancets, and the Accu-Answer Uric Acid Test Strips ready. Wash and dry your hands thoroughly.



Image 4.1: Visual guide for the testing process with an Accu-Answer meter and test strips.

Step-by-Step Guide:

1. **Insert a Test Strip:** Insert a new Accu-Answer Uric Acid Test Strip into the test strip port of your Accu-Answer Multi-Monitoring device. This action will automatically turn on the meter.
2. **Match Code Number:** Compare the code number displayed on the meter with the code number printed on the test strip vial. If they do not match, adjust the meter's code manually according to your device's instructions.
3. **Ready for Testing:** Once the code numbers match and the meter is ready, a blood drop symbol will appear on the display, indicating the meter is prepared for testing.

5. Operating Instructions

After completing the setup steps, proceed with the blood sample collection and testing:

4. **Apply Blood Sample:** Use a lancing device to obtain a blood sample. Touch the tip of the test strip to the blood drop. The strip will draw blood into the confirmation window. Ensure the confirmation window is completely filled.
5. **View Results:** The meter will begin counting down. After approximately 5 seconds, your uric acid test result will appear on the display. The result is automatically stored in the meter's memory.
6. **Eject Used Strip:** Press the ejection button on your meter to safely dispose of the used test strip. The meter will automatically turn off after approximately 2 minutes of inactivity.

6. Storage and Handling

Proper storage and handling of your test strips are crucial for accurate results and product longevity.



Image 6.1: A single vial of Accu-Answer Uric Acid Test Strips with three individual strips laid out, showing their design.

- Store test strips in their original vial with the cap tightly closed.
- Keep the vial in a cool, dry place, away from direct sunlight and heat.
- Do not refrigerate or freeze test strips.
- Use test strips within the expiration date printed on the vial and packaging.
- Once the vial is opened, use the strips within the specified period (typically 3 months, refer to vial label).
- Do not bend, cut, or alter the test strips in any way.

7. Product Specifications

Specification	Detail
Brand	Accu-Answer
Model	UA01 Uric Acid Test Strip
Number of Pieces	25 strips per vial
Special Feature	Rapid Test
Specific Product Uses	Uric Acid Monitoring
Compatibility	Accu-Answer Multi-Monitoring devices only

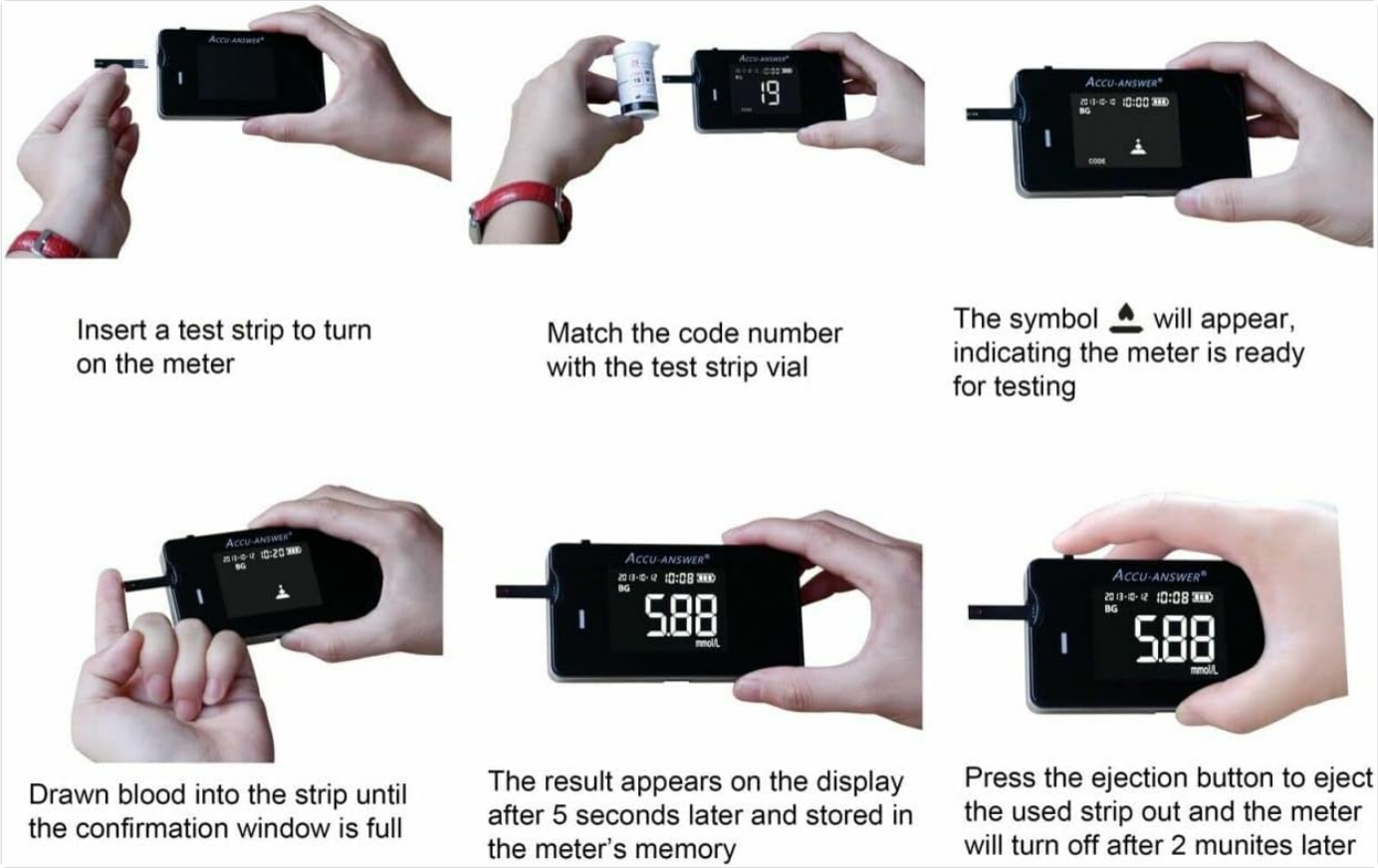


Image 7.1: An Accu-Answer Multi-Monitoring device displayed with different types of test strips, including uric acid strips, highlighting its versatility.

8. Troubleshooting

If you encounter issues while using the Accu-Answer Uric Acid Test Strips, consider the following:

- **Inaccurate Results:** Ensure the test strip is not expired, the vial cap was tightly closed, and the code on the meter matches the vial. Verify proper blood sample application.
- **Meter Error Message:** Refer to your Accu-Answer Multi-Monitoring device's user manual for specific error codes and troubleshooting steps.
- **No Result Displayed:** Check if the test strip was inserted correctly and if enough blood was applied to the strip's confirmation window.

If problems persist, contact Accu-Answer customer support or consult a healthcare professional.

9. Warranty and Support

For information regarding the warranty of your Accu-Answer Multi-Monitoring device or the test strips, please refer to the documentation included with your monitoring system or visit the official Accu-Answer website. For technical support or inquiries, please contact Accu-Answer customer service through the contact information provided in your device's manual.

© 2023 Accu-Answer. All rights reserved.

Ask a question about this manual

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. Privacy details.

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question

Open the live page

https://manuals.plus/asin/B08SHHDZFB

More documents and user manuals: **https://manuals.plus/**