

Manuals.plus /
› BEBIRD /
› BEBIRD Note3 Pro Max User Manual

BEBIRD Note3 Pro Max

BEBIRD Note3 Pro Max User Manual

Ultimate Version Smart Visual Ear Cleaning Tweezer

1. INTRODUCTION

The BEBIRD Note3 Pro Max is an advanced smart ear wax removal tool designed for safe and effective ear cleaning. It features a high-definition camera, dual operating modes (tweezer and ear picker), and seamless integration with a mobile application for real-time visual guidance. This manual provides comprehensive instructions for setting up, operating, maintaining, and troubleshooting your device to ensure optimal performance and user safety.



bebird® | Note3 Pro Max

Ultimate Version



Image: BEBIRD Note3 Pro Max device alongside a smartphone displaying the live ear canal view.

2. SETUP GUIDE

2.1. App Download and Installation

To begin, download the official BEBIRD application on your smartphone or tablet. The app is compatible with both Android and iOS systems.

- Scan the QR code provided in the product packaging or on the device itself.
- Alternatively, visit the official BEBIRD app download page: <http://www.bebird.cn/app/en.html>
- Follow the on-screen instructions to install the application.

1 Scan QR code to download app



2 Start ear-cleaning rod



3 Connect WiFi Bebird-XXXX



4 Return the app Click on "start"



About Note3 Pro Max



So far, only doctors do advanced ear cleaning. With BEBIRD, the story becomes simple and effective. The Bebird note3 pro is the top model of the BEBIRD family, bringing a fresh design and the all-important self-care you do often, making it easier than ever

20,000,000+	100+	200+
Use and download	Countries and regions	Hard core patent








Image: Visual guide for app download and initial connection steps.

2.2. Device Connection

Once the app is installed, connect your BEBIRD Note3 Pro Max to your device:

1. Power on the BEBIRD Note3 Pro Max. The device will emit a Wi-Fi signal.
2. Go to your phone's Wi-Fi settings and connect to the network named "bebird-XXXX" (where XXXX is a unique identifier for your device).
3. Open the BEBIRD app. The device should automatically connect. If not, select the "Note3 Pro Max" model from the device list within the app.
4. Tap "Start" to begin the visual ear cleaning process.

Your browser does not support the video tag.

Video: Demonstrates the setup and connection process for the BEBIRD Note3 Pro Max.

3. OPERATING INSTRUCTIONS


3.1. Understanding Dual Modes: Tweezer and Ear Picker

The BEBIRD Note3 Pro Max offers two primary modes for ear wax removal, catering to different needs and types of earwax:


- **Tweezer Mode:** Utilizes a high-precision micro-mechanical arm with rubber-covered tips. This mode is ideal for grasping and removing complex or larger foreign bodies and stubborn earwax. The open and close degree protection system ensures safe operation.
- **Ear Picker Mode:** Features a spade-shaped food-grade silicone ear scoop. This mode is suitable for general ear cleaning, effective for both oily and dry earwax.

Simulated Use Scenarios

Earpick Mode

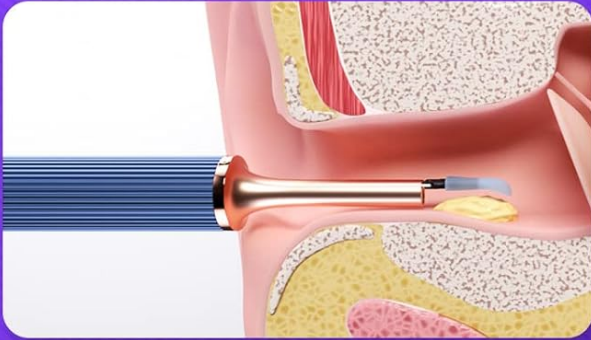

Tweezers mode

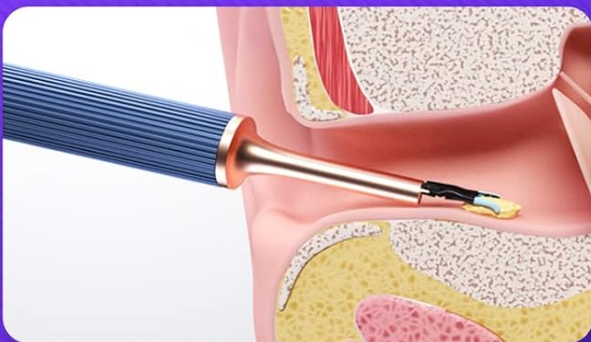

Tweezers Mode & Earpick Mode

Send real-time in-ear video to APP


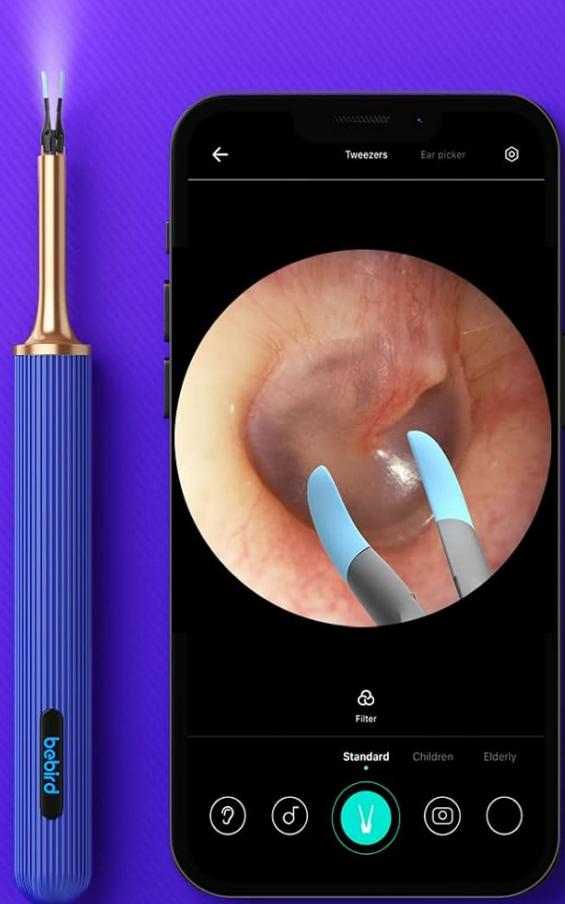


Image: Simulated use scenarios illustrating the ear pick and tweezer modes within the ear canal.



Enjoy a God-Like ear picking experience

The latest technology ensures you handle the buildup accurately without damaging the ear canal and eardrum

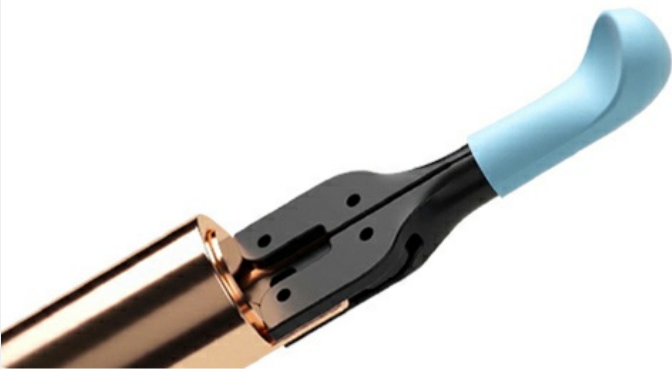

100%
Soft Smooth

30mm
Depth Focus

Six
Gyroscope

5 LED
Highlighting

Integrates proven micro-robotic arm technology

Clip tweezer mode with smaller tweezer tip for narrow ear canals and stubborn ear wax



80g
Force

10.2°
Degree

0.02mm
Precision

Safe/On-off



Image: Detailed view of the micro-robotic arm technology and ear spoon attachment mechanism.

3.2. Visual Cleaning and Safety Features

The device integrates a 10-megapixel HD endoscope with a 360° wide-angle lens and 6-axis gyroscope for precise and safe operation. The latest Wi-Fi chip ensures high-quality, real-time image transmission to your smartphone.

- **Real-time View:** Observe the ear canal clearly on your phone screen, allowing for accurate and controlled cleaning.
- **6-Axis Gyroscope:** Provides stable and accurate measurement, ensuring precise movement and preventing accidental contact with sensitive areas.
- **LED Shadowless Light:** Illuminates the ear canal effectively, eliminating shadows for a clear view.
- **Non-Sensing Temperature:** The device maintains a temperature close to the human body, ensuring comfort and safety during use, even for children.
- **App Features:** The BEBIRD app supports left and right ear modes, allows you to take pictures and record videos of the cleaning process, and can be used for daily checks of ear, nose, and mouth health.

State-of-the-art Technology

2-IN-1

Tweezer and Rod

5-LED

Astral Lamp

6-Axis

Gyroscope

High-Quality

Anti-slip Material

Super-WiFi

Fast WiFi Chip

Type-C

Charging Base

1000W HD

Professional Otoscope

77F°

Otoscope Non-sensing
temperature

85F°

Otoscope Non-sensing
temperature

#1
Winner

BRAND NEW TECHNOLOGY

bebird

Image: Overview of the state-of-the-art technology integrated into the BEBIRD Note3 Pro Max.

3.3. Accessories and Family Use

The BEBIRD Note3 Pro Max comes with 12 different accessories in 4 distinct colors, making it suitable for the entire family. The food-grade silicone ear spoon heads and robotic arm are round and skin-friendly, ensuring a safe and comfortable experience for all ages.



Image: Contents of the BEBIRD Note3 Pro Max package, highlighting the various ear spoons and cleaning tools.

4. MAINTENANCE

4.1. Cleaning the Device

Regular cleaning ensures hygiene and extends the lifespan of your device:

- After each use, remove the ear spoon or tweezer head.
- Use the provided cleaning brush to remove any visible earwax or debris from the tip and camera lens.
- For thorough cleaning, use the alcohol cotton swabs to gently wipe the camera lens and the accessory tips. Ensure the lens is completely dry before next use.
- Do not submerge the main device body in water.

4.2. Charging

The BEBIRD Note3 Pro Max features a convenient Type-C magnetic charging base.

- Place the device onto the magnetic charging base. Ensure it is properly seated for charging to begin.
- Connect the Type-C cable to the charging base and a suitable power source (e.g., USB wall adapter, computer USB port).
- The device's indicator light will show charging status.

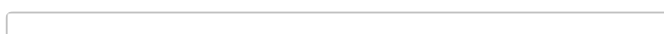


Image: The magnetic charging base and accessory storage for the BEBIRD Note3 Pro Max.

5. TROUBLESHOOTING

If you encounter issues with your BEBIRD Note3 Pro Max, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
No camera signal / Black screen in app	◦ Device not connected to Wi-Fi. ◦ App not recognizing device. ◦ Camera malfunction.	◦ Ensure your phone is connected to the "bebird-XXXX" Wi-Fi network. ◦ Restart the BEBIRD app and the device. ◦ Check if the device is fully charged. ◦ If the issue persists, contact customer support.
Video latency or hiccups	◦ Weak Wi-Fi signal. ◦ Interference from other devices. ◦ Outdated app version.	◦ Ensure you are close to the device for a strong Wi-Fi connection. ◦ Close other applications running in the background on your phone. ◦ Update the BEBIRD app to the latest version.
Device not charging	◦ Improper placement on charging base. ◦ Faulty charging cable or adapter.	◦ Ensure the device is correctly seated on the magnetic charging base. ◦ Try a different Type-C cable and power adapter.
Tweezer tips or ear spoons break easily	◦ Excessive force during use. ◦ Wear and tear.	◦ Use gentle pressure when operating the device. ◦ Replacement accessories can be purchased from authorized retailers.

6. SPECIFICATIONS

Feature	Detail
Model	Note3 Pro Max
Camera Resolution	10 Megapixel HD
Lens Length	3.7 mm
Gyroscope	6-Axis Smart Gyroscope
Connectivity	High-Speed Wi-Fi Chip
Battery Capacity	300mAh (approx.)
Charging Interface	Type-C Magnetic Charging Base
Accessories Included	12 different ear spoons (3 types, 4 colors), cleaning brush, alcohol cotton swabs, liquid cotton swabs.
Package Dimensions	8.19 x 3.54 x 2.6 inches
Weight	12.31 ounces
Manufacturer	BEBIRD
Date First Available	April 1, 2021

		Ultimate bebird Note3 Pro Max			
C3	x17			Tvifm	Removal Tools
					
3Million	3Million	Camera	10Million Pixe	0.92Million	N/A
15-20mm	15-20mm	Depth Focus	10-30mm	15-20mm	N/A
3.9mm	3.9mm	Camera Size	3.7mm	3.9mm	N/A
3-Axis	3-Axis	Gyroscope	6-Axis	3-Axis	N/A
250mAh	300mAh	Battery	300mAh	200mAh	N/A
N/A	N/A	Head Structure	Rod/Tweezers	N/A	N/A
N/A	Type C	Fast Charge	Type C	N/A	N/A
		 Note3 pro max has upgraded various core components and processed it with the world's top production line, which is better than ever			

Image: Detailed comparison of BEBIRD models, highlighting the Note3 Pro Max's features.

7. WARRANTY AND SUPPORT

For warranty information, technical support, or any inquiries regarding your BEBIRD Note3 Pro Max, please refer to the official BEBIRD website or contact their customer service directly. Keep your purchase receipt as proof of purchase for warranty claims.

- **Manufacturer:** BEBIRD
- **Official Website:** Please visit the BEBIRD official website for the most up-to-date support information and contact details.
- **App Support:** For app-related issues, ensure your app is updated to the latest version.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question

Open the live page

<https://manuals.plus/asin/B091J7363G>

More documents and user manuals: <https://manuals.plus/>