

GROWNSY BC023

GROWNSY Electric Baby Nasal Aspirator Instruction Manual

Model: BC023

INTRODUCTION

Thank you for choosing the GROWNSY Electric Baby Nasal Aspirator. This device is designed to provide gentle and effective relief for babies and toddlers experiencing nasal congestion. Featuring adjustable suction levels, soothing music, and a built-in light, it aims to make the nose-clearing process comfortable for your child. Please read this manual thoroughly before use and keep it for future reference.

WHAT'S INCLUDED

Your GROWNSY Electric Baby Nasal Aspirator package includes the following items:

- 1 x Electric Nasal Aspirator Unit
- 3 x Food-Grade Silicone Tips (various sizes for different age ranges)
- 1 x USB-C Charging Cable
- 1 x Storage Case
- 1 x User Manual

Certified for Safety & Newborn Care

3 Food-grade Silicone Tips



Image: All components of the GROWNSY Electric Baby Nasal Aspirator, including the main unit, three silicone tips, USB-C charging cable, and storage case.

SETUP

1. **Initial Charge:** Before first use, fully charge the device using the provided USB-C cable. The battery indicator on the display will show charging status. A full charge can last up to 30 days with typical use.
2. **Assemble Suction Chamber:** Attach the transparent mucus collection chamber to the main aspirator unit. Ensure it is twisted clockwise until it locks securely into place.
3. **Attach Silicone Tip:** Select the appropriate food-grade silicone tip for your child's age and nasal passage size. Gently push the chosen tip onto the top of the suction chamber until it is firmly seated.

Essential Item for A Baby Registry



✓ Gets snot out quickly


✓ Gentle on baby

✓ multifunctional

✓ Easy to clean

BEST ELECTRIC NASAL ASPIRATOR

what to expect.

 babycentre.

 Babygaga




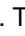
 babylist

 Best Products Reviews

Parents the **BUMP**

Image: The GROWNSY Electric Baby Nasal Aspirator with a silicone tip attached, ready for use.

OPERATING INSTRUCTIONS

- Turn On:** Press and hold the power button  for 3 seconds to turn on the device. The display will illuminate.
- Adjust Suction Level:** Press the fan icon button  to cycle through 3 adjustable suction levels (SP-1: Weak, SP-2: Moderate, SP-3: Strong). Choose the level appropriate for the type and consistency of nasal mucus.
- Activate Light & Music:** Optionally, press the light bulb icon button  to turn on the built-in soothing light, and the music note icon button  to play calming lullabies. These features can help distract and comfort your baby during use.

4. **Perform Aspiration:** Gently insert the silicone tip into your baby's nostril. Ensure a snug fit to create effective suction. The device will automatically clear mucus. Repeat for the other nostril as needed.
5. **Turn Off:** Press and hold the power button for 3 seconds to turn off the device after use.

Patent Design

- ✔ Instant Relief
- ✔ Noninvasive and Safe
- ✔ Automatic Suction



Image: A baby being comforted by the aspirator's light and music functions during use.

Video: Official GROWNSY video demonstrating the use of the Baby Nasal Aspirator, highlighting its features and ease of operation.


MAINTENANCE

1. **Disassembly:** After each use, remove the silicone tip and twist the mucus collection chamber

counter-clockwise to detach it from the main unit.

2. **Cleaning Tips & Chamber:** Rinse the silicone tips and mucus collection chamber thoroughly with warm soapy water. The silicone nozzle can be sterilized in boiling water.
3. **Cleaning Main Unit:** The main aspirator unit is not waterproof. Do not immerse it in water. Wipe the exterior with a damp cloth.
4. **Drying:** Ensure all detachable parts are completely dry before reassembling or storing the device.
5. **Motor Maintenance:** If liquid accidentally enters the motor, remove the storage chamber and run the machine for 5-10 minutes to drain the liquid. Wipe away any stains. This helps prevent damage and prolongs device life.
6. **Storage:** Store the clean and dry aspirator and its components in the provided storage case, away from direct sunlight and extreme temperatures.

USB-C Rechargeable Convenience

 Enjoy 120Min of continuous use



No filters required, battery-included:
hassle-free nasal aspirator for life

Image: The GROWNSY Electric Baby Nasal Aspirator disassembled into its main components for easy cleaning.

TROUBLESHOOTING

- **Device not turning on:** Ensure the device is fully charged. The unit will not power on while charging. If the battery is low, connect it to a USB power source using the provided cable.
- **Weak suction:** Check that the suction chamber and silicone tip are securely attached and properly sealed. Ensure there are no blockages in the tip or chamber.
- **Power button unresponsive:** Some users have reported needing to press the power button firmly for 3 seconds to activate.
- **Liquid in motor:** Follow the motor maintenance steps outlined in the Maintenance section.

SPECIFICATIONS

Feature	Detail
Brand	GROWNSY
Model Name	Electric Nasal Aspirator
Model Number	BC023
Age Range	0-12 Years Old
Material Type	Plastic, Food-Grade Silicone (PVC Free, Latex Free, Lead Free, Nitrosamine Free, Phthalate Free, Bisphenol-A (BPA) Free)
Power Source	Battery Powered (USB-C Rechargeable)
Suction Levels	3 Adjustable Levels
Special Features	Music and Light Soothing Function, Anti-backflow Technology, Larger Suction Chamber
Water Resistance Level	Not waterproof (main unit), Silicone nozzle sterilizable

WARRANTY AND SUPPORT

For warranty information or customer support, please refer to the contact details provided with your product packaging or visit the official GROWNSY website. Ensure to keep your purchase receipt for warranty claims.

Important Notes:

- The device is designed for operation by adults. Keep it out of reach of children to prevent accidental damage.
- Always follow the user manual for correct assembly and operation to avoid damaging the device.

