

McKesson 5059

McKesson Quintet AC Blood Glucose Test Strips - Instruction Manual

Model: 5059

1. INTRODUCTION

The McKesson Quintet AC Blood Glucose Test Strips are designed for use with the McKesson QUINTET AC Blood Glucose Monitoring System and QUINTET AC Blood Glucose Control Solution. These strips are intended for measuring glucose in fresh capillary whole blood samples.

They are suitable for self-testing at home and for multiple-patient use in professional healthcare settings, such as long-term care facilities and medical offices. The auto-coding feature ensures quick and accurate testing, aiding in the effective monitoring of diabetes.

Key Features:

- Accurate results and quick response time.
- Designed for use with McKesson QUINTET AC Blood Glucose Monitoring System.
- Features a unique glove-grip paddle for caregiver's convenience.
- Requires a minimal sample size of 1 μ L.
- Auto-coding for simplified and precise testing.
- Single-use strips for hygiene and accuracy.

2. PRODUCT COMPONENTS

Each package of McKesson Quintet AC Blood Glucose Test Strips contains the following:




An overview of the McKesson Quintet AC Blood Glucose Test Strips packaging, including the box, a sealed vial of strips, and a single test strip.

QUINTET AC

McKES

■ Blood Glucose Test Strips

25

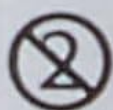

PER VIAL

DO NOT REUSE

MFR # 5

Read the insert before using. Use with McKesson Quintet AC Blood Glucose Meter. For *in vitro* diagnostic use only. For testing glucose in capillary whole blood use only. For testing glucose in capillary whole blood use only. Don't store strips outside of vial. Use within 3 months of opening the vial. Re-cap vial immediately after use. Store between temperatures 39-86°F (4-30°C). **DO NOT FREEZE.**

101-165260-051
ENUS0

Technical Support? Call 1-888-481-8485
Distributed by McKesson Medical-Surgical Inc.
Richmond, VA 23233 PVN A0915 Made in Taiwan

A close-up view of the McKesson Quintet AC Blood Glucose Test Strips vial, showing the label with product information and usage instructions.



A detailed view of a single McKesson Quintet AC Blood Glucose Test Strip, highlighting its design and the area for blood application.

3. SETUP

Before using the test strips, ensure you have the McKesson QUINTET AC Blood Glucose Monitoring System and, if needed, the control solution. Always read the owner's manual for your blood glucose meter for complete instructions. Upon opening a new vial of test strips, note the discard date. The strips should be used within 3 months of opening the vial. Do not use strips past their expiration date printed on the box and vial.



The back of the test strip box provides important information regarding in vitro diagnostic use, sample sources, and storage conditions.

4. OPERATING INSTRUCTIONS

Follow these general steps for blood glucose testing. Always refer to your McKesson QUINTET AC Blood Glucose Meter's owner's manual for specific instructions on meter operation.

1. **Prepare:** Wash your hands thoroughly with soap and warm water. Dry them completely.
2. **Obtain Blood Sample:** Use a lancing device to obtain a blood sample from your fingertip, palm, or forearm, as instructed by your healthcare professional and meter manual.
3. **Apply Blood to Strip:** Gently touch the tip of the test strip to the blood sample. The strip will draw the blood into the reaction area. Ensure the blood sample is sufficient to fill the strip's confirmation window.
4. **Insert Strip into Meter:** Insert the test strip into the McKesson QUINTET AC Blood Glucose Meter. The meter will automatically code and begin the test.
5. **Read Results:** The blood glucose result will appear on the meter's display within a few seconds.

Your browser does not support the video tag.

This video demonstrates the use of McKesson QUINTET AC Blood Glucose Test Strips with the monitoring system, showing the process of obtaining a blood sample and applying it to the strip.



A hand is shown holding a McKesson Quintet AC test strip, with a small blood sample visible on the strip's application area, ready for testing.

5. MAINTENANCE AND STORAGE

Proper handling and storage of your test strips are crucial for accurate results and product longevity.

- Always keep test strips in their original vial with the cap tightly sealed immediately after removing a strip.
- Never transfer strips from one vial to another.
- Never reuse strips; they are for single-use only. Discard used test strips properly after each test.
- Store test strips between temperatures of 39-86°F (4-30°C). Do not freeze.
- Avoid storing strips outside of the vial for extended periods.
- This product is not made with natural rubber latex.



A hand is shown carefully inserting a used test strip back into its vial, emphasizing the importance of proper disposal and re-capping the vial.

6. TROUBLESHOOTING

If you encounter issues during testing, consider the following common troubleshooting tips:

- **Error Codes:** Refer to your McKesson QUINTET AC Blood Glucose Meter's owner's manual for explanations of specific error codes and recommended actions.
- **Insufficient Blood Sample:** If the meter indicates an insufficient blood sample, discard the used strip and repeat the test with a new strip and a fresh blood sample. Ensure the blood drop is large enough to be fully absorbed by the strip.
- **Expired or Damaged Strips:** Check the expiration date on the vial and box. Do not use expired or visibly damaged strips. Ensure the vial cap was tightly sealed.
- **Storage Conditions:** Verify that the strips have been stored within the recommended temperature range and away from direct sunlight or extreme humidity.

- **Control Solution Test:** If you suspect inaccurate results, perform a control solution test as described in your meter's manual to verify the meter and strip functionality.

For persistent issues, contact customer support.

7. SPECIFICATIONS

Feature	Detail
Product Dimensions	1.4 x 3.2 x 4.8 inches
Item Weight	6.7 ounces
Manufacturer	McKesson Brand
ASIN	B09X6DK6KS
Item Model Number	5059
Date First Available	April 4, 2022
Number of Pieces (per pack)	150
Special Feature	Accurate
Test Type	Blood Glucose Test
Specific Uses For Product	Blood, Blood Glucose Test, Glucose Test

8. WARRANTY AND SUPPORT

For general questions or technical support regarding McKesson Quintet AC Blood Glucose Test Strips, please use the contact information below:

General Questions: Call 1-800-777-4908

Technical Support: Call 1-888-481-8485

McKesson Medical-Surgical Inc. distributes this product. For detailed warranty information, please refer to the documentation included with your McKesson QUINTET AC Blood Glucose Monitoring System or contact McKesson directly.

© 2023 McKesson. All rights reserved. Information subject to change without notice.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question