



[Manuals.plus](#) /

› [IBEET](#) /

› [IBEET P10 Replacement Kit: Tubing with Swivel and Frame Instruction Manual](#)

IBEET P10

IBEET P10 Replacement Kit: Tubing with Swivel and Frame Instruction Manual

Model: P10

[Introduction](#) [Product Overview](#) [Setup & Specifications](#) [Support](#)
[Installation](#) [Operation](#) [Maintenance](#) [Troubleshooting](#)

1. INTRODUCTION

This manual provides instructions for the IBEET P10 Replacement Kit, which includes tubing with a swivel and a frame. This kit is designed to replace original components for the P10 system, ensuring proper function and comfort. Please read these instructions carefully before installation and use.

2. PRODUCT OVERVIEW

2.1 Components

- Tubing with 360-degree swivel
- Frame (headgear strap not included)

2.2 Key Features

- **Upgraded Vent Holes:** Designed to vent air from both sides, similar to original parts, to reduce noise and condensation.
- **Secure Fit:** Engineered to fit snugly with other system components, creating a reliable seal and minimizing hissing sounds.
- **Durable Tubing:** The 18.5-inch short tube is designed for enhanced sturdiness and durability, reducing the likelihood of punctures.
- **Comfort and Flexibility:** Maintains the original angle for secure placement and features a 360-degree swivel on the short tube for hose pivot.

DIRECT REPLACEMENT

Our product is designed to be more similar to the original structure, allowing you to sleep more quietly and comfortably.



Gentle exhaust



Snug fits without leaks



Ultra-flexible material



Figure 1: Overview of the IBEET P10 Replacement Kit components, highlighting the frame, tubing, and gentle exhaust design.

360° ROTATION

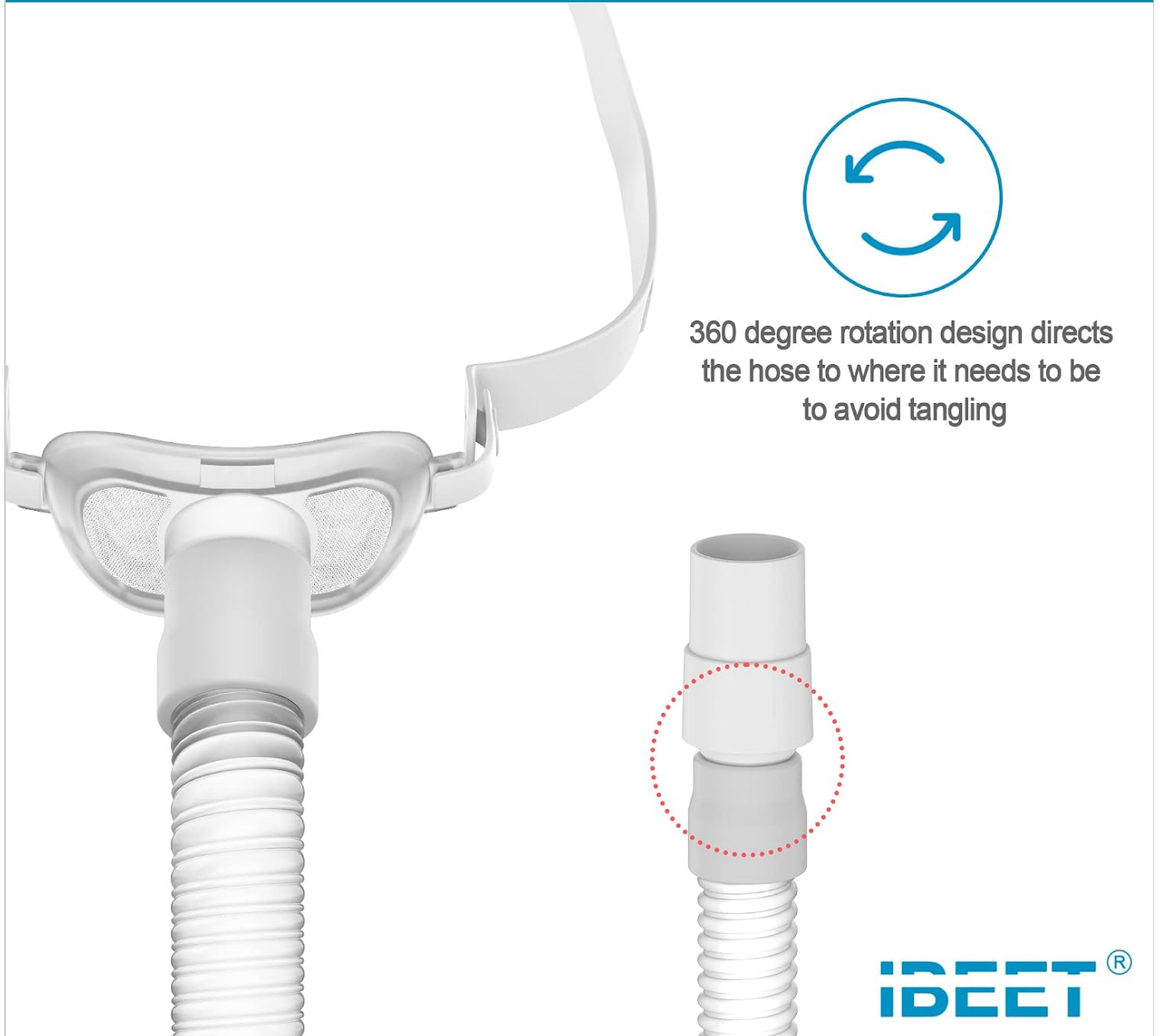


Figure 2: Detail of the 360-degree rotation feature on the tubing, designed to prevent tangling and allow flexible hose positioning.

3. SETUP & INSTALLATION

Follow these steps to properly install your IBEET P10 Replacement Kit.

3.1 Connecting the Tubing

1. Ensure the short tubing is connected to a standard 22mm hose. This connection must be made first before attaching to the device.
2. Once the short tubing is securely connected to the 22mm hose, connect the other end of the 22mm hose to your P10 device.

Suitable for various 22 mm diameter hoses

IDEET®

NOTE:
The short tube should connect with 22mm hose first.



Figure 3: Illustration of connecting the short tube to a 22mm hose, a crucial first step in the setup process.

Correct Use for Reliable Performance

Should connect to 22mm hose before connecting to the device.



Figure 4: Visual instructions for the correct connection sequence: short tube to 22mm hose, then 22mm hose to the device.

3.2 Attaching the Frame

The frame is designed to integrate with your existing P10 nasal pillow and headgear strap (not included in this kit).

1. Align the clips of your nasal pillow with the corresponding slots on the replacement frame and press firmly until they click into place.
2. Secure your headgear strap at the frame nodes. Ensure the headgear slots are properly secured for a stable fit.

Perfect for P10 Nasal Pillow & Strap (Not Included)



Fit clips into corresponding slots.



Secure strap at frame nodes.

Figure 5: Detailed steps for assembling the nasal pillow and headgear strap onto the replacement frame.

3.3 Installation Video Guide

Your browser does not support the video tag.

Video 1: This video demonstrates the installation process for the P10 replacement supplies, including how to connect the components.

4. OPERATION

Once installed, the IBEET P10 Replacement Kit operates in conjunction with your existing P10 system. Ensure all connections are secure before initiating therapy. The 360-degree swivel on the short tube allows for natural movement and helps prevent hose tangling during use.

5. MAINTENANCE

Regular cleaning and inspection of your replacement kit components are essential for hygiene and optimal performance.

- **Daily Cleaning:** Disconnect the tubing and frame from your device and nasal pillow. Wash the tubing and frame with mild soap and warm water. Rinse thoroughly and allow to air dry completely away from direct sunlight.
- **Weekly Inspection:** Inspect the tubing and frame for any signs of wear, tears, or damage. Replace components if any damage is observed to maintain an effective seal and prevent air leaks.
- **Avoid Harsh Chemicals:** Do not use bleach, alcohol, or other harsh cleaning agents, as these can damage the materials.

6. TROUBLESHOOTING

Issue	Possible Cause	Solution
-------	----------------	----------

Issue	Possible Cause	Solution
Air Leakage / Hissing Sound	Loose connection between tubing and hose, or hose and device. Damage to tubing or frame.	Ensure all connections are secure. Inspect tubing and frame for cracks or tears; replace if damaged.
Tubing detaches easily	Improper connection or worn connectors.	Verify the tubing is fully inserted into the 22mm hose. If the issue persists, the connectors may be worn and require replacement.
Excessive Noise or Condensation	Vent holes are obstructed or damaged.	Check vent holes for blockages and ensure they are clear. If vent holes are damaged, replace the frame.

7. SPECIFICATIONS

- **Brand:** IBEET
- **Model:** P10 Replacement Kit
- **Material:** Transparent
- **Product Dimensions:** 18.5 inches (Length) x 18.5 inches (Width) -*Note: This likely refers to the overall packaging or a flat measurement, the tubing length is 18.5 inches.*
- **Item Weight:** 0.07 Kilograms
- **Outside Diameter (Tubing Connection):** 22 Millimeters
- **Components Included:** Tubing with Swivel, Frame (Headgear strap and nasal pillow not included)

8. WARRANTY AND SUPPORT

IBEET products are designed for quality and reliability. While specific warranty details are not provided in this manual, if you encounter any concerns or issues with your IBEET P10 Replacement Kit after use, please contact the seller or manufacturer for assistance. Your satisfaction is important, and support is available to address product-related inquiries. For further assistance, please refer to the contact information provided with your purchase or visit the official IBEET store online.

