BlessReach WB-2106

BlessReach Foldable Rollator Walker with Seat (Model WB-2106)

User Manual

1. Introduction

Thank you for choosing the BlessReach Foldable Rollator Walker with Seat, Model WB-2106. This heavy-duty aluminum rolling walker is designed to provide reliable mobility support with enhanced safety and comfort. It features an impressive 450 lb weight capacity, large 8-inch wheels for various terrains, an 18-inch wide seat, and a 6-position adjustable handle height to suit users from 5.5' to 6.5' tall. Please read this manual thoroughly before assembly and use to ensure safe and proper operation.

2. Safety Information

- **Weight Capacity:** Do not exceed the maximum weight capacity of 450 lbs (204 kg).
- **Brakes:** Always engage the parking brakes before sitting on the seat or when the walker is stationary. Regularly check brake functionality.
- **Stability:** Ensure all parts are securely locked in place before each use. Do not use on uneven or slippery surfaces without extreme caution.
- **Adjustments:** Ensure all height adjustment mechanisms are securely locked after making changes.
- **Folding:** Exercise caution when folding or unfolding the walker to prevent pinching.
- **Supervision:** Users with cognitive impairments or balance issues should be supervised.

3. Product Components

Familiarize yourself with the main components of your BlessReach Rollator Walker:



Figure 3.1: Key Components of the BlessReach Rollator Walker

1. **Handle:** Ergonomic grips for comfortable pushing.
2. **Hand Brake:** Controls speed and locks the wheels.
3. **Back Strap:** Provides back support when seated.
4. **Height Adjustment Lever:** For adjusting handle height.
5. **Seat:** Padded seat for resting.
6. **Folding Release Button:** Mechanism for folding the walker.
7. **Storage Bag:** Under-seat storage for personal items.
8. **Swivel Wheel:** Front wheels for maneuverability.
9. **Brake Adjustment Screw:** For fine-tuning brake tension.

4. Setup and Assembly

The BlessReach Rollator Walker requires minimal assembly. Follow these steps:

1. **Unfold the Frame:** Carefully pull the sides of the walker apart until the frame locks into its open position. Ensure the folding mechanism is fully engaged.
2. **Install Handles:** Insert the handle tubes into the main frame. Align the holes and ensure the locking pins snap into place. Adjust to your desired height (see Section 5.3).
3. **Attach Backrest:** Secure the backrest to the designated slots on the handles.
4. **Check Brakes:** Verify that the brakes are functioning correctly after assembly.

For a visual guide, please watch the assembly video:

Your browser does not support the video tag.

Video 4.1: BlessReach Rollator Walker Assembly Instructions. This video demonstrates the step-by-step process of assembling the rollator walker, including unfolding the frame, installing handles, and attaching the backrest.



Figure 4.1: Foldable and Adjustable Features with Dimensions

5. Operating Instructions

5.1 Using the Walker for Mobility

- Stand behind the walker with both hands on the handles.
- Push the walker forward a comfortable distance.
- Step forward with one foot, then the other, maintaining balance.
- Keep your body upright and centered between the handles.

5.2 Brake Operation

The rollator features a reliable braking system:



Figure 5.1: Reliable Braking System

- **Lift-to-Slow (Resistance Brake):** Gently pull the brake levers upwards (labeled 'R' in the diagram) to apply resistance and control your speed, especially when going downhill.
- **Press-to-Lock (Parking Brake):** Push the brake levers downwards (labeled 'P' in the diagram) until they click into place to engage the parking brake. This locks the wheels, preventing the walker from rolling when you sit or stand. To release, pull the levers upwards.

5.3 Height Adjustment

The handles can be adjusted to 6 different positions to accommodate users from 5.5' to 6.5' tall.

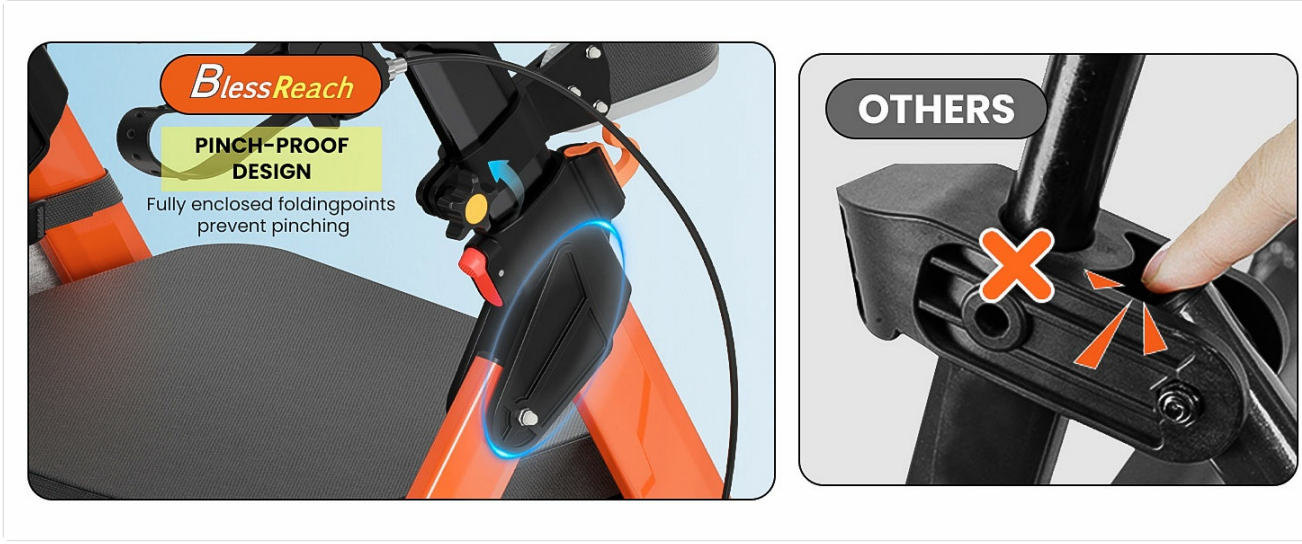


Figure 5.2: Dual-Adjustable Height

- Locate the height adjustment lever on each handle support.
- Press the lever or button to release the locking pin.
- Slide the handle up or down to the desired height. Ensure your wrists are straight and elbows are slightly bent when holding the handles.
- Release the lever/button and ensure the locking pin fully engages in one of the adjustment holes. Repeat for the other handle.

5.4 Using the Seat

The rollator features an 18-inch wide padded seat and a breathable backrest for comfortable resting.



Figure 5.3: Optimal Comfort Features

- Ensure the walker is on a flat, stable surface.
- Engage both parking brakes (press down on the brake levers).
- Carefully turn around and sit on the seat. Use the handles for support.
- The back strap provides support while seated.
- To stand up, push off the handles and stand slowly. Release the parking brakes before moving.

5.5 Folding for Storage and Transport

The rollator features an effortless 3-second folding mechanism for portability.



Figure 5.4: Convenient Folding and Storage

- Ensure the storage bag is empty.
- Lift the seat upwards using the designated strap or handle.
- Locate the folding release button (refer to Figure 3.1). Press the button and push the sides of the walker together until it is fully folded.
- The compact unit can be stored in car trunks or small spaces.

Watch this video for a demonstration of the rollator in use:

Your browser does not support the video tag.

Video 5.1: Heavy Duty Rollator Walker in Use. This video showcases the BlessReach rollator walker's features and ease of use for seniors in various environments.

6. Maintenance

6.1 Cleaning

- Wipe down the frame with a damp cloth and mild detergent.
- Do not use abrasive cleaners or solvents.
- Ensure all parts are dry before storage or next use.

6.2 Brake Adjustment

If the brakes feel too loose or too tight, they can be adjusted.

Your browser does not support the video tag.

Video 6.1: Brake Adjustment Instructions. This video provides detailed instructions on how to adjust the brakes on your rollator walker for optimal performance and safety.

- Locate the brake adjustment screw near the wheel.
- Turn the small nut to the right to tighten the brake cable, increasing brake tension.
- Turn the small nut to the left to loosen the brake cable, decreasing brake tension.
- Test the brakes after adjustment to ensure they engage and release properly.

7. Troubleshooting

- **Wheels not rolling smoothly:** Check for debris in the wheels. Ensure all wheel bolts are secure.
- **Brakes not engaging/releasing:** Refer to Section 6.2 for brake adjustment. Ensure brake cables are not kinked or damaged.

- **Walker feels unstable:** Verify that all locking pins for height adjustment and folding mechanisms are fully engaged. Ensure the frame is completely unfolded and locked.
- **Difficulty folding:** Ensure the seat is lifted and the folding release mechanism is fully activated.

8. Specifications

Feature	Detail
Model Number	WB-2106
Weight Capacity	450 lbs (204 kg)
Wheel Size	8 inches (puncture-proof)
Seat Dimensions	18 inches wide
Adjustable Handle Height	6 positions, for users 5.5' to 6.5'
Product Dimensions (Unfolded)	27.5"D x 23.6"W x 41"H
Minimum Height	33 inches
Maximum Height	40 inches
Product Weight	23.3 lbs (10.5 kg)
Material	Heavy-Duty Aluminum Frame
Color	Grizzled
Foldable	Yes



Figure 8.1: Enhanced Durability Features

9. Warranty and Support

BlessReach stands behind the quality of its products. For specific warranty information, please refer to the warranty card included with your purchase or visit the official BlessReach website. If you have any questions, require assistance with assembly, operation, or maintenance, or need to report an issue, please contact BlessReach

customer support. Contact details can typically be found on the product packaging or the BlessReach website.

Ask a question about this manual

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. Privacy details.

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question