

**Manuals+**

[Q & A](#) | [Deep Search](#) | [Upload](#)

**Manuals.plus /**

- › **Osim /**
- › OSIM uSqueeze 2 Smart Leg Massager Instruction Manual (Model OS-393S)

**Osim OS-393S**

# OSIM uSqueeze 2 Smart Leg Massager Instruction Manual

Model: OS-393S

## 1. INTRODUCTION

---

This manual provides essential information for the safe and effective use of your OSIM uSqueeze 2 Smart Leg Massager. Please read it thoroughly before operating the device and retain it for future reference. The OSIM uSqueeze 2 Smart is designed to provide a comprehensive leg massage experience, targeting feet, calves, and thighs with various massage techniques.



*Image 1.1: The OSIM uSqueeze 2 Smart Leg Massager, showcasing its ergonomic design and dark grey finish with purple accents.*

## **2. SAFETY INFORMATION**

To ensure safe operation and prevent injury or damage, please observe the following safety precautions:

- Always unplug the massager from the power outlet immediately after use and before cleaning.
- Do not use the massager near water or in damp environments.
- Keep the power cord away from heated surfaces.
- Do not operate the massager if it has a damaged cord or plug, or if it is not working properly.
- Consult a physician before using if you have medical conditions such as heart disease, pregnancy, diabetes, or if you use a pacemaker.
- Do not use the massager on swollen, inflamed, or injured areas.
- This appliance is not intended for use by persons with reduced physical, sensory, or mental capabilities unless supervised or instructed.
- Ensure the voltage specified on the product matches your local power supply.

### 3. PRODUCT FEATURES

The OSIM uSqueez 2 Smart Leg Massager incorporates advanced features for a comprehensive massage experience:

#### 3.1 Massage Techniques

- **Kneading Plates:** Emulate the hands of a professional masseuse for deep tissue relief.
- **Acupressure Massage Nodes:** Target specific acupressure points on the feet, calves, and thighs.
- **Reflex Rollers:** Provide stimulating reflexology massage to the soles of the feet.
- **Warm Air:** Integrated heat therapy to enhance relaxation and circulation.
- **Vibration:** Gentle vibrations to soothe muscles and improve comfort.



Image 3.1: Diagram illustrating the five core massage techniques: Kneading Plates, Acupressure Massage Nodes, Reflex Rollers, Warm Air, and Vibration.

#### 3.2 Auto Massage Programs

The device offers 6 pre-programmed massage routines tailored for different needs:

- **Reflexology:** Focuses on stimulating reflex zones in the feet.
- **Calves:** Targets calf muscles for relief and relaxation.
- **Thighs:** Provides massage to the thigh area.
- **Sleep:** A gentle program designed to promote relaxation before sleep.
- **Circulation:** A program aimed at improving blood flow.
- **Deep Comfort:** A comprehensive program for overall leg comfort.



Image 3.2: The control panel of the uSqueeze 2 Smart, displaying options for vibration, warm air, kneading, and the six auto massage programs.

### 3.3 Adjustable Design

- **158° Adjustment:** The kneading plates can be adjusted up to 158 degrees to target feet, calves, and thighs effectively.
- **Wide & Narrow Modes:** The kneading plates can be adjusted to accommodate different leg sizes and allow for a stronger massage experience in Narrow mode.

## 5 positions for full leg coverage



### Up to 158° adjustment to target calves and legs

uSqueeze 2 Smart features adjustable kneading plates with 5 positions that expand the massage coverage up to 158°. It can be adjusted between wide and narrow to fit different leg sizes for the most optimal massage enjoyment.

Image 3.3: Visual representation of the 5 adjustable positions for full leg coverage and the Wide/Narrow modes for the kneading plates.

### 3.4 Additional Features

- **Sleek Design:** Encased with textured leatherette, the uSqueeze 2 features a modern aesthetic.
- **Washable Cloth:** The fabric covering the base massager can be easily removed and washed for hygiene.

## Comfortable leg massage everywhere



Office



Living Room



Bedroom

### Additional Features



#### Sleek design

Encased with textured leatherette, the uSqueeze 2 features a sleek design and silhouette that will fit well in modern homes.



#### Washable cloth for hygiene

The fabric covering the base massager can be easily removed and washed for better hygiene.

Image 3.4: The OSIM uSqueeze 2 Smart shown in different environments (office, living room, bedroom) and close-ups of its sleek design and washable cloth feature.

## 4. SETUP

Follow these steps to set up your OSIM uSqueeze 2 Smart Leg Massager:

1. **Unpacking:** Carefully remove the massager from its packaging. Ensure all components, including the user manual, are present.
2. **Placement:** Place the massager on a flat, stable surface. Ensure there is enough space around the unit for comfortable use and proper ventilation.
3. **Power Connection:** Connect the power cord to the massager's power inlet, then plug it into a suitable 220-240V ~ 50-60Hz electrical outlet.
4. **Initial Inspection:** Before first use, inspect the massager for any visible damage.

## 5. OPERATING INSTRUCTIONS

To begin using your OSIM uSqueez 2 Smart Leg Massager:

1. **Positioning:** Sit comfortably in a chair or sofa. Place your feet, calves, or thighs into the massage compartments. Adjust the angle of the massager to target the desired area using the adjustable plates.
2. **Power On:** Press the power button on the control panel to turn on the massager.
3. **Select Program:** Choose one of the 6 auto massage programs (Reflexology, Calves, Thighs, Sleep, Circulation, Deep Comfort) by pressing the corresponding button.
4. **Adjust Settings:**
  - **Intensity:** Adjust the massage intensity as desired.
  - **Heat Therapy:** Activate or deactivate the warm air function.
  - **Kneading Width:** Select between Wide or Narrow mode for the kneading plates based on your leg size and preferred massage strength.
  - **Vibration:** Turn the vibration function on or off.
5. **Duration:** The massager has an auto-timer of 15 minutes. It will automatically shut off after this period. You can manually stop the massage at any time by pressing the power button.
6. **Power Off:** After your massage session, press the power button to turn off the unit and unplug it from the power outlet.



## uSqueez 2

The strongest reflexology massage

FOR DETOX & GOOD HEALTH



**AsiaNo.1**

Image 5.1: A user comfortably seated on a sofa, utilizing the OSIM uSqueez 2 Smart Leg Massager for a full leg massage.

## 6. MAINTENANCE

---

Proper maintenance ensures the longevity and hygiene of your massager:

- **Cleaning:** Always unplug the massager before cleaning. Wipe the exterior with a soft, damp cloth. Do not use abrasive cleaners or solvents.
- **Washable Cloth:** The fabric covers inside the massage compartments can be removed and hand-washed with mild detergent. Allow them to air dry completely before reattaching.
- **Storage:** Store the massager in a cool, dry place, away from direct sunlight and extreme temperatures. Avoid placing heavy objects on top of the unit.
- **Cord Care:** Do not wrap the power cord tightly around the unit. Avoid twisting or bending the cord excessively.

## 7. TROUBLESHOOTING

---

If you encounter any issues with your OSIM uSqueez 2 Smart Leg Massager, refer to the following common problems and solutions:

- **Massager does not turn on:**
  - Ensure the power cord is securely plugged into both the massager and a working electrical outlet.
  - Check if the power outlet is functional by plugging in another appliance.
- **Massage feels weak or inconsistent:**
  - Adjust the kneading plates to 'Narrow' mode for a stronger massage.
  - Ensure your legs are properly positioned within the compartments.
  - Select a different auto program or adjust intensity settings.
- **Unusual noises during operation:**
  - Turn off the massager and check for any obstructions in the massage mechanisms.
  - If the noise persists, discontinue use and contact customer support.
- **Heat function not working:**
  - Ensure the warm air function is activated on the control panel.
  - Allow a few minutes for the heat to build up.

If the problem persists after attempting these solutions, please contact OSIM customer support.


## 8. SPECIFICATIONS

---

Feature	Detail
Model Name	uSqueez 2 Smart
Model Number	OS-393S

<b>Feature</b>	<b>Detail</b>
Voltage	220-240V ~ 50-60Hz
Power Consumption	60W
Auto Timer	15 minutes
Vibration Speed	Approx. 1750 - 2600 rpm
Kneading Speed	Approx. 4 - 8 rpm
Warm Air Temperature	Approx. 38 - 45°C
Product Dimensions (LxWxH)	51.5 x 44.5 x 47 cm
Item Weight	10 kg
Material	Acrylonitrile Butadiene Styrene
Included Components	uSqueez2 Smart Foot Massager, User Manual

## PRODUCT SPECIFICATIONS

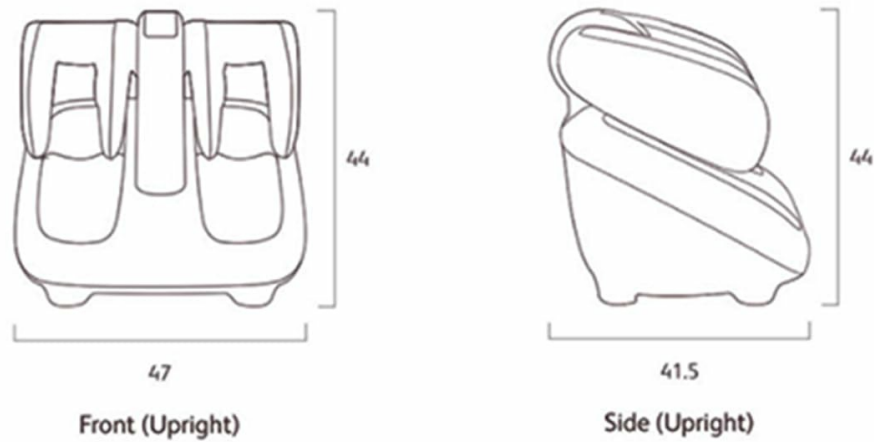

<b>Model</b> uSqueez 2 Smart	<b>Power consumption</b> 60W	<b>Kneading speed</b> Approx. 4 - 8 rpm	<b>Warm air temperature</b> Approx. 38 - 45°C
<b>Model number</b> OS-393S	<b>Auto timer</b> 15 minutes	<b>Vibration speed</b> Approx. 1750 - 2600 rpm	<b>Weight</b> Gross weight - Approx. 11.5kg Net weight - Approx. 10kg
<b>Voltage</b> 220-240V - 50/60Hz			

Image 8.1: Diagram illustrating the dimensions and key specifications of the OSIM uSqueez 2 Smart Leg Massager.

## 9. WARRANTY AND SUPPORT

Your OSIM uSqueez 2 Smart Leg Massager comes with a **1-Year Warranty** from the OSIM brand.

For warranty claims, technical support, or any inquiries, please contact:

- **Manufacturer:** Osim International - Health care Company, Headquarters: 65 Ubi Avenue 1, 408939, Singapore
- **Importer (India):** Osim India (A Division Of Paramount Surgimed Limited), No. 1, L. S. C., Okhla Industrial Area, Phase- 2, New Delhi – 110020, Delhi
- **Contact Number:** 09311660866 (for India inquiries)

Please retain your proof of purchase for warranty validation.

