

[Manuals.plus](#) /

› [yuwell](#) /

› yuwell Blood Glucose Test Strips and Lancets User Manual

yuwell test strip set 50pcs(model 710&582 only)

yuwell Blood Glucose Test Strips and Lancets User Manual

For Model 710 & 582 Blood Glucose Meters

PRODUCT OVERVIEW

This manual provides essential information for the proper use and care of yuwell Blood Glucose Test Strips and Lancets, designed exclusively for use with yuwell Blood Glucose Meters Models 710 and 582. These components are crucial for accurate blood glucose monitoring.



Image: yuwell Blood Glucose Test Strips and Lancets.

COMPATIBILITY

The yuwell Blood Glucose Test Strips and Lancets are specifically engineered to work only with yuwell Blood Glucose Meters Model 710 and Model 582. Using these strips and lancets with other branded meters or different yuwell models may lead to inaccurate results or device malfunction.

Only Compatible With: yuwell blood glucose meter 710 & 582



yuwell blood glucose meter only accepts yuwell branded test strips
yuwell test strips will NOT work with other branded blood glucose meters
The blood glucose meter only accepts test strips of the same model from the yuwell brand

Image: Compatibility of yuwell Test Strips with specific yuwell Blood Glucose Meters (Models 710 & 582).

SETUP AND PREPARATION

Before performing a blood glucose test, ensure you have all necessary components ready and follow these steps:

1. Wash your hands thoroughly with soap and warm water. Dry them completely.
2. Gather your yuwell Blood Glucose Meter (Model 710 or 582), a yuwell test strip, and a yuwell lancet.
3. Prepare the lancing device according to its specific instructions, inserting a fresh lancet.
4. Ensure the test strip bottle is properly sealed when not in use to maintain strip integrity.

OPERATING INSTRUCTIONS

Follow these steps for accurate blood glucose measurement:

Instructions



① Insert a test strip into the meter



② Use yuwell lancing device to obtain a drop of blood



③ Touch the drop of blood to the the test strip



④ The result will appear on the display

TIPS FOR USE (In order to make the test results more accurate)

1. After inserting the test strip, the monitor will turn on automatically
2. Do not drop blood on the test strip before inserting it, otherwise code 306 will appear and measurement will not be possible
3. After the test strip bottle is opened, the test strip inside needs to be used up within 6 months

Image: Step-by-step instructions for performing a blood glucose test.

1. **Insert Test Strip:** Insert a fresh yuwell test strip into the test strip port of your yuwell blood glucose meter. The meter will automatically turn on.
2. **Obtain Blood Sample:** Use the yuwell lancing device to obtain a small drop of blood from your fingertip.
3. **Apply Blood to Strip:** Gently touch the tip of the test strip to the drop of blood. The strip will draw the blood in automatically.
4. **Read Result:** The blood glucose result will appear on the meter's display within a few seconds.

Important Tips for Accurate Results:

- Do not drop blood on the test strip before inserting it into the meter. This may cause an error code (e.g., "CODE 306") and prevent measurement.
- Once a test strip bottle is opened, the test strips inside should be used within 6 months to ensure accuracy.
- The yuwell blood glucose monitor recognizes the encrypted batch code on each test strip, eliminating the need for manual code entry.

- Calibration liquid is not required for yuwell blood glucose kits.

Some tips for using yuwell glucometer 710



After inserting the test strip, the monitor will turn on automatically. The screen will show "CODE 306", and then the symbol "BLOOD" will flashing which means the monitor is ready to test.



The monitor could recognize the encrypted batch code on each test strip, so you don't need to manually insert any codes for each new batch. To save power, yuwell blood glucose monitor will keep closed if there is no operation.



The liquid for calibration is not a necessity for a blood glucose kit. Without it will not affect the normal use of the meter.

Image: Additional tips for using the yuwell glucometer 710 for optimal performance.

MAINTENANCE AND STORAGE

- Store test strips and lancets in their original sealed containers in a cool, dry place, away from direct sunlight and heat.
- Do not store test strips outside their original bottle.
- After opening, use test strips within the validity period indicated on the bottle, typically 6 months.
- Dispose of used test strips and lancets properly in a sharps container or according to local regulations for biohazardous waste.
- Keep the meter clean by wiping it with a soft, damp cloth. Do not immerse the meter in water or use harsh cleaning agents.

TROUBLESHOOTING

If you encounter issues during testing, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Meter displays "CODE 306" or similar error.	Blood applied to strip before insertion into meter, or strip is damaged/expired.	Ensure strip is inserted first. Use a new, unexpired strip.
No blood drawn into strip.	Insufficient blood sample, or strip not touched correctly.	Ensure a sufficient blood drop is formed. Touch the tip of the strip to the blood drop.
Inaccurate readings.	Expired strips, improper storage, or incorrect testing procedure.	Check strip expiry date. Store strips correctly. Review operating instructions.

If problems persist, contact yuwell customer support.

SPECIFICATIONS

- **Product Name:** yuwell Blood Glucose Test Strips and Lancets
- **Compatible Meters:** yuwell Blood Glucose Meter Model 710, yuwell Blood Glucose Meter Model 582
- **Test Strip Technology:** Updated Five-electrode Technology
- **Package Contents:** 50pcs Test Strips and Lancets
- **Item Model Number:** test strip set 50pcs(model 710&582 only)
- **Package Dimensions:** 4.88 x 2.6 x 1.73 inches
- **Package Weight:** 2.82 ounces

Updated Five-electrode Technology

Double corrected electrodes,
multiple anti-interference, more accurate measurements

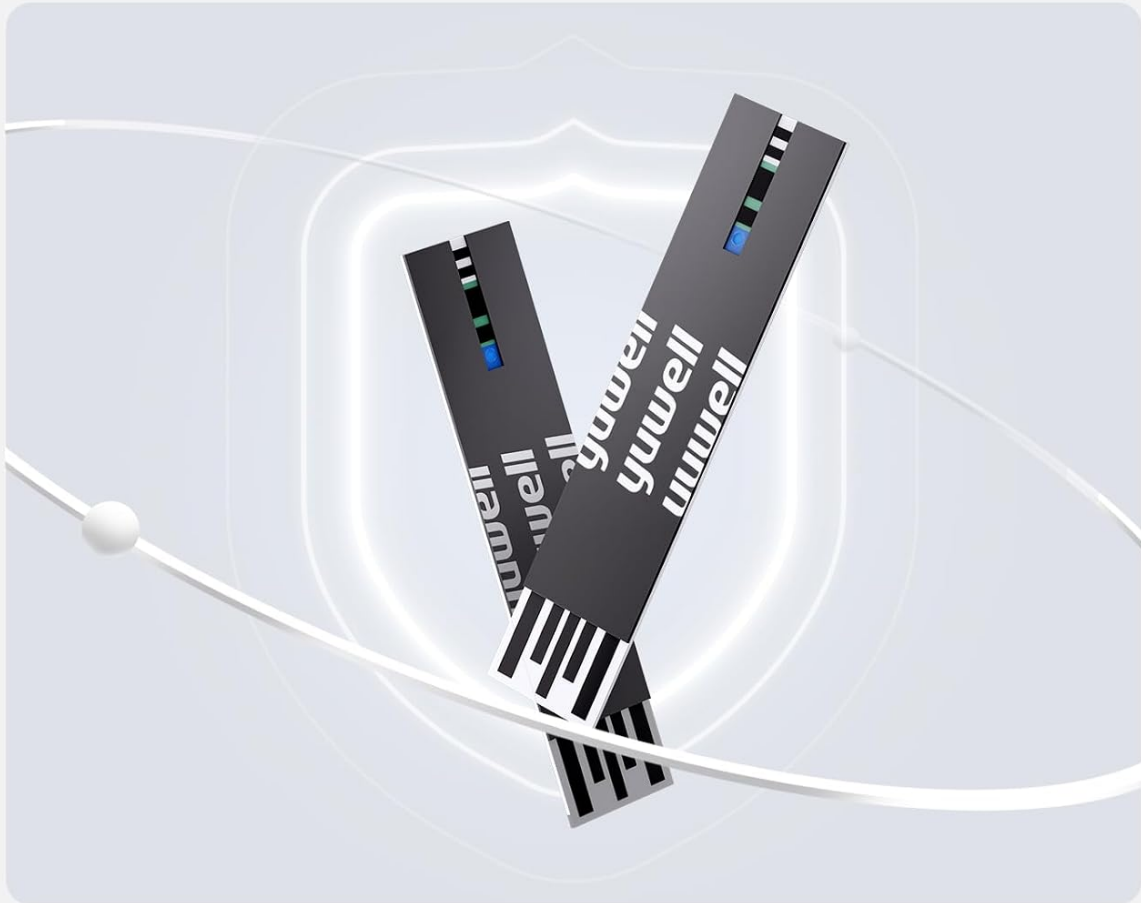


Image: Illustration of the Updated Five-electrode Technology in yuwell test strips.

IMPORTANT SAFETY INFORMATION

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

Always consult with a healthcare professional for diagnosis and treatment of any health condition.

WARRANTY AND SUPPORT

For warranty information or technical support regarding your yuwell Blood Glucose Test Strips, Lancets, or Meter, please refer to the documentation included with your yuwell Blood Glucose Meter or visit the official yuwell website.

You can also visit the [yuwell Store on Amazon](#) for more product information and support resources.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question

Open the live page

<https://manuals.plus/asin/B0C13R2ZN5>

More documents and user manuals: <https://manuals.plus/>