

[Manuals.plus](#) /

› [Sawgmore](#) /

› Sawgmore TH2D4 Water Flosser Replacement Tips User Manual

Sawgmore TH2D4

Sawgmore TH2D4 Water Flosser Replacement Tips User Manual

Model: TH2D4 Compatible Tips

INTRODUCTION

This manual provides instructions for the proper use and maintenance of your Sawgmore TH2D4 Water Flosser Replacement Tips. These tips are designed to be compatible with the Sawgmore TH2D4 Oral Irrigator, enhancing your oral hygiene routine with specialized cleaning functions.



Image: Overview of the Sawgmore TH2D4 compatible replacement tips.

SAFETY INFORMATION

Please read all instructions before using the replacement tips. Improper use may lead to injury or damage to the product. These

tips are for oral hygiene use only. Consult your dentist before use if you have any oral health concerns.

- Do not share tips with others for hygiene reasons.
- Replace tips every 3 months, or sooner if they show signs of wear.
- Keep out of reach of children.

PRODUCT OVERVIEW: TIP TYPES AND FUNCTIONS

The Sawgmore TH2D4 replacement tip set includes various specialized tips designed for comprehensive oral care.

6 NOZZLES With Different Functions



CLASSIC JET TIP*2
General cleaning




TONGUE CLEANER
Removes tongue coating buildup for fresher breath



ORTHODONTIC TIP
Remove food debris in the brace



MASSAGE TIP
Covers more area, massage gum and prevent gum bleeding

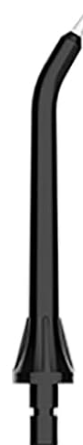

PERIODONTAL TIP
Rinses and massage deep below the gumline



Image: Detailed view of each tip type and its specific function.

- **Classic Jet Tip (x2):** Designed for general cleaning, effectively removing plaque and food debris between teeth and below the gumline.
- **Tongue Cleaner:** Features a spoon-like shape to gently remove bacteria and food particles from the tongue surface, promoting fresher breath.
- **Orthodontic Tip:** Specifically designed with a tapered brush to clean around braces, wires, and other orthodontic appliances.
- **Massage Tip:** Provides a gentle massage to the gums, which can help improve gum health and circulation.
- **Periodontal Tip:** Features a soft rubber tip for gentle, low-pressure delivery of water or therapeutic rinses into periodontal pockets and below the gumline.

SETUP AND INSTALLATION

1. **Select a Tip:** Choose the appropriate replacement tip for your desired cleaning function.
2. **Insert the Tip:** With the water flosser unit turned off, insert the stem of the chosen tip into the handle's opening until it clicks securely into place. Ensure it is fully seated to prevent leaks.
3. **Remove the Tip:** To remove a tip, press the tip eject button (if available on your flosser model) or firmly pull the tip straight out from the handle.

OPERATING INSTRUCTIONS

Once the desired tip is installed, fill your Sawgmore TH2D4 Oral Irrigator with water and select your preferred pressure setting. For optimal results, follow these general guidelines:

USING PULSING WATER FLOWS TO CLEAN THE TEETH

Use 1400 times per minutes frequency conversion motor, and the maximum water pressure can reach 90PSI.

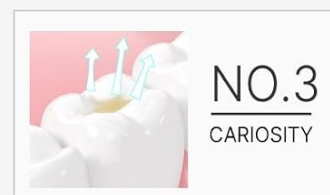
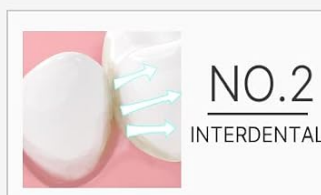
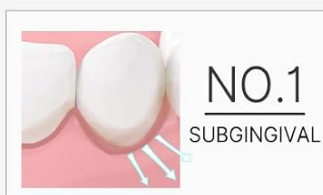


Image: Demonstrating the pulsating water flow for effective cleaning.

- **General Use:** Lean over a sink and place the tip in your mouth. Turn on the device. Direct the water stream at a 90-degree angle to your gumline, moving along the teeth.

- **Classic Jet Tip:** Use for daily cleaning, focusing on the spaces between teeth and the gumline.
- **Tongue Cleaner:** Set the pressure to low. Place the tip on the back of your tongue and move it forward, cleaning the entire surface.
- **Orthodontic Tip:** Carefully guide the tip along the gumline and around orthodontic brackets and wires.
- **Massage Tip:** Use a lower pressure setting and gently massage the gumline.
- **Periodontal Tip:** Use at the lowest pressure setting. Place the soft tip against the tooth at a 45-degree angle and gently insert it into the gum pocket.

MAINTENANCE AND CLEANING

Regular cleaning of your replacement tips is essential for hygiene and longevity.

1. **Daily Cleaning:** After each use, rinse the tip thoroughly under running water to remove any debris.
2. **Deep Cleaning:** Periodically, soak the tips in a solution of 2 parts water to 1 part white vinegar for 5-10 minutes to remove mineral deposits and sanitize. Rinse thoroughly afterward.
3. **Storage:** Store tips in a clean, dry place, away from direct sunlight.

TROUBLESHOOTING

Issue	Possible Cause	Solution
Tip not fitting securely	Incorrect tip for model, or not fully inserted.	Ensure the tip is compatible with Sawgmore TH2D4. Push firmly until it clicks.
Reduced water flow from tip	Mineral buildup or blockage.	Soak the tip in a vinegar solution as described in Maintenance.
Tip appears worn or damaged	Normal wear and tear.	Replace the tip with a new one.

SPECIFICATIONS

These replacement tips are designed for the Sawgmore TH2D4 Portable Oral Irrigator.

- **Brand:** Sawgmore
- **Compatible Model:** TH2D4 Oral Irrigator
- **Tip Types:** Classic Jet, Tongue Cleaner, Orthodontic, Massage, Periodontal
- **UPC:** 744243953135
- **Manufacturer Reference:** SA-122 (for the irrigator unit)

Sawgmore®

Portable Oral Irrigator

Model:TH2D4



WARRANTY INFORMATION

Specific warranty details for these replacement tips are not provided in the product information. Please refer to the warranty information included with your Sawgmore TH2D4 Oral Irrigator unit or contact the manufacturer directly for details regarding accessories.

CUSTOMER SUPPORT

For further assistance or inquiries regarding your Sawgmore TH2D4 Water Flosser Replacement Tips, please refer to the contact information provided with your original Sawgmore Oral Irrigator purchase or visit the official Sawgmore website.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Email

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Details

Submit question