

[Manuals.plus](#) /

› [yeghhbm](#) /

› yeghhbm M607 Rechargeable Hearing Aids User Manual

## yeghhbm M607

# yeghhbm M607 Rechargeable Hearing Aids User Manual

Model: M607 | Brand: yeghhbm

## 1. INTRODUCTION

This manual provides essential information for the proper use and maintenance of your yeghhbm M607 Rechargeable Hearing Aids. These devices are designed to be nearly invisible, offering an ergonomic fit and incorporating layered noise reduction technology for enhanced sound clarity. The portable charging case ensures convenience and extended use.

## 2. PACKAGE CONTENTS

Please verify that all items listed below are included in your package:

- Left and Right Hearing Aids
- Charging Case with LED Power Display
- USB-A to USB-C Power Cord
- Multiple Sizes of Earplugs (Ear Domes)
- Cleaning Brushes
- User Manual

# The Package Includes

Complete Accessories, Matching Models, The Best Choice Of Similar Products



Earplugs



Hearing Aid




Cleaning Brush



Power Cord



Charging Case

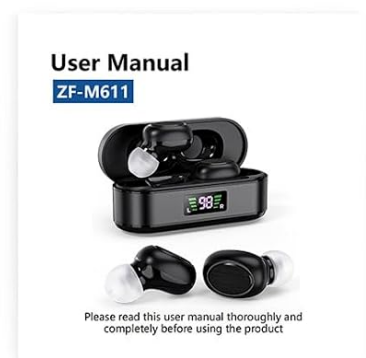

Instructions


**Image 2.1:** All items included in the yeghhbm M607 package: earplugs of various sizes, the left and right hearing aids, cleaning brushes, a power cord, the charging case, and the user manual.

## 3. PRODUCT OVERVIEW

Familiarize yourself with the components and design of your hearing aids.

# Structure & Function

Ergonomic curved body, better fit the ear canal

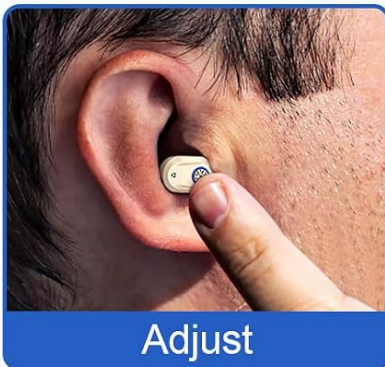

Adjust



Charge



Easy to carry

**Image 3.1:** This diagram illustrates the key components of the yeghhbm M607 hearing aid, including the ear dome, volume adjustment wheel, microphone, charging points, and speaker. It also shows how to adjust, charge, and carry the device.

## 3.1 Key Components

- **Ear Dome:** Soft silicone tip for comfortable fit in the ear canal.
- **Volume Adjustment:** Wheel for adjusting sound amplification levels.
- **Microphone:** Captures surrounding sounds.
- **Charging Point:** Connects to the charging case for power.
- **Speaker:** Delivers amplified sound into the ear.

## 3.2 Design and Dimensions

The M607 hearing aids feature an ergonomic curved body designed for a snug and comfortable fit within the ear canal, making them discreet and nearly invisible. Each hearing aid weighs approximately 2.2 grams.

## Unique Design — & Nearly Invisible

Ideal size to feel comfortable in the cana while large enough for easier operation



Natural Fitting  
**Comfortable** in the ear



**Weight**  
2.2g

**Size**  
0.49"\*0.95\*0.44

**Image 3.2:** This image highlights the compact and nearly invisible design of the yeghhbm M607 hearing aid, showing its natural fit in the ear canal. Dimensions for both the hearing aid (0.49"x0.95"x0.44", 2.2g) and the charging case (2.36"x2.36"/6cm) are provided.

## 4. SETUP

### 4.1 Initial Charging

Before first use, fully charge both the hearing aids and the charging case. Connect the provided USB-C cable to the charging port on the case and a suitable USB power source. The LED display on the case will indicate the charging status and battery percentage. A full charge typically takes 2-3 hours, providing over 20 hours of use for the hearing aids.

# Fast Charge & Larger capacity

Fully Charged

≤3hrs

For Hearing Aid  
And Charging Box

Working Time

≥20hrs

For Hearing Aid



Backup Power

≥100hrs

5 Times For Both Units 10 Times For Single Unit



**Image 4.1:** This graphic details the charging capabilities and battery life of the yeghhbm M607 hearing aids. It shows that the hearing aids and charging box fully charge in under 3 hours, providing over 20 hours of working time for the hearing aids and over 100 hours of backup power from the case.

## 4.2 Selecting and Fitting Ear Domes

Choose the ear dome size that provides the most comfortable and secure fit in your ear canal. A proper fit is crucial for optimal sound quality and to prevent feedback (whistling). Gently insert the hearing aid into your ear until it feels snug.

## 5. OPERATING INSTRUCTIONS

### 5.1 Power On/Off

- **To Power On:** Simply remove the hearing aids from the charging case. They will automatically power on.
- **To Power Off:** Place the hearing aids back into the charging case. They will automatically power off and begin charging.

## 5.2 Volume Adjustment

Use the small volume adjustment wheel located on the hearing aid to increase or decrease the sound amplification to your preferred level. Adjust gradually to avoid sudden loud sounds.

## 5.3 Noise Reduction

The M607 utilizes layered noise reduction and directional feedback canceling technology to process sound more naturally. While this technology aims to reduce background noise, it does not eliminate all ambient sounds, allowing for a more natural hearing experience.



# Simple & Effective

Before After

Noise Reduction

Reduce the background noises and process the sound as naturally as possible.  
Make perfect gifts for elderly family members with hearing loss!

**Image 5.1:** The image displays the yeghhbm M607 hearing aids in their charging case, alongside a graphic illustrating the effect of noise reduction technology. The text emphasizes the reduction of background noise for clearer sound.

## 6. CHARGING

The portable charging case is designed for convenient charging of your hearing aids. It is recommended to charge the

hearing aids and case once a day for optimal performance.

- **Charging the Hearing Aids:** Place the hearing aids into their respective slots in the charging case. They will automatically begin to charge. The LED display on the case will show the individual charge level for each hearing aid (L for left, R for right).
- **Charging the Case:** Connect the USB-C cable to the charging port on the bottom of the case and a power source. The central LED display will show the battery percentage of the charging case.



**Image 6.1:** The yeghhbm M607 hearing aids are shown with their portable charging case. Two hearing aids are placed inside the open case, while two others are displayed in front, highlighting their compact design.

## 7. MAINTENANCE

Regular cleaning and proper care will extend the lifespan and ensure optimal performance of your hearing aids.

- **Cleaning:** Use the provided cleaning brushes to gently remove any earwax or debris from the ear domes, microphone, and speaker openings. Do not use liquids or harsh chemicals.

- **Storage:** When not in use, always store the hearing aids in their charging case to protect them and keep them charged.
- **Avoid Moisture:** Keep the hearing aids and charging case dry. Avoid exposure to water, humidity, and extreme temperatures.

## 8. TROUBLESHOOTING

If you encounter issues with your yeghhbm M607 hearing aids, please refer to the following common solutions:

| Problem                           | Possible Cause                                                              | Solution                                                                                                                                                                                                                    |
|-----------------------------------|-----------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| No sound or low sound             | Low battery, improper fit, volume too low, ear dome blocked.                | Ensure hearing aids are charged. Check ear dome for proper fit and blockages. Adjust volume.                                                                                                                                |
| Whistling or feedback sound       | Improper ear dome size, hearing aid not inserted correctly, earwax buildup. | Try a different size ear dome for a tighter seal. Reinsert the hearing aid properly. Clean the ear dome and speaker opening. Turn off the hearing aid before inserting, then turn it on once in place to minimize feedback. |
| Charging case not charging        | Cable not connected properly, faulty cable or power source.                 | Ensure the USB-C cable is securely connected to both the case and the power source. Try a different cable or power adapter.                                                                                                 |
| Hearing aids not charging in case | Charging points obstructed, case battery low.                               | Clean the charging points on both the hearing aids and the case. Ensure the charging case itself is charged.                                                                                                                |

## 9. SPECIFICATIONS

- **Model Name:** M607
- **Brand:** yeghhbm
- **Power Source:** Rechargeable (2 Lithium Ion batteries included)
- **Special Features:** Adjustable, Lightweight, Noise Cancellation, Rechargeable
- **Item Model Number:** M607Brown
- **Package Dimensions:** 3.86 x 3.82 x 2.32 inches
- **Item Weight:** 5.61 ounces (package)
- **Manufacturer:** yeghhbm
- **Date First Available:** August 10, 2023

## 10. WARRANTY AND SUPPORT

If you encounter any problems when using your yeghhbm M607 hearing aids, please contact our customer support. We are committed to providing a satisfactory answer within 24 hours.  
For support, please refer to the contact information provided with your purchase or on the official yeghhbm website.

## ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

### Name

Anonymous

### Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

### Question

Example: How do I reset this device or fix this error?

### Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question