

INIA HAUTE

INIA Laser Hair Removal Device User Manual

Model: HAUTE (3023281259)

Brand: INIA

1. INTRODUCTION

The INIA HAUTE IPL Hair Removal Device is designed for safe and effective hair removal in the comfort of your home. Utilizing advanced Intense Pulsed Light (IPL) technology with a unique Double-Pulse feature and integrated ice-cooling, it aims to provide long-lasting hair reduction results.

This manual provides essential information on how to safely and effectively use your INIA HAUTE device, including setup, operation, maintenance, and troubleshooting.



Figure 1.1: Comparison of INIA HAUTE at-home treatment versus traditional beauty salon laser hair removal, highlighting convenience and cost-effectiveness.

2. SAFETY INFORMATION

Before using the INIA HAUTE device, please read all safety instructions carefully. This device is intended for hair removal. Do not use it for any other purpose.

2.1 Skin Tone and Hair Color Compatibility

The effectiveness and safety of IPL hair removal depend on your skin tone and hair color. This device is most effective on lighter skin tones and darker hair colors. It is not suitable for very dark skin tones or very light hair colors (white, grey, red, or very light blonde).

INIA

Skin Tone & Body Hair Type Chart

Please review the charts below to ensure your eligibility for safe & effective IPL use.

Suitable Skin Tones



Suitable Body Hair Colors



Figure 2.1: Skin Tone & Body Hair Type Chart. Please review this chart to ensure your eligibility for safe and effective IPL use. Darker skin tones and very light hair colors are generally not suitable.

2.2 General Safety Precautions

- Always wear the provided protective eyewear during treatment to protect your eyes from the intense light flashes.
- Do not use on tattooed areas, moles, birthmarks, or areas with skin conditions like eczema, psoriasis, or open wounds.
- Avoid direct exposure to the light window.
- Keep the device away from water and other liquids.
- Do not use if you are pregnant, breastfeeding, or have certain medical conditions (consult a doctor if unsure).
- Keep out of reach of children.

3. PACKAGE CONTENTS

Upon opening the package, verify that all items are present and in good condition:

- INIA HAUTE IPL Hair Removal Device
- Power Adapter
- Protective Eyewear
- Razor
- Storage Bag
- User Manual (this document)



Figure 3.1: INIA HAUTE device and standard accessories.

4. GETTING STARTED (SETUP)

Follow these steps to prepare your device for its first use:

1. **Prepare the Skin:** Shave the area you intend to treat. Ensure the skin is clean, dry, and free of any lotions, creams, or deodorants.
2. **Connect Power:** Plug the power adapter into the device and then into a power outlet. The device will power on automatically or indicate readiness.
3. **Wear Protective Eyewear:** Always put on the provided protective eyewear before beginning treatment to shield your eyes from the intense light.

5. OPERATING INSTRUCTIONS

The INIA HAUTE device offers efficient and comfortable hair removal with its advanced features.

5.1 Power On/Off and Intensity Levels

- **Power On:** The device typically powers on when connected to the power adapter.
- **Adjust Intensity:** Press the power button briefly to cycle through the 5 intensity levels. Start with the lowest intensity and gradually increase it based on your comfort and skin reaction.
- **Power Off:** Unplug the device from the power outlet to turn it off.

5.2 Treatment Modes

The device features two primary modes for effective treatment:

- **Double-Pulse Mode:** This mode emits two pulses within milliseconds, maximizing energy absorption by hair follicles for enhanced effectiveness. This is the default and recommended mode for optimal results.
- **Auto-Glide Mode:** For faster treatment of larger areas, activate Auto-Glide mode. In this mode, the device will automatically flash when in continuous contact with the skin, allowing for smooth gliding across the treatment area. The device can achieve up to 5 flashes per second in this mode.

INIA

Industry Record 5 Flashes per Second

Lightning Fast
and Incredibly Effective



**8 Minutes Full Body
Hair Removal**

Others

2 Flashes per Second	1.2 Flashes per Second	1.4 Flashes per Second	0.6 Flashes per Second

Figure 5.1: The device's rapid flashing capability for efficient full-body treatment.

5.3 Ultra Sapphire Ice-Cooling Function


The INIA HAUTE device incorporates an advanced ice-cooling feature to maintain a comfortable skin temperature during treatment, typically around 50°F (10°C). This helps minimize discomfort and protect the skin from heat injury.

- The ice-cooling function activates automatically when the device is in use.
- Ensure the treatment window is in full contact with the skin for optimal cooling effect.

INIA

Ultra Sapphire Ice-Cooling

Keep the temperature below 60°F for up to 20 minutes.

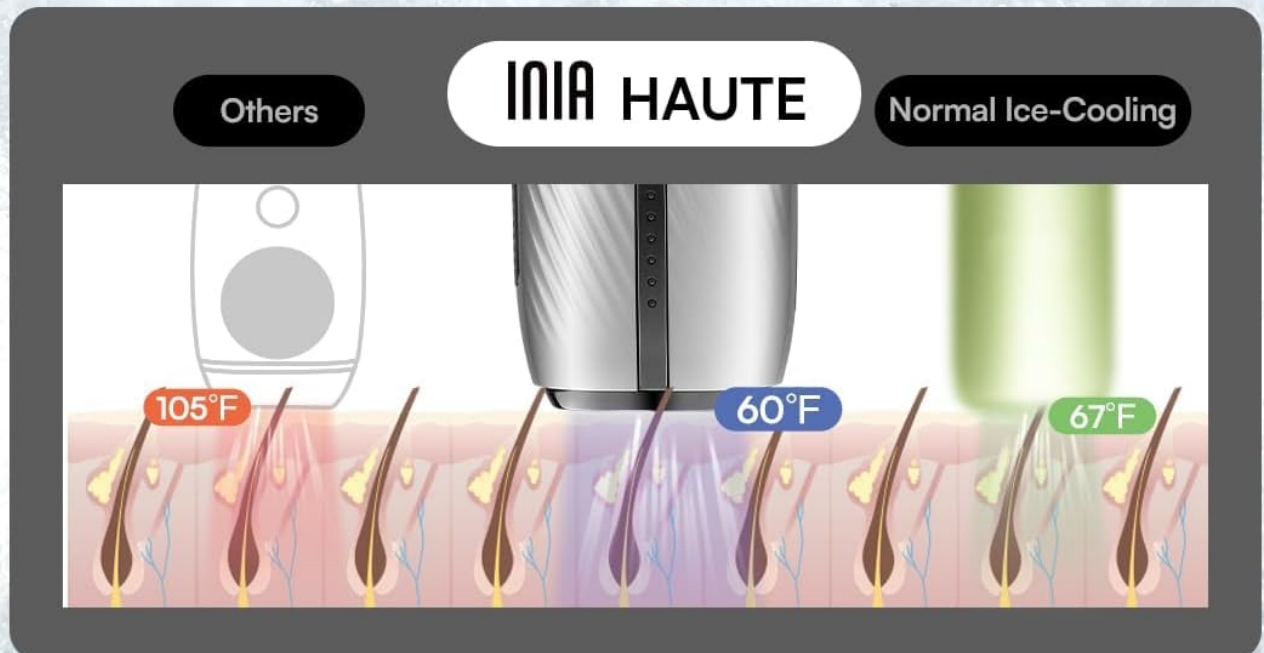

*High temperature burns skin, Low temperature cares for skin.

Figure 5.2: Illustration of the Ultra Sapphire Ice-Cooling feature, designed to keep the skin cool during treatment.

5.4 Application Technique

1. Place the treatment window firmly against the skin, ensuring full contact.
2. Press the flash button (or glide in Auto-Glide mode) to deliver a pulse.
3. Move the device to the next adjacent area, ensuring slight overlap to avoid missed spots.
4. Treat each area once per session.

INIA

At Home vs Beauty Salon

Freedom and Savings



	INIA HAUTE	Salon laser hair removal
Treatment Outcomes	Long lasting hair reduction	Long lasting hair reduction
Body Areas Safety	Head to toe	Head to toe
On-Site Requirement	FDA-cleared equipment	FDA-cleared equipment
Appointment Required	NO Need	Need
Costs for Tips	NO Need	Need
Cost per Body Area	NO Need	Need
Cost per Treatment	NO Need	Need
TOTAL COST	NO Need	Need
	\$199	\$4,000+

Figure 5.3: Proper handling and application of the device during treatment.

5.5 Treatment Schedule

For optimal results, follow a consistent treatment schedule. Initial results may be visible within 3 weeks, with significant hair reduction over time.

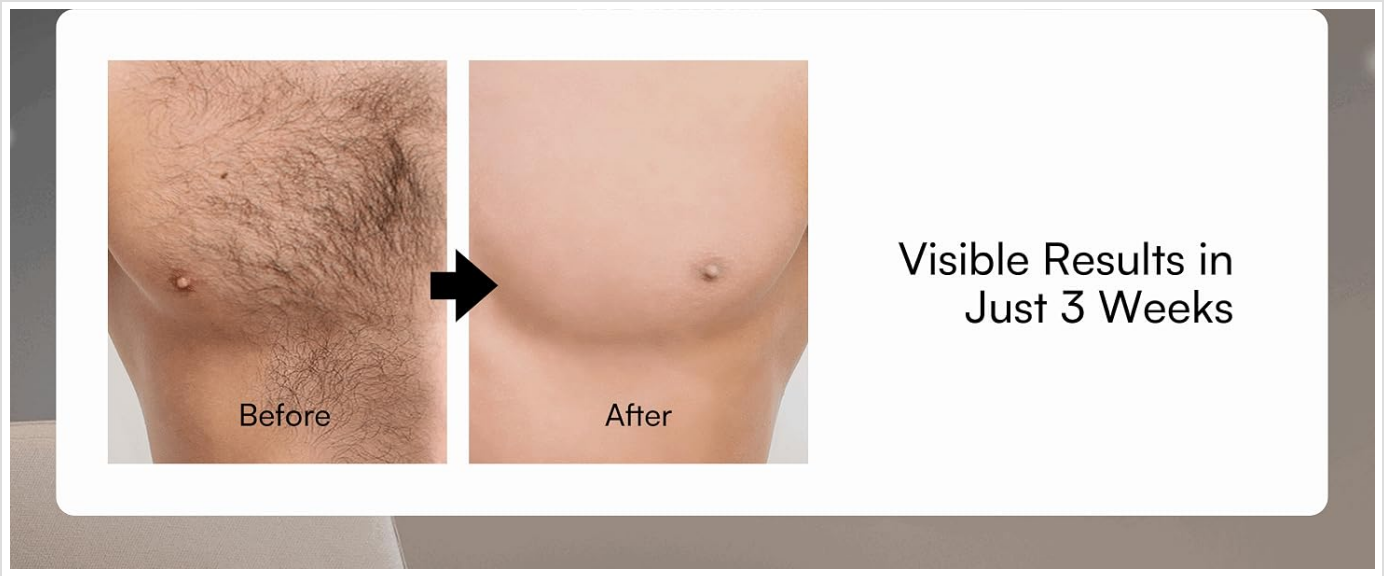


Figure 5.4: Example of visible results after consistent use.

6. CARE AND MAINTENANCE

Proper care and maintenance will ensure the longevity and optimal performance of your INIA HAUTE device.

6.1 Cleaning the Device

- Always unplug the device before cleaning.
- Use a soft, dry cloth to wipe the device body.
- For the light window, use a cotton swab or a soft cloth slightly dampened with alcohol to gently clean any residue. Ensure the window is completely dry before next use.
- Do not use abrasive cleaners or immerse the device in water.

6.2 Storage

- Store the device in its original packaging or the provided storage bag.
- Keep it in a cool, dry place, away from direct sunlight and moisture.
- Ensure the device is unplugged and cooled down before storing.

7. TROUBLESHOOTING

If you encounter any issues with your INIA HAUTE device, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Device does not power on.	Not properly plugged in; Power outlet issue.	Ensure the power adapter is securely connected to the device and the outlet. Try a different outlet.
Device does not flash.	Light window not in full contact with skin; Incorrect mode selected; Device malfunction.	Ensure the light window is fully pressed against the skin. Check if the device is in the correct operating mode. If the issue persists, contact customer support.

Problem	Possible Cause	Solution
Skin feels too hot or uncomfortable.	Intensity level too high; Improper skin preparation.	Reduce the intensity level. Ensure skin is clean and dry. Verify skin tone compatibility.
No hair reduction after several sessions.	Inconsistent use; Incorrect intensity; Incompatible skin/hair type.	Ensure consistent treatment schedule. Increase intensity if comfortable. Re-check skin tone and hair color compatibility chart.

If the problem persists after trying these solutions, please contact INIA customer support for further assistance.

8. SPECIFICATIONS

Feature	Detail
Model Name	HAUTE
Item Model Number	3023281259
Brand	INIA
Color	Black
Power Source	Corded Electric
Special Features	Double-Pulse, Ice-cooling, Rechargeable (Note: 'Rechargeable' might be a mislabel, as 'Corded Electric' is also listed. Device operates via cord.)
Flashes per Second	Up to 5
Energy Output	Up to 21J (Joules)
Ice-Cooling Temperature	Below 60°F (approx. 50°F)
Product Dimensions	10.43 x 6.88 x 0.1 inches
Item Weight	2.54 Pounds
Date First Available	December 18, 2023

9. WARRANTY AND SUPPORT

9.1 Warranty Information

The INIA HAUTE device comes with a 1-year warranty. This warranty can be extended to 2 years upon product registration. Free replacements are provided within the warranty period for manufacturing defects or if the device runs out of flashes.

9.2 Customer Support

For any questions, technical assistance, or warranty claims, please contact INIA customer support. Refer to the contact

information provided in your product packaging or visit the official INIA website for support details.

Seller: [INIA Official US](#)

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question