

[manuals.plus](#) /

› [Bitvae](#) /

› [Bitvae Water Dental Flossers for Teeth, Cordless Flosser with 3 Modes & 5 Intensities, Water Teeth Cleaner Picks for Cleaning, IPX7 Waterproof & Rechargeable, Portable for Travel, Pink User Manual](#)

Bitvae D2

Bitvae Water Dental Flosser User Manual

Model: D2

PRODUCT OVERVIEW

The Bitvae Cordless Water Dental Flosser is designed to provide effective and convenient oral hygiene. It features multiple cleaning modes and intensity settings to cater to various dental needs, from gentle cleaning to deep plaque removal. Its portable design, long battery life, and waterproof construction make it suitable for home use and travel.



The Bitvae Water Dental Flosser, including its packaging and the six interchangeable tips provided for comprehensive oral care.

Safety Information

Please read all safety instructions before using the Bitvae Water Dental Flosser to ensure safe operation and prevent injury or damage.

- Do not submerge the device completely in water for extended periods, despite its IPX7 waterproof rating.
- Keep the device out of reach of children.
- Consult a dental professional before use if you have certain medical conditions or have recently undergone oral surgery.
- Always follow proper charging and battery handling instructions to prevent electrical hazards or battery damage.

- Use only clean water or dentist-recommended mouthwash in the reservoir. Do not use abrasive or corrosive liquids.

COMPONENTS AND WHAT'S IN THE BOX

Your Bitvae Water Dental Flosser package includes the following items:

- Bitvae Cordless Water Dental Flosser Unit
- 5 Standard Jet Tips
- 1 Tongue Cleaner Tip
- DC-USB Charging Cable
- User Manual

SETUP INSTRUCTIONS

Initial Charging

Before first use, fully charge the device. Connect the DC-USB cable to the charging port on the flosser and plug the USB end into a 5V1A adapter or a compatible USB port (e.g., power bank, laptop). A full charge takes approximately 4 hours and provides up to 40 days of use.

Attaching a Nozzle

Select your desired tip (standard or tongue cleaner). Insert the tip into the hole on the top of the flosser handle until it clicks into place. To remove, press the pop-up knob on the side of the handle and pull the tip out.

Filling the Water Tank

The flosser features a 260ML wide-mouth water tank. To fill, detach the water tank from the main unit by twisting or pulling it down (refer to the diagram below). Fill the tank with lukewarm water. Reattach the tank securely to the main unit.

9 Easy to Clean OZ (260ML)



The wide-mouth design of the water tank allows for easy filling and cleaning.

OPERATING INSTRUCTIONS

The Bitvae Water Dental Flosser offers 3 modes and 5 intensity settings for a total of 15 cleaning combinations.

Selecting Mode and Intensity

Press the 'Mode' button to cycle through the available modes: Clean, Jet, and Massage. The selected mode will be indicated by an LED light. Use the 'Intensity' button to adjust the water pressure from 1 (gentle) to 5 (powerful), ranging from 16-136 psi.

Using the Flosser

1. Fill the water tank and attach a nozzle.
2. Select your desired mode and intensity.

3. Lean over a sink and place the tip into your mouth, directing it at a 90-degree angle to your gumline.
4. Press the 'Power' button to turn on the device. Allow the water to flow from your mouth into the sink.
5. Move the tip along your gumline, pausing briefly between teeth. Ensure you clean both the front and back surfaces of your teeth.
6. Continue for approximately 60 seconds or until the water tank is empty.
7. Press the 'Power' button again to turn off the device.



Visual guide for the five main steps of operating the water flosser.

Deep Cleaning and Improves Gum Health



The water flosser effectively cleans between teeth and around braces, promoting improved gum health.

Deeply Removes Food & Stains

Enjoy a thorough, comfortable teeth-cleaning experience

Bitvae


1400

Pulsations per minute

0.6mm

Water column

20-120psi

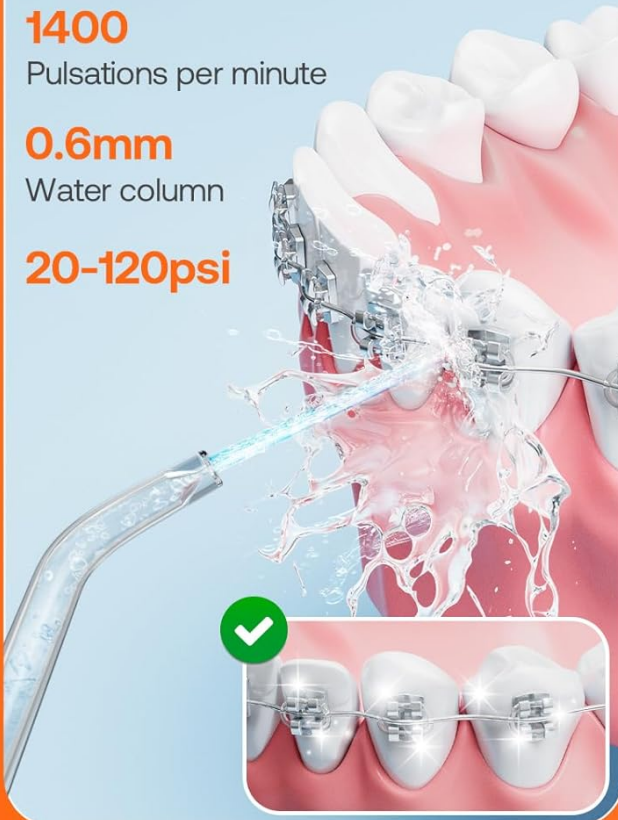

Other

1200

Pulsations per minute

0.66mm

Water column

30-105psi



The Bitvae flosser's high pulsation rate and fine water column effectively remove food particles and stains.

MAINTENANCE

Cleaning the Water Tank

The wide-mouth design of the 260ML water tank allows for easy and thorough cleaning. After each use, detach the tank and rinse it with clean water. For deeper cleaning, you can use a small brush to remove any residue. Ensure the tank is completely dry before reattaching or storing.

9 Easy to Clean OZ (260ML)



260_{ML}

The detachable wide-mouth water tank is designed for easy cleaning, preventing water fouling.

Cleaning the Nozzles

Rinse the tips under running water after each use. For a more thorough cleaning, you can soak them in a solution of 2 parts water to 1 part white vinegar for 5-10 minutes, then rinse thoroughly.

Battery Care

To maintain optimal battery life, avoid completely draining the battery before recharging. If the device will not be used for an extended period, charge it fully before storage and recharge every few months.

Long Battery Life

Up to 40 Days

40 Days



14 Days



Bitvae Others

The Bitvae Water Dental Flosser boasts a long battery life of up to 40 days on a single charge.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on or stops working suddenly.	Low battery; mineral buildup from hard water.	Charge the device fully. If using hard water, try using distilled water. For mineral buildup, add a few drops of vinegar to the water reservoir and run the device briefly, then rinse thoroughly.

Problem	Possible Cause	Solution
Weak water pressure.	Low battery; clogged nozzle; incorrect intensity setting.	Recharge the device. Clean the nozzle thoroughly. Increase the intensity setting.
Water tank leaks.	Water tank not properly attached.	Ensure the water tank is securely attached and properly sealed to the main unit.

SPECIFICATIONS

- **Model Number:** D2
- **Dimensions:** 6.6 x 9.4 x 22.3 cm
- **Weight:** 380 g
- **Water Tank Capacity:** 260 ML (9 oz)
- **Pressure Range:** 16-136 psi (5 intensities)
- **Modes:** 3 (Clean, Jet, Massage)
- **Battery Life:** Up to 40 days (approx. 4 hours charge time)
- **Waterproof Rating:** IPX7
- **Power Source:** Battery Powered (1 Lithium Ion battery included)
- **Country of Origin:** China

WARRANTY AND SUPPORT

Bitvae offers a **24-month quality service** for this product. For any issues, questions, or support needs, please contact Bitvae customer service through the retailer where you purchased the product or visit the official Bitvae website for contact information.

For more information, you may refer to the official Bitvae brand page [Bitvae Brand Store](#)

ADDITIONAL FEATURES

IPX7 Waterproof Design

The IPX7 waterproof rating allows for safe use in the shower, providing convenience and flexibility for your oral hygiene routine.



IPX7 Waterproof

Allows use in the shower

The IPX7 waterproof design ensures the flosser can be safely used in wet environments like the shower.

Portable for Travel

Its compact and cordless design makes the Bitvae Water Dental Flosser ideal for travel, ensuring you can maintain your oral hygiene routine anywhere.

Portable for Travel

At Home, in The Office, or on The Go



The compact size makes the flosser easy to pack and carry for travel.

Family Sharing Tips

The inclusion of multiple standard tips and a tongue cleaner tip allows for individual use by different family members, promoting shared oral hygiene with distinct tips.

2-Colored Tips for Family Sharing




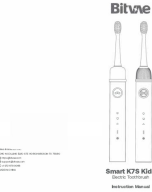
- 5 Normal Nozzles
- 1 Tongue Scraping Nozzle



The flosser comes with multiple tips, allowing for family sharing while maintaining hygiene.

© 2025 Bitvae. All rights reserved.

Related Documents - D2

	<p>Bitvae Oral Irrigator User Manual and Guide</p> <p>Comprehensive user manual for the Bitvae Oral Irrigator, covering important safeguards, product design, usage instructions, specifications, and warranty information. Learn how to charge, operate, and maintain your oral irrigator for optimal dental hygiene.</p>
	<p>Bitvae S3 Sonic Electric Toothbrush User Manual</p> <p>Comprehensive user manual for the Bitvae S3 Sonic Electric Toothbrush, covering features, operation, troubleshooting, and safety information.</p>
	<p>Bitvae Daily R2 Rotary Electric Toothbrush Instruction Manual</p> <p>Comprehensive instruction manual for the Bitvae Daily R2 Rotary Electric Toothbrush, covering features, specifications, usage guidelines, and safety information.</p>
	<p>Bitvae Smart K7S Kids Electric Toothbrush: Instruction Manual</p> <p>Comprehensive instruction manual for the Bitvae Smart K7S Kids Electric Toothbrush, covering features, operation, safety guidelines, and after-sales service. Learn how to use the Bluetooth app and maintain your child's oral hygiene.</p>

