

DermRays B0D1DRYKSM

DermRays LED Light Therapy Mask Instruction Manual

Model: B0D1DRYKSM

1. INTRODUCTION

This manual provides comprehensive instructions for the safe and effective use of your DermRays LED Light Therapy Mask. Please read this manual thoroughly before operating the device and retain it for future reference.

The DermRays LED Light Therapy Mask is designed for at-home skin care, utilizing various light wavelengths to address specific skin concerns such as fine lines, wrinkles, redness, and blemishes.

2. SAFETY INFORMATION

No specific safety information was provided in the product data. Always consult a healthcare professional before starting any new skin treatment, especially if you have pre-existing skin conditions, are pregnant, or have light sensitivity. Avoid direct eye exposure to the LED lights. Do not use on broken or irritated skin. Discontinue use if irritation occurs.

3. PRODUCT COMPONENTS

Your DermRays LED Light Therapy Mask package includes the following items:

- DermRays LED Light Therapy Mask (Face and Neck piece)
- Remote Control
- USB Charging Cable
- Adjustable Straps (4 pieces)
- Instruction Manual (this document)



Image: The DermRays LED Light Therapy Mask and its retail packaging, showcasing the mask, remote, and other accessories.

Your browser does not support the video tag.

Video: An unboxing demonstration of the DermRays LED Light Therapy Mask, showing the mask, remote, straps, and charging cable.

4. SETUP

1. **Prepare Your Skin:** Ensure your face and neck are clean and dry before each session. Remove all makeup and skincare products.
2. **Attach Straps:** Secure the adjustable straps to the mask for a comfortable and snug fit around your head and neck.
3. **Connect Remote:** Plug the mask's USB cable into the remote control unit.
4. **Position the Mask:** Place the mask gently over your face, ensuring the eye and mouth openings are aligned. Adjust the straps so the mask sits comfortably and securely.

HOW TO USE

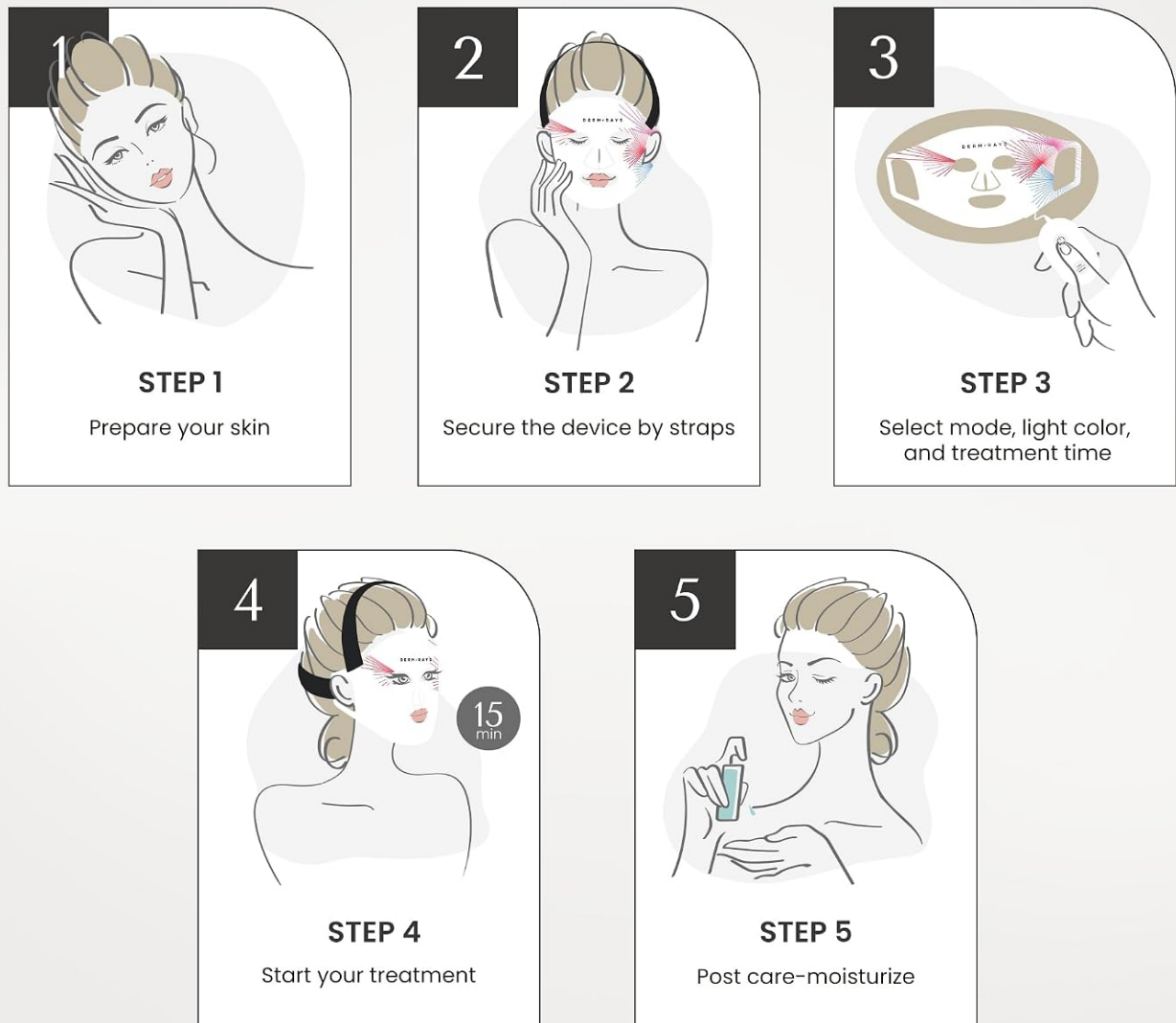


Image: Visual guide illustrating the five steps for using the DermRays LED Light Therapy Mask, from skin preparation to post-treatment moisturizing.

Your browser does not support the video tag.

Video: A demonstration on how to properly wear and use the LED Light Therapy Mask, including strap adjustment and remote control operation.

5. OPERATING INSTRUCTIONS

The DermRays LED Light Therapy Mask features four distinct light modes, each targeting specific skin concerns. Each session is automatically set for 15 minutes.

Light Modes:

- **Red Light (630nm):** Stimulates collagen production, improves skin elasticity, and reduces the appearance of fine lines and wrinkles.
- **Blue Light (460nm):** Targets acne-causing bacteria, helps control oil production, and reduces blemishes.

- **Infrared Light (850nm):** Promotes skin healing, reduces inflammation, and calms the skin.
- **Mixed Light:** Combines the benefits of all modes for a comprehensive skin treatment.

Using the Remote Control:

The remote control allows you to select your desired light mode. Refer to the markings on the remote for mode selection (M1, M2, M3, M4 for different light combinations or automatic/manual modes).



4 Lights Therapy



15 Mins Treatment



100 LEDs

Image: The remote control for the DermRays LED Light Therapy Mask, showing buttons for automatic and manual mode selection, and indicating specific light durations for each mode.

Your browser does not support the video tag.

Video: A demonstration of the four color light modes of the LED mask and their respective benefits for skin treatment.

Treatment Duration and Frequency:

Each session lasts 15 minutes, after which the device will automatically turn off. For optimal results, it is recommended to use the mask 3-4 times per week.

15 Mins for Face and Neck



Image: The LED Light Therapy Mask emitting red light, highlighting its use for 15-minute face and neck treatments to reduce fine lines and wrinkles.

HOW TO USE




STEP 1

Prepare your skin



STEP 2

Secure the device
by straps

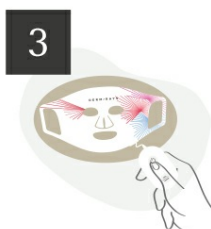

STEP 3

Select mode, light color,
and treatment time



STEP 4

Start your treatment



STEP 5

Post care-moisturize

Image: A cross-section diagram of skin layers (epidermis, dermis, fat, hypodermis) illustrating how blue, red, and near-infrared (NIR) light penetrate to different depths.

6. MAINTENANCE

- **Cleaning:** After each use, gently wipe the mask with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the mask in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective pouch to prevent damage.
- **Charging:** Charge the remote control unit as needed using the provided USB cable. Ensure the device is fully charged before use for uninterrupted sessions.

7. TROUBLESHOOTING

- **Device Not Turning On:**
 - Ensure the remote control is fully charged.
 - Check that the USB cable is securely connected between the mask and the remote.
- **Lights Not Activating:**
 - Verify that a light mode has been selected on the remote control.
 - Confirm the device is powered on.
- **Discomfort During Use:**
 - Adjust the straps to ensure the mask is not too tight.
 - Ensure the mask is positioned correctly on your face.
 - If discomfort persists, discontinue use and consult a healthcare professional.

8. SPECIFICATIONS

Feature	Detail
Brand	DermRays

Model	B0D1DRYKSM
Light Colors	4 (Red, Blue, Infrared, Mixed)
Treatment Duration	15 minutes (automatic shut-off)
Power Source	USB rechargeable remote control
Material	Soft Silicone
Package Dimensions	15.8 x 8.75 x 2.37 inches
Weight	2.2 Pounds
Manufacturer	Guangdong Newdermo Biotech Co.,Ltd

9. WARRANTY AND SUPPORT

Please refer to the product packaging or contact DermRays customer support for detailed warranty information and assistance. The product data indicates a 1-year warranty, though some reviews mention a 2-year warranty. For the most accurate information, consult the official product documentation or manufacturer directly.

© 2023 DermRays. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question

