

[Manuals.plus](#) /

› [Rayovac](#) /

› Rayovac Extra Size 312 Zinc-Air Hearing Aid Batteries Instruction Manual

Rayovac 312

Rayovac Extra Size 312 Zinc-Air Hearing Aid Batteries Instruction Manual

Model: 312 | Brand: Rayovac

INTRODUCTION

This manual provides essential instructions for the safe and effective use of Rayovac Extra Size 312 Zinc-Air Hearing Aid Batteries. These batteries are designed to provide long-lasting, mercury-free power for compatible hearing devices. Please read this manual thoroughly before use.

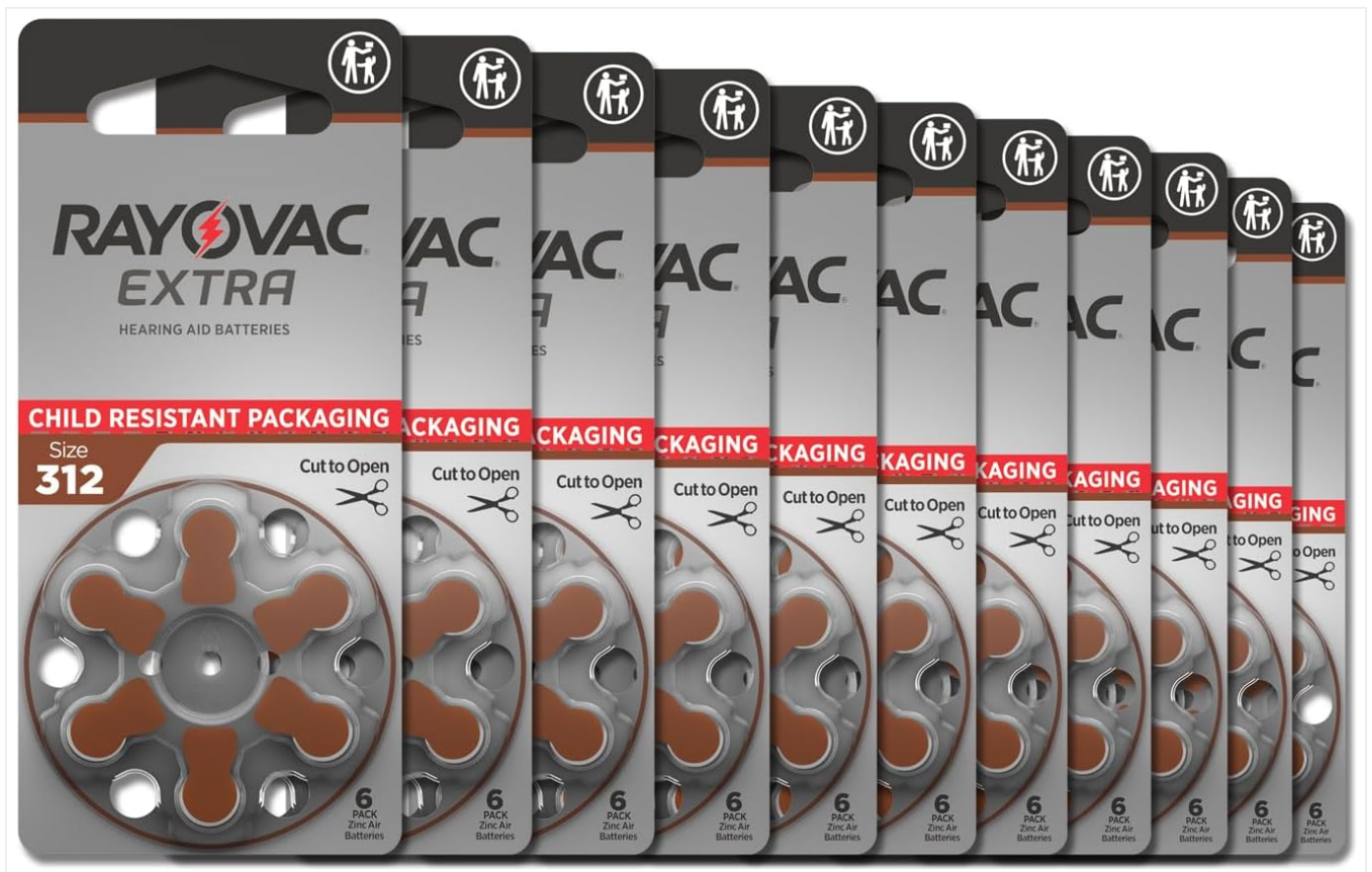


Image: Rayovac Extra Size 312 Hearing Aid Batteries, 60-count packaging.

SETUP AND INSTALLATION

Rayovac Extra Size 312 batteries feature child-resistant packaging to enhance safety. Follow these steps to safely access and install your batteries:

1. **Opening the Packaging:** The battery packaging is designed to be child-resistant and requires scissors to open. Carefully cut along the indicated line on the package to access the battery wheel. It is recommended to perform this over a table to prevent batteries from falling.

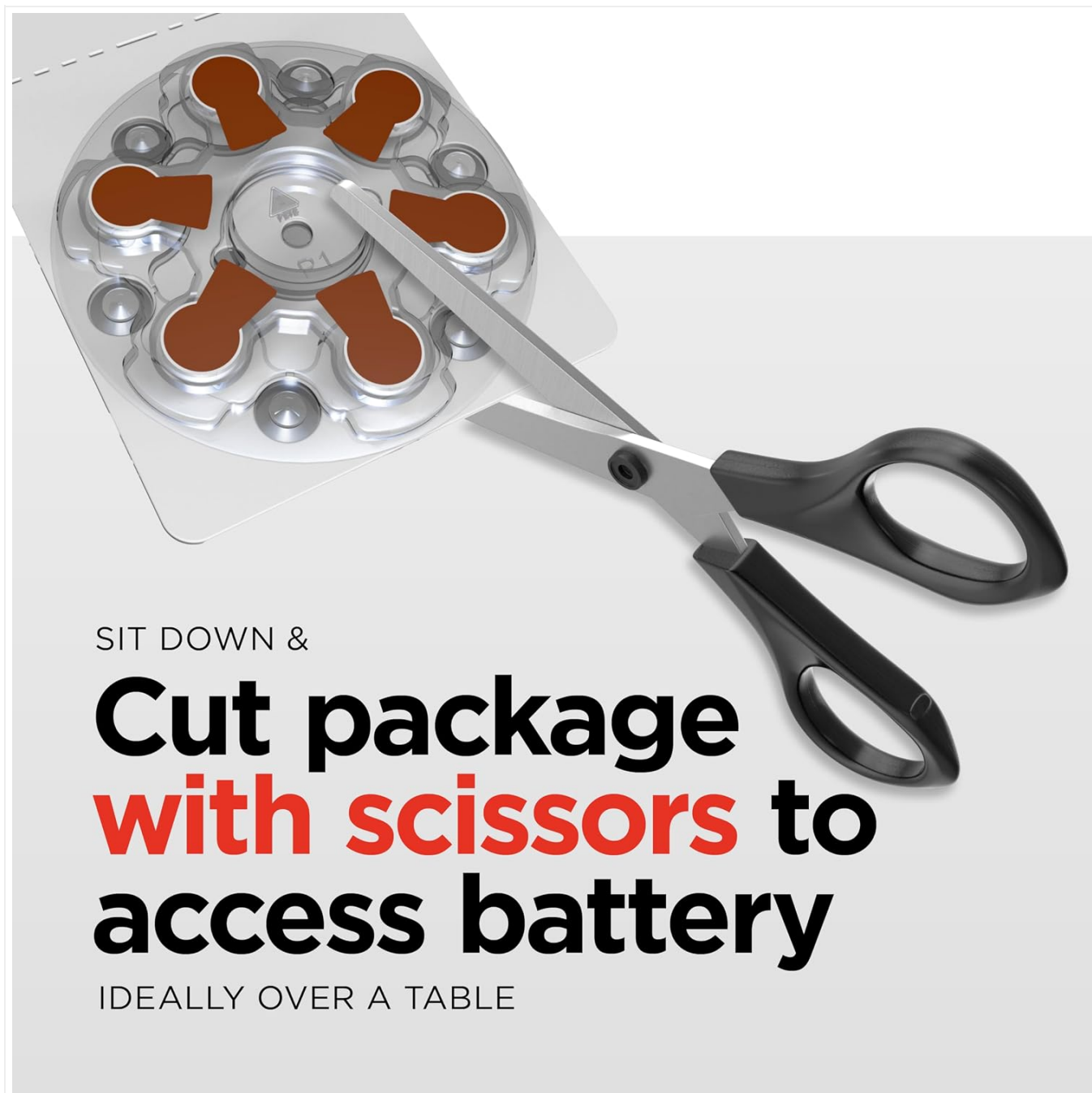


Image: Instructions for cutting open the child-resistant battery packaging.

2. **Removing a Battery:** Once the package is open, carefully remove one battery from the wheel.
3. **Activating the Battery:** Rayovac Extra batteries are zinc-air batteries. They are activated by air. Remove the protective tab from the back of the battery. For optimal performance, allow the battery to sit un-tabbed in the open air for approximately one minute before inserting it into your hearing device. This allows air to fully activate the zinc-air chemistry.
4. **Inserting into Device:** Refer to your hearing aid's instruction manual for specific guidance on how to correctly insert the battery. Ensure the battery is oriented with the correct polarity (+/-) as indicated by your device.

OPERATING INSTRUCTIONS

Once properly installed and activated, your Rayovac Extra Size 312 batteries will power your hearing aid. The lifespan of the battery will vary depending on the hearing aid model, usage patterns, and environmental conditions.

- **Power Indicator:** Most modern hearing aids have a low battery indicator. Replace the battery when this indicator appears or when you notice a decrease in sound quality or volume.
- **Consistent Performance:** Rayovac Extra batteries are designed for consistent power output throughout their life cycle.

MAINTENANCE AND STORAGE

- **Storage:** Store unused batteries in their original, unopened packaging at room temperature, away from direct sunlight and extreme temperatures. Do not refrigerate.
- **Handling:** Always handle batteries with clean, dry hands. Avoid touching the activated side of the battery (where the tab was removed) with your fingers, as oils and dirt can affect performance.
- **Device Care:** When not using your hearing aid, open the battery door to conserve battery life and allow moisture to escape, which can help prevent corrosion.

TROUBLESHOOTING

- **Hearing Aid Not Working:**
 - Ensure the battery tab has been removed and the battery has been allowed to air activate for one minute.
 - Check that the battery is inserted correctly with the proper polarity.
 - Verify the battery door is fully closed.
 - Try a new battery. If the issue persists, consult your hearing aid's manual or a hearing care professional.
- **Short Battery Life:**
 - Ensure the battery tab was removed just prior to use and not prematurely.
 - Check for excessive moisture or extreme temperatures, which can reduce battery life.
 - High-drain features on your hearing aid (e.g., streaming, Bluetooth) can consume more power.
 - If consistently experiencing short life with new batteries, consult your hearing care professional to check your device.

SPECIFICATIONS

Battery Type	Zinc-Air
Size	312 (Brown Tab)
Voltage	1.45V
Chemistry	Mercury-Free
Number of Batteries per Pack	6
Total Batteries (as packaged)	60
Product Barcode (UPC/EAN)	10012800539670
Package Dimensions	4.61 x 2.28 x 2.24 inches (for 60-count package)
Item Weight	3.2 ounces (for 60-count package)
Manufacturer	Energizer Trading Limited

Date First Available	June 22, 2024
----------------------	---------------

SAFETY WARNINGS AND DISPOSAL

WARNING:

- **Keep out of reach of children.** These batteries are small and can be swallowed. If swallowed, promptly seek medical attention and call the National Battery Ingestion Hotline at 800-498-8666.
- **Do not dispose of in fire or charge.** Batteries may explode or leak, causing injury.
- **Immediately dispose of used batteries.** Place them in a secure container out of reach of children and pets until proper disposal can occur.
- **Avoid short-circuiting.** Do not carry loose batteries in pockets or purses with metal objects like coins or keys.

Disposal: Dispose of used batteries according to local regulations. Many communities have recycling programs for batteries. Do not throw batteries into household trash.

Size 312 6 PACK
Zinc Air
Batteries

312CR-6XEUS
1.45V

HEARING AID BATTERIES



For best results, let a new battery sit un-tabbed for one minute before inserting into device.

Warning: (1) Keep out of reach of children. Keep in original package until use. If swallowed, promptly see doctor; have doctor phone 800-498-8666. (2) do not dispose of in fire or charge – may explode or leak causing injury. Immediately dispose of used batteries.



REPLACES ALL SIZE 312

©2024 Energizer. Trademarks and certain designs are owned by Energizer. Made in the United Kingdom for Energizer Brands, LLC, St. Louis, MO 63141



Cut to Open



RAYOVAC
EXTRA

BEST BEFORE



Mercury Free
UPN-175837

WARRANTY AND SUPPORT

For specific warranty information or product support, please refer to the official Rayovac website or contact Energizer Trading Limited, the manufacturer. Contact details are typically available on the manufacturer's official channels.

© 2024 Rayovac. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question