

[Manuals.plus](#) /

> [UPBRIGHT](#) /

> UpBright AC/DC Adapter User Manual

UPBRIGHT DA60-220S24-W1, JE01-40-006HM

UpBright AC/DC Adapter User Manual

Brand: UPBRIGHT

Models: DA60-220S24-W1, JE01-40-006HM

1. INTRODUCTION

This user manual provides essential information for the safe and effective use of your UpBright AC/DC Adapter. This power supply is designed to be compatible with various Resvent iBreeze APAP and CPAP machine series, ensuring reliable power delivery for your medical device.

Please read this manual thoroughly before using the adapter and retain it for future reference.

2. SAFETY INFORMATION

Always observe the following safety precautions to prevent injury or damage to the adapter or connected devices:

- Use the adapter only with compatible devices as specified in this manual.
- Do not expose the adapter to water, moisture, or extreme temperatures.
- Avoid dropping or subjecting the adapter to severe impact.
- Do not attempt to open, modify, or repair the adapter. Refer all servicing to qualified personnel.
- Ensure the power outlet is easily accessible and unplug the adapter during electrical storms or when unused for extended periods.
- Keep the adapter away from heat sources and ensure adequate ventilation during use.
- This adapter is for **indoor use only**.

3. PACKAGE CONTENTS


Verify that all items are present in the package:

- 1 x UpBright AC/DC Adapter
- 1 x AC Power Cord (suitable for your region)



Image 3.1: The UpBright AC/DC Adapter shown with its accompanying AC power cord.

Come With A Suitable AC Power Cord (Only One)
to Fit Your Country/Area

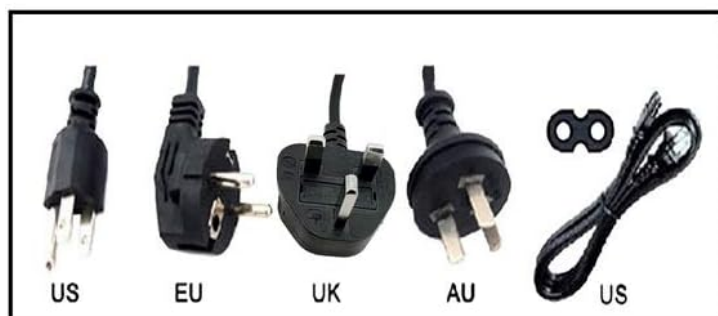

THIS PRODUCT COMES WITH A UPBRIGHT®

LIMITED WARRANTY CERTIFICATE
+ UPBRIGHT® LIMITED WARRANTY



Image 3.2: Examples of different AC power cord plug types (US, EU, UK, AU) that may be included depending on your region.

4. SETUP INSTRUCTIONS

Follow these steps to set up your UpBright AC/DC Adapter:

1. **Connect the AC Power Cord:** Insert the AC power cord firmly into the input port of the UpBright AC/DC Adapter.
2. **Connect to Device:** Plug the DC output connector of the adapter into the power input port of your compatible Resvent iBreeze APAP/CPAP machine. Ensure a secure connection.
3. **Connect to Power Outlet:** Insert the AC power cord plug into a standard wall power outlet (100-240VAC, 50/60Hz).
4. **Verify Power:** Once connected, your Resvent iBreeze machine should receive power. Refer to your machine's manual for specific power indicators.



Image 4.1: The UpBright AC/DC Adapter connected to a wall power outlet.

UpBright Power Supply Adapter Compatible with:



Image 4.2: An example of a Resvent iBreeze CPAP machine, which is compatible with this power adapter.

5. OPERATING INSTRUCTIONS

The UpBright AC/DC Adapter provides stable power to your Resvent iBreeze machine. It is designed for continuous operation within its specified voltage and current limits.

5.1. Compatibility

This adapter is compatible with the following Resvent iBreeze APAP/CPAP Machine Series models:

- iBreeze 20A
- iBreeze 20A Pro
- iBreeze 20C
- iBreeze 20C Pro
- iBreeze 25S
- iBreeze 25A
- iBreeze 25ST
- iBreeze 25STA
- iBreeze 30ST
- iBreeze 30STA
- Resvent Model: DA60-220S24-W1 (Power Supply)
- SHANGHAI BIAOJUN ELECTRONICS TECHNOLOGY Model: JE01-40-006HM (Power Supply)

Always ensure your device's power requirements match the adapter's output specifications to prevent damage.

6. MAINTENANCE AND CARE

Proper care will extend the life of your UpBright AC/DC Adapter:

- **Cleaning:** Disconnect the adapter from all power sources and devices before cleaning. Use a soft, dry cloth to wipe the surface. Do not use liquid cleaners or abrasive materials.
- **Storage:** Store the adapter in a cool, dry place away from direct sunlight and extreme temperatures when not in use. Avoid tightly winding the cables to prevent damage.
- **Handling:** Do not pull on the cord to unplug the adapter from the wall outlet or the device. Always grasp the plug itself.

7. TROUBLESHOOTING

If you encounter issues with your UpBright AC/DC Adapter, try the following steps:

- **No Power to Device:**
 - Ensure the adapter is securely plugged into both the wall outlet and the Resvent iBreeze machine.
 - Test the wall outlet with another device to confirm it is functioning.
 - Check for any visible damage to the adapter or cables. If damaged, discontinue use.
- **Intermittent Power:**
 - Ensure all connections are firm and not loose.
 - Avoid bending or crimping the power cords.
- **Adapter is Hot:**
 - It is normal for the adapter to feel warm during operation. However, if it is excessively hot to the touch, unplug it immediately.
 - Ensure adequate ventilation around the adapter and that it is not covered.
 - Verify that the connected device is compatible and not drawing excessive power.

If the issue persists after attempting these steps, please contact UpBright customer support.

8. SPECIFICATIONS

Feature	Specification
Input Voltage	100-240VAC, 50/60Hz (World Wide Input Voltage)
Output Wattage	24 watts
Protection Features	OVP (Over Voltage output Protection), OCP (Over Current output Protection), SCP (Short Circuit output Protection)
Product Dimensions	6 x 3 x 2 inches
Item Weight	10 ounces
Connector Type	Barrel Connector
Color	Black

9. WARRANTY AND SUPPORT

UpBright stands behind the quality of its products.

- **Refund Policy:** 30 days refund period from the date of purchase.
- **Exchange Policy:** 24 months exchange period for defective units.

For technical support, warranty claims, or any inquiries, please visit the official UpBright store or contact their customer service through the platform where the product was purchased.

[Visit the Official UPBRIGHT Store on Amazon](#)

© 2024 UPBRIGHT. All rights reserved.

This manual is for informational purposes only. Specifications are subject to change without notice.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question

Open the live page

<https://manuals.plus/asin/B0DCNQW19J>

More documents and user manuals: <https://manuals.plus/>