

[Manuals.plus](#) /

> [WELLLEO](#) /

> Water Flosser for Teeth, 300ML Cordless Oral Irrigator Rechargeable Dental Flosser Picks with 5 Modes 5 Jet Tips, IPX7 Waterproof Teeth Cleaner for Braces Gum Care, Home Travel Use White User Manual

WELLLEO HF-P12

WELLLEO Water Flosser for Teeth

USER MANUAL

Brand: WELLLEO | Model: HF-P12

INTRODUCTION AND OVERVIEW

The WELLLEO Water Flosser is engineered to provide an effective and gentle solution for comprehensive oral hygiene. Utilizing a strong pulsating water flow, it targets stubborn debris and food particles, promoting healthier teeth and gums. Its versatile design includes five adjustable cleaning modes to cater to various needs, from sensitive gums to deep periodontal cleaning. The device is cordless, rechargeable, and features a large water reservoir for convenience at home or while traveling.

WHAT'S IN THE BOX

Upon unboxing your WELLLEO Water Flosser, you will find the following components:

- 1 x Water Flosser Unit
- 5 x Jet Tips (various types for different needs)
- 1 x USB-C Charging Cable (adapter not included)
- 1 x Cleaning Brush
- 1 x User Manual
- 1 x Reservoir Cleaning Tablet (10pcs)



Image: All components included in the WELLLEO Water Flosser package.

SETUP AND GETTING STARTED

Follow these steps to set up your water flosser for first use:

1. **Charge the Device:** Before first use, fully charge the water flosser using the provided USB-C cable. A full charge can last up to 30 days.
2. **Attach a Jet Tip:** Select the desired jet tip and firmly insert it into the opening at the top of the flosser unit until it clicks into place.
3. **Fill the Reservoir:** Detach the 300ML water reservoir from the main unit. Fill it with warm water for a


more comfortable experience, especially if you have sensitive teeth. Reattach the reservoir securely.

4. **Select Mode:** Press the 'M' button to cycle through the 5 cleaning modes: Strong (140Psi), High (100Psi), Normal (80Psi), Soft (60Psi), and Low (40Psi). The device remembers your last used mode.
5. **Position for Use:** Lean over a sink and place the jet tip into your mouth, directing it at a 90-degree angle towards your gum line. Slightly close your mouth to prevent splashing.
6. **Start Cleaning:** Press the power button to begin flossing. Move the tip along your teeth and gums, pausing briefly on each tooth. The 2-minute timer ensures optimal cleaning duration.

Reservoir Cleaning Tablets



Add 1Pcs cleaning tablet and wait 30 minutes

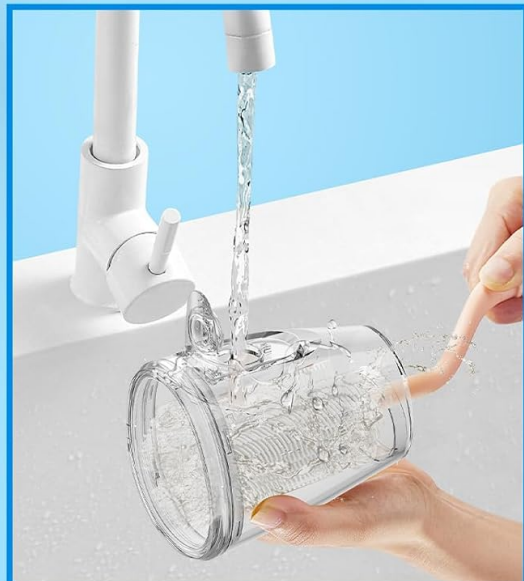

Clean with a brush, then rinse twice with clean water

Note: clean the water tank once a week

Image: Visual guide for initial setup and usage.

Video: Official demonstration of the WellLeo Water Flosser's features and usage.

OPERATING INSTRUCTIONS

The WELLLEO Water Flosser offers various modes to customize your oral care routine:

- **Strong Mode (140Psi):** For deep cleaning and effective removal of stubborn debris.
- **High Mode (100Psi):** Ideal for gum massage and stimulating blood circulation.
- **Normal Mode (80Psi):** Suitable for daily use and regular cleaning.
- **Soft Mode (60Psi):** Provides gentle cleaning for sensitive gums.

- **Low Mode (40Psi):** Recommended for first-time users and those with very sensitive teeth.

The smart memory mode automatically recalls your last used water flow intensity and cleaning mode, streamlining your routine. A built-in 2-minute timer ensures you clean for the dentist-recommended duration.



Image: Detailed breakdown of the 5 cleaning modes.

Strong, High-pulsed Stream



Image: The 300ML large capacity reservoir for extended use.

CLEANING AND MAINTENANCE

Proper cleaning and maintenance will ensure the longevity and hygiene of your water flosser:

1. **Daily Cleaning:** After each use, remove the jet tip and rinse it under clean water. Empty any remaining water from the reservoir and allow all parts to air dry.
2. **Weekly Deep Cleaning:** For thorough cleaning, detach the reservoir. Add 1 cleaning tablet (provided) to the reservoir filled with warm water and let it dissolve for 30 minutes. Use the cleaning brush to scrub the inside of the reservoir. Rinse the reservoir twice with clean water.
3. **Waterproof Design:** The device is IPX7 waterproof, allowing safe use in the bathroom or shower. Always ensure the waterproof cover on the USB-C charging port is securely closed during use to prevent water entry.
4. **Tip Replacement:** It is recommended to replace jet tips every 3 months for optimal performance and hygiene.

IPX7 Waterproof



Image: The IPX7 waterproof design for safe use.

Video: Official demonstration of the WellLeo Water Flosser's features and cleaning process.

TROUBLESHOOTING

If you encounter any issues with your WELLLEO Water Flosser, please refer to the following common solutions:

- **No Water Flow:** Ensure the reservoir is properly filled and securely attached. Check if the jet tip is correctly inserted and not clogged.
- **Weak Water Pressure:** Verify the selected mode is appropriate for your desired pressure. Ensure the battery is sufficiently charged.
- **Device Not Turning On:** Confirm the device is fully charged. If the device overheats, the overheat protection system will cut off power; allow it to cool down before attempting to use it again.
- **Charging Issues:** Ensure the USB-C cable is properly connected to both the device and a 5V 1A adapter (not included). Check the charging port for any debris.

For persistent issues, please contact customer support.

SPECIFICATIONS

Feature	Detail
Product Dimensions	2.68 x 3.35 x 8.86 inches
Item Weight	13.4 ounces
Model Number	HF-P12
Battery	2500mAh Rechargeable (1 AAA battery included for initial function)
Water Tank Capacity	300ML
Water Pressure Range	40-140 PSI (5 Modes)
Pulses Per Minute	1300
Waterproof Rating	IPX7
Power Source	Battery Powered

WARRANTY AND SUPPORT

The WELLLEO Water Flosser comes with a **12-month warranty** to ensure your satisfaction. For any questions, concerns, or support needs, please utilize the **24/7 email support** provided by WELLLEO. Refer to the contact information in the included user manual for details.