

COMFIER B0DMVV8YZD

COMFIER Shiatsu Back Massager with Heat, Tapping, Kneading, and Compression (Model B0DMVV8YZD)

User Instruction Manual

1. INTRODUCTION

Thank you for choosing the COMFIER Shiatsu Back Massager. This device is designed to provide a comprehensive massage experience, combining shiatsu kneading, tapping, compression, and heat therapy. It is suitable for use on most chairs, including office chairs, sofas, and recliners, offering relaxation for your neck, back, and hips. Please read this manual thoroughly before use to ensure proper operation and safety.

2. IMPORTANT SAFETY INFORMATION

To reduce the risk of electric shock, fire, and injury, please observe the following precautions:

- Always unplug the massager from the electrical outlet immediately after using and before cleaning.
- Do not use while bathing or in a shower.
- Do not place or store the unit where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Never use pins or other metallic fasteners with this appliance.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- Keep the cord away from heated surfaces.
- Never use while sleeping or drowsy.
- Never drop or insert any object into any opening.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use on infants, invalids, or on a sleeping or unconscious person.
- Do not use on insensitive skin or on a person with poor blood circulation.
- Do not crush. Avoid sharp folds.

3. PACKAGE CONTENTS

- COMFIER Shiatsu Back Massager Cushion
- Remote Control
- Power Adapter
- User Manual

4. PRODUCT OVERVIEW



Figure 1: COMFIER Shiatsu Back Massager with its remote control and illustrated massage areas.

The COMFIER Shiatsu Back Massager is a versatile cushion designed to fit most chairs. It features:

- **Adjustable Neck Massage Nodes:** For targeted neck and shoulder relief.
- **Shiatsu & Tapping Massage Nodes:** Located in the backrest for deep tissue massage.
- **Compression Airbags:** On the sides and seat for a soothing compression massage.
- **Vibration Motors:** In the seat cushion for added relaxation.
- **Heat Function:** For both the back and neck massage areas.
- **Integrated Remote Control:** For easy selection of massage modes and intensities.

5. SETUP INSTRUCTIONS

1. **Unpack the Massager:** Carefully remove the massager from its packaging.
2. **Position the Massager:** Place the massager cushion on your desired chair (e.g., office chair, sofa, recliner). Ensure the chair has a sturdy backrest to support the massager.
3. **Secure the Massager:** Use the integrated straps on the back of the massager to secure it firmly to your chair. This prevents it from slipping during use.
4. **Connect Power:** Plug the power adapter into the massager's power input port, then plug the adapter into a standard electrical outlet.
5. **Remote Control:** The remote control is tethered to the massager and can be stored in the side pocket when not in use.

6. OPERATING INSTRUCTIONS

Compression Massage

Massage on waist & hips ,Adjustable 3 intensity settings are available.



Figure 2: Remote Control Layout

The massager is operated using the attached remote control. Familiarize yourself with the buttons:

- **Power Button:** Turns the massager ON/OFF.
- **Vibration (1, 2, 3):** Selects one of three vibration intensities for the seat cushion.
- **Demo:** Activates a pre-programmed demonstration mode that cycles through all massage functions.
- **Neck Massage (V/A, Circle Left/Right):** Adjusts the height of the neck massage nodes (V for down, A for up) and changes the rotation direction.
- **2D/3D Shiatsu:** Toggles between 2D and 3D shiatsu massage for the back.
- **Rolling:** Activates a rolling massage function for the back.
- **Spot (V/A):** Allows you to pinpoint the massage to a specific area of your back (V for down, A for up).
- **Air Back / Air Intensity / Air Seat:** Controls the compression airbags for the back and seat, and adjusts

their intensity.

- **Full, Upper, Lower Back:** Selects the massage area for the back (full, upper, or lower).
- **Heat:** Turns the heat function ON/OFF for the back and neck massage nodes.

Massage Functions:

1. **Tapping Massage:** The massager offers 4 tapping massage modes for full back, upper back, and lower back. It has 3 intensity levels.

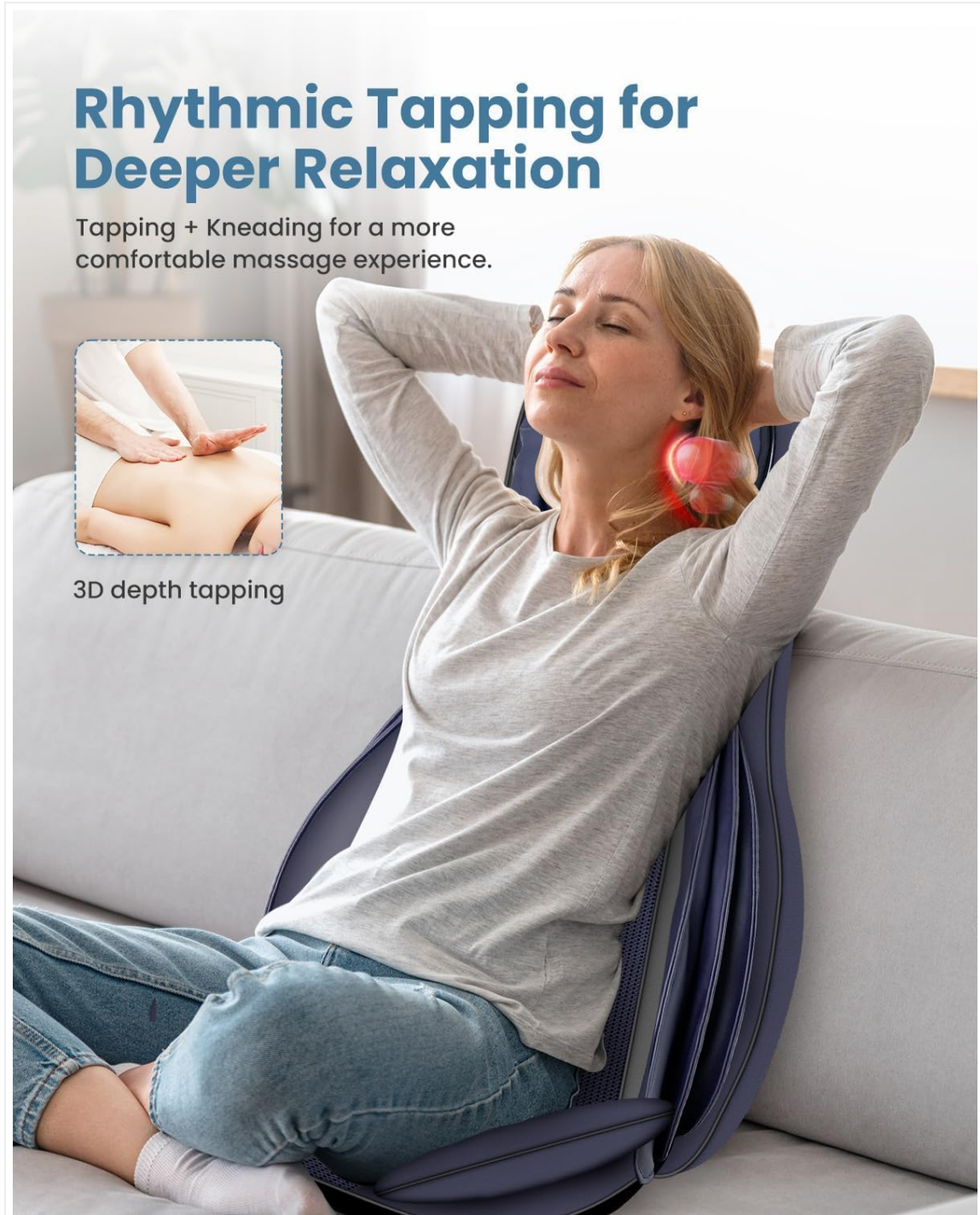


Figure 3: Rhythmic Tapping for Deeper Relaxation

2. **Adjustable Neck Shiatsu:** The neck massage nodes are height-adjustable to suit users between 5' to 5'9". The neck cover is removable and washable for hygiene.

Spot Massage

SPOT massage function allows you to concentrate massage on one area for pin-point relaxation

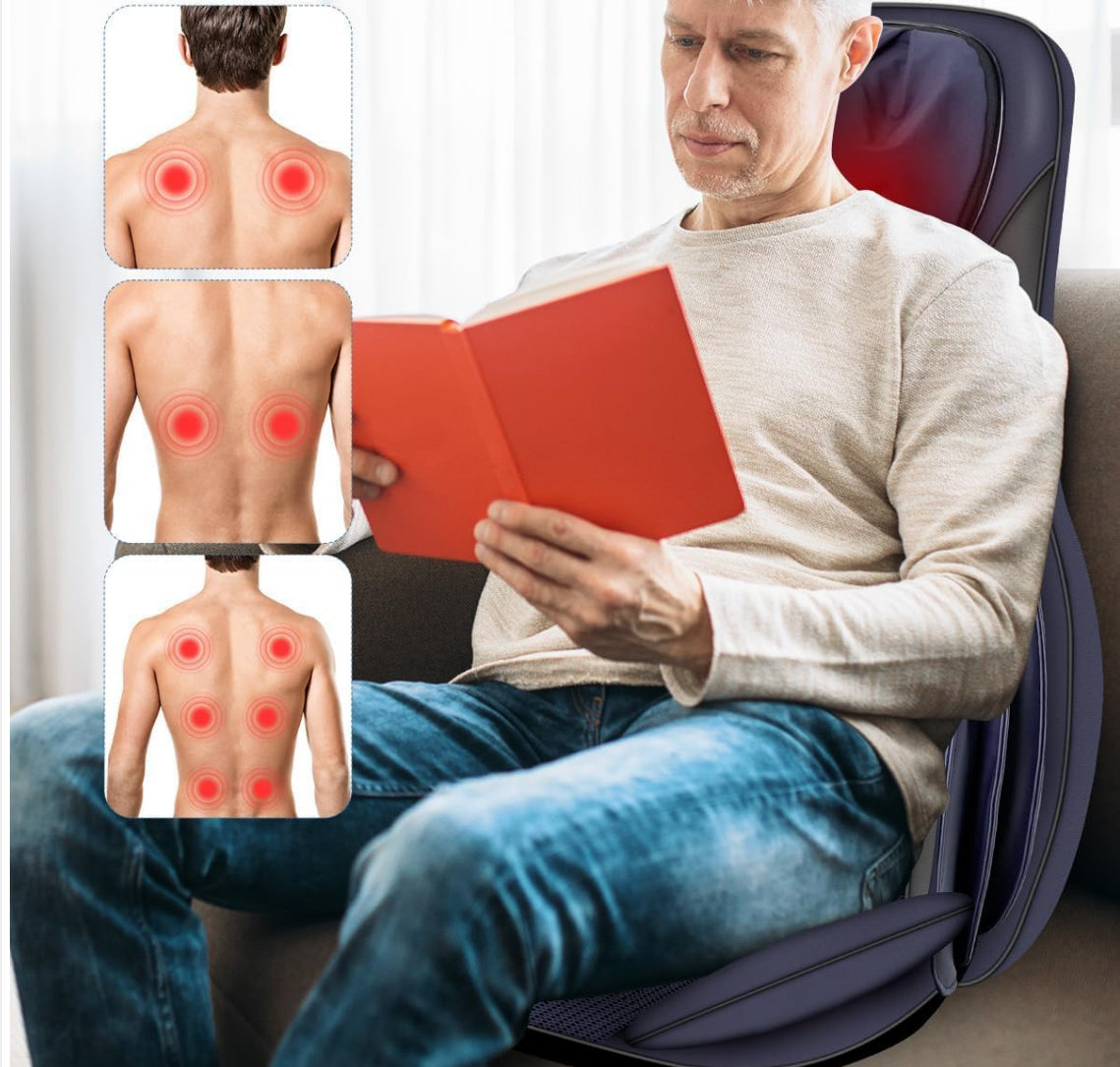


Figure 4: Adjustable Neck Shiatsu Massager

3. **SPOT Massage:** Use the SPOT function to concentrate the massage on a specific area of your back. You can also choose full back, upper back, or lower back massage for the target area.

Timing Settings & Auto off

10/20/30 Mins auto shut off timer



Figure 5: SPOT Massage Function

4. **Compression Massage:** Features 3 compression intensities for the lower back and hip area, providing a multi-layered soothing experience.



Figure 6: Compression Massage on Waist & Hips


5. **Vibration Massage:** The seat pad offers 3 levels of vibration massage.

6. **Heat Function:** Provides soothing heat to the back and neck areas.

Auto Shut-Off Timer:


The massager includes an auto shut-off timer with settings for 10, 20, or 30 minutes to prevent overuse and ensure safety.

Adjustable Neck Shiatsu Massager

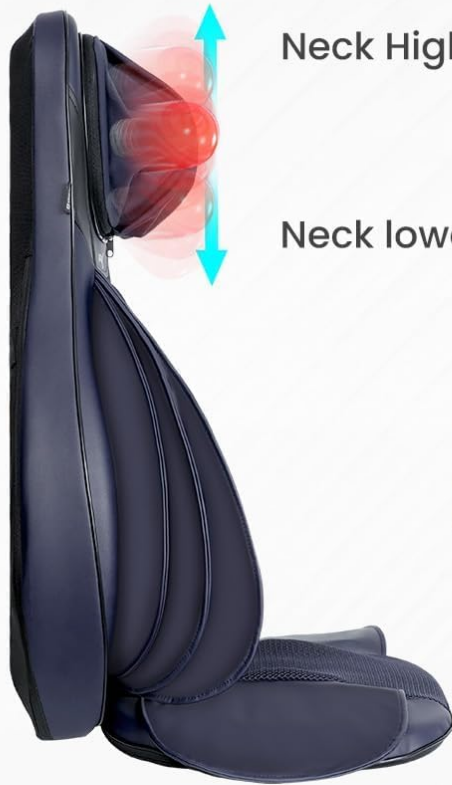


Neck Highest: 5'9

Neck lowest: 5'

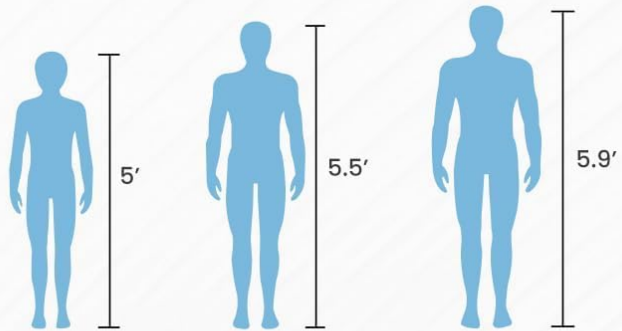

Suitable for most users
(height 5' to 5'9")

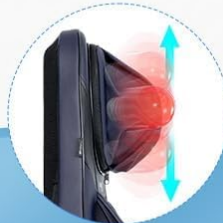

Neck Shiatsu
Auto-Walk Function



Removable and
washable fabric cover

Figure 7: Timing Settings & Auto Off

7. MAINTENANCE AND CARE

- **Cleaning:** Always unplug the massager before cleaning. Wipe the surface with a soft, damp cloth. Do not use abrasive cleaners or harsh chemicals. The neck cover is removable and washable.
- **Storage:** Store the massager in a cool, dry place when not in use. Avoid placing heavy objects on it.
- **Avoid Damage:** Do not attempt to repair the massager yourself. Contact customer support for any issues.

8. TROUBLESHOOTING

Fault	Solutions
-------	-----------

Fault	Solutions
Device switches off suddenly.	Device may have been in use for 20 mins. Auto-shut-off function activated. Continue using by pressing the "power" button after resting the unit for 15 mins. Over-heat protection function activated. Device is switched off for your safety. Please rest device for 15-30 mins before restarting.
Device fails to start.	Check if adapter is plugged into socket. Check if adapter output is connected to device power inlet wire.
How the auto-timer works and how to use this function well.	The auto-timer function starts and defaults to a timer of 20 minutes when device turned on. When you press the "timer" button on the controller once, the auto-timer becomes 30 minutes. When you press it again, the auto-timer becomes 10 minutes. You can choose the timing setting you want, and the device will stop work after timer runs out.
If your height is between 5-5.9ft, but the shiatsu massager roller can't reach your neck after adjusting the roller height by the up/down button.	You can place a seat cushion or a pillow on the sitting part. By this way, the using experience is better.






9. SPECIFICATIONS

- **Product Dimensions:** 20.47 x 9.65 x 30.71 inches
- **Item Weight:** 17.8 Pounds
- **Brand:** COMFIER
- **Model:** B0DMVV8YZD
- **Color:** Navy Blue
- **Manufacturer:** COMFIER

10. WARRANTY AND SUPPORT

This product comes with a manufacturer's warranty. For detailed warranty information or customer support, please refer to the warranty card included in your package or contact COMFIER customer service directly. Statements regarding this product have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

Related Documents - B0DMVV8YZD

<div><div><div>COMFIER</div><div>Shiatsu massagecushion with heat:en:en: 02 Shiatsu massagecushion with heat:fr:fr: 03 Cushion de massage Shiatsu avec chaleur:fr:fr: 04 Cushion per massage Shiatsu riscaldato:it:it: 05 Cojín masajeador Shiatsu con calefacción:es:es: 06</div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-2606MC</div></div></div>	<p>COMFIER CF-2606MC Shiatsu Massage Cushion with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2606MC Shiatsu massage cushion with heat. Includes safety instructions, operating guide, features, maintenance, and warranty information.</p>
<div><div><div>COMFIER</div><div>Shiatsu Neck and Back Massager with Heat</div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-2303-KR</div></div></div>	<p>COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat, covering safety instructions, features, technical data, setup, operation, controller functions, care, cleaning, troubleshooting, and warranty information.</p>
<div><div><div>COMFIER</div><div>Shiatsu Back Support Massage Cushion:en:en: 01 Shiatsu Back Support Massage Cushion:fr:fr: 02 Cushion de massage Shiatsu pour le dos:fr:fr: 03 Cushion massaggiante Shiatsu per il sostegno della schiena:it:it: 04 Cojín de masaje Shiatsu para la espalda:es:es: 05 Shiatsu rugstuen massagecushion:nl:nl: 06</div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-1503S</div></div></div>	<p>Comfier CF-1503S Shiatsu Back Support Massage Cushion User Manual</p> <p>User manual for the Comfier CF-1503S Shiatsu Back Support Massage Cushion, detailing features, operation, maintenance, safety instructions, and warranty information.</p>
<div><div><div>COMFIER</div><div>Shiatsu massage pillow with heat</div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-1228</div></div></div>	<p>Comfier CF-1228 Shiatsu Massage Pillow with Heat User Manual</p> <p>User manual for the Comfier CF-1228 Shiatsu massage pillow with heat, detailing safety instructions, operating procedures, features, maintenance, disposal, technical specifications, and warranty information.</p>
<div><div><div>COMFIER</div><div>Shiatsu Neck and Back Massager</div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-2309A</div></div></div>	<p>Comfier Shiatsu Neck and Back Massager CF-2309A User Manual</p> <p>User manual for the Comfier Shiatsu Neck and Back Massager, model CF-2309A. Includes safety instructions, features, technical data, setup and operation, controller instructions, care and cleaning, troubleshooting, and warranty information.</p>



[Comfier Shiatsu & Percussion Massage Cushion User Manual \(CF-1903\)](#)

User manual for the Comfier Shiatsu & Percussion Massage Cushion (Model CF-1903), detailing safety instructions, features, setup, operation, care, troubleshooting, and warranty information.