

[Manuals.plus](#) /

> [gelpal](#) /

> gelpal 7-1 Colors Red Light Therapy Facial Skin Care Mask (Model MZ-01) User Manual

gelpal MZ-01

gelpal 7-1 Colors Red Light Therapy Facial Skin Care Mask (Model MZ-01) User Manual

Model: MZ-01 | Brand: gelpal

1. INTRODUCTION

This manual provides essential information for the safe and effective use of your gelpal 7-1 Colors Red Light Therapy Facial Skin Care Mask. Please read these instructions thoroughly before operating the device and retain for future reference.

The LED Facial Mask is an over-the-counter, facemask-shaped device that emits light energy in the red, blue, and amber spectrum for the treatment of wrinkles or mild to moderate inflammatory conditions on the face.

2. SAFETY INFORMATION

Important Safeguards:

- Do not use if you are pregnant or breastfeeding.
- Consult a physician before use if you have any pre-existing skin conditions, light sensitivity, or are taking medications that may cause photosensitivity.
- Avoid direct eye exposure to the LED lights. Use the provided eye pads for protection.
- Keep out of reach of children.
- Do not use on broken, irritated, or infected skin.
- Discontinue use if irritation or discomfort occurs.
- Do not immerse the device in water or other liquids.
- Use only the provided charging cable.

3. PRODUCT FEATURES

The gelpal 7-1 Colors LED Facial Skin Care Mask is designed with advanced technology for comprehensive skin

care:

- **7+1 Color LED Light Modes:** This advanced LED treatment mask features 7 colors of light (Red, Blue, Yellow, Green, Purple, Cyan, White) providing a comprehensive skin solution.
- **Cordless & Portable Design:** Enjoy light therapy anytime, anywhere. The lightweight and portable design allows for use on the go. A full charge provides power for 3-4 treatment courses (approximately 2 weeks of typical use).
- **Automatic Timing:** The mask is equipped with automatic timing and timed shutdown functions, automatically closing after a 20-minute treatment period.
- **High-Performance LEDs:** Features 90 high-performance LED beads that precisely target over 60 facial zones for rejuvenation in just 15-minute sessions.
- **Universal Skin Type Compatibility:** Gentle yet effective, designed for all skin types—dry, oily, sensitive, or combination—to address various concerns while maintaining comfort and balance.



Figure 1: Key components of the gelpal LED Facial Mask.

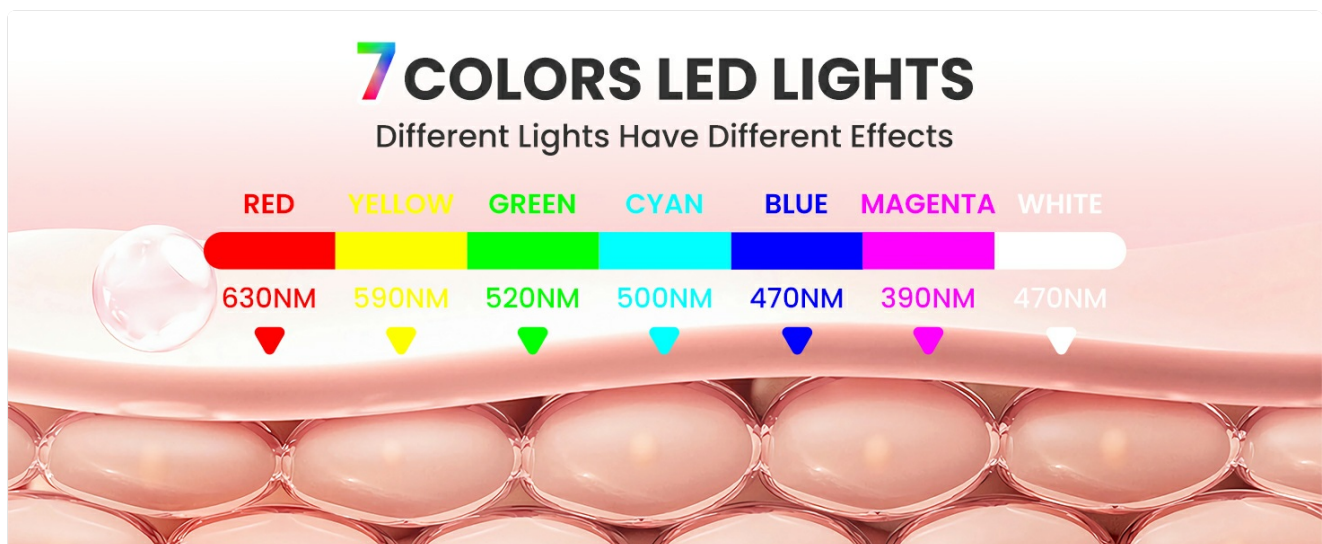


Figure 2: Portable and rechargeable design details.

4. SETUP

1. **Unpack:** Carefully remove the LED Facial Mask and all accessories from the packaging.
2. **Charge:** Before first use, fully charge the device using the provided USB Type-C charging cable. The charging

port is located on the chin area of the mask. A full charge typically takes 1-2 hours. The product will emit a beep sound when charging is complete.

3. **Attach Strap:** Secure the adjustable head strap to the sides of the mask for a comfortable and snug fit.
4. **Insert Eye Pads:** Place the silicone eye pads into the designated areas around the eye openings to protect your eyes from direct light exposure.

5. OPERATING INSTRUCTIONS

1. **Prepare Skin:** Cleanse and dry your face thoroughly before each session.
2. **Wear Mask:** Place the mask on your face, ensuring the eye pads are correctly positioned and the strap provides a comfortable fit.
3. **Power On:** Press and hold the power button (indicated by a hand icon) on the forehead of the mask for approximately 3 seconds until the lights illuminate.
4. **Select Light Mode:** Briefly press the power button to cycle through the different LED light colors. Each color offers specific benefits:
 - **Red Light (630nm):** Deep repair, wrinkle reduction, skin rejuvenation.
 - **Blue Light (470nm):** Acne care, oil control, calming properties.
 - **Yellow Light (590nm):** Brightens skin, detoxifies, improves texture.
 - **Green Light (520nm):** Soothes, balances, and improves complexion.
 - **Purple Light (390nm):** Targets acne scars and promotes healing (combination of red and blue light benefits).
 - **Cyan Light (500nm):** Calms inflamed skin, cooling and soothing effect, good for stressed skin.
 - **White Light (470nm):** Supports skin tightening, helps with fine lines, overall skin rejuvenation.
5. **Treatment Session:** Each session is automatically set for 15-20 minutes. The device will automatically shut off after the session is complete.
6. **Frequency:** For optimal results, use 3-4 times per week.
7. **Power Off:** To manually turn off the device, press and hold the power button until the lights extinguish.



Figure 3: LED Light Color Modes and Benefits.

Your browser does not support the video tag.

Video 1: A brief demonstration of the gelpal 7-1 Colors LED Facial Skin Care Mask, showing how to power it on and cycle through the different light colors.

6. MAINTENANCE

- **Cleaning:** After each use, gently wipe the mask with a soft, damp cloth. Do not use abrasive cleaners or harsh chemicals. Ensure the device is completely dry before storing.
- **Storage:** Store the mask in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Battery Care:** To prolong battery life, avoid fully discharging the device frequently. Charge it when the battery indicator is low.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or improper power button press.	Ensure the device is fully charged. Press and hold the power button for 3 seconds.
Lights are dim or flickering.	Low battery.	Recharge the device.
Discomfort during use.	Mask not properly positioned or light intensity too high (if adjustable, though this model doesn't specify).	Adjust the mask for a comfortable fit. Ensure eye pads are correctly placed. Discontinue use if discomfort persists.

8. SPECIFICATIONS

- **Model Number:** MZ-01
- **Brand:** gelpal
- **LED Colors:** 7 (Red, Blue, Yellow, Green, Purple, Cyan, White)
- **LED Beads:** 90
- **Power Source:** Rechargeable (1 x 9V battery included)
- **Charging Time:** 1-2 hours
- **Usage per Charge:** 3-4 treatment courses
- **Session Duration:** 15-20 minutes (automatic shutdown)
- **Dimensions:** 9.29 x 8.39 x 4.88 inches
- **Weight:** 15.2 ounces
- **UPC:** 721545087492

9. WARRANTY AND SUPPORT

For warranty information or customer support, please refer to the contact details provided with your product packaging or visit the official gelpal website. Keep your purchase receipt for warranty claims.