

[Manuals.plus](#) /

› [CWTUL](#) /

› CWTUL Laser Hair Growth Cap TW280 User Manual

CWTUL TW280

CWTUL Laser Hair Growth Cap TW280 User Manual

PRODUCT OVERVIEW

The CWTUL Laser Hair Growth Cap TW280 is a dome-shaped Low-Level Laser Therapy (LLLT) device. It is designed to promote hair growth in both women and men by exposing the entire scalp to photobiostimulation. The cap features 80 laser diodes, each emitting at 650nm with 41mW power.

The device consists of an outer plastic cover and a protective inner liner that houses the electronics and laser array. It is powered by an included battery pack, ensuring portability and ease of use.



Image: The CWTUL Laser Hair Growth Cap, displaying its external appearance and the internal array of red laser diodes.

TARGET USERS

SPECIALLY DESIGNED FOR
MEN AND WOMEN
WITH THINNING HAIR
STOP HAIR LOSS EARLY



Image: A detailed diagram explaining the mechanism of Low-Level Laser Photobiomodulation (LLLT) on hair follicles, highlighting its effects

SAFETY INFORMATION

Please read and understand all safety warnings before using the device.

- **For external use only:** Do not ingest or allow the product to come into contact with eyes, mouth, or other mucous membranes.
- **Allergic reactions:** Discontinue use if you experience scalp irritation, redness, or rash.
- **Pregnancy and breastfeeding:** Not recommended for use by pregnant or nursing women.
- **Children:** Keep out of reach of children. Not for use by individuals under 18 years old.
- **Scalp conditions:** Do not use on damaged, sunburned, or irritated scalp.
- **Drug interactions:** Consult a healthcare professional if you're using other medications, especially those applied to the scalp.
- **Cardiovascular concerns:** Some treatments may affect blood pressure. Consult a doctor if you have heart disease.
- **Absorption:** Avoid using on other body parts, as this may lead to unwanted hair growth.
- **Sun exposure:** Some treatments may increase sensitivity to sunlight.
- **Proper application:** Follow instructions carefully regarding dosage and frequency of application.
- **Discontinuation effects:** Stopping treatment may lead to reversal of hair regrowth.
- **Flammability:** Some products may be flammable; keep away from fire or flame.
- **Storage:** Store at room temperature and away from direct sunlight.
- **Expiration:** Do not use after the expiration date.
- **Medical conditions:** Consult a healthcare provider before use if you have any underlying medical conditions.

SETUP INSTRUCTIONS

1. **Charge the Battery Pack:** Before initial use, fully charge the included battery pack using the provided charging cable and a compatible power adapter.
2. **Connect the Battery Pack:** Connect the fully charged battery pack to the designated port on the Laser Hair Growth Cap. Ensure a secure connection.
3. **Position the Cap:** Place the Laser Hair Growth Cap on your head. Adjust it to ensure it fits comfortably and covers your entire scalp evenly.

Over 300,000 Reverse Hair Loss

Expect to see noticeable results in 3-6 months

INDIVIDUAL RESULTS MAY VARY

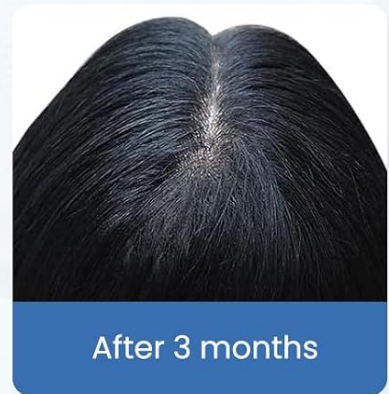
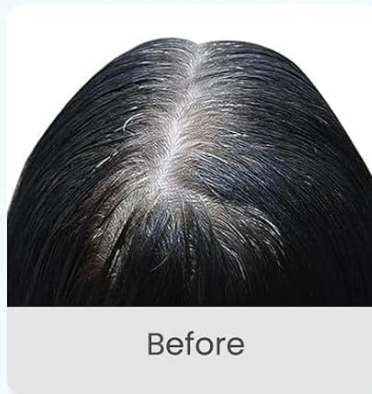
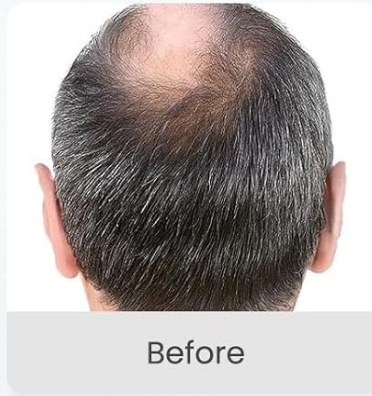


Image: The CWTUL Laser Hair Growth Cap with dimensions, highlighting its adjustable Velcro strap for a comfortable fit.

OPERATING INSTRUCTIONS

1. **Power On:** Once the cap is comfortably positioned and the battery pack is connected, turn on the cap using the power switch.
2. **Select Mode:** The cap features multiple modes for hair health and stress relief. Select the desired mode using the control panel. Available modes typically include:
 - M1: Nourish and promote hair
 - M2: Protect hair
 - M3: Control oil production and nourish hair
 - M4: Promote hair density and growth
3. **Treatment Duration:** Wear the cap for the recommended duration of 10-30 minutes per session. It can be used daily or every other day. The device typically features a 30-minute auto shut-off function.

4. **Power Off and Store:** After the session, turn off the cap and remove it. Store the cap and battery pack in a clean, dry place away from direct sunlight.

SUITABLE FOR MOST PEOPLE TO WEAR



TIPS: there may be 1 inch error due to manual measurement

Image: A user wearing the cap, illustrating the control panel and the four distinct modes for hair health and stress relief.

10-30 Minute Sessions, Daily or Every Other Day



0-1 Month

No noticeable changes or slight shedding of weaker hairs.



2-3 Months

Visible improvement in hair growth.



4-5 Months

Initial reduction in hair loss.



6 Months After

Significant results with continued use.

Image: Examples of multiple use scenarios, demonstrating how the cap can be worn during various daily activities.

Suitable for Men & Women



Alopecia



Pattern baldness



Thinning



Image: A visual timeline detailing expected hair growth progress over 6 months of consistent use.

Your browser does not support the video tag.

Video: An official product video demonstrating the CWTUL Laser Hair Growth Cap and its features.

MAINTENANCE

- **Cleaning:** Gently wipe the cap's exterior with a soft, damp cloth. Do not use harsh chemicals or abrasive cleaners. Ensure the device is completely dry before storage or next use.
- **Storage:** Store the cap in its original packaging or a protective case in a cool, dry place. Avoid extreme temperatures and direct sunlight.
- **Battery Care:** To prolong battery life, avoid fully discharging the battery frequently. Charge it regularly, even if not in constant use.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on	Battery not charged or not connected properly.	Ensure battery pack is fully charged and securely connected to the cap.
Lights are dim or flickering	Low battery charge.	Recharge the battery pack.
Cap feels uncomfortable	Improper adjustment.	Adjust the Velcro strap for a more comfortable and secure fit.

SPECIFICATIONS

- **Model Number:** TW280
- **Laser Diodes:** 80
- **Wavelength:** 650nm
- **Power per Diode:** 41mW
- **Treatment Modes:** 4 (M1-M4)
- **Auto Shut-off:** 30 minutes
- **Power Source:** Rechargeable Battery Pack
- **Package Dimensions:** 9.76 x 8.39 x 2.48 inches
- **Item Weight:** 10.55 ounces
- **UPC:** 701042567574
- **Manufacturer:** CWTUL

WARRANTY INFORMATION

Warranty information for the CWTUL Laser Hair Growth Cap TW280 is not provided in the available product data. Please refer to the product packaging or contact the manufacturer directly for warranty details.

CUSTOMER SUPPORT

For technical assistance, product inquiries, or customer service, please contact CWTUL customer support through their official website or the contact information provided with your purchase.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question

Open the live page

<https://manuals.plus/asin/B0F1N5S5D1>

More documents and user manuals: <https://manuals.plus/>