

[Manuals.plus](#) /

› [NCVI](#) /

› NCVI Electric Breast Pump 8132 User Manual

NCVI 8132

NCVI Electric Breast Pump 8132 User Manual

Model: 8132

1. INTRODUCTION

Thank you for choosing the NCVI Electric Breast Pump 8132. This manual provides essential information for the safe and effective use of your breast pump. Please read all instructions carefully before first use and retain this manual for future reference.

The NCVI Electric Breast Pump 8132 is a portable, hospital-grade double breast pump designed for efficient milk expression. It features 4 pumping modes, 9 suction levels, 5 speed settings, a smart display with a night light, and a 4400mAh rechargeable battery for convenient use at home or on the go. It includes multi-size flanges (17-24mm) to ensure a comfortable fit.

2. IMPORTANT SAFETY INFORMATION

Please observe the following safety precautions to prevent injury or damage to the device:

- Always inspect all pump components before each use for damage or wear. Replace any damaged parts immediately.
- Ensure all parts that come into contact with breast milk are cleaned and sterilized according to the 'Cleaning and Maintenance' section before first use and after each subsequent use.
- Do not immerse the pump motor unit in water or other liquids. Clean the motor unit with a damp cloth only.
- Keep the device out of reach of children.
- Use only the original power adapter and USB cable provided with the pump.
- Do not use the pump while bathing or showering.
- If you experience pain or discomfort during pumping, reduce the suction level or stop pumping. Consult a healthcare professional if discomfort persists.
- This device is for personal use only. Sharing the pump may pose health risks.

3. PACKAGE CONTENTS

Verify that all items are present in your package:

- Pump Motor Unit (x1)
- Flanges (24mm) (x2)

- Pump Body (x2)
- Bottles (x2)
- Replacement Valves (x2)
- Inserts (17mm, 19mm, 21mm) (x2 sets)
- Nipple & Cap (x2)
- Air Tube (x2)
- USB Type-C Cable (x1)
- Nursing Pads (x6)
- Milk Storage Bags (x10)
- User Manual (x1)
- Nipple Ruler (x1)

Find the Pump That Fits Your Lifestyle

Feature	NCVI 8132	NCVI 8102	NCVI 8782
			
Pumping Modes	4	4	4
Suction Levels	9	9	9
Night Light	✓	✓	×
Speed Settings	5 Speeds	—	—
Battery Capacity	4400mAh	4400mAh	2200mAh
Micro-Vibration Comfort Tech	✓	×	×

4. PRODUCT COMPONENTS

Familiarize yourself with the main parts of your NCVI Electric Breast Pump 8132:

- **Pump Motor Unit:** The central unit housing the controls, display, and battery.
- **Flange:** The part that cups the breast and creates a seal for suction. Available in multiple sizes for optimal fit.
- **Pump Body:** Connects the flange to the bottle and contains the valve.
- **Valve:** Regulates milk flow into the bottle and prevents backflow.
- **Bottle:** Collects expressed milk.
- **Air Tube:** Connects the pump body to the motor unit, providing suction.

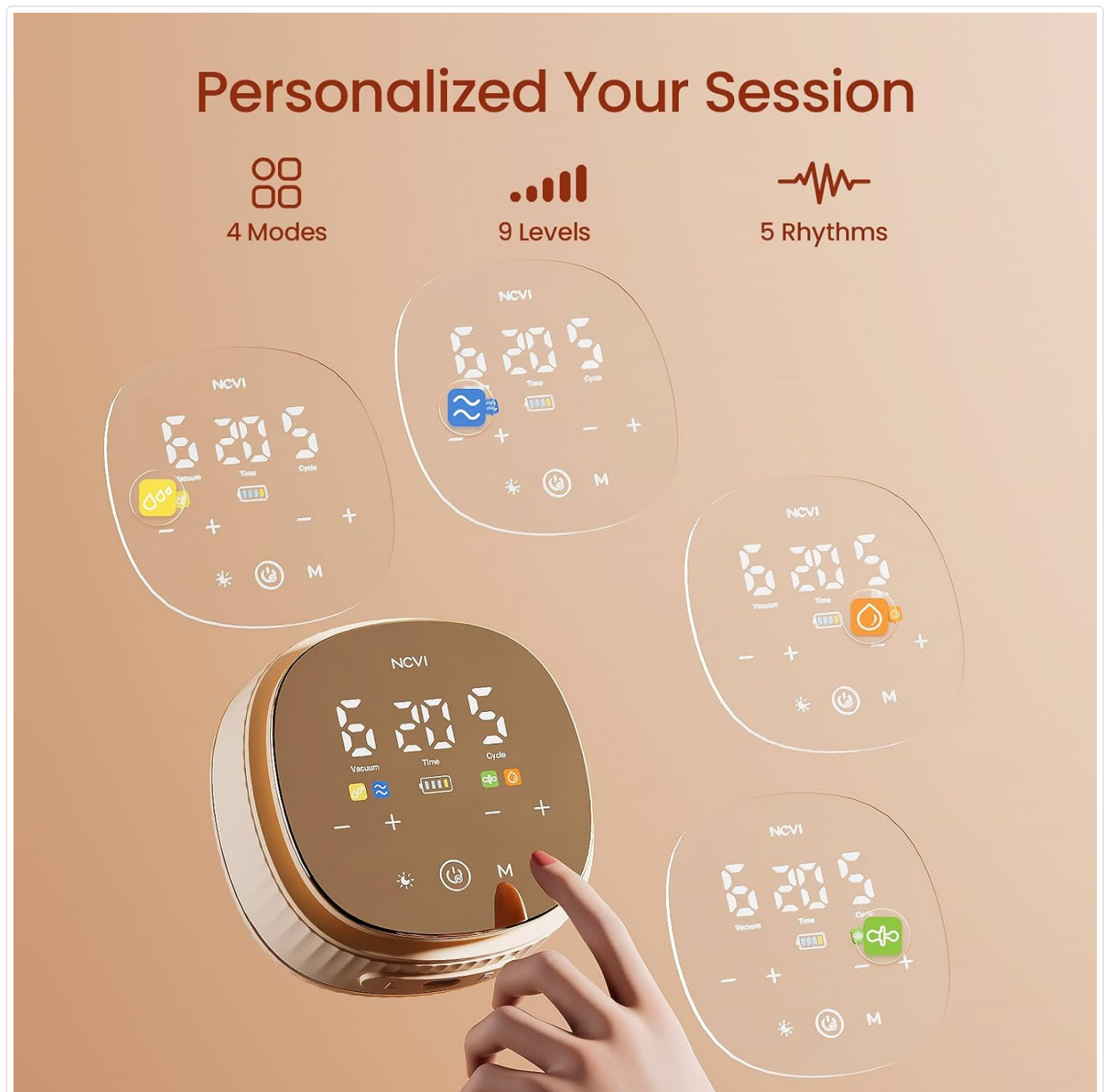


Figure 4.1: Flange sizes and comfortable fit.

5. ASSEMBLY AND SETUP

Follow these steps to assemble your breast pump:

1. **Clean and Sterilize:** Before first use, disassemble all parts that come into contact with breast milk (flanges, pump bodies, valves, bottles, nipples, caps). Wash them with warm soapy water, rinse thoroughly, and sterilize according to manufacturer guidelines (e.g., boiling, steam sterilizer).
2. **Assemble the Pump Body:** Insert the valve into the bottom of the pump body. Ensure it is securely seated.
3. **Attach Flange:** Securely attach the appropriate size flange to the pump body. Use the nipple ruler to determine your correct flange size for optimal comfort and efficiency.
4. **Connect Bottle:** Screw the milk collection bottle onto the bottom of the pump body.
5. **Connect Air Tube:** Attach one end of the air tube to the pump body and the other end to the designated port on the pump motor unit. If using as a double pump, connect both tubes.
6. **Charge Battery:** Before initial use, fully charge the pump motor unit using the provided USB Type-C cable and a compatible power adapter. The battery indicator on the display will show charging status.

Your browser does not support the video tag.

Video 5.1: Instructions on how to install the NCVI Breast Pump 8132 components.

Your browser does not support the video tag.

Video 5.2: An unboxing and first pumping experience guide for the NCVI Breast Pump 8132.

6. OPERATING INSTRUCTIONS

The NCVI Electric Breast Pump 8132 offers customizable settings for a personalized pumping experience.

Power On/Off

Press and hold the power button (O) on the motor unit to turn the pump on or off. The smart display will illuminate.

Modes and Levels

The pump features 4 modes and 9 suction levels, along with 5 speed settings:

- **Multitronic Mode:** Variable frequency to optimize milk expression.
- **Massage Mode:** Gentle stimulation to initiate milk flow.
- **Suck Simulate Mode:** Variable frequency to encourage more milk.
- **Express Mode:** Slow and deep suction to empty the breast efficiently.

Press the 'M' button to cycle through the different modes. Use the '+' and '-' buttons to adjust the suction level (1-9) and speed settings (1-5) within each mode.

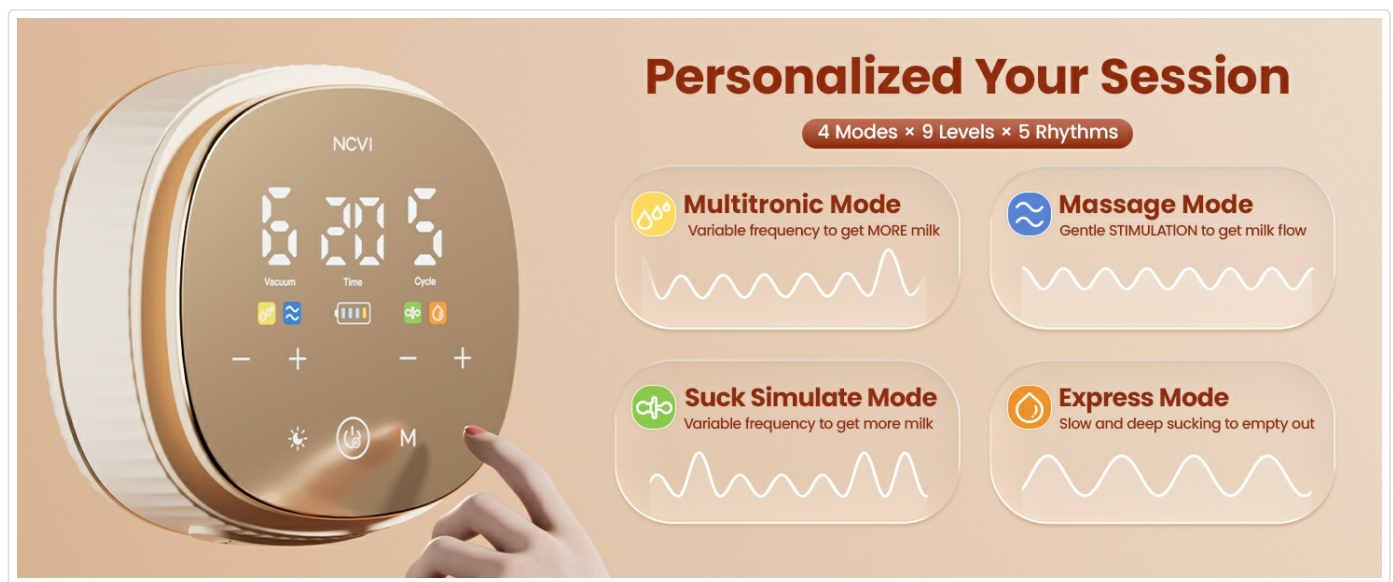


Figure 6.1: Customizable pumping modes and settings.

Night Light Function

The integrated night light provides soft illumination for discreet pumping in low-light conditions. Press the light bulb icon button to turn the night light on or off.

Comfortable Fit



4 Flange Sizes



Soft Silicone



360°Innovative Design

Figure 6.2: Smart night light for convenient use in the dark.

Battery Life

The 4400mAh rechargeable battery provides extended wireless use. A full charge supports approximately 8 pumping sessions or 210 minutes of operation. The battery indicator on the display shows the remaining charge.

Effortless Night Pumping

Built-in warm-glow night light

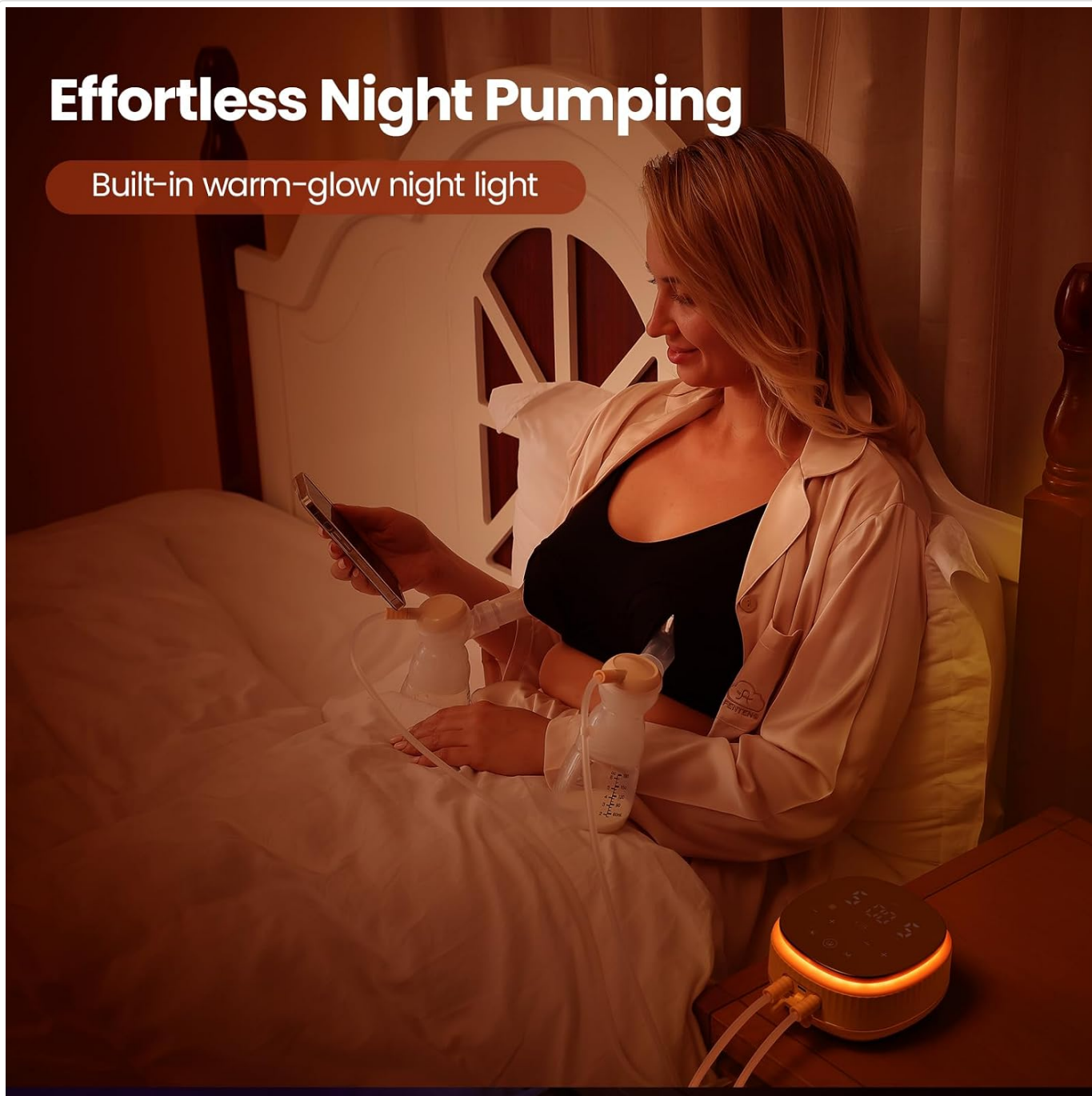


Figure 6.3: Long-lasting battery for portable use.

Your browser does not support the video tag.

Video 6.4: Information on the customizable features and micro-vibration technology of the NCVI Electric Pump 8132.

7. CLEANING AND MAINTENANCE

Proper cleaning and maintenance are crucial for hygiene and the longevity of your breast pump.

After Each Use

1. **Disassemble:** Separate all parts that came into contact with breast milk (flanges, pump bodies, valves, bottles, nipples, caps).
2. **Wash:** Rinse parts with cold water to remove milk residue. Then, wash thoroughly with warm soapy water using a dedicated brush.
3. **Rinse:** Rinse all parts under running water to remove soap residue.

4. **Sterilize:** Sterilize parts by boiling for 5 minutes, using a steam sterilizer, or following your preferred sterilization method.
5. **Dry:** Air dry all components on a clean drying rack or use a clean cloth. Ensure parts are completely dry before storing to prevent mold growth.

Motor Unit Care

The pump motor unit should not be washed or immersed in water. Wipe the exterior of the motor unit with a clean, damp cloth. Do not use abrasive cleaners or solvents.



Figure 7.1: Easy cleaning process for pump parts.

8. TROUBLESHOOTING

If you encounter issues with your NCVI Electric Breast Pump 8132, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
No suction or weak suction	<ul style="list-style-type: none">• Improper assembly• Damaged valve• Air leak in tubing or connections• Low battery	<ul style="list-style-type: none">• Ensure all parts are correctly assembled and sealed.• Check the valve for tears or damage; replace if necessary.• Verify air tubes are securely connected and free from kinks or damage.• Charge the pump motor unit.
Pump not turning on	<ul style="list-style-type: none">• Low or depleted battery• Power button not pressed correctly	<ul style="list-style-type: none">• Charge the pump motor unit.• Press and hold the power button firmly for a few seconds.
Pain or discomfort during pumping	<ul style="list-style-type: none">• Incorrect flange size• Suction level too high	<ul style="list-style-type: none">• Use the nipple ruler to ensure the correct flange size.• Reduce the suction level. If pain persists, stop pumping and consult a healthcare professional.
Milk backflow into tubing	<ul style="list-style-type: none">• Damaged or improperly seated valve• Overfilled bottle	<ul style="list-style-type: none">• Check the valve for damage and ensure it is securely in place. Replace if damaged.• Do not overfill bottles; ensure milk does not reach the valve.

If the problem persists after trying these solutions, please contact NCVI customer support.

9. TECHNICAL SPECIFICATIONS

Feature	Specification
Model Number	8132 (AM8132)
Pumping Modes	4 (Multitronic, Massage, Suck Simulate, Express)
Suction Levels	9
Speed Settings	5
Battery Capacity	4400mAh
Power Source	Battery Powered
Noise Level	40 Decibels
Flange Sizes Included	17mm, 19mm, 21mm, 24mm
Material Type	Polypropylene, BPA Free

Feature	Specification
Capacity (Bottles)	6 ounces
Package Dimensions	9.72 x 6.65 x 6.02 inches
Item Weight	2.98 Pounds

10. WARRANTY AND SUPPORT

The NCVI Electric Breast Pump 8132 is covered by a manufacturer's warranty. Please refer to the warranty card included in your package for specific terms and conditions. For technical support, troubleshooting assistance, or to inquire about replacement parts, please contact NCVI customer service through the contact information provided on the official NCVI website or your purchase platform.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question

Open the live page

<https://manuals.plus/asin/B0F1YJ7VBB>

More documents and user manuals: <https://manuals.plus/>