

Ni-SHEN WD-JY-B002

Ni-SHEN Dog Ultrasonic Teeth Cleaner Instruction Manual

Model: WD-JY-B002

INTRODUCTION

This manual provides instructions for the safe and effective use of your Ni-SHEN Dog Ultrasonic Teeth Cleaner. This device is designed to assist in removing plaque, tartar, and calculus from your pet's teeth using high-frequency vibration technology. Please read this manual thoroughly before first use and retain it for future reference.

IMPORTANT SAFETY INFORMATION

- **Do not use on humans.** This device is specifically designed for pet dental care.
- Keep away from young children to prevent choking hazards from small parts.
- Always supervise the device when charging to prevent electrical hazards.
- Replace the device immediately if the battery becomes damaged or swollen.
- Consult a veterinarian if your pet has existing dental conditions or if you are unsure about using this device.

PACKAGE CONTENTS

Verify that all items are present in your package:

- 1 x Ni-SHEN Ultrasonic Teeth Cleaner Device
- 2 x Replaceable Working Tips (Sharp Head, Flat Head)
- 1 x Dental Probe (double-ended scraper)
- 1 x Mouth Mirror
- 1 x Dental Tweezers
- 1 x USB Charging Cable

DURABLE & IPX8 WATERPROOF



Image: The complete Ni-SHEN Dog Ultrasonic Teeth Cleaner kit, showing the main device, two interchangeable cleaning tips, a mouth mirror, dental tweezers, and a dental probe.

PRODUCT FEATURES

High-Frequency Vibration Technology

The device utilizes advanced high-frequency vibration, reaching up to 44KHZ or 2 million vibrations per minute, to effectively dislodge plaque, tartar, and calculus from pet teeth.



Image: Illustration of the device's high-frequency vibration technology, demonstrating its effectiveness in removing plaque and tartar from a dog's teeth.

Intelligent Recognition Function

The cleaner is designed to vibrate only when its cleaning tip makes contact with hard surfaces like teeth or plaque. It automatically stops vibrating when touching soft tissues such as gums or skin, ensuring safety and preventing discomfort for your pet.

Sharp Head

For plaque and stains on the teeth.



Image: Visual explanation of the intelligent recognition feature, showing the device activating on teeth and deactivating upon contact with gums.

5 Cleaning Modes

The device offers five distinct cleaning modes to suit various needs and pet sensitivities:

- **Soft:** Recommended for first-time users and sensitive pets.
- **Comfort:** Gentle cleaning.
- **Standard:** Regular cleaning.
- **Strong:** For more stubborn plaque.
- **Super Strong:** Maximum power for significant tartar buildup.

5 CLEANING MODES

It is recommend to start with the "Soft" mode for the first use

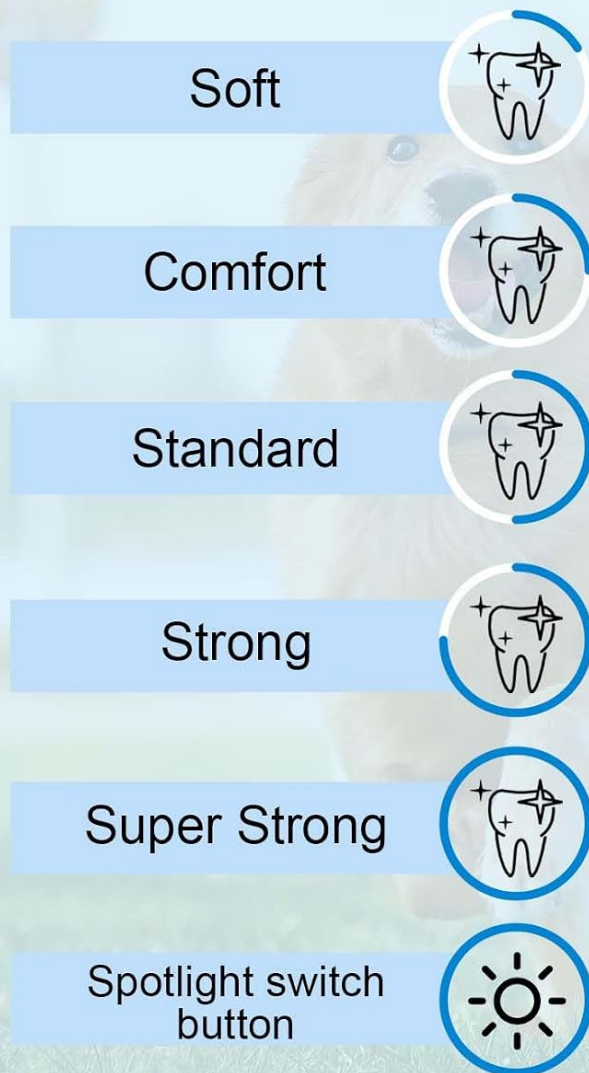


Image: The device displaying its five cleaning modes and the spotlight switch button.

LED Ring Light

Equipped with 5 bright LED ring lights, the device illuminates hard-to-see areas in your pet's mouth, allowing for more thorough cleaning. The lights can be turned on or off independently.

Flat Head

For large tartar on the surface of teeth.



Image: The device's tip with the LED light activated, highlighting its ability to illuminate the cleaning area.

IPX8 Waterproof and USB Rechargeable

The device features an IPX8 waterproof rating, allowing it to be rinsed directly with water for easy cleaning. It is powered by a robust lithium battery that fully charges in 2-3 hours via USB, providing extended working time.



6 SIMPLE STEPS TO USE



1.Tighten the cleaning head before using.



2.Rinse the cleaning head with water to clean the cleaning head.



3.Short press to turn on the device to start working,we recommend using soft mode for the first use.



4.Place the cleaning head on the teeth and slowly move it to clean.



5.Short press the button to switch to the next mode according to your cleaning needs.



6.Long press for 2-3s to turn off the device. Rinse thoroughly with water after use.

Image: The device shown under running water, illustrating its IPX8 waterproof design.

Automatic Shut-off

For safety and to prevent overheating, the device automatically powers off after 10 minutes of continuous operation.

Tartar Scrapor

Reaching into every corner of the mouth to quickly remove tartar.



Image: A visual representation of the device's 10-minute automatic shut-off feature.

COMPONENTS OVERVIEW

The kit includes specialized tools for comprehensive dental care:

- **Sharp Head Tip:** Designed for removing plaque and stains from tooth surfaces.
- **Flat Head Tip:** Ideal for addressing larger tartar deposits on the surface of teeth.
- **Dental Probe:** Helps to remove hard-to-reach food residues between teeth.
- **Mouth Mirror:** Allows visibility of inner and outer areas of teeth that are otherwise difficult to see.
- **Dental Tweezers:** For handling small items or assisting with debris removal.

Dental Probe

Say Goodbye to Pet Oral Disease



Loss of Appetite



Bad Breath



Periodontal Disease



Gingival Hyperemia

Image: A dental probe shown alongside its application for cleaning between a dog's teeth.

Mouth Mirror

Dental Probe

Help to get rid of hard to-reach food residues between the teeth.



Image: A mouth mirror being used to examine a dog's oral cavity.

Sharp Head Tip

Self-contained Safe Mode

Auto power off after 10 minutes of operation to prevent overheating

Open

Auto-off

Image: The sharp head tip shown in use for targeting plaque and stains on a dog's tooth.

Flat Head Tip

IPX8



Image: The flat head tip demonstrated for removing larger tartar deposits from a dog's tooth surface.

SETUP

1. **Initial Charge:** Before first use, fully charge the device using the provided USB cable. A full charge typically takes 2-3 hours.
2. **Attach Cleaning Tip:** Select the appropriate cleaning tip (Sharp Head or Flat Head) and securely attach it to the device. Ensure it is tightened properly.

OPERATING INSTRUCTIONS

Follow these steps for effective and safe use:

1. **Prepare the Device:** Ensure the selected cleaning head is securely tightened onto the device.
2. **Rinse Head:** Briefly rinse the cleaning head with water before use.
3. **Power On and Select Mode:** Short press the power button to turn on the device. For the first use, it is recommended to start with the "Soft" mode. You can cycle through modes by short pressing the mode button.
4. **Begin Cleaning:** Gently place the cleaning head on your pet's teeth where plaque or tartar is present. Move it slowly and carefully. The device will only vibrate when it touches hard surfaces (plaque/tartar), not soft gums.
5. **Adjust Mode (if needed):** If your pet tolerates it, you can short press the mode button to switch to a stronger cleaning mode for more stubborn areas.
6. **Power Off and Clean:** After cleaning, long press the power button for 2-3 seconds to turn off the device. Rinse the cleaning head thoroughly with water to remove any debris.



Image: Step-by-step visual instructions for operating the ultrasonic teeth cleaner.

MAINTENANCE

- **Cleaning the Device:** After each use, rinse the cleaning head and the device body under running water due to its IPX8 waterproof rating. Ensure all debris is removed.
- **Storage:** Store the device and its accessories in a clean, dry place, away from direct sunlight and extreme temperatures.
- **Battery Care:** Recharge the device when the battery indicator is low. Avoid completely draining the battery frequently to prolong its lifespan.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device not vibrating	Cleaning tip not touching a hard surface (teeth/plaque) or device is off.	Ensure the device is powered on. The intelligent recognition function activates vibration only upon contact with hard surfaces. Apply the tip directly to plaque or tartar.
Pet shows discomfort/fear	Sensitivity to sound/vibration or unfamiliarity with the device.	Start with the "Soft" mode. Introduce the device gradually with short sessions. Reward your pet for cooperation. Some pets may require more time to adapt.
Ineffective plaque/tartar removal	Incorrect technique or insufficient power for stubborn buildup.	Ensure the cleaning tip makes direct contact with the plaque. For stubborn areas, try a higher cleaning mode after your pet has adapted to the device. Multiple sessions may be required.
Battery not holding charge	Battery degradation or damage.	If the battery is damaged or swollen, discontinue use and replace the device. Ensure you are using the correct charging cable and power source.

SPECIFICATIONS

Feature	Detail
Model Number	WD-JY-B002
Brand	Ni-SHEN
Power Source	Battery Powered (Lithium Battery)
Charging Time	2-3 hours for full charge
Operation Modes	5 (Soft, Comfort, Standard, Strong, Super Strong)
Waterproof Rating	IPX8
Automatic Shut-off	After 10 minutes
Item Weight	0.36 Kilograms
Handle Material	Plastic

WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the contact details provided with your purchase or visit the official Ni-SHEN website. Keep your purchase receipt as proof of purchase.

© 2026 Ni-SHEN. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question

and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question