

[Manuals.plus](#) /

› [Facemoon](#) /

› Facemoon Cordless Heated Eye Mask Instruction Manual (Model B0F2LQYLCS)

Facemoon B0F2LQYLCS

Facemoon Cordless Heated Eye Mask Instruction Manual

Model: B0F2LQYLCS

PRODUCT OVERVIEW

The Facemoon Cordless Heated Eye Mask is designed to provide gentle, consistent warmth for dry and tired eyes. It features a rechargeable, cordless design for convenient use at home, work, or while traveling. The mask includes a flaxseed-filled nose bridge for improved light blocking and heat retention, three adjustable heat levels, and four timer settings with an automatic shut-off function for safety and comfort. The soft microfiber cover and adjustable strap ensure a comfortable fit.

FACEMOON

Cordless Warm Compress for Dry, Tired Eyes

Soothing warmth without cords,
microwaves, or messy reheating

- ✓ Dry Eyes
- ✓ Screen Fatigue
- ✓ Rechargeable Comfort
- ✓ No Microwave Needed

FSA / HSA
Eligible




- 3 Heat Levels
- 4 Timer Settings
- Auto Shut-Off



Flaxseed Nose Bridge
for Better Fit &
Focused Warmth

Image: The Facemoon Cordless Heated Eye Mask highlighting its key features like cordless design, heat levels, and timer options.

Steady Warm Compress, Made for Daily Eye Care

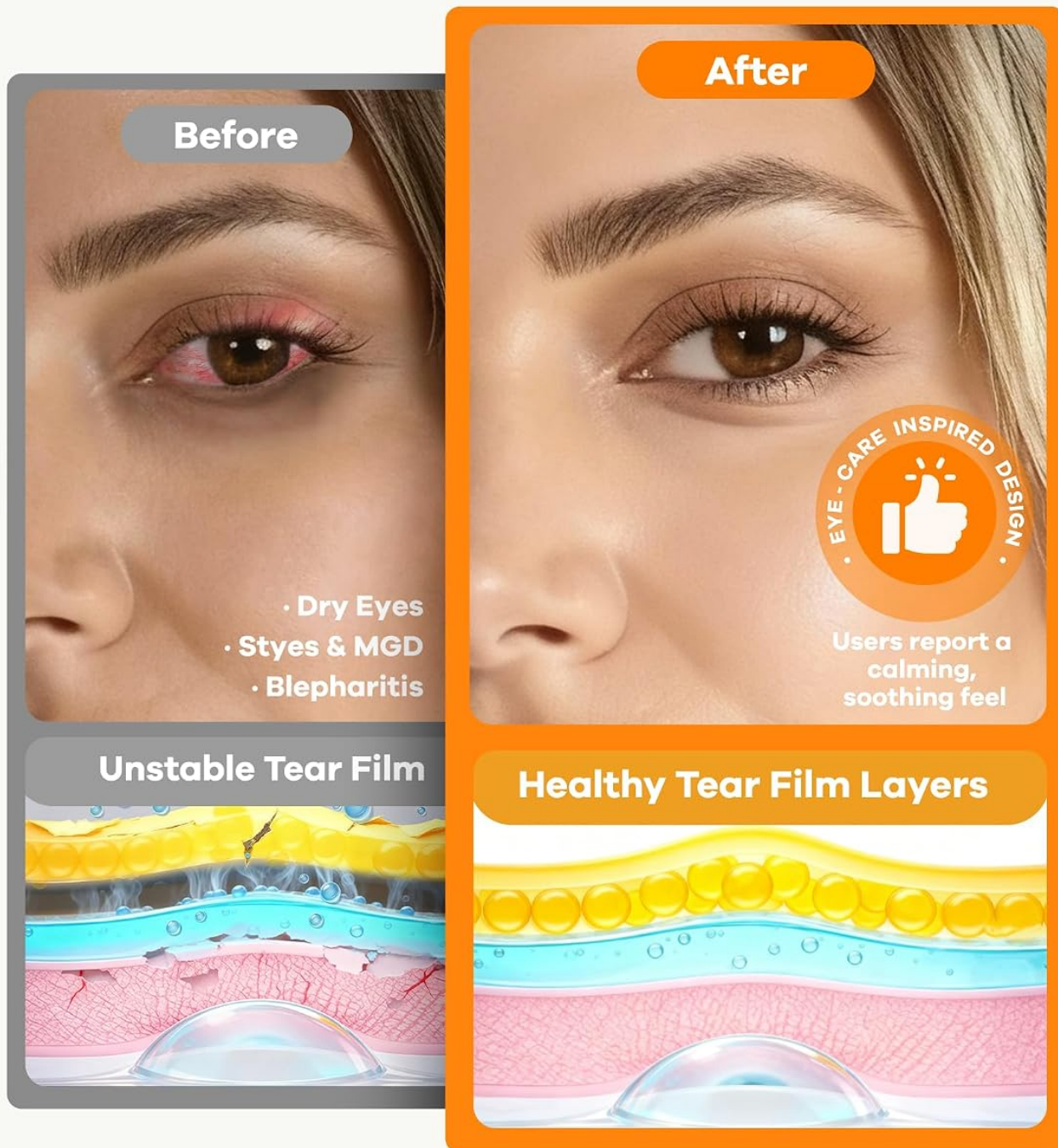

Gentle, consistent warmth helps support a more comfortable routine for dry, tired eyes


Image: Visual representation of the eye mask's effect on dry eyes, showing 'before' and 'after' states and tear film layers.

SETUP INSTRUCTIONS

1. The battery may lose some charge during shipping. Please use the original Type-C cable and a USB adapter under 20W to charge the device. The indicator light will turn on during charging.
2. Unplug the charging cable before turning on the device.

How to Charge

Please Charge Fully before the First Use

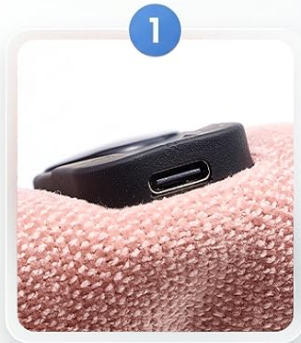

Insert the Type-C Connector.



Connect to a USB Power Supply.



LCD Screen shows charging status.



Fully charged and ready for use.

Image: A visual guide demonstrating the steps to charge the Facemoon Heated Eye Mask using a Type-C cable and USB power supply.

OPERATING INSTRUCTIONS

1. Press and hold the right button for 2 seconds to turn the device on or off.
2. Use the right button to adjust the heat level and the left button to set the timer.

FACEMOON

Choose Your Heat Set Your Time

3 heat levels and 4 timer options help you
create a routine that feels right for you



Auto Shut-Off
for Peace of Mind



FSA / HSA
Eligible



3 Heat Levels



4 Timer Options



Start on low and adjust to your comfort level

Image: The eye mask's control panel showing options for 3 heat levels (104°F, 113°F, 140°F) and 4 timer options (15, 30, 45, 60 minutes).

Video: Official seller video demonstrating the features and use of the heated eye mask for dry eyes.

Video: Official seller video showcasing the newest heated eye mask model and its benefits for dry eyes.

MAINTENANCE

The Facemoon Heated Eye Mask features a detachable cover for easy cleaning. The cover is machine-washable. Always

remove the electronic heating core before washing the cover. Hand wash is recommended for longevity.



Image: The eye mask cover is shown being removed, indicating its washable design for hygiene.

Cordless Comfort, Ready When You Need It



Image: A person demonstrating how to remove the washable cover from the eye mask to keep it fresh after each use.

TROUBLESHOOTING

- **Is the heat strong enough to actually help with dry eyes?**

The warmth is designed to be comfortably warm, not hot. It is sufficient to help loosen oil glands and relieve dryness, as recommended by eye care professionals. Consistent gentle warmth is more effective than intense heat.

- **How to get a 'moist heat' effect with this mask?**

- Why doesn't the battery last long?

Q Is the heat strong enough to actually help with dry eyes?

How to get the 'moist heat' effect with this mask?

Why doesn't the battery last long?

Image: A section detailing Frequently Asked Questions (FQA) regarding the heated eye mask's performance and usage.

Model Number (ASIN)	B0F2LQYLCS
Brand	Facemoon
Package Dimensions	10.24 x 5 x 1.46 inches
Weight	10.86 ounces
Batteries	1 Lithium Ion battery (included)
Heat Levels	3 (Low, Medium, High)
Timer Settings	4 (with auto shut-off)
Material Type	Foam, Hypoallergenic Microfiber
Care Instructions	Hand Wash Only (detachable cover)
Special Features	Cordless, Electric, USB Type-C, Portable, Stable Temperature, Adjustable Strap, Flaxseed Nose Bridge

1. Use the mask for a minimum of 10 minutes to initiate its therapeutic effect. A duration of 20–30 minutes per session is sufficient. It is recommended to use the mask once or twice daily. Prolonged use is not advised.
2. The heat generated by the mask helps to liquefy and release meibomian gland secretions, thereby reducing

inflammation. To avoid the risk of burns, ensure that the temperature is within a comfortable and safe range.

3. For best results, wear the mask while seated or in a slightly reclined position. Lying flat may apply unnecessary pressure to the eyes and is not recommended.
4. Temporary blurred vision may occur after use due to increased oil secretion stimulated by the heat. This is a normal and short-term response and should resolve on its own.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

WARRANTY INFORMATION

Facemoon offers a **1-Year Hassle-Free Replacement** for this heated eye mask. If any issue arises with materials or workmanship within 12 months of purchase, please contact customer support for a replacement or refund.

CUSTOMER SUPPORT

Our friendly support team is available to assist with any questions or concerns regarding your Facemoon heated eye mask. Please refer to the product packaging or the official Facemoon website for contact details.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Email

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Details

Submit question

Open the live page

<https://manuals.plus/asin/B0F2LQYLCS>

More documents and user manuals: <https://manuals.plus/>