

[Manuals.plus](#) /

› [WAFLOO](#) /

› WAFLOO Y826 Portable Water Flosser Instruction Manual

WAFLOO Y826

WAFLOO Y826 Portable Water Flosser Instruction Manual

Model: Y826 | Brand: WAFLOO

INTRODUCTION

The WAFLOO Y826 Portable Water Flosser is designed for effective and convenient oral hygiene, especially for travel. This ultra-compact and lightweight device is engineered to provide dentist-level cleaning, removing plaque and improving gum health. Its telescopic water tank and multiple cleaning modes make it a versatile tool for daily oral care, orthodontic users, and those on the go.

WHAT'S IN THE BOX

- 1 x WAFLOO Y826 Portable Water Flosser
- 4 x Multi-Functional Nozzles (2 Standard, 1 Orthodontic, 1 Tongue Cleaner)
- 1 x USB-C Charging Cable
- 1 x Travel Pouch
- 1 x User Manual



Image: The WAFLOO Y826 Portable Water Flosser, along with its four interchangeable nozzles, USB-C charging cable, and travel pouch.

SETUP

1. **Charge the Device:** Before first use, fully charge the water flosser using the provided USB-C cable. A full charge

takes approximately 3 hours and provides up to 30 days of use. The intelligent pressure compensation ensures steady water pressure even when the battery is low.

2. **Attach Nozzle:** Select the desired nozzle (2 standard, 1 orthodontic, or 1 tongue cleaner). Insert the nozzle into the opening at the top of the device until it clicks into place. The nozzle can rotate 360° for comprehensive cleaning.
3. **Fill Water Tank:** Pull down the telescopic water tank to extend it. Open the wide-mouth cap on the tank and fill it with lukewarm water. The detachable design makes refilling and cleaning easy.

Your browser does not support the video tag.

Video: A demonstration of unboxing, assembling, and initial setup of the WAFLOO Travel Water Flosser.



Image: Visual steps showing how to insert the tip, pull down the tank, add water, choose a mode, start flossing, and store the tip.

OPERATING INSTRUCTIONS

1. **Power On/Off:** Press the power button (usually indicated by a power symbol) to turn the device on or off.
2. **Select Cleaning Mode:** Press the 'M' button to cycle through the 5 available cleaning modes: **Child** (30 PSI), **Soft** (60 PSI), **Normal** (106 PSI), **Pulse** (80–120 PSI), and **Strong** (130 PSI). The smart memory function will remember your last used mode.
3. **Begin Flossing:** Lean over a sink and place the nozzle tip into your mouth, aiming it at your gumline. Turn on the device and gently guide the water stream along your gumline and between your teeth.
4. **Battery Level Display:** Long-press the 'M' button to check the battery level. 5 mode lights indicate 100% battery, 4 lights for 80%, and so on. A red light indicates low battery.

5 Most Popular Cleaning Modes



Strong – 130 psi

Deep cleaning for a truly refreshed mouth.



Normal – 106 psi

Daily dental care made easy and effective.



Soft – 60 psi

Designed for sensitive teeth and a softer touch.



Pulse – 80 to 120 psi

Dynamic pressure to massage gums and protect enamel.



Child – 30 psi

Safe and gentle for children and first-time users.



Image: The WAFLOO Water Flosser displaying its 5 cleaning modes: Strong, Normal, Soft, Pulse, and Child, with corresponding PSI levels.

Battery Level Display

Press & hold M button to check power

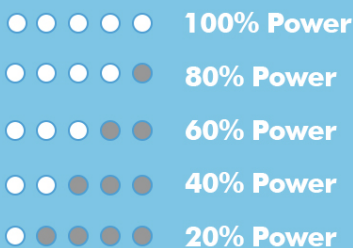


Image: The WAFLOO Water Flosser illustrating its battery level display, showing power percentages via illuminated mode lights.

MAINTENANCE

- **Cleaning the Tank:** The detachable wide-mouth tank is easy to clean. After each use, rinse the tank with water and allow it to air dry completely to prevent bacterial growth.
- **Cleaning Nozzles:** Rinse the nozzles thoroughly after each use.
- **Nozzle Replacement:** It is suggested to replace the nozzle every three months for optimal hygiene and performance.
- **Upgraded Filter:** The device features an upgraded filter that efficiently removes impurities, protecting both the flosser and your mouth.
- **Nozzle Storage:** The built-in storage at the top keeps standard nozzles organized, clean, and safe during travel.

Upgraded Water Filter

Cleaner Water, Safer Flossing



Efficiently filters
impurities



Protects device
and nozzles

Image: Close-up of the WAFLOO Water Flosser's upgraded water filter, highlighting its function in providing cleaner water for flossing.

Convenient Nozzle Storage

Keeps your tips clean,
organized, and ready to use

Note:
Storage compartment is compatible with
standard nozzles only

Image: The WAFLOO Water Flosser demonstrating its convenient nozzle storage feature, keeping tips clean and organized.

SAFETY PRECAUTIONS

- Do not use while sleeping or if you have a medical condition that makes you prone to falling.
- Keep away from children.
- Do not submerge the charging port in water.
- Always unplug the device and allow it to cool before charging.
- Start from the lower gear and gently press the nozzle against the tooth surface to avoid excessive distance.
- Be sure to drain the water and air dry after each use to prevent bacterial growth.
- Suggest using purified water with the dental irrigator.
- Do not direct the water flow towards the tongue, ears, nose, or other sensitive areas.

SPECIFICATIONS

Model Number	Y826
-----------------	------

Brand	WAFLOO
Package Dimensions	6.46 x 3.31 x 3.23 inches
Item Weight	10.86 ounces
Battery	1 A battery required (included), 30-day life on full charge
Charging Time	3 hours via Type-C port
Waterproof Rating	IPX7
Cleaning Modes	5 (Child, Soft, Normal, Pulse, Strong)
Water Pressure	30 PSI (Child) to 130 PSI (Strong)
Pulse Frequency	1400-1800 times/min
Material Type	Acrylonitrile Butadiene Styrene (ABS)
Additional Features	Easy Clean Water Tank, Travel Bag Included, 2-minute Smart Timer, Memory Function, Power Check Function, Pocket Size Design, Ultra-Portable, Leakproof, Simplified Controls, Low Noiseless, 360° Rotating Nozzle, Removable Water Tank, Quick Change Nozzle

WARRANTY & SUPPORT

The WAFLOO Y826 Portable Water Flosser comes with a**2-year manufacturer's warranty**. We also offer**24-hour after-sales service** and a **free 60-day refund/replacement** policy. For any technical support or customer service inquiries, please refer to the contact information provided in your user manual or visit the WAFLOO brand store.



Image: WAFLOO's commitment to customer satisfaction, detailing 24-hour customer support, 2-year warranty, and lifetime technical support.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question

Open the live page

<https://manuals.plus/asin/B0F7R814SX>

More documents and user manuals: <https://manuals.plus/>