**Saved from manuals.plus.**
**Open the live page for links to other documents and user manuals.**



### Oricle JH-D26-HA

# Oricle TrueFit Rechargeable Behind-The-Ear Hearing Aids (Model JH-D26-HA) Instruction Manual

## Product Overview

Experience clearer hearing with effortless comfort. Designed for those with mild-to-moderate hearing loss, the Oricle TrueFit hearing aids enhance everyday conversations without overwhelming amplification. Featuring advanced 16-channel DSP, 16 EQ band adjustments, and 4 environment-specific modes, these devices deliver crystal-clear, natural sound tailored to your surroundings. They are simple to use, powering on automatically when removed from the charging case. The snug, barely-there fit and LED-charging case make them a smart, comfortable, and convenient choice for daily life.

## Key Features

- **Designed for Mild-to-Moderate Hearing Loss:** Ideal for boosting everyday conversations and enhancing clarity through noise without excessive amplification. Intelligently filters and fine-tunes sound.
- **Clear, Natural Sound:** Equipped with advanced 16-channel DSP digital noise reduction technology, 16 EQ band adjustments, and 4 listening modes tailored for different environments. Reduces background noise and enhances speech clarity.
- **Easy to Use:** Automatic on/off when placing in or removing from the charging case. Manual power on by holding the “–” button for 5 seconds and off for 3 seconds. The charging case features a clear LED display showing the exact battery percentage.
- **Comfortable & Easy to Operate:** Includes 3 ear dome sizes for a secure and comfortable fit, intuitive button controls with 11 volume levels, and 4 adaptive modes.
- **Discreet & Stylish:** Sleek, low-profile design that blends in naturally.

## Package Contents

Each Oricle TrueFit set includes:

- Oricle TrueFit Hearing Aids (Left and Right)
- Compact Charging Case with LED Power Display
- Tube Cleaning Tool
- USB-C Charging Cable
- 3 Pairs of Ear Domes (different sizes for personalized fit)

## Product Components

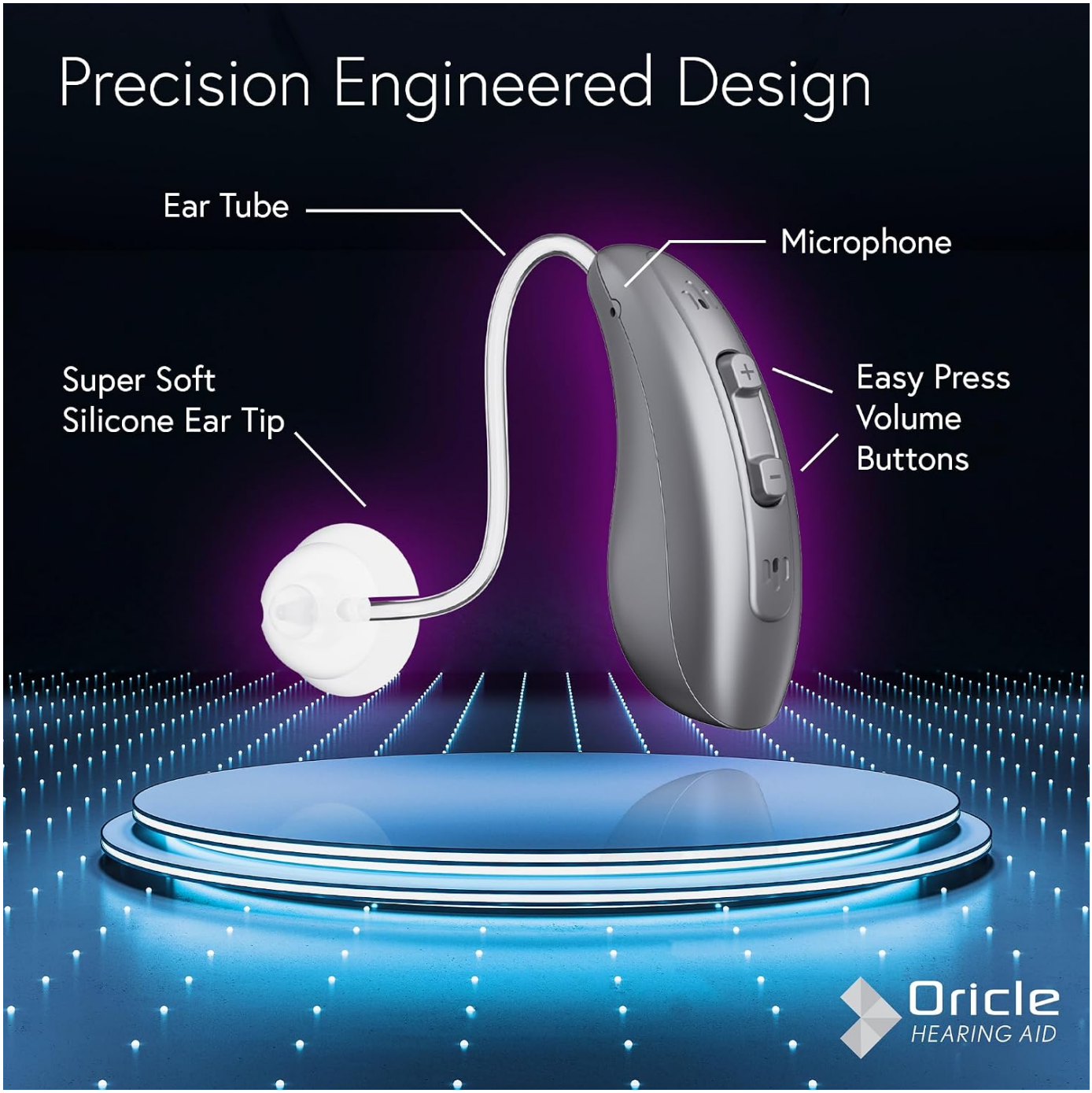



Image: A detailed diagram highlighting the ear tube, microphone, and easy press volume buttons on the Oricle TrueFit hearing aid.

## Setup and First Use

1. **Unboxing:** Carefully remove the hearing aids and charging case from the packaging.
2. **Initial Charge:** Before first use, ensure the hearing aids are fully charged. Place them in the charging case. The charging case features a clear LED display showing the exact battery percentage. Use the provided USB-C cable to charge the case.
3. **Ear Dome Selection:** Select the appropriate ear dome size for a secure and comfortable fit. Three sizes are included.
4. **Power On/Off:** The hearing aids power on automatically when removed from the charging case and power off when placed back in. To manually power on, hold the "–" button for 5 seconds. To manually power off, hold the "–" button for 3 seconds.

## Operating Instructions

### Volume Control

The Oricle TrueFit hearing aids feature intuitive button controls with 11 volume levels. Use the '+' and '-' buttons on the device to adjust the volume to your preference.

### Listening Modes

The device is equipped with 4 adaptive listening modes tailored for different environments. These modes adjust to reduce background noise and enhance speech clarity. Refer to the detailed manual for specific instructions on switching between modes.




Image: Visual examples of environments where Oricle TrueFit hearing aids improve sound, including noisy restaurants, nature walks, watching television, and phone conversations.

## Battery and Charging

The Oricle TrueFit hearing aids utilize a rechargeable lithium-ion battery. The compact charging case provides over 30 hours of usage time and can be fully charged in approximately 1 hour. The LED display on the case indicates the current

battery percentage.




Image: An illustration of the Oricle TrueFit charging case displaying 100% battery, alongside a large battery icon indicating long usage time.

## Maintenance and Care

- **Cleaning:** Use the provided tube cleaning tool to remove any wax buildup or debris from the ear domes and sound tubes. Regular cleaning ensures optimal sound quality.
- **Storage:** Store the hearing aids in their charging case when not in use to protect them and keep them charged.
- **Avoid Moisture:** Keep the hearing aids and charging case dry. Avoid exposure to water or excessive humidity.

## Troubleshooting

If you experience issues with your Oricle TrueFit hearing aids, consider the following:

- **No Sound/Low Sound:** Ensure the devices are fully charged. Check for any blockages in the ear domes or sound tubes and clean them if necessary. Verify the volume level is adjusted appropriately.

- **Feedback/Whistling:** This can occur if the ear dome is not properly sealed in the ear canal. Try a different size ear dome for a better fit. Ensure the device is inserted correctly.
- **Intermittent Sound:** Check battery level. Ensure the devices are clean and free from debris.

For further assistance, please contact Oricle customer support.

## Specifications

| Feature | Detail |
|---|---|
| Model Number | JH-D26-HA |
| Technology | Advanced 16-channel DSP Digital Noise Reduction |
| EQ Band Adjustments | 16 |
| Listening Modes | 4 Environment-Specific Modes |
| Volume Levels | 11 |
| Battery Type | Lithium Ion (included) |
| Charging Port | USB-C |
| Usage Time (per charge) | 30+ Hours |
| Charging Case Dimensions | 4.72 x 4.69 x 3.23 inches |
| Item Weight | 11.68 ounces |

## Warranty and Support

For information regarding product warranty, returns, or technical support, please refer to the official Oricle website or contact Oricle customer service directly. Contact details are typically found in the product packaging or on the manufacturer's website.

## Important Safety Information

- This hearing aid is for adults with signs of mild to moderate hearing loss.
- Do not use if you are younger than 18.
- Consult a doctor if you have any of the problems listed in the product packaging, such as fluid coming out of your ear, sudden hearing loss, or pain.
- Do not use if you have a birth defect or an unusual shape of your ear.
- Keep out of reach of children.


### Ask a question about this manual

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

**Name**

Anonymous

**Email**

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. Privacy details.

**Question**

Example: How do I reset this device or fix this error?

**Details**

Model number, symptoms, what you tried, or the section of the manual you are using.

**Submit question**

---

**Open the live page**

**https://manuals.plus/asin/B0FB15DHMW**

More documents and user manuals: **https://manuals.plus/**