

MLAY T11

MLAY T11 IPL Laser Hair Remover User Manual

Model: T11 | Brand: MLAY

1. INTRODUCTION

The MLAY T11 IPL Laser Hair Remover utilizes advanced Intense Pulsed Light (IPL) technology to provide effective and long-lasting hair reduction for home use. Designed for comfort and ease of use, this device helps achieve up to 96% visible hair reduction after consistent use, offering a convenient alternative to salon treatments.

The MLAY T11 features a maximum energy output of 14.4 Joules (J) and a wavelength range of 510–1200 nm. Its wide treatment area of 3.5 cm² allows for efficient coverage of larger body parts. The device offers 3 adjustable energy levels and 2 operating modes (automatic and manual) to suit various body areas and individual sensitivities.



The MLAY T11 IPL Laser Hair Remover.

Frequência de Tratamento Recomendada



Semanas 1-3
3 vezes por semana



Semanas 4-6
2 vezes por semana



Após 6 semanas
1-2 vezes por mês,
para manutenção






Illustration of IPL technology targeting hair follicles for long-lasting hair reduction.

2. IMPORTANT SAFETY INFORMATION

Please read all safety instructions carefully before using the MLAY T11 device. Failure to follow these instructions may result in injury or device damage.

2.1 Skin Tone and Hair Color Suitability

IPL technology is most effective on lighter skin tones with darker hair. It is not suitable for all skin tones and hair colors. Consult the suitability chart below before use.

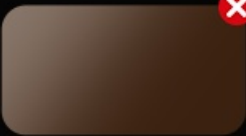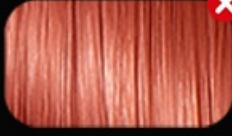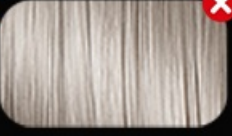	 	 	
	Castanho	Castanho Escuro	
 	 	 	 
Branca	Marfim	Bege	Castanho Claro
 	 	 	 
Preto	Castanho Escuro	Castanho	Dourado
	 	 	 
	Ruivo	Branco	Cinza

2.2 Contraindications

Do not use the MLAY T11 device on the following areas or conditions:

- Tattoos, permanent makeup, dark spots, birthmarks, or dermal fillers.
- Sunburned skin or skin with active allergic reactions.
- Not recommended for pregnant women.
- Do not use near eyes without proper eye protection.



Always wear the provided protective glasses during operation to shield your eyes from the intense pulsed light. The device is not water resistant; do not immerse it in water or use it in wet environments.

3. PACKAGE CONTENTS

Your MLAY T11 IPL Laser Hair Remover package includes:

- MLAY T11A IPL Laser Hair Remover (AC 100-240V)
- Power Adapter
- User Manual
- Protective Glasses
- Shaving Blade
- Cotton Swabs

4. SETUP

Before your first treatment, ensure your skin is clean, dry, and freshly shaved. Do not wax, pluck, or use depilatory creams, as these methods remove the hair follicle, which is necessary for IPL treatment to be effective.

1. **Prepare the Area:** Shave the desired treatment area.
2. **Connect Power:** Plug the power adapter into the device and then into an electrical outlet.
3. **Wear Protection:** Put on the protective glasses provided.



Preparation steps: Shave, power on, and apply treatment.

5. OPERATING INSTRUCTIONS

5.1 Powering On/Off


To turn the device on or off, press and hold the Power button for more than 2 seconds.

5.2 Adjusting Energy Levels


Briefly press the Power button to cycle through the 3 intensity levels. Start with Level 1 and gradually increase as your skin tolerates it. The energy levels are approximately: Level 1 (7.2J), Level 2 (10.8J), Level 3 (14.4J).

Tratamento Seguro para Todo o Corpo




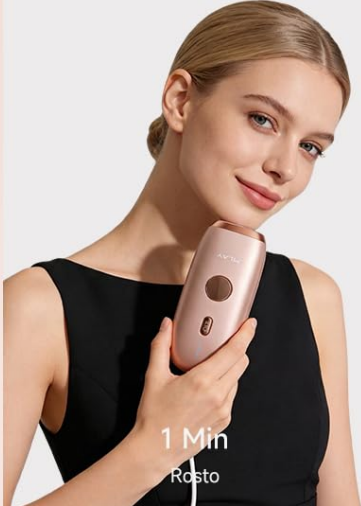


1 Min
Rosto

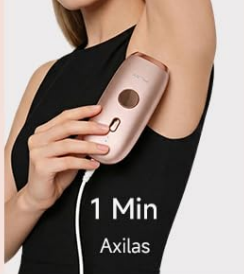

1 Min
Axilas

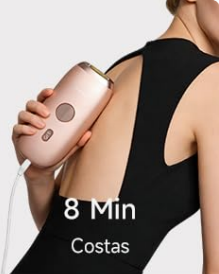

8 Min
Costas

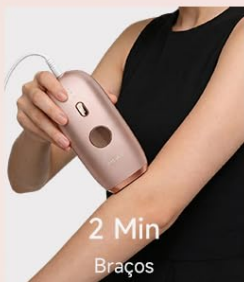

2 Min
Braços



6 Min
Pernas

Recommended treatment frequency and energy levels.

5.3 Selecting Operating Modes

The MLAY T11 offers two modes: Automatic and Manual. Press and hold the flash button for 3 seconds to switch between these modes.

- **Manual Mode:** Ideal for small, sensitive areas like the face, armpits, and bikini line. Press the flash button for each pulse.

- **Automatic Mode:** Suitable for larger areas like legs, arms, and back. Hold the flash button, and the device will automatically emit pulses as you glide it over the skin.

Contraindicações

- Tatuagens, pintas e sardas
- Queimaduras solares
- pele com alergia ativa



Automatic and Manual modes with intensity levels.

5.4 Application

Place the device's treatment window firmly against the skin, ensuring full contact at a 90-degree angle. The device will only emit a flash when it detects full skin contact, indicated by the indicator lights.



Tips for proper device application.

5.5 Recommended Treatment Areas and Times

The MLAY T11 is suitable for use on various body areas. Adjust the energy level according to the area's sensitivity.

- **Face:** Level 1 (approx. 1 minute)
- **Armpits:** Level 3 (approx. 1 minute per armpit)
- **Arms:** Level 2 (approx. 2 minutes per arm)
- **Bikini Line:** Level 3
- **Legs:** Level 2 (approx. 6 minutes per leg)
- **Back:** Level 3 (approx. 8 minutes)



Início



Após 1-3 semanas
2 a 3 vezes por semana




Após 4-6 semanas
1 a 2 vezes por semana



Após 7 semanas
1 vez por semana

*Os resultados podem variar de pessoa para pessoa. Você pode ajustar a frequência de acordo com a sua necessidade e a tolerância da sua pele.

Personalized use for different body areas.

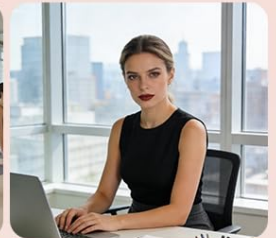

Estimated treatment times for safe whole-body treatment.

5.6 Treatment Frequency

For optimal results, follow the recommended treatment schedule:

- **Weeks 1-3:** 3 times per week.
- **Weeks 4-6:** 2 times per week.
- **After 6 weeks:** 1-2 times per month for maintenance.

Dicas para o uso do dispositivo IPL MLAY

1

Certifique-se de que a combinação entre o tom de pele e a cor dos pelos é adequada para o uso do IPL.

2

Não utilize o aparelho sobre tatuagens, maquiagem permanente, manchas escuras, marcas de nascença, sinais, verrugas ou preenchimentos dérmicos.

3

Durante o uso, mantenha o aparelho em bom contato com a pele, posicionado a 90°, perpendicular à superfície da pele.

Visible hair reduction progression over 4 and 8 weeks.

6. CARE AND MAINTENANCE

Proper care ensures the longevity and performance of your MLAY T11 device:

- **Cleaning:** After each use, wipe the treatment window and the device body with a dry, soft cloth or the provided cotton swabs.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture.
- **Avoid Liquids:** Do not immerse the device in water or use abrasive cleaners, alcohol, or harsh chemicals.

7. TROUBLESHOOTING

If you encounter any issues with your MLAY T11, refer to the following common troubleshooting steps:

- **Device Not Turning On:** Ensure the power adapter is securely connected to both the device and a working electrical outlet.
- **No Flash Emitted:** Verify that the treatment window is in full and firm contact with your skin. The device has a built-in safety sensor that prevents flashing if contact is not complete. Also, confirm that your skin tone and hair color are suitable for IPL treatment as per the suitability chart.
- **Discomfort During Use:** If you experience discomfort or pain, reduce the energy level. If discomfort persists, discontinue use and consult a healthcare professional.
- **Skin Redness/Irritation:** Mild redness is normal after treatment and should subside within a few hours. If redness or irritation is severe or persists, apply a soothing lotion and avoid further treatment on the affected area until it heals.

For issues not covered here, please contact MLAY customer support.

8. TECHNICAL SPECIFICATIONS

- **Model Number:** T11
- **Brand:** MLAY
- **Dimensions:** 21.2 x 12.2 x 0.1 cm
- **Weight:** 860 g

- **Power Source:** AC/DC
- **Voltage:** 100-240 Volts
- **Energy Levels:** 3 (approx. 7.2J, 10.8J, 14.4J)
- **Wavelength:** 510-1200 nm
- **Treatment Area:** 3.5 cm²
- **Water Resistance:** Not water resistant
- **Color:** Orange (as per specifications, actual product may vary)

9. WARRANTY AND CUSTOMER SUPPORT

The MLAY T11 IPL Laser Hair Remover comes with a **12-month warranty** from the date of purchase, covering manufacturing defects.

For any questions, concerns, or technical assistance, please contact MLAY customer support. Our professional support team is available to help you with your device.

© 2026 MLAY. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question