

COMEBY J090102

COMEBY Body Sculpting Machine Model J090102 User Manual

Professional Cordless Cellulite Massage Tool

1. INTRODUCTION

Thank you for choosing the COMEBY Body Sculpting Machine Model J090102. This device is designed to assist in body contouring and skin revitalization through a combination of massage modes, intensity levels, and suction settings. Please read this manual thoroughly before use to ensure safe and effective operation.

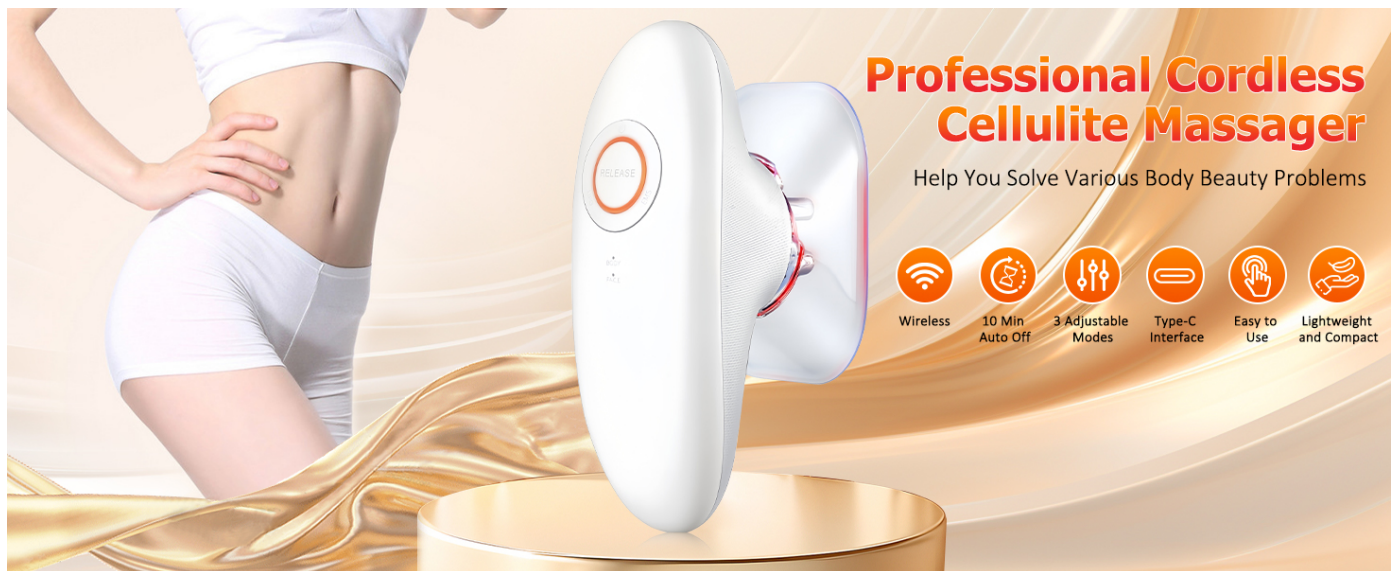


Image 1.1: Overview of the COMEBY Body Sculpting Machine and its intended applications.

2. SAFETY INFORMATION

- Do not use the device on broken, irritated, or infected skin.
- Consult a physician before use if you have any medical conditions, including but not limited to heart conditions, epilepsy, or if you are pregnant.

- Keep out of reach of children.
- Do not immerse the device in water.
- Use only the provided charging cable.
- Discontinue use if you experience discomfort or adverse reactions.
- The device is for external use only.

3. PACKAGE CONTENTS

Verify that all items are present in the package:

- Main Device
- Large Suction Cup (for body)
- Small Suction Cup (for face)
- Type-C Charging Cable
- Filter Elements (10 pieces)
- User Manual

Perfect Body Shape

Works for all skin types



Image 3.1: Included components of the COMEBY Body Sculpting Machine.

4. PRODUCT OVERVIEW

The COMEBY Body Sculpting Machine is a cordless, handheld device featuring multiple functions for skin care. It includes interchangeable suction cups for different body areas.



Image 4.1: Key components and controls of the device.

Device Functions:

- **Vacuum Therapy:** Utilizes suction to stimulate skin.
- **Sculpting Therapy:** Helps to contour and firm the skin.
- **Red-Light Therapy:** Emits red light to support skin revitalization.



Image 4.2: Explanation of the 3-in-1 functions.

5. SETUP

5.1 Charging the Device

Before first use, fully charge the device. Connect the provided Type-C charging cable to the device and a suitable USB power adapter. The indicator light will show charging status and turn off when fully charged.

5.2 Attaching Suction Cups

Select the appropriate suction cup for your desired treatment area. Align the suction cup with the device head and gently twist until it clicks into place. Ensure it is securely attached before use.

Your browser does not support the video tag.

6. OPERATING INSTRUCTIONS

6.1 Preparation

1. Ensure your skin is clean and dry.
2. Apply a thin layer of your preferred body gel or lotion to the treatment area to ensure smooth gliding of the device. Aloe vera gel is recommended.



Image 6.1: Step-by-step guide for device preparation and use.

6.2 Powering On and Mode Selection

1. Press and hold the power button to turn on the device.
2. Press the mode button to cycle through the 2 available massage modes.
3. Adjust the intensity level (5 levels) and suction setting (3 settings) as desired using the respective controls.
4. The device features a 10-minute auto-shutoff for convenience and safety.

6.3 Application Areas and Techniques

The device is suitable for various body parts. Choose the appropriate suction cup:

- **Small Suction Cup:** Ideal for face, arms, and calves.
- **Large Suction Cup:** Recommended for abdomen, waist, thighs, and glutes.

Ensure the suction cup fully covers the skin to be effective.



Image 6.2: Recommended application areas for the device.

6.4 Specific Mode Usage

- **Face Mode:** Designed to increase skin elasticity. Use on each area for 2-3 minutes. No suction cup needed for face mode.
- **Body Mode:** Promotes local circulation. Use the large suction cup for legs, abdomen, and buttocks. Regular use leads to better outcomes.

How does it work

Deeply activates and exercises muscle movement, and stimulates muscle contraction to achieve the effect of firming skin and improving contour

20-40KPa



Image 6.3: Device in use on the face.

SHAPING A GOOD BODY

Reshape Your Body for 10 Minutes Every Day



Stretch Marks
Removal



Body Sculpting



Body Relaxation



Tightening Skin



Firm and Tone Skin



Face & Neck Lifting



Image 6.4: Device in use on the thigh.

HUMANIZED DESIGN

Small and lightweight, easy to carry



Image 6.5: Device in use on the arm.

3 in 1 Body Sculpting Machine You Deserve


Lights up when turned on 3 in 1



**Vacuum
Therapy**



**Sculpting
Therapy**

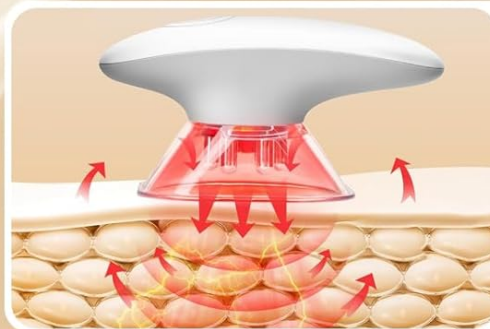

**Red-Light
Therapy**

Image 6.6: Device in use on the abdomen.

7. MAINTENANCE

7.1 Cleaning the Device

After each use, turn off the device and disconnect it from the charger. Wipe the device body and suction cups with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.

7.2 Replacing Filter Elements

The device uses filter elements to maintain hygiene and performance. Replace the filter element regularly, or if you notice reduced suction power. To replace, carefully remove the used filter and insert a new one.

7.3 Storage

Store the device in a cool, dry place away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case when not in use.

8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on	Low battery; Device not charged	Charge the device fully using the provided Type-C cable.
Weak or no suction	Suction cup not properly attached; Filter element clogged; Insufficient gel/lotion on skin	Ensure suction cup is securely attached. Replace the filter element. Apply more body gel or lotion to the skin.
Device stops during use	10-minute auto-shutoff activated; Low battery	Restart the device if desired. Charge the device if the battery is low.
Skin irritation after use	Excessive intensity/suction; Allergic reaction to gel/lotion; Prolonged use on one area	Reduce intensity/suction. Try a different gel/lotion. Ensure even movement and do not stay on one area for too long. Discontinue use if irritation persists.

9. SPECIFICATIONS

Feature	Detail
Model Number	J090102
Product Dimensions	7.8 x 6.5 x 3.1 inches
Item Weight	15.52 ounces
Power Source	Cordless, Type-C Fast Charging
Massage Modes	2
Intensity Levels	5
Suction Settings	3
Auto-Shutoff	10 minutes
Material Type	Plastic
Color	White

10. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the product packaging or contact COMEBY customer service directly. Keep your purchase receipt as proof of purchase for any warranty claims.

For further assistance, please visit the official COMEBY website or contact their support team.

