

BBLOVE AES-U632

BBLOVE Blood Pressure Monitor AES-U632 User Manual

Model: AES-U632

1. IMPORTANT SAFETY INFORMATION

Please read this manual thoroughly before using the BBLOVE Blood Pressure Monitor. Retain this manual for future reference.

- Consult your physician before using this device if you have any medical conditions.
- Do not self-diagnose or self-treat based on measurements from this device. Always consult a healthcare professional.
- Keep the device out of reach of children.
- Do not use the device in a moving vehicle or on an airplane.
- Avoid using the device near strong electromagnetic fields.
- This device is intended for home use and is not a substitute for professional medical advice or treatment.

2. PRODUCT OVERVIEW

The BBLOVE Blood Pressure Monitor AES-U632 is designed for accurate and convenient blood pressure measurement at home. It features a large LED display, voice announcement, and dual-user memory.



Medical-Grade Home Blood Pressure Monitor

Best Solution for At-home Self-monitoring



Medical-Grade Accuracy



One-Click Automatic Measurement



Easy To Use



Image: Overview of BBLOVE Blood Pressure Monitor features compared to other products.

Key Features:

- One-Touch Automatic Measurement
- Dual User Mode with 90 readings storage per user
- Clear LED Display with three-color warning indicator
- Voice Announcement for easy reading
- Rechargeable via Type-C (200+ measurements per charge)
- Adjustable Universal Cuff (8.7-16.5 inches)
- Irregular Heartbeat Detection

One-Touch operation

Easy to use, and all functions are one-touch operations.



One-touch measurement



One-touch user switching



One-touch memory browsing



Image: Display and icon indicators on the BBLOVE Blood Pressure Monitor.

3. SETUP

3.1 What's in the Box

Upon unpacking, ensure all components are present:

Innovative blood pressure monitor



Image: Contents of the product packaging, including the monitor, cuff, and accessories.

- 1 x Blood Pressure Monitor
- 1 x Wide Range Cuff
- 1 x Storage Bag
- 1 x Instruction Manual
- 1 x Charging Cable (Type-C)

3.2 Charging the Device

The monitor is equipped with a rechargeable battery. Use the provided Type-C cable to charge the device. A full charge supports over 200 measurements. The device will indicate low battery when charging is required.

4. OPERATING INSTRUCTIONS

4.1 Preparing for Measurement

For accurate readings, follow these guidelines:

- Sit quietly for at least 5 minutes before measurement.
- Avoid talking or moving during measurement.
- Measure before taking breakfast or medication in the morning.
- Wait at least 1-1.5 minutes between measurements.
- Ensure at least 30 minutes have passed since physical activity.



Image: Visual tips for preparing to take a blood pressure measurement.

4.2 Applying the Cuff

1. Place the cuff directly on the bare skin of your upper left arm. The bottom edge of the cuff should be 2-3 cm above your elbow.
2. Ensure the arm cuff is at the same level as your heart.

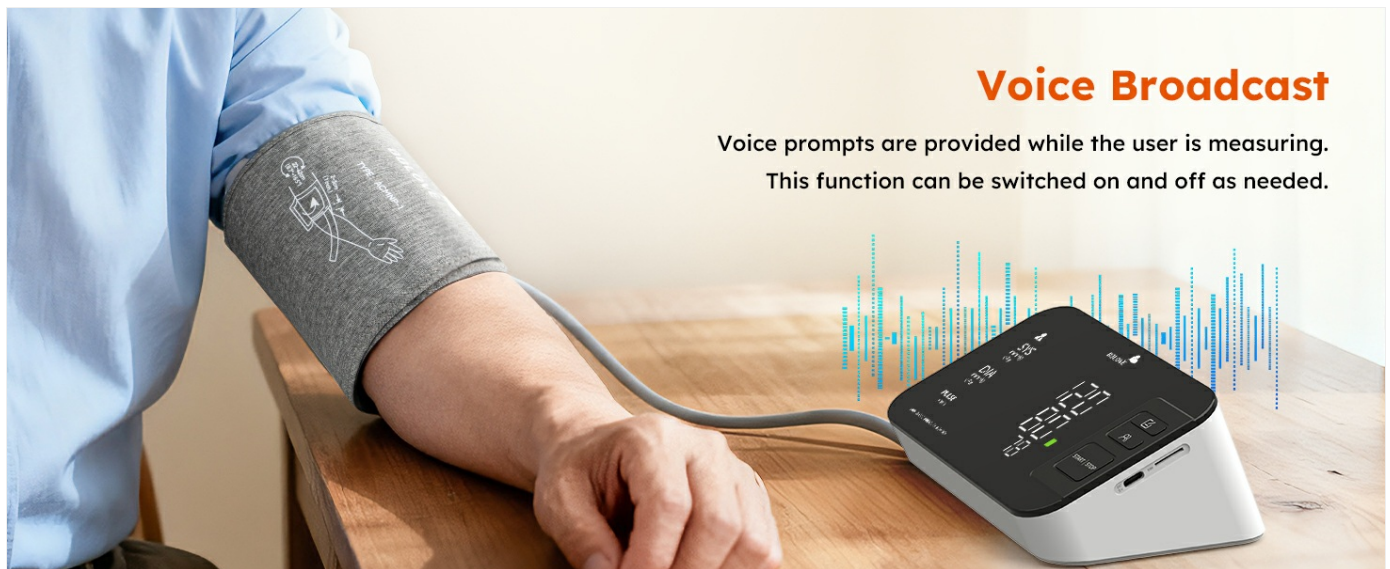


Image: Step-by-step guide on how to correctly position the arm cuff for measurement.

4.3 Taking a Measurement

1. With the cuff correctly applied, press the **START/STOP** button to begin the measurement. Keep your arm in a relaxed position.
2. The cuff will inflate, and the measurement will proceed automatically.
3. Once the device finishes measuring, the results will be displayed on the LED screen.

4.4 Dual User Mode

The device supports two users (User A and User B) and stores 90 readings separately for each. To switch between users, press the user switch button (often indicated by a person icon) before taking a measurement.



Image: Display showing separate readings for User A and User B.

4.5 Silent Mode

To enable or disable the voice announcement:

1. In the power-off state, press and hold the user switch button for three seconds. The display will show "SP ON" (Voice On) or "SP OFF" (Voice Off).
2. Press the user switch button again to toggle between "SP ON" and "SP OFF".
3. Press **START/STOP** to confirm and restart in the selected mode.



Irregular Heart Rate Detection

Monitoring your blood pressure at home is crucial for managing hypertension. Our blood pressure monitor also detects irregular heart rate, which can help prevent heart attacks and strokes. Simple at-home monitoring aids in early diagnosis and follow-up care.



Image: Visual guide for activating or deactivating silent mode.

4.6 Understanding Your Readings

The monitor displays Systolic (SYS), Diastolic (DIA), and Pulse readings. The three-color indicator provides a quick visual reference for your blood pressure category:



Image: Chart explaining the multi-color blood pressure indicator and corresponding ranges.

The device also features Irregular Heartbeat Detection. If an irregular heartbeat is detected, a corresponding icon will appear on the display. Consult your physician if this occurs frequently.



Image: Illustration of the irregular heart rate detection feature.

5. MAINTENANCE

5.1 Cleaning

- Wipe the monitor and cuff with a soft, dry cloth.
- If necessary, use a cloth lightly dampened with water or a mild detergent, then wipe dry.
- Do not use abrasive cleaners or immerse the device in water.

5.2 Storage

- Store the device in the provided storage bag when not in use.
- Keep the monitor in a cool, dry place, away from direct sunlight and extreme temperatures.
- Avoid dropping the device or subjecting it to strong impacts.

6. TROUBLESHOOTING

If you encounter issues with your BBLOVE Blood Pressure Monitor, refer to the following common problems and

solutions:

Problem	Possible Cause	Solution
No display when powered on	Low battery or device not charged	Charge the device using the Type-C cable.
Inaccurate readings	Incorrect cuff placement, movement during measurement, or talking.	Ensure correct cuff placement, remain still and silent during measurement. Refer to section 4.1 and 4.2.
Cuff does not inflate or inflates too little	Air hose disconnected or leak in cuff.	Check that the air hose is securely connected to the monitor. Inspect the cuff for any damage.
Irregular Heartbeat icon appears frequently	Detection of irregular heart rhythm.	This indicates an irregular heartbeat. Consult a healthcare professional for advice.

If the problem persists after attempting these solutions, please contact customer support.

7. SPECIFICATIONS

Feature	Specification
Model Number	AES-U632
Product Dimensions	5.12 x 5.12 x 3.15 inches
Weight	1.34 Pounds
Cuff Size	Universal (8.7-16.5 inches)
Power Source	Rechargeable Battery (Type-C)
Display Type	LED
Memory	2 Users, 90 readings each
Manufacturer	Alicn Medical Shenzhen Inc
UPC	634052119054, 739768508441, 634052123372, 634052122528

8. WARRANTY AND SUPPORT

8.1 Warranty Information

Please refer to the warranty card included with your product for specific warranty terms and conditions. Typically, BBLOVE products come with a limited warranty covering manufacturing defects from the date of purchase.

8.2 Customer Support

For technical assistance, troubleshooting, or warranty claims, please contact BBLOVE customer support through the retailer where the product was purchased or visit the official BBLOVE website for contact information.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question