

[Manuals.plus](#) /

› [Ulike](#) /

› Ulike Air 10 Deluxe IPL Laser Hair Removal Device Instruction Manual

Ulike UI20 DB

Ulike Air 10 Deluxe IPL Laser Hair Removal Device Instruction Manual

Model: UI20 DB

INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your Ulike Air 10 Deluxe IPL Laser Hair Removal Device. Please read all instructions carefully before first use and retain this manual for future reference. The Ulike Air 10 Deluxe is designed for at-home hair reduction, featuring advanced IPL technology with sapphire ice-cooling for a comfortable experience. It includes dual lights and SHR technology to target stubborn hair and a smart sensor for skin tone adaptation.

IMPORTANT SAFETY INFORMATION

For external use only.

Skin Safety Test: IPL technology requires a contrast between the melanin in your hair and the melanin in your skin. It is not safe for all hair and skin types. A Skin Safety Test must be completed before your first treatment session. Refer to the compatibility chart below.

Skin and Hair Color Compatibility

The Ulike Air 10 Deluxe is suitable for specific skin tones and hair colors. Consult the chart to determine if this device is appropriate for your skin and hair type.



Image: A chart illustrating compatible skin tones (Very Fair, Fair, Medium, Olive Light Brown) and hair colors (Black, Dark Brown, Brown, Light Brown) for the Ulike Air 10 Deluxe IPL device, alongside non-recommended types (Brown, Dark Brown skin; Red, Light Blonde, White Gray hair).

Do not use if: Your skin tone or hair color falls under the "Not Recommend" category in the chart. Using the device on incompatible skin or hair types may result in adverse reactions or ineffective treatment.

SETUP

1. **Unpack Contents:** Carefully remove the Ulike Air 10 Deluxe device and all accessories from the deluxe storage bag. The package includes the Air 10 device, a deluxe bag, a razor, and protective glasses.



Image: The Ulike Air 10 Deluxe IPL device being placed into its pink deluxe storage case, alongside other accessories, illustrating the complete kit.

2. **Prepare Skin:** Before each treatment, ensure the area to be treated is clean, dry, and free of lotions or creams. Use the included razor to shave any visible hair in the treatment area. Do not wax or pluck, as IPL targets the hair follicle.

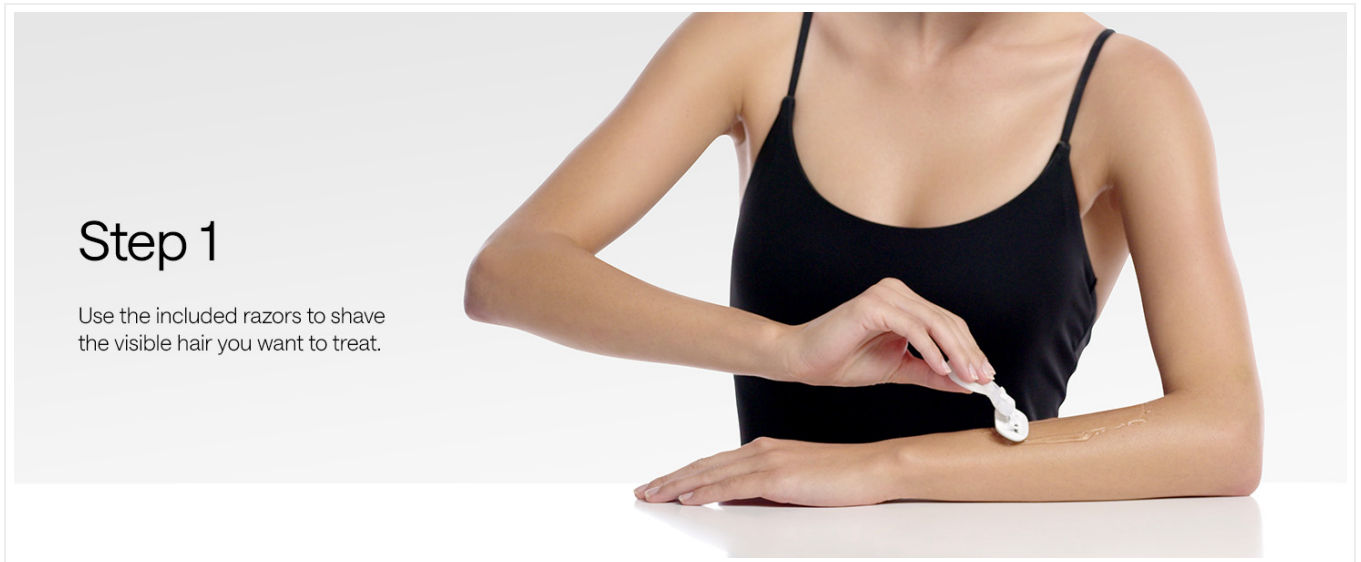


Image: A person shaving visible hair on their arm with a razor, illustrating the first step of preparation for IPL treatment.

3. **Wear Protective Glasses:** Always wear the provided protective glasses during operation to shield your eyes from the intense light pulses.
4. **Connect Power:** Plug the device into a power outlet using the provided adapter. The device is corded electric and requires continuous power during use.

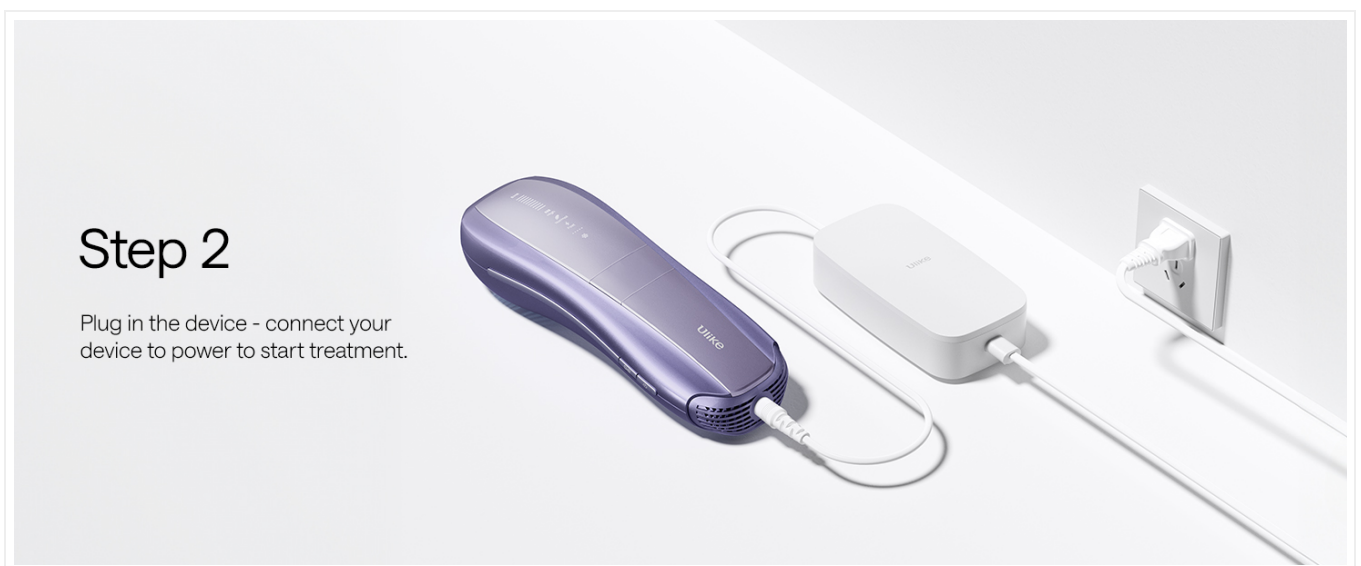


Image: The Ulike Air 10 Deluxe IPL device connected to its power adapter, which is plugged into a wall outlet, demonstrating the power connection step.

OPERATING INSTRUCTIONS

Power On and Mode Selection

- **Power On:** Press and hold the power button to turn on the device.
- **Skin Sensor:** The intelligent skin sensor will automatically detect your skin tone and adjust the energy intensity to a safe and effective level. The device features 10 intensity levels.

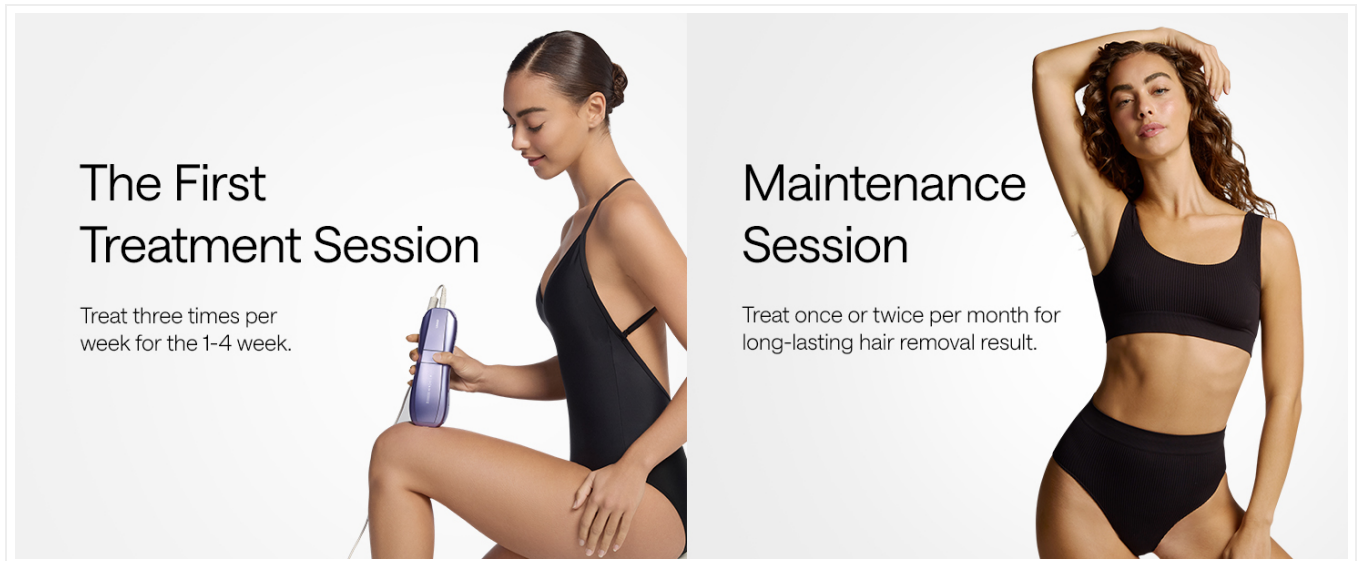


Image: The Ulike Air 10 Deluxe IPL device being used on a person's leg, with an overlay showing skin tone detection and 10 intensity levels, highlighting the SkinSensor feature.

- **Select Mode:** Choose the appropriate treatment mode based on the area and hair type:
 - **Fast Mode:** Delivers double low-energy flashes, suitable for small touch-up areas like the face.
 - **Normal Mode:** Provides a single high-energy flash for finer to normal hair, ideal for arms and legs.
 - **High Mode:** Emits two high-energy flashes for coarser, thicker hair, such as bikini areas.
 - **SHR Mode:** Features dual lamps for double coverage and 26J SHR mode, specifically designed for stubborn, coarse hair.



Image: The Ulike Air 10 Deluxe IPL device on a person's leg, highlighting the SHR mode with indicators for 26 Joules, 4 flashes per session, and 2x power output, designed for stubborn hair.

Treatment Application

- **Position Device:** Place the treatment window firmly against your skin. The device will only flash when full contact is made.
- **Ice-Cooling Technology:** The sapphire ice-cooling feature will activate automatically, cooling the skin to approximately 65°F (18°C) to ensure a comfortable and nearly painless treatment experience.



Image: A close-up of the Ulike Air 10 Deluxe IPL device's treatment window surrounded by ice, illustrating its ice-cooling technology.

- **Flash Application:**
 - **Auto-Glide Mode:** For continuous treatment, hold the flash button down and glide the device slowly over the skin. It will automatically emit light pulses every 0.25 seconds, allowing for quick full-body treatment.
 - **Manual Mode:** For precise treatment of smaller areas, press the flash button once for each pulse.
- **Treatment Schedule:**

- **Initial Treatment:** Treat three times per week for the first 1-4 weeks.
- **Maintenance:** After the initial phase, treat once or twice per month to maintain long-lasting hair reduction results.



Image: A before and after comparison showing significant hair density reduction in the underarm area after 2 weeks of using the Ulike Air 10 Deluxe IPL device.

MAINTENANCE

- **Cleaning:** After each use, ensure the device is unplugged and cool. Gently wipe the treatment window with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the Ulike Air 10 Deluxe device and its accessories in the provided deluxe pink leather case. This keeps the device protected and organized, extending its lifespan. Store in a cool, dry place away from direct sunlight and out of reach of children.

TROUBLESHOOTING

If you encounter any issues with your Ulike Air 10 Deluxe device, please refer to the following common solutions:

- **Device Not Turning On:**
 - Ensure the power adapter is securely plugged into both the device and a working electrical outlet.
 - Check if the power button is pressed and held for a few seconds to activate the device.
- **No Flash Emitted:**
 - Verify that the treatment window is in full and firm contact with your skin. The device has a safety mechanism that prevents flashing if not properly placed.
 - Confirm that your skin tone is compatible with the device (refer to the Skin and Hair Color Compatibility section). The device will not flash on very dark skin tones for safety reasons.
- **Reduced or Slow Results:**
 - Ensure you are following the recommended treatment schedule (three times per week for initial phase). Consistency is key for effective hair reduction.
 - Confirm that the hair in the treatment area was shaved, not waxed or plucked, before the session. IPL targets the hair follicle, which must be present.
 - Adjust the intensity level if comfortable. While the smart sensor adjusts automatically, ensuring proper contact and consistent application can improve results.
 - Hair growth cycles vary. It may take several weeks to see significant reduction. Continue with the recommended treatment plan.
- **Discomfort or Pain During Treatment:**
 - Although the device features ice-cooling technology, some sensitivity may occur. Ensure the device is making full contact to allow the cooling plate to work effectively.
 - If discomfort persists, reduce the intensity level.
 - Ensure your skin is clean and dry. Avoid using on irritated or sensitive skin.

If the issue persists after attempting these troubleshooting steps, please contact Ulike customer support for further assistance.

SPECIFICATIONS

Model Name	UI20 DB
Item Model Number	UI20 DB
Manufacturer Part Number	UI06PR New-Gift
Brand	Ulike
ASIN	B0G39WBP1S
UPC	850067658894
Package Dimensions	11.14 x 6.97 x 6.18 inches
Item Weight	3.81 pounds
Power Source	Corded Electric
Head Type	Dual-Head
Color	Purple
Recommended Uses	Armpit, Face, Leg, Lip
Additional Features	Ergonomic Design, Sapphire Ice-Cooling, Smart Skin Sensor, Dual Lights, SHR Technology

WARRANTY AND SUPPORT

- **Warranty:** Ulike offers a 100-day money-back guarantee and 5 years of sales service for this product. Please retain your proof of purchase for warranty claims.
- **Customer Support:** For any questions, concerns, or technical assistance, please contact Ulike customer support. Refer to the official Ulike website or your product packaging for contact details.

