

Braun PL7432

Braun Skin i-Expert IPL Device PL7432 User Manual

Model: PL7432

1. INTRODUCTION

The Braun Skin i-Expert IPL device is designed for effective and gentle permanent hair reduction at home. This smart IPL device learns and adapts to your skin and hair type, offering a personalized treatment experience. It is clinically tested and dermatologically accredited for safe use.

This manual provides essential information for the safe and efficient use of your device, including setup, operation, and maintenance. Please read it thoroughly before first use.

2. WHAT'S IN THE BOX

Your Braun Skin i-Expert IPL device package includes:

- Braun Skin i-Expert IPL Device with Smart Flex Head
- 3 additional precision heads for various body areas
- Storage pouch
- Mini facial hair remover
- Gillette Venus razor (model may vary)

#BraunLucePulsataSmart

L'alternativa
intelligente al laser.

Il meglio di Braun.
Goditi la libertà di
2 anni di pelle liscia.¹

VOGUE

Miglior dispositivo per
l'epilazione

¹Seguendo il regime raccomandato, i risultati individuali possono variare.

Image: The Braun Skin i-Expert IPL device shown with its various accessories, including different treatment heads, a mini facial hair remover, and a Gillette Venus razor, all neatly arranged.


3. SAFETY INFORMATION AND SKIN SUITABILITY

The Braun Skin i-Expert IPL device is designed for use on suitable skin tones and hair colors. It is crucial to check your compatibility before starting treatment.


Skin Tone and Hair Color Compatibility

The device is effective on light to medium skin tones and natural blonde, brown, or black hair. It is not suitable for very dark skin tones or very light, red, grey, or white hair, as these do not contain enough melanin to absorb the light energy.

Inquadra per scoprire
se l'epilatore a Luce
Pulsata di Braun è
adatto a te:

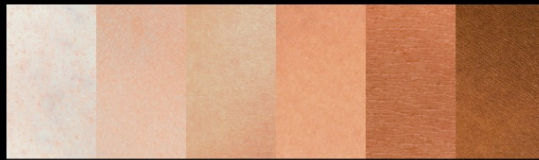

Tonalità della pelle idonee

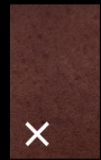

Non idoneo



Colori dei peli idonei



Non idoneo

Image: A chart illustrating suitable skin tones and hair colors for IPL treatment, along with a QR code to check personal compatibility. The chart shows a range of light to medium skin tones and dark hair colors as suitable, while very dark skin and very light hair colors are marked as unsuitable.

For a personalized compatibility check, scan the QR code in the image above or visit tipcheck.pg.com.

Important Safety Guidelines

- Ensure the device is in full contact with the skin for proper function and safety.
- Do not use on tattoos, permanent makeup, dark brown or black spots, freckles, moles, birthmarks, or warts.
- Avoid use on skin with dermal fillers.
- No protective eyewear is necessary during use.

Consigli per l'uso dell'epilatore a Luce pulsata

1

Assicurati che la combinazione di tonalità della pelle e colore dei peli sia adatta all'epilatore a Luce Pulsata.

2

Assicurati che il dispositivo sia a contatto con la pelle poiché il dispositivo funziona solo così, grazie alle nostre caratteristiche di sicurezza.

3

La luce pulsata non si può usare su tatuaggi o trucco permanente, macchie marrone scuro o nero, lentiggini, voglie, nei, verruche e pelle con filler dermici.

→ Non sono necessari occhiali di sicurezza.

Image: A graphic displaying three key usage tips for the IPL device, emphasizing skin tone and hair color compatibility, full skin contact for safety, and areas to avoid such as tattoos and moles.

4. SETUP AND OPERATION

Follow these steps for optimal use of your Braun Skin i-Expert IPL device:

1. **Prepare the Skin:** Ensure the skin is clean, dry, and shaved. Use the included Gillette Venus razor for shaving if needed.
2. **Connect to the App:** Download the free Braun Skin i-Expert app. Connect your device to the app via Bluetooth or Wi-Fi for personalized guidance and treatment tracking.

3. **Select Treatment Area and Mode:** Within the app, select the body area you wish to treat (e.g., lower legs, underarms, face, bikini line). Choose a treatment mode based on your goals: Fastest Results, Balanced, or Skin Comfort.
4. **Skin Compatibility Check:** The app will guide you through a skin test to ensure compatibility. The device's Smart SkinProtect sensor automatically reads your skin tone 80 times per second and adjusts the light intensity for safe and effective treatment.
5. **Attach the Correct Head:** The device comes with a Smart Flex head and additional precision heads. The app will recommend the appropriate head for each body part. Attach the recommended head firmly to the device.
6. **Begin Treatment:** Place the device head flat against your skin. The device will emit light pulses. Glide the device evenly over the skin, ensuring full contact. The app provides real-time coverage guidance to prevent missed spots.
7. **Moisturize:** After each session, moisturize the treated skin.

Using the Smart Flex Head and Precision Heads

The Smart Flex head adapts to your body's contours for comprehensive coverage, even in difficult areas. Use the smaller precision heads for sensitive or hard-to-reach areas like the face, bikini line, or underarms.


Testina flessibile Smart Flex

Copertura al

100%

grazie alle ali flessibili.



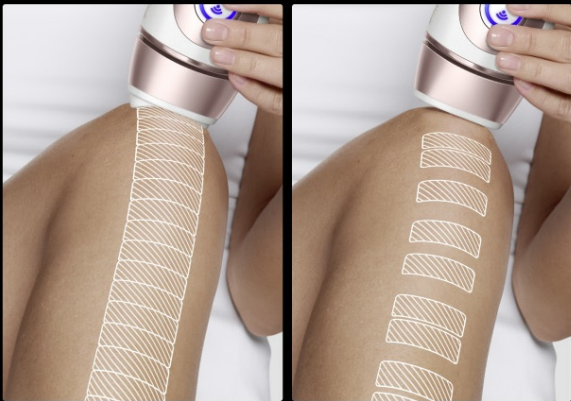

Testina flessibile Smart Flex

Testina precedente

Image: A close-up of the Smart Flex head demonstrating its flexible wings, designed to adapt to body curves for 100% coverage. An illustration shows how the flexible head covers the skin more effectively than a rigid head.

Adatto anche alla zona intima:

 Zona bikini: Monte di Venere e linea bikini.	 Zona situata tra le gambe: Grandi labbra e perineo. Non utilizzare sulle membrane mucose (come le piccole labbra o la vagina).	 Zona situata tra i glutei: Attorno all'ano. Non utilizzare sull'ano.	 Area intima maschile: Pube e zona pantaloncini. Non utilizzare sull'asta del pene e sullo scroto.
---	--	--	---

→ Ginecologicamente testato.

Image: A series of images showing the IPL device with different attachment heads for various body parts, including legs/arms, face/small areas, underarms, chest/back, and intimate areas.

App Features

The Braun Skin i-Expert app acts as your personal guide, offering:

- Personalized treatment plans based on body area, hair density, and desired results.
- Real-time progress tracking and session management.
- Smart rescheduling if you miss a session.
- Intelligent head recognition to recommend the best head for each body part.

Your browser does not support the video tag.

Video: A demonstration of the Braun Skin i-Expert IPL device and its accompanying app. The video shows a user selecting a treatment area, performing a skin test, attaching different heads, and tracking treatment progress in real-time on the app. It highlights the device's ability to adapt to different skin tones and hair types for safe and effective hair removal across various body parts.

5. MAINTENANCE

To ensure the longevity and optimal performance of your Braun Skin i-Expert IPL device, follow these maintenance guidelines:

- **Cleaning:** After each use, ensure the device is unplugged and cool. Wipe the treatment window and the exterior of the device with a dry, lint-free cloth. Do not use abrasive cleaners or solvents.
- **Storage:** Store the device and its accessories in the provided pouch in a cool, dry place, away from direct sunlight and moisture.
- **Inspection:** Regularly inspect the device and power cord for any signs of damage. Do not use the device if it appears damaged.

6. TROUBLESHOOTING

If you encounter issues with your device, refer to the following common troubleshooting tips:

Problem	Possible Cause	Solution
Device not flashing	Not in full contact with skin; skin tone too dark; hair color too light.	Ensure full skin contact. Check skin and hair compatibility using the app or iplcheck.pg.com .
App not connecting to device	Bluetooth/Wi-Fi off; device not paired; app outdated.	Ensure Bluetooth/Wi-Fi is enabled. Follow in-app pairing instructions. Update the app to the latest version.
Skin irritation after use	Incorrect power setting; sensitive skin.	Ensure the device's sensor is correctly reading your skin tone. Use a lower intensity setting or the Skin Comfort mode. Consult a dermatologist if irritation persists.

If the problem persists, please contact Braun customer support.

7. SPECIFICATIONS

- **Model Number:** PL7432
- **Power Source:** Corded Electric

- **Voltage:** 240 Volts
- **Batteries:** 1 AA battery required (included)
- **Water Resistance Level:** Not Water Resistant
- **Country of Origin:** Germany
- **Recommended Uses:** Underarms, Legs (lower and upper), Bikini line
- **Dimensions:** 28 x 19.7 x 9.8 cm; 1.59 kg

8. WARRANTY AND SUPPORT

Your Braun Skin i-Expert IPL device comes with a manufacturer's warranty. Please refer to the warranty card included in your packaging for specific terms and conditions. For further assistance, technical support, or warranty claims, please contact Braun customer service.

Braun also offers an exclusive 'White Gloves' service for personalized assistance and immediate device support, including logistics and appointment scheduling with service centers.

Esclusivo Servizio White Gloves.

Assistenza Personalizzata: Goditi il supporto dedicato del nostro team.

Supporto Immediato al Dispositivo: In caso di problemi, gestiremo tutta la logistica, compresa la programmazione degli appuntamenti con il centro di assistenza.



Image: A promotional image for Braun's 'Exclusive White Gloves Service', detailing personalized assistance and immediate device support. A person is shown using the IPL device on their leg while holding a smartphone, implying connected support.

© 2024 Braun. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for

conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question