



ATIK CAMERA 4120EX Series 4 Color Camera User Manual

[Home](#) » [ATIK CAMERA](#) » ATIK CAMERA 4120EX Series 4 Color Camera User Manual 

Contents

- 1 [ATIK CAMERA 4120EX Series 4 Color Camera](#)
- 2 [Introduction](#)
 - 2.1 [Further information](#)
- 3 [Pack Contents](#)
- 4 [Getting to know your camera](#)
- 5 [Technical Information Summary](#)
- 6 [Declaration of Conformity](#)
- 7 [Servicing and Repairs](#)
- 8 [Warranty](#)
- 9 [Documents / Resources](#)
- 10 [Related Posts](#)



ATIK CAMERA 4120EX Series 4 Color Camera



Introduction

Congratulations on your purchase of an Atik Series 4 Camera. This manual will help you get the most out of your Atik camera so please take the time to read it thoroughly and you'll soon be ready to discover new worlds.

Atik Cameras provide exceptional value for money, superior performance and unparalleled ease of use. They are the result of extensive research and development, each one having been designed and built with the requirements of the most demanding astro imager in mind. Your Atik camera incorporates state-of-the-art design and materials, and will be your trusted astrophotography companion for a long time to come.

Further information

For installation instructions and other useful information please refer to the Quickstart guide which was shipped with the camera. Information on the software is given in the ArtemisCapture guide, a PDF file of which is copied to your computer hard drive by the software installer. Further information, including a number of 'Getting Started' tutorial videos, is available on our website at <http://www.atik-cameras.com>

Pack Contents

This pack includes:



1. Atik Series 4 camera
2. USB cable
3. Car lighter type power cable
4. CD with software and manuals
5. Quickstart Guide

Getting to know your camera

Camera Parts



1. Mini USB port
2. 2.1mm centre-positive 12V DC input
3. Desiccant port

Sensor

The sensors available in the Atik series 4 cameras are listed below:

	Monochrome	Colour
Atik 420	SONY ICX274AL	SONY ICX274AQ
Atik 450	SONY ICX655ALA	SONY ICX655AQA
Atik 414EX	SONY ICX825ALA	SONY ICX825AQA
Atik 428EX	SONY ICX674ALG	SONY ICX674AQG
Atik 460EX	SONY ICX694ALG	SONY ICX694AQG
Atik 490EX	SONY ICX814ALG	SONY ICX814AQG
Atik 4120EX	N/A	SONY ICX834AQG

Due to the very clean nature of these CCDs, dark frames are unnecessary in most cases. This is due to the low noise which stands at an amazing 5 electrons RMS typical.

Optical Window

The optical window used in front of the CCD is a BK7 with BBAR coatings on both sides ensuring that no reflection will appear in your image.

Analog to Digital Converter (ADC)

The Analog to Digital Converter (ADC) is a 16-bit ADC. This means that your Atik camera will allow you to record subtle levels of gray, providing you with enhanced dynamic range when capturing an image.

Power Consumption

Your Atik camera was designed to have low power consumption: please see the table in section 4 for details

specific to your camera model. The 2.1 mm centre-positive DC input is compatible with a wide range of 12V DC supplies.

WARNING: If you have purchased the optional mains power adaptor, please note that it is for indoor or observatory use only. There is a risk of electric shock if the adaptor is used in damp environments or outside. If in doubt, do not use the adaptor and consult a trained electrician.

USB Port

The Atik Series 4 cameras use a USB 2.0 high-speed interface, allowing full-frame image downloads in only a few seconds. A higher speed “preview” mode is also present, which provides roughly twice the speed.

Cooling

All Series 4 cameras have a robust, regulated cooling system, meaning that the CCD is kept at the desired temperature throughout the imaging session. The very low readout noise, combined with efficient cooling, means that dark frames are not necessary in most situations, allowing you to spend more time imaging and less time taking calibration frames. If very long exposures are to be used or high precision data is required (exo-planet searches, for example), dark frame libraries will be easy to maintain.

The cooling delta varies between models and full specifications are available in the Technical Information Summary on page 7. Please note: the cooling delta is how far below the ambient temperature the camera can achieve, rather than a fixed temperature it will cool to.

After switching on the camera, it is advisable to allow up to 5 minutes before taking images in order to allow the temperature to stabilize.

Replacing the desiccant

Your camera includes a high-performance molecular sieve desiccant tablet which is used to avoid condensation in the CCD chamber. Although this desiccant will last a long time it may eventually need replacement or recharging, in which case please follow the instructions below:

- To replace the desiccant, open the desiccant port with the supplied plastic tool and a screwdriver and take out the used tablet.
- To recharge the desiccant, place the tablet in an electric oven at 200°C for 2 hours. Take the tablet out the oven (CAUTION: tablet will be very hot) and allow it to cool down.
- Place the recharged tablet back in the camera, replace the port cover and tighten (no need to over-tighten). Wait 24 hours before connecting the camera again.

You will notice that there is a filter inside the desiccant chamber. The purpose of this is to avoid contamination of the chamber when the desiccant is being replaced. The filter is very fragile and should not be touched.

Technical Information Summary

	420	450	414EX	428EX	460EX	490EX	4120EX
Image sensor	Sony ICX274	Sony ICX655	Sony ICX825	Sony ICX674	Sony ICX694	Sony ICX814/5	Sony ICX834A QG
Resolution	1620×1220	2448×2050	1392×1040	1932×1452	2750×2200	3380×2704	4242×2830

The following table gives the angular resolution per pixel with certain focal distances. The formula to calculate any other focal length is:

$(\text{Pixel Size } (\mu\text{m}) / \text{Focal Distance (mm)}) * 206.3 = \text{angular resolution (arcseconds/pixel)}$

Focal Length	Atik 420	Atik 450	Atik 414EX	Atik 428EX	Atik 460EX	Atik 490EX	Atik 4120EX
250	3,63	2,85	5,32	3,75	3,75	3,04	2,56
300	3,03	2,37	4,44	3,12	3,12	2,54	2,13
350	2,59	2,03	3,80	2,68	2,68	2,17	1,83
400	2,27	1,78	3,32	2,34	2,34	1,90	1,60
450	2,02	1,58	2,99	2,08	2,08	1,69	1,42
500	1,82	1,42	2,66	1,87	1,87	1,52	1,28
550	1,65	1,29	2,42	1,70	1,70	1,38	1,16
600	1,51	1,19	2,22	1,56	1,56	1,27	1,07
650	1,40	1,09	2,05	1,44	1,44	1,17	0,98
700	1,30	1,02	1,90	1,34	1,34	1,09	0,91
750	1,21	0,95	1,77	1,25	1,25	1,01	0,85
800	1,13	0,89	1,66	1,17	1,17	0,95	0,80
850	1,07	0,84	1,57	1,10	1,10	0,90	0,75
900	1,01	0,79	1,48	1,04	1,04	0,85	0,71

Declaration of Conformity

EU Declaration of Conformity



This product carries the CE Mark in accordance with the related European Directive. CE Marking is the responsibility of:

Perseu, SA


R. Dr. Agostinho Neto, 1D
2690-576 Sta Iria da Azoia
Portugal

Critical Applications

This product is not designed for any “critical applications”. “Critical applications” means life support systems, medical applications, connections to medical devices, commercial transportations, nuclear facilities or systems or any other applications where product failure could lead to injury to persons or loss of life or catastrophic property damage.

This product is not a toy.

This is a class A product. In a domestic environment this product may cause radio interference in which case the user may be required to take adequate measures.

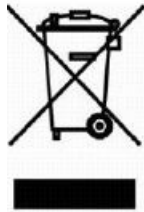

Disposal of the camera

When no longer required do not dispose of this electronic device with general household waste. To minimise pollution and protect the environment the camera should be recycled. Local recycling drop off points available under the Waste from Electrical and Electronic Equipment (WEEE) regulations which will accept the camera. For further information contact Perseu SA at the above address, or the shop from which the camera was bought.

Servicing and Repairs

Repairs, servicing and upgrades are available through your local dealer or by emailing support@atik-cameras.com

Please note that modifications to the camera and/or accessories which are undertaken without the manufacturer's written permission will void the warranty.

Warranty


The equipment is guaranteed against defective design, manufacture or materials for a period of one year from the date of purchase.

This means that Atik Cameras will repair or replace the equipment at its sole option, at no charge to the purchaser for parts or for labour, if the fault is reported within the guarantee period, provided however that Atik Cameras is able to duplicate the defect or problem at its facilities. This warranty does not apply to damage that occurred as a result of abuse or misuse, abnormal service or handling, damage which may have been caused either directly or indirectly by another product, or if the equipment has been altered or modified in any way, or if the damage was caused by repairs or service provided or attempted by anyone other than Atik Cameras. This warranty does not include or provide for incidental or consequential damages.

To exercise your rights under this warranty, you must return the equipment to the dealer from whom it was purchased together with proof of purchase and a clear description of the fault. If it's not possible to return the equipment to your dealer, you should contact Atik Cameras. Equipment returned to Atik Cameras must be sent in appropriate packaging and at your expense (insurance is recommended), together with proof of purchase, a return address and a clear description of the fault.

This does not affect your statutory rights.

Documents / Resources

 <small>Atik Series 4 User Manual Version 1.1 - 10/2019</small>	ATIK CAMERA 4120EX Series 4 Color Camera [pdf] User Manual 4120EX Series 4 Color Camera, 4120EX, Series 4 Color Camera
---	---