



DEXIS Imaging Suite Software Manager Instructions

[Home](#) » [DEXIS](#) » DEXIS Imaging Suite Software Manager Instructions 

Contents

- [1 DEXIS Imaging Suite Software Manage](#)
- [2 DSM – DEXIS Software Manager](#)
- [3 Main User Interface](#)
- [4 Download](#)
- [5 Install](#)
- [6 Download Server](#)
- [7 Events](#)
- [8 Settings](#)
- [9 Documents / Resources](#)
 - [9.1 References](#)
- [10 Related Posts](#)

DEXIS

DEXIS Imaging Suite Software Manage



DSM – DEXIS Software Manager


DEXIS Software Manager is a software product from DEXIS Imaging that allows the dental practice to keep up to date with DEXIS's latest offerings. DSM utilizes cloud technology to provide automated download and installation of software and updates without the need for physical media such as a DVD or USB drives. It is currently packaged with the DEXIS IO Sensor driver download on tech.dexis.com, and with the USB drive that ships with DEXIS Titanium / DEXIS IXS sensors.

Once installed, DEXIS Software Manager examines the current DEXIS software installed on the computer and presents a customized list of software updates and links to other products of interest. Then pick and choose which updates to download and install. DEXIS Software Manager can be configured to run daily so you are aware of new products the same day they are released.

DEXIS Software Manager runs as an application in the system tray. Double click the system tray icon to open the main user interface.

Main User Interface

The left column displays two available functions; Installed and Add Feature.

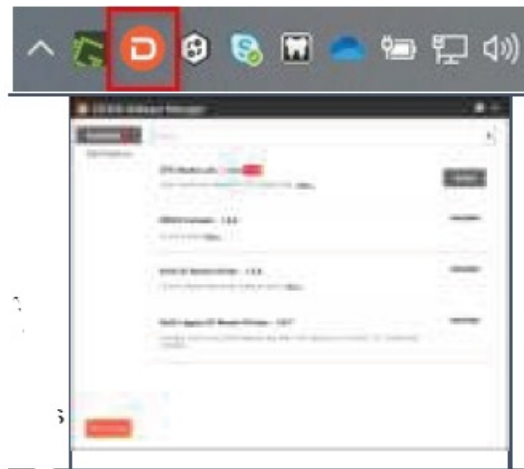

- Installed – lists software that is installed on the system and may have available updates.
- Add Feature – lists new software that is available for download and installation. The middle column displays the name, version and detailed information of the software installers which are available for Installed/Add Features. The right side displays whether the software is Installed or a Modify button is provided.

A Check Now button on the bottom left is for a query of new updates. When a new update is available, the Updates tab highlights with the number of new available installers.

Download



Download progress is displayed as well as the option to Cancel the download.
A notification is displayed when the new software download is complete and ready to install.

Install



When the software is downloaded, the button name will change to Install.
Click the Install button to initiate the installation.

Download Server

The above scenarios require that the DEXIS Software Manager be installed on computers where the software is to be installed and that machines have internet access.

If this is not the desired setup for the office, the office may pick one machine with internet access to install DEXIS Software Manager, and manually install the software to other computers.

Events



DEXIS Software Manager captures logs for events that happen on the machine.
These logs are helpful in determining the health of the system and providing technical support with valuable information. Events can be viewed by selecting the Settings icon Events.
The user interface allows for filtering of these events by various categories.


Settings



- Automatically check for updates – detect and download updates and then copy the installer from the specified download location and manually install on the client machines. DEXIS Software Manager, by default, is configured to run on a schedule to automatically detect, download, and install updates. Daily checks and installation updates can run at a specific time of their choosing.
- Show all product lines – receive updates regarding all DEXIS product lines or filter to only show updates relevant to the current product line installed.
- Automatically download updates – Downloads may be configured as a daily scheduled task. Once an installer is downloaded, it can be run by clicking the Install button on the Main screen. Downloads automatically go to C:\ProgramData\DEXIS Imaging\Software Manager\downloads.

Dental Imaging Technologies Corporation 450 Commerce Dr Quakertown, PA 18951 USA Customer Service: 1-888-883-3947
Technical Support: 1-800-769-2909

Documents / Resources

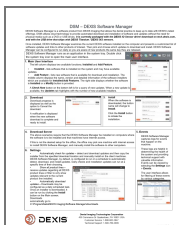

[DEXIS Imaging Suite Software Manager](#) [pdf] Instructions
Imaging Suite Software Manager, Software Manager, Imaging Suite, Manager, Software

References

- [Log in | tech.dexis.com](https://tech.dexis.com)