



KERN PBS Precision Balance User Manual

[Home](#) » [KERN](#) » KERN PBS Precision Balance User Manual 

KERN PBS Precision Balance



Contents

- [1 INTRODUCTION](#)
- [2 Features](#)
- [3 PRODUCT OVERVIEW](#)
- [4 Accessories](#)
- [5 Technical data](#)
- [6 Pictograms](#)
- [7 Documents / Resources](#)
- [8 Related Posts](#)

INTRODUCTION

Multifunctional laboratory balance with single-cell weighing system, also with EC type approval [M]

Features

- KERN PBJ: Internal adjustment in the case of a change in temperature and time-controlled at defined intervals, guarantees high degree of accuracy and makes the balance independent of its location of use
- KERN PBS: Adjusting program CAL for quick setting of the balance accuracy using an external test weight
- Metal housing: robust and sturdy
- Dosage aid
- Weighing with tolerance range (checkweighing): a visual signal helps with portioning, dispensing or grading
- Totalising of weights
- Identification number: 4 digits, printed on calibration protocol freely programmable
- Automatic data output to the PC/printer each time the balance is steady
- 1 Draught shield standard for models with [d] = 0,001 g, weighing space W×D×H 180×193×87 mm
- Protective working cover included with delivery

PRODUCT OVERVIEW



Accessories

- Protective working cover, scope of delivery: 5 items, for models with weighing plate size A KERN PBS-A01S05

B KERN PBS-A02S05

- 2 Set for density determination of liquids and solids with density ≥ 1 , for models with weighing plate size

A KERN PBS-A04

B KERN PBS-A03

- RS-232/Ethernet adapter for connection to an IP-based Ethernet network, KERN YKI-01
- Minimum weight of sample, smallest weight to be weighed, depending on the required process accuracy, only in combination with a DAkkS calibration certificate, KERN 969-103
- Equipment qualification: compliant qualification concept which includes the following validation services: Installation Qualification (IQ), Operating Qualification (OQ), Further details see 208
- Further details, plenty of further accessories and suitable printers see Accessories


3 Single-cell advanced technology:

- Fully automatic manufactured weighing cell from one piece of material
- Stable temperature behaviour
- Short stabilisation time: steady weight values within approx. 3 s under laboratory conditions
- Shock proof construction
- High corner load performance

Technical data

- Large backlit LCD display, digit height 14 mm
- Dimensions weighing surface, Stainless Steel A W×D 108×105 mm B W×D 170×180 mm, see larger picture
- Overall dimensions W×D×H 209×322×78 mm (without draught shield)
- Net weight approx. 3,2 kg
- Permissible ambient temperature 10 °C/30 °C

STANDARD

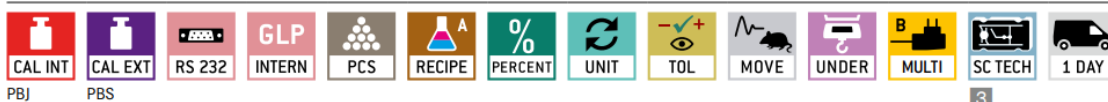

OPTION




FACTORY



PBj

Model KERN	Weighing capacity [Max] g	Readability [d] g	Verification value [e] g	Minimal load [Min] g	Linearity g	Weighing plate		Option			
								Verification		DAkkS Calibr. Certificate	
								KERN		KERN	
PBS 620-3 M	620	0,001	–	–	$\pm 0,002$	A		–		963-103	
PBS 4200-2 M	4200	0,01	–	–	$\pm 0,02$	B		–		963-127	
PBS 6200-2 M	6200	0,01	–	–	$\pm 0,02$	B		–		963-104	

Note: For applications that require verification, please order verification at the same time, initial verification at a later date is not possible.



Verification at the factory, we need to know the full address of the location of use.

PBJ 620-3M	620	0,001	0,01	0,1	$\pm 0,002$	A		965-201		963-103	
PBJ 4200-2 M	4200	0,01	0,1	0,5	$\pm 0,02$	B		965-216		963-127	
PBJ 6200-2 M	6200	0,01	0,1	1	$\pm 0,02$	B		965-202		963-104	
PBJ 8200-1 M	8200	0,1	1	5	$\pm 0,2$	B		965-217		963-128	

Pictograms



• Internal adjusting:

Quick setting up of the balance's accuracy with internal adjusting weight (motordriven)



• Adjusting program CAL:

For quick setting up of the balance's accuracy.

External adjusting weight required



- **Easy Touch:**

Suitable for the connection, data transmission and control through PC or tablet.



- **Memory:**

Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.



- **Alibi memory:**

Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.



- **Data interface RS-232:**

To connect the balance to a printer, PC or network



- **RS-485 data interface:**


To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances.

Network in bus topology is possible



- **USB data interface:**

To connect the balance to a printer, PC or other peripherals

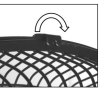

- **Bluetooth* data interface:**

To transfer data from the balance to a printer, PC or other peripherals



- **WiFi data interface:**

To transfer data from the balance to a printer, PC or other peripherals



- **Control outputs (optocoupler, digital I/O):**

To connect relays, signal lamps, valves, etc.



- **Analogue interface:**

to connect a suitable peripheral device for analogue processing of the measurements



- **Interface for second balance:**

For direct connection of a second balance



- **Network interface:**

For connecting the scale to an Ethernet network



- **KERN Communication Protocol (KCP):**

It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems



- **GLP/ISO log:**

The balance displays serial number, user ID, weight, date and time, regardless of a printer connection



- **GLP/ISO log:**

With weight, date and time. Only with KERN printers.



- **Piece counting:**

Reference quantities selectable. Display can be switched from piece to weight



- **Recipe level A:**

The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out



- **Recipe level B:**

Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display



- **Totalising level A:**

The weights of similar items can be added together and the total can be printed out



- **Percentage determination:**

Determining the deviation in % from the target value (100 %)



- **Weighing units:**

Can be switched to e.g. nonmetric units. See balance model. Please refer to KERN's website for more details



- **Weighing with tolerance range:**

(Checkweighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model



- **Hold function:**

(Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value



- **Protection against dust and water splashes IPxx:**

The type of protection is shown in the pictogram.



- **Suspended weighing:**

Load support with hook on the underside of the balance



- **Battery operation:**

Ready for battery operation. The battery type is specified for each device



- **Rechargeable battery pack:**

Rechargeable set



- **Universal plug-in power supply:**

with universal input and optional input socket adapters for A) EU, CH, GB; B) EU, CH, GB, USA; C) EU, CH, GB, USA, AUS



- **Plug-in power supply:**

230V/50Hz in standard version for EU, CH.

On request GB, USA or AUS version available



- **Integrated power supply unit:**

Integrated in balance. 230V/50Hz standard EU.

More standards e.g. GB, USA or AUS on request



- **Weighing principle: Strain gauges:**

Electrical resistor on an elastic deforming body



- **Weighing principle: Tuning fork:**

A resonating body is electromagnetically excited, causing it to oscillate



- **Weighing principle: Electromagnetic force compensation:**

Coil inside a permanent magnet. For the most accurate weighings



- **Weighing principle: Single cell technology:**

Advanced version of the force compensation principle with the highest level of precision



- **Verification possible:**

The time required for verification is specified in the pictogram



- **DAkkS calibration possible (DKD):**

The time required for DAkkS calibration is shown in days in the pictogram



- **Factory calibration (ISO):**

The time required for Factory calibration is shown in days in the pictogram



- **Package shipment:**

The time required for internal shipping preparations is shown in days in the pictogram



- **Pallet shipment:**

The time required for internal shipping preparations is shown in days in the pictogram

KERN – Precision is our business Your KERN specialist dealer:

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg – 2500 kg.

In combination with a DAkkS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAkkS calibration laboratory today is one of the most modern and bestequipped DAkkS calibration laboratories for balances, test weights and force measurement in Europe.

Thanks to the high level of automation, we can carry out DAkkS calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

Range of services:

- DAkkS calibration of balances with a maximum load of up to 50 t
- DAkkS calibration of weights in the range of 1 mg – 2500 kg
- Volume determination and measuring of magnetic susceptibility (magnetic characteristics) for test weights
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAkkS calibration certificates in the following languages DE, EN, FR, IT, ES, NL, PL
- Conformity evaluation and reverification of balances and test weights

Documents / Resources



[KERN PBS Precision Balance \[pdf\] User Manual](#)
PBS Precision Balance, PBS, Precision Balance, Balance

