



Contents [[hide](#)]

- [1 Linshang LS177 Spectrophotometer](#)
- [2 Product Introduction](#)
- [3 Parameters](#)
- [4 Features](#)
- [5 Operation](#)
- [6 Measurement](#)
- [7 Menu bar](#)
- [8 Bluetooth Connection](#)
- [9 Device connection](#)
- [10 PC Software](#)
- [11 Precaution](#)
- [12 Plastic cuvettes](#)
- [13 Packing List](#)
- [14 Documents / Resources](#)
 - [14.1 References](#)



Linshang LS177 Spectrophotometer



- QR code for app installation <https://www.linshangtech.cn/a/170/>
- Please read this manual carefully before using and reserve it for reference.

Product Introduction

The instrument is a multifunctional large aperture spectrophotometer. It has a 7-inch IPS color display and capacity touch panel, which has an excellent operating experience. The instrument integrates multiple color charts and can quickly match the measured color to the closest color code. It also has color difference comparison function, color difference threshold can be set and a variety of color difference formulas can be selected to achieve rapid QC testing. The instrument can also connect to APP for color measurement and sharing. It is also equipped with a powerful PC software to meet various needs of customers. The instrument meets the requirements about grade 1 colorimeter in the Chinese JJG 595-2002 Colorimeters and Color Difference Meters in all directions employing full-spectrum LED light source and advanced spectral sensor.

Standards for the product

- JJG 595-2002 Colorimeters and Color Difference Meters
- GB/T 3978-2008 Standard illuminants and geometric conditions
- GBT 7921-2008 Uniform color space and color difference formula
- GB/T 11186.1-1989 Methods for measuring the colour of paint films—Part 1: Principles
- GB/T 11186.2-1989 Methods for measuring the colour of paint films—Part 2: Colour

measurement

- GB/T 11186.3-1989 Methods for measuring the colour of paint films—Part 3: Calculation of colour differences
- GB/T 39822-2021 Plastics—Determination of yellowness index and change in yellowness index
- GB/T 17749-2008 Methods of whiteness specification
- ASTM E313-98 Standard Practice for Calculating Yellowness and Whiteness Indices from Instrumentally Measured Color Coordinates

Parameters

Illumination geometry	D/8, specular component include (SCI)
Illumination light source	Full spectrum LED light source
Spectral range	400-700nm
Spectral interval	10nm
Measuring aperture	50mm
Measurement conditions	Light source D65, field of view 10°
Measuring time	About 1.5s

Color space	CIE Lab, Luv, LCh, Yxy, CMYK, RGB, WI-98, WI-Gauz, WI-Hunter, YI-98, WI-R457, Reflectance
Color difference formula	ΔE^*ab , ΔE^*uv , ΔE^*94 , $\Delta E^*cmc(2:1)$, $\Delta E^*cmc(1:1)$, $\Delta E^*cmc(1.4:1)$, ΔE^*00
Repeatability	Standard deviation ΔE^*ab is within 0.03 (Measurement condition: the average value of 30 measurements on the whiteboard at an interval of 3s after calibration)
Inter-instrument agreement	$\Delta E^*ab \leq 0.4$, based on avg. of 12 BCRA series II tiles
Dimension	21.1 * 26.4 * 15.5cm (L*W*H)
Weight	About 3.29kg
Display	7-inch 1024*600 dot matrix IPS color screen
Language	Simplified Chinese, English
Data transmission	USB, Bluetooth
Operating temperature range	0~45°C, 0~85%RH (no condensation)
Storage temperature range	-25~55°C, 0~85%RH (no condensation)
Supply Voltage	AC100~277V 50/60HZ
Operating Current	0.4A
Operating Power Consumption	80W

Features

1. The measuring aperture is up to 50mm, which is especially suitable for measuring uneven surfaces such as masterbatches, powders and coarse cloths.
2. The instrument is equipped with special plastic cuvettes for the measurement of masterbatches, powder and opaque liquid.
3. The instrument is fashioned with a full-spectrum LED light source and spectral sensor, and it meets the requirements of the national metrology grade 1 machine.
4. Adopt a 7-inch capacity touch panel with human-machine interaction interface for convenient operation.
5. Contains multiple electronic color charts, which can quickly match the closest color code with matching accuracy rate greater than 90%.
6. QC testing function: the color difference threshold can be set and a variety of color

difference formulas are optional.

7. A variety of color spaces are available, suitable for various needs of color measurement.
8. It can store 1000 colors and 1000 color comparison records with large storage space.
9. It can connect to APP for color sharing.
10. Equip with PC software which has stronger color difference statistical analysis function, spectral reflectance measurement and curve data export EXCEL function.
11. The instrument has a built-in compensation optical path with stable values, allowing long-period calibration.


Operation

Power on/off

- **Power on:** Plug in the power supply, short press the right side button to turn on the instrument.
- **Power off:** Long press the right side button to turn off or press “Power off” in the menu bar.

Calibration

After entering “calibration” interface, you can perform calibration operation as the calibration animation, or skip the calibration. It is recommended to calibrate after a long time of no usage.

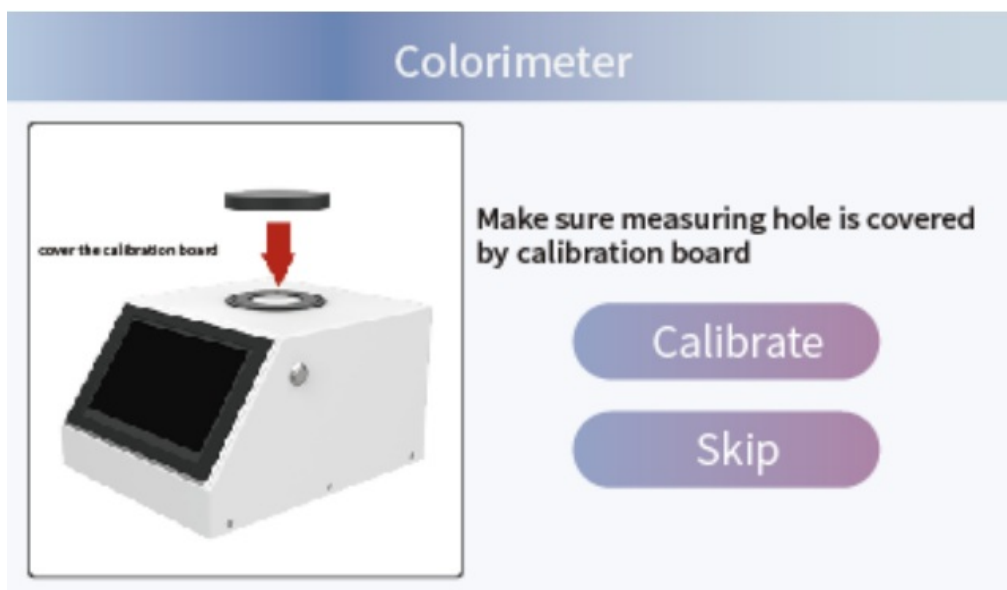

For failed calibration, the possible reasons are as below:


- The standard tile is not closed properly;
- The standard tile needs to clean;
- The attenuation of the light source leads to failure of normal use and it needs to be returned to the factory for inspection and repair.

Measurement


The default measurement parameter of the instrument is “Lab”. Users can click the “Param” at the top of the screen or the “Param Settings” in menu bar to enter the setting interface to select the parameters to be measured (CIE Lab, CIE Luv, CIE LCh, Yxy, CMYK, RGB, WI-98, WI-Gauz, WI-Hunter, YI-98, WI-R457, Spectra).



Entering the measurement interface, you can click the ” Compare” on the top of the screen to switch the compare measurement mode.

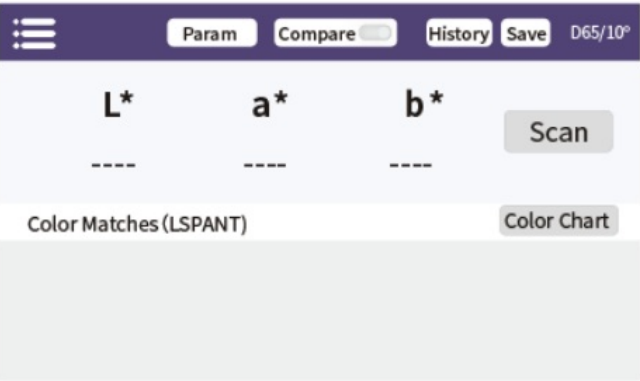

Compare off

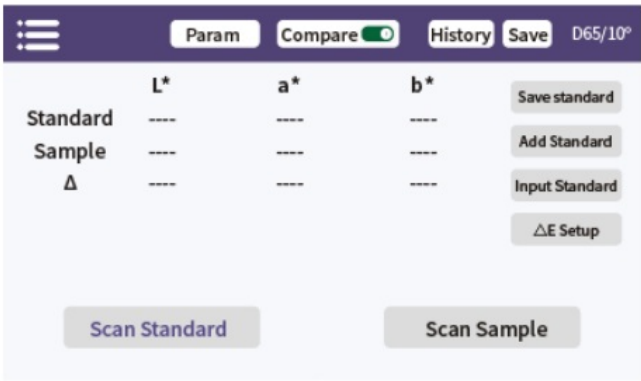

Compare on

Compare on

The screenshot shows a color measurement interface with a purple header bar containing a menu icon, 'Param', 'Compare' (with a green toggle), 'History', 'Save', and 'D65/10°'. Below the header is a table with columns L*, a*, and b*. The 'Standard' row shows values 38.84, 2.05, and 2.21. The 'Sample' row shows values 37.68, 2.37, and 3.19. The difference row (Δ) shows values -1.16, 0.32, and 0.98. Below the table, there are color difference indicators: 'Slight black' for ΔL, 'Slight red' for Δa, and 'Slight yellow' for Δb. A green bar displays 'ΔE*ab=0.01', 'PASS', and '(ΔE≤5.0)'. To the right of the table are buttons for 'Save standard', 'Add Standard', 'Input Standard', and 'ΔE Setup'. At the bottom are two large buttons: 'Scan Standard' and 'Scan Sample'.

	L*	a*	b*
Standard	38.84	2.05	2.21
Sample	37.68	2.37	3.19
Δ	-1.16	0.32	0.98

Slight black Slight red Slight yellow

ΔE*ab=0.01 PASS (ΔE≤5.0)

Scan Standard Scan Sample

Compare on

In the measurement interface, press the “Scan Standard” button to set current value as standard. Press the “Scan Sample” button to use current value as sample. If short press the button on the right side, another measurement will also be taken and the current value will replace the previous one. When measuring or adding a standard and measuring a sample, the color difference between the two measured values will be calculated with the following results:

- **The difference value**

The difference Δ is the sample data minus the standard data;

- **Color Bias**

When the measurement parameter is selected as Lab, the color bias will be judged by ΔL , Δa , and Δb ;

- **Δ E value**

When the measurement parameter is selected as Lab LCh Luv Yxy CMYK or RGB, ΔE is calculated using the ΔE^*ab color difference formula;

- **QC test**

Suitable for color difference inspection in production process. If ΔE is higher than the set threshold, it will read “NG” and shown in red; if ΔE is less than or equal to the set threshold, it will read “Pass” and display in green.

	L*	a*	b*
Standard	74.26	-0.21	8.60
Sample	74.31	-0.18	0.31
Δ	0.05	0.03	0.31
Normal	Normal	Normal	More yellow
ΔE*ab=0.32 PASS (ΔE≤5.0)			

Pass

	L*	a*	b*
Standard	58.22	25.95	18.72
Sample	58.52	27.22	20.32
Δ	0.30	1.27	1.60
Slight white	Slight Red	Slight Yellow	
ΔE*ab=2.06 NG (ΔE≤2.0)			

NG

The following steps can be performed after completing compare measurement:

Save

To save the comparison result, press “Save” on the upper right corner. You can choose the default or name the record manually.

History

If users need to view the saved results, click “History” to view the records, and records can be paged, deleted, renamed and searched.

Save standard

After the user has finished measuring the standard value, if needed to save it for later calling, click the ” Save standard” to save it, and the saved value can be queried in the add standard interface.

Add standard

After clicking “Add standard”, the interface will jump to “Add standard”, click the value you want to add, and the selected value will be used as the standard, and will automatically return to the measurement interface to compare with the current sample or the upcoming sample.

Input standard


The user can manually input the Lab, Luv, LCh, Yxy, CMYK, RGB, WI-98, WI-Gauz, WI-Hunter, YI-98, and WI-R457 values as standard.

ΔE Setup

Setting the ΔE formula and threshold for QC inspection.

Compare off

In the measurement interface, click on “Scan” or short-press the button on the right side of the instrument to measure.

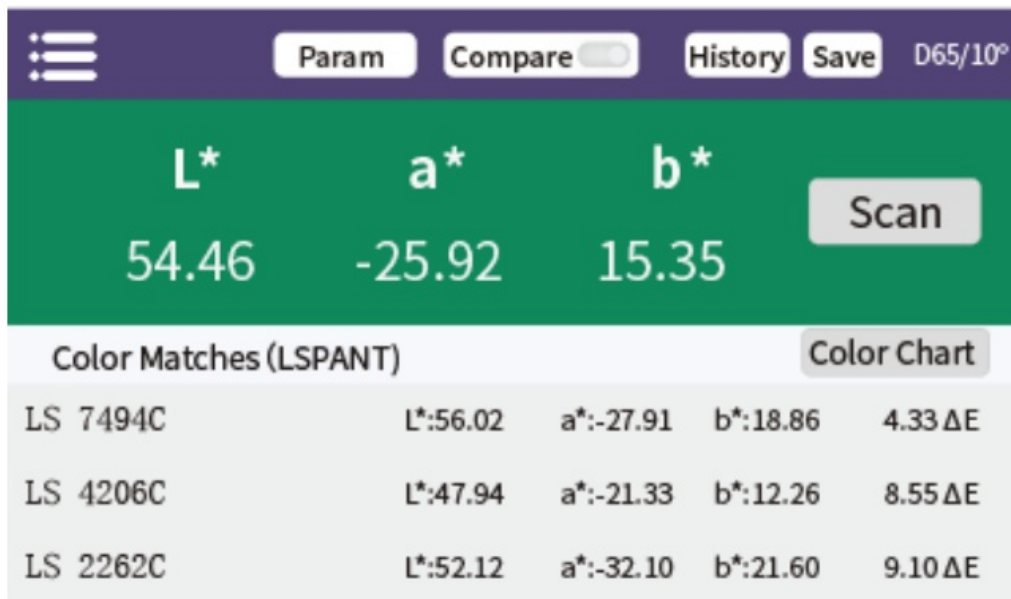

The screenshot shows a measurement interface with a purple header bar containing a menu icon, 'Param', 'Compare' (with a toggle switch), 'History', 'Save', and 'D65/10°'. Below the header, a green section displays color data: L* (54.46), a* (-25.92), b* (15.35), and a 'Scan' button. Below this is a table titled 'Color Matches (LSPANT)' with a 'Color Chart' button. The table lists three matches: LS 7494C, LS 4206C, and LS 2262C, each with L*, a*, b* values and a ΔE value.

Color Matches (LSPANT)				Color Chart
LS 7494C	L*:56.02	a*:-27.91	b*:18.86	4.33 ΔE
LS 4206C	L*:47.94	a*:-21.33	b*:12.26	8.55 ΔE
LS 2262C	L*:52.12	a*:-32.10	b*:21.60	9.10 ΔE

Compare off

Matching colors: According to the selected color chart by the user, the closest three color chart numbers are matched for display. The following steps can be performed after completing measurement:

Save

To save the comparison result, press “Save” on the upper right corner. You can choose default or name the record manually.


History

If users need to view the saved results, click “History” to view the records, and records can be paged, deleted, renamed and searched.

Color chart

Select one set of color charts for matching, then samples will be only compared within this chart. Three closest color codes will be matched.

Menu bar

In the measurement interface, click the “” icon on the upper left corner to pop up the menu bar, with the following options: Language, Calibration, Param Settings, Device info, Reset to defaults and Power off.

Calibration

Can enter the calibration interface to perform the calibration operation of the instrument.

Language

Support Simplified Chinese and English

Param Settings

Select the parameter desired to be measured and displayed. Shortcut keys for the same functions are available at the top of the measurement interface

Device Info

“Device Info” displays the basic information of the instrument, including: Serial Number, Firmware Version, and Device Status.

Reset to defaults

To restore parameters to the factory settings and clear the data in the measurement interface

Power Off

Press “Power off” to power off the instrument.

Bluetooth Connection

App Installation

The colorimeter APP supports Android operating system 7.0 and above as well as iOS system 7.0 and above. Choose one of the following methods to install:

1. Scan the QR code on the cover with a mobile browser or WeChat, follow the prompts to download and install L scolor.
2. Search for “LScolor” from the iOS App Store or Google Play App market and follow the prompts to download and install LScolor.



App icon

Precautions

Open the App during the installation process or after the first installation, the mobile phone will prompt the permission of setting, the customer needs to allow all of them, otherwise the App cannot search for the device, and App cannot be used.

Device connection

Turn on the colorimeter, open the App, after the display of the LOGO. Enter the “Connect Device” interface and automatically start searching for Bluetooth devices. The searched devices will be displayed at the bottom of the interface. Click the device that matches the SN of the colorimeter and wait for the device to connect. Once connected, it will jump to “Calibration interface.”

Precautions

1. If there is a Bluetooth connection error or failure, you can restart the App and the colorimeter, or turn off/start Bluetooth in the settings interface and reconnect.
2. Click the “Browse” button. App is not connected to the colorimeter, you can enter the App to view historical data.

APP operation

The APP has the same function of color scan and color compare as the instrument with the same operation. The only difference is that the APP has color sharing function, and can import and export colors.

Color sharing

There are two ways to share colors:


- Scan a color in the “Color Scan” interface, press the “Share” icon to share the scanned color in image format or data format (the sharing methods are: WhatsApp,

Facebook, email, etc.).

- In the “My Colors” folder, select the colors to share and click “Share” icon to share the colors in data format.


Color import

Users can import colors shared by others into their own App, but the colors must be scanned by the same model of colorimeter. Both single color and batch of colors can be imported as follows: Take WeChat sharing and import on Android as an example, find the file to import with a name of “import_mycolor” and long press the file to be imported.

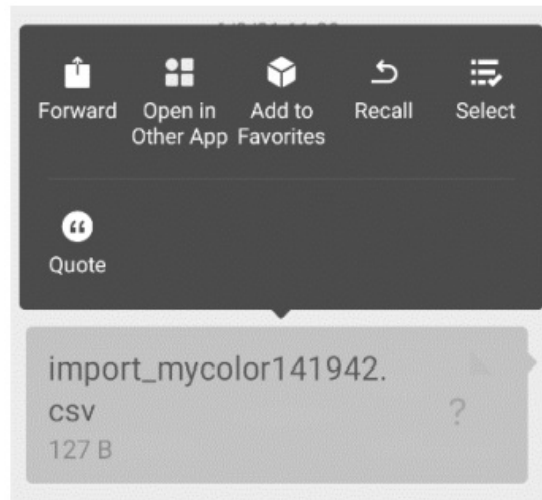

If use the iOS system, click the file and download it first, and then select “Open with other app” to realize data sharing.


- If the above interface appears, click “Open with other app” and select the App “LScolor”, after confirmation, jump to the import interface and select the file to complete the import.

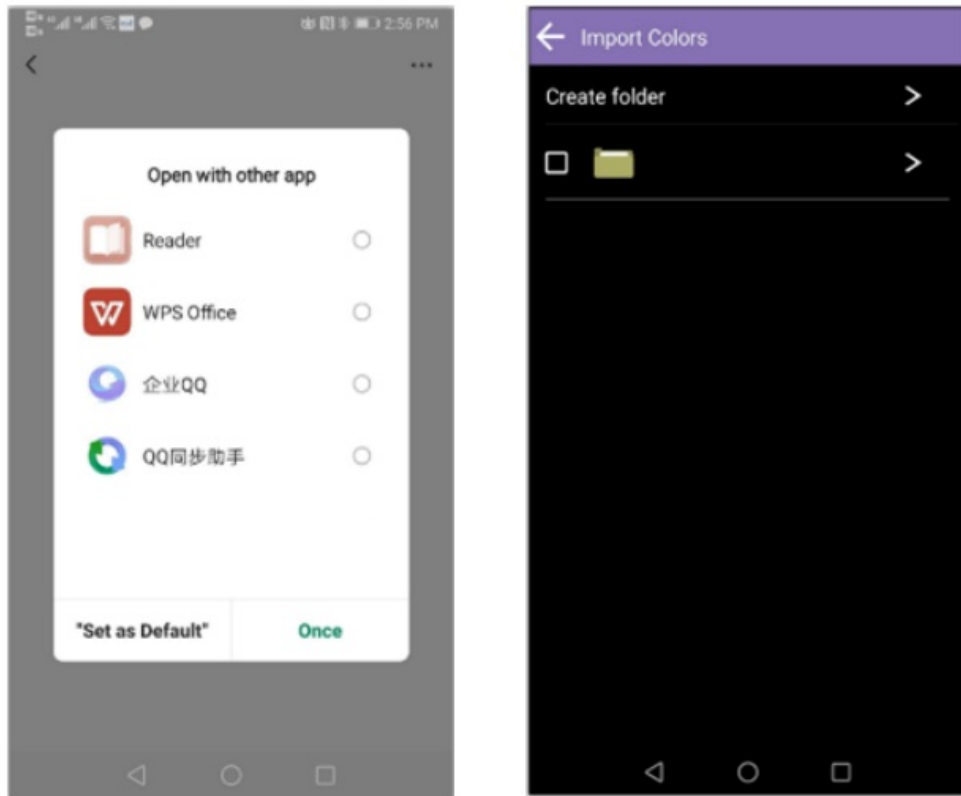

- If the import has completed, it will prompt “Import successful”. After pressing “Confirm”, the imported color will be displayed. After the import is completed, the import can be checked as the first record under the folder.
- If the product model does not match, it will prompt “Unmatched device model “. Please connect to the correct model before importing.

PC Software

You can use the PC software to connect the computer via USB. The software has functions of color compare measurement, reading compare record, export compare data to Excel, spectral reflectance measurement, qualified number, unqualified number, total number statistics, report generation and printing, etc. For specific operation, please refer to the “Color Haze Meter Software Operation Manual”.

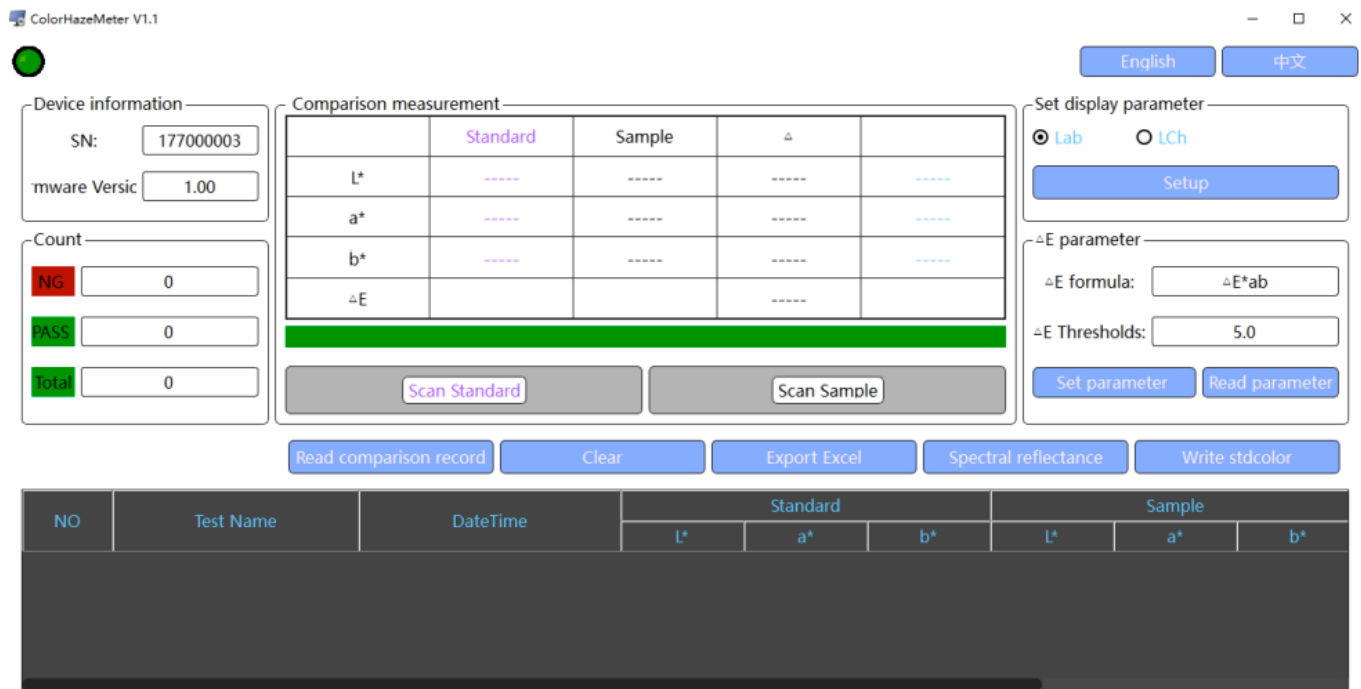

Precaution

1. When the colorimeter has not been used for a long time, it is recommended to perform calibration before using.
2. Please ensure that the sample is evenly colored with a flat and clean surface; otherwise, it will affect the measurement accuracy.
3. The calibration cover should be attached to the iron shell after it is removed from the main unit to avoid fouling of the standard plate.
4. Do not apply any objects to the glass of the measuring port, as this will damage the instrument and affect the measurement accuracy and operational safety.
5. When the instrument is not in use, put on the calibration cover to prevent dust from entering or prolonged humidity, which may affect the measurement accuracy.
6. It's recommended to calibrate the instrument once each year. And we offer calibration service.
7. Due to color differences in the LCD display, the color displayed on the instrument screen is for reference only.

Plastic cuvettes

In order to meet the measurement requirements of masterbatches, powders, and opaque liquids, the instrument is equipped with two plastic cuvettes.



Packing List

No.	Description	Quantity	Unit
1	Colorimeter	1	Set
2	USB cable	1	pcs
3	Cleaning cloth	1	pcs
4	Power cable	1	pcs
5	User manual	1	pcs
6	Certificate / Warranty card	1	pcs
7	Calibration Report	1	pcs
8	Plastic Cuvettes	2	pcs


Service

1. The colorimeter has one-year warranty. If the meter works abnormally, please send the whole instrument to our company for maintenance.
2. Provide users with spare parts and lifelong maintenance services.
3. Provide the users with the gauge calibration service.
4. Free technical support for long term.

- **Manufacturer:** Shenzhen Linshang Technology Co., Ltd.
- **Website:** www.linshangtech.com
- **Service hotline:** +86-755-86263411

- Email: sales21@linshangtech.com

Documents / Resources

	Linshang LS177 Spectrophotometer [pdf] User Manual LS177, LS177 Spectrophotometer, Spectrophotometer
-----------------------------------------------------------------------------------	-------------------------------------------------------------------------------------------------------------------------

References

- [User Manual](#)

Linshang, LS177, LS177 Spectrophotometer,

Linshang Spectrophotometer

Leave a comment

Your email address will not be published. Required fields are marked *

Comment *

Name

Email

Website

☐ Save my name, email, and website in this browser for the next time I comment.

Post Comment

Search:

e.g. whirlpool wrf535swhz

Search

[Manuals+](#) | [Upload](#) | [Deep Search](#) | [Privacy Policy](#) | [@manuals.plus](#) | [YouTube](#)

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.