

Fast and easy maintenance- the Inogen benefit.

INOGEN ONE[®] G5/Inogen Rove 6 Column Change Guide



Your Equipment Provider:

Your Physician:

 
Inogen, Inc.
859 Ward Drive Suite 200.
Goleta, CA 93111
1-877-4-INOGEN (Toll Free)
1-877-446-6436 (Toll Free)
1-805-562-0515 (International)
info@inogen.net
www.inogen.com

96-09240-00-01 D



Column Change Guide

INOGEN ONE® G5/ Inogen Rove 6


To change the columns on your Inogen One® G5/Rove 6, first:

- 1 Turn off the Inogen One® G5/Rove 6 concentrator by pressing the power button to shut down the device.
- 2 Remove the Inogen One® G5/Rove 6 concentrator from the carrying case (if attached).




Image of device may vary.


NOTE Column change instructions are only to be used when maintenance is required and are not intended for practice purpose.

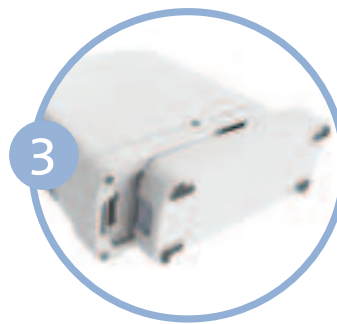


Remove the battery from the Inogen One® G5/Rove 6 concentrator.

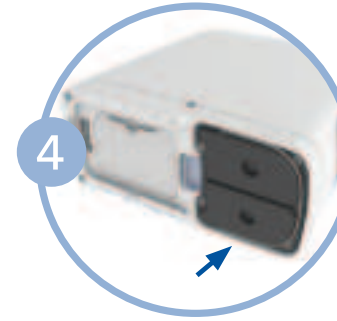


Place the Inogen One® G5/Rove 6 concentrator on its side so that the underside is visible. The metal column assembly can be seen on one side of the device.

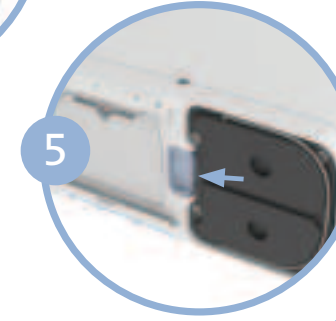


Unlock the column assembly by pushing the latch button away from the columns.

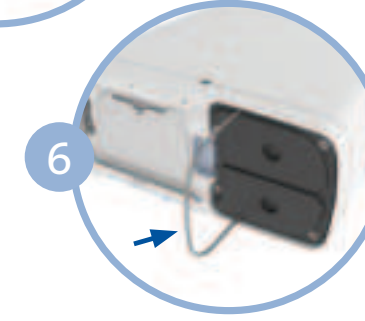


While holding the button open, slide the column assembly out of the device by pulling on the column pull handle.

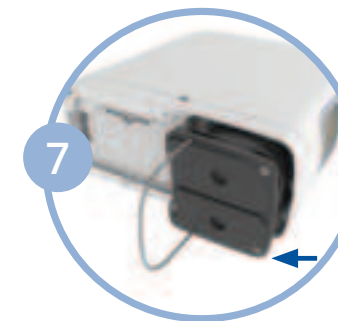

Remove the columns completely from the Inogen One® G5/Rove 6. Both columns are removed as one piece.

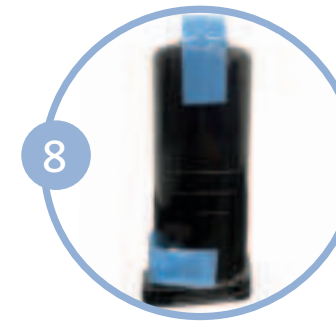


If tape is applied to column assembly, remove tape

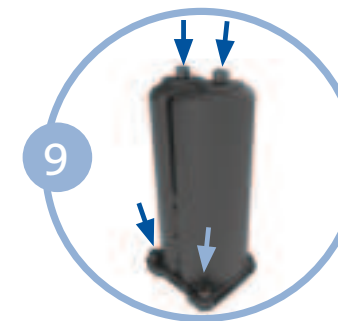


Remove dust caps of new coulumn assembly. Make sure there is no dust or debris where the dust caps were located.



Insert column assembly into the Inogen One® G5/Rove 6 concentrator. Do not leave the column ends exposed; column assembly should be inserted into the Inogen One® G5/Rove 6 as soon as the dust caps have been removed.

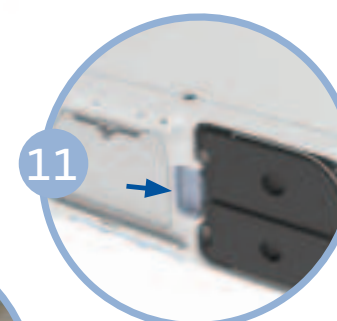

Push the column assembly into the device such that the columns are fully seated into the Inogen One® G5/Rove 6 concentrator. The spring loaded latch button should fully return to the closed position.

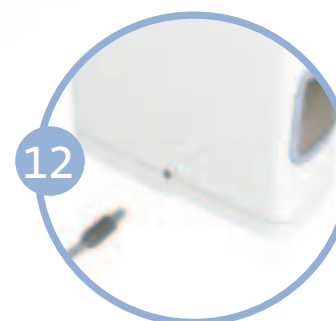


Connect the AC power supply cord to the Inogen One® G5/Rove 6 and plug the power supply AC cord into an electrical outlet. Do Not Power on the Inogen One® G5/Rove 6 concentrator.

Continues on back →



Press the alert/bell button once and screen will display "sieve reset success." Press the power button to turn on the Inogen One® G5/Rove 6, and use normally.

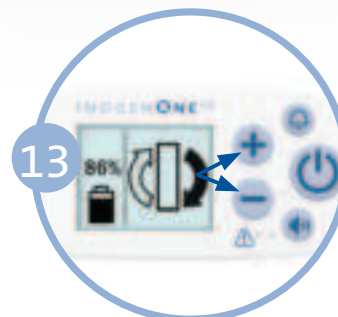

Press and hold the plus (+) and (-) minus buttons for 5 seconds. The screen will display "sieve reset" message. Release buttons once message is displayed on screen.