

For further assistance, please
visit [Tenovi.com/support](https://tenovi.com/support).

tenovi

© 2025 Tenovi Co.

All rights reserved. Republication, duplication or modification
in whole or in part is prohibited without prior written permission
from Tenovi. All trademarks are owned by or licensed to
Tenovi, Co. unless otherwise noted.

Revision 05.01.2024


Power Your Gateway

Plug your Gateway into a power outlet with the provided power cord and wait until the LED on the Gateway turns **RED**.



Insert Batteries



Insert batteries into the back of the glucometer. **Stand within 10 feet** of the gateway and **press the top button** on the glucometer. The gateway will beep indicating that it has connected successfully to the glucometer.

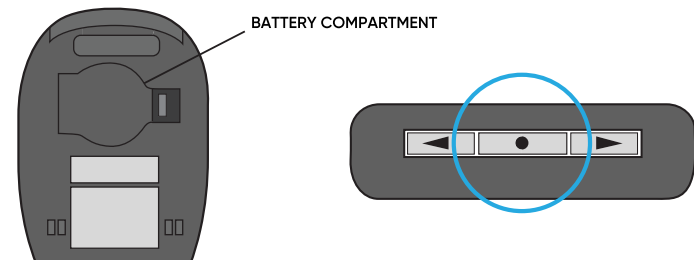

Preparation

3


Insert a test strip with the contact bars facing upwards into the meter's test strip port. **Insert a new disposable lancet firmly into the lancet holder.**

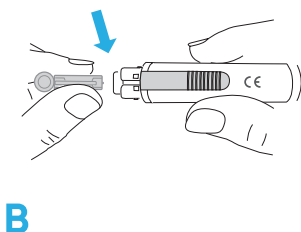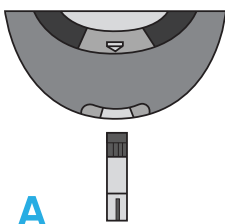

Transmit Readings

5



Remove the test strip and the screen will turn off. Stand **within 10 feet** of the gateway and **press the power button** on top of the glucometer to transmit the reading to the Gateway. You can take readings while away from the Gateway. **Press this button** the next time you are close to your Gateway to transmit all of your readings.

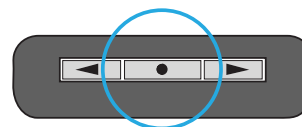

Take A Sample

4

Use prepped lancet to take a blood sample. Quickly apply the blood sample to the end of the test strip. When the measurement concludes, the results are shown on the display.

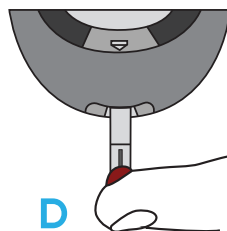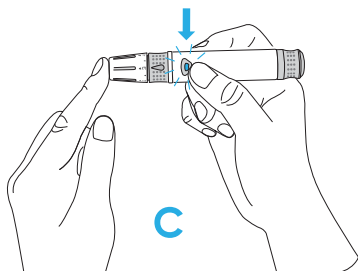

Daily Measurements

6

Take one measurement each day unless advised otherwise by your healthcare provider. The LED will turn **RED** each morning and then turn **GREEN** once a measurement has been successfully transmitted.



Measurement Needed



Measurement Successful