




NATIONAL INSTRUMENTS PXIe-4322 Analog Output Module User Guide

[Home](#) » [NATIONAL INSTRUMENTS](#) » NATIONAL INSTRUMENTS PXIe-4322 Analog Output Module User Guide



NATIONAL INSTRUMENTS PXIe-4322 Analog Output Module

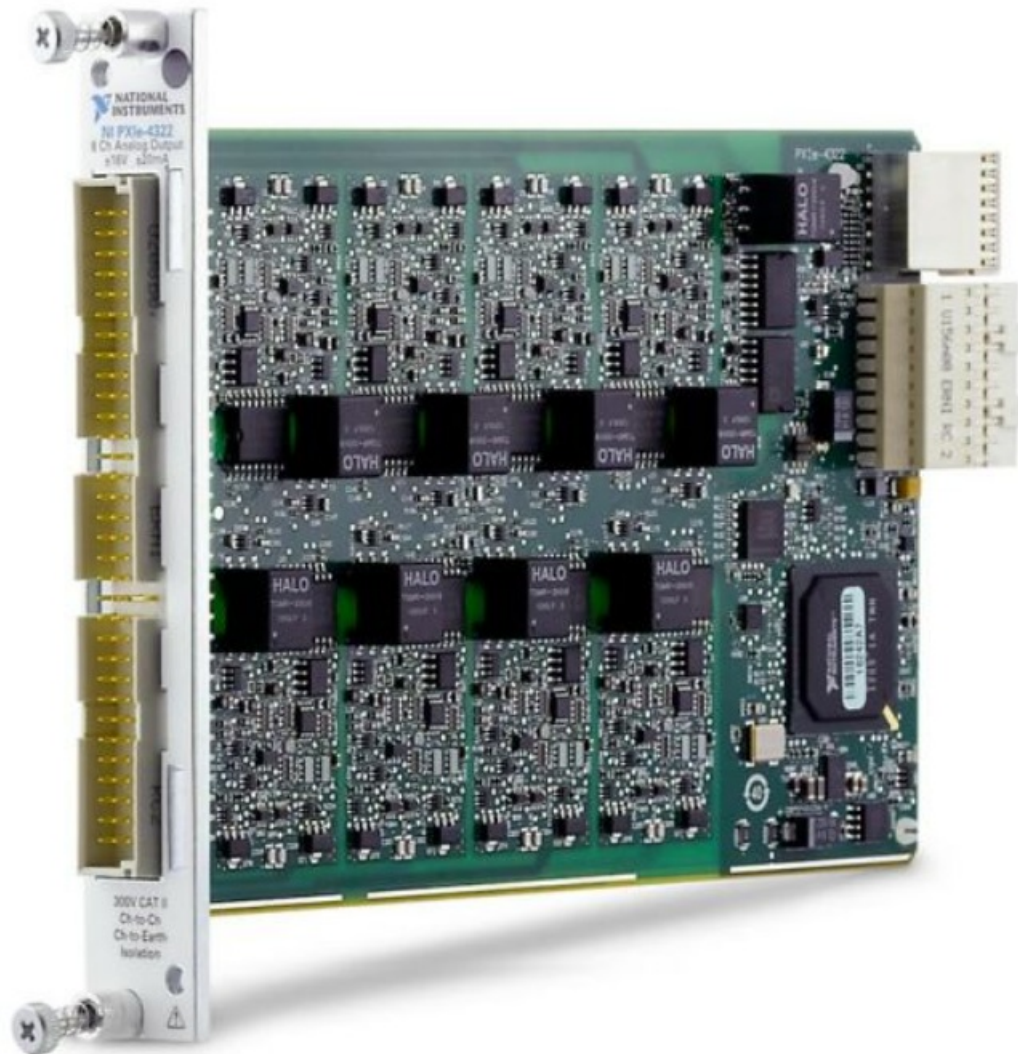

Contents

- [1 COMPREHENSIVE SERVICE](#)
- [2 SELL YOUR SURPLUS](#)
- [3 OBSOLETE NI HARDWARE IN STOCK & READY TO SHIP](#)
- [4 Manufacturer](#)
- [5 Board Assembly Part Numbers](#)
- [6 Volatile Memory](#)
- [7 Non-Volatile Memory \(incl. Media Storage\)](#)
- [8 Procedures](#)
- [9 Terms and Definitions](#)
- [10 Customer Support](#)
- [11 Documents / Resources](#)
 - [11.1 References](#)
- [12 Related Posts](#)

COMPREHENSIVE SERVICE

We offer competitive repair and calibration service, as well as easily accessible documentation and free downloadable resources.

SELL YOUR SURPLUS

We buy new, decommissioned, and surplus parts from every NI series. We work out the best solution to suit your individual needs.



Sell For Cash



Get Credit Receive a



Trade-In Deal

OBSOLETE NI HARDWARE IN STOCK & READY TO SHIP

We stock New, New Surplus, Refurbished, and Reconditioned NI Hardware.

Bridging the gap between the manufacturer and your legacy test system.

Request a Quote  [CLICK HERE PXIe-4322](#)

Manufacturer

National Instruments

Board Assembly Part Numbers

(Refer to Procedure 1 for identification procedure):

Part Number and Revision	Description
153917C-01L or later	PXI-4322

Volatile Memory

N/A

Non-Volatile Memory (incl. Media Storage)

Target Data	Type	Size	Battery Backup	User1 Accessible	System Accessible	Sanitization Procedure
Channel configuration	CPLD	Lattice LCMXO2-1200	No	No	No	None
Device configuration <ul style="list-style-type: none"> • Device configuration information • Calibration metadata • Calibration data2 	Flash	1 MB	No	No Yes No	Yes Yes Yes	None Procedure 2 None

1 Refer to Terms and Definitions section for clarification of User and System Accessible

2 Calibration constants that are stored on the device include information for the device's full operating range. Any implications resulting from partial self-calibration can be eliminated by running the full self-calibration procedure.

Procedures

Procedure 1 – Board Assembly Part Number identification:

To determine the Board Assembly Part Number and Revision, refer to the label applied to the surface of your product. The Assembly Part Number should be formatted as "P/N: #####a-##L.

Procedure 2 – Device Configuration Flash (Calibration Metadata):

The user-accessible areas of the Device Configuration Flash are exposed through a calibration Applications Programming Interface (API) in LabVIEW. To clear the calibration meta-data area, complete the following steps:

1. To clear the calibration password, overwrite the existing Calibration Password twice with non sensitive password using the DAQmx Change External Calibration Password function in the NIDAQmx Calibration API.
2. To clear the user-defined information, overwrite the existing User-Defined Information twice with non-sensitive information using the User-Defined Information property in the NI-DAQmx Calibration API.

Terms and Definitions

Cycle Power:

The process of completely removing power from the device and its components and allowing for adequate discharge. This process includes a complete shutdown of the PC and/or chassis containing the device; a reboot is not sufficient for the completion of this process.

Volatile Memory:

Requires power to maintain the stored information. When power is removed from this memory, its contents are

lost. This type of memory typically contains application specific data such as capture waveforms.

Non-Volatile Memory:

Power is not required to maintain the stored information. Device retains its contents when power is removed. This type of memory typically contains information necessary to boot, configure, or calibrate the product or may include device power up states.

User Accessible:

The component is read and/or write addressable such that a user can store arbitrary information to the component from the host using a publicly distributed NI tool, such as a Driver API, the System Configuration API, or MAX.

System Accessible:

The component is read and/or write addressable from the host without the need to physically alter the product.

Clearing:

Per NIST Special Publication 800-88 Revision 1, “clearing” is a logical technique to sanitize data in all User Accessible storage locations for protection against simple non-invasive data recovery techniques using the same interface available to the user; typically applied through the standard read and write commands to the storage device.

Sanitization:

Per NIST Special Publication 800-88 Revision 1, “sanitization” is a process to render access to “Target Data” on the media infeasible for a given level of effort. In this document, clearing is the degree of sanitization described.

Notice: This document is subject to change without notice. For the most recent version, visit ni.com/manuals

Customer Support

Contact: 866-275-6964

support@ni.com



1-800-915-6216



www.apexwaves.com



sales@apexwaves.com

All trademarks, brands, and names are the property of their respective owners.



Documents / Resources



[NATIONAL INSTRUMENTS PXIe-4322 Analog Output Module](#) [pdf] User Guide
PXIe-4322, PXIe-4322 Analog Output Module, Analog Output Module, Output Module, Module

References

- [NI Product Documentation - NI](#)
- [AW PXIe-4322 National Instruments PXI Analog Output Module | Apex Waves](#)

Manuals+.