

Palintest Lumiso Expert Benchtop Kit User Guide

Home » Palintest » Palintest Lumiso Expert Benchtop Kit User Guide 🖫





Lumiso Expert Quick Start Guide



Contents

- 1 Registration and technical support
- 2 Getting started
- 3 Achieving accurate results
- 4 Documents / Resources
 - 4.1 References
- **5 Related Posts**

Registration and technical support

Thank you for choosing a Palintest Lumiso

Expert Photometer. Please register your instrument now at www.palintest.com/products/lumiso-expert or scan the QR Code.




There, you will also find invaluable support for your new instrument, including:

- User Manual
- · Supporting documents for all test methods
- · Firmware updates
- · Palintest Connect cloud based data storage and management

Getting started

Onscreen instructions are available for all test methods.

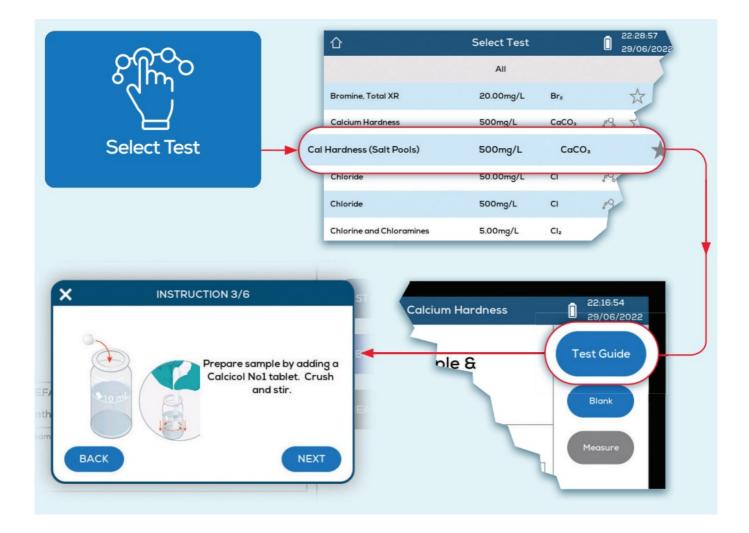

Achieving accurate results

Blanking


- Blanking is critical and sets the zero level for each test.
- Usually, the blank is a sample without any reagents added.
- The onscreen instructions show blanking as the first few steps, however, the "blank" can be prepared in a separate tube and blanking carried out just prior to the measurement step.

Best practice when testing

Rinse equipment thoroughly with deionised water or the water being tested.




Measure volumes accurately as shown (or use a pipettor or syringe).

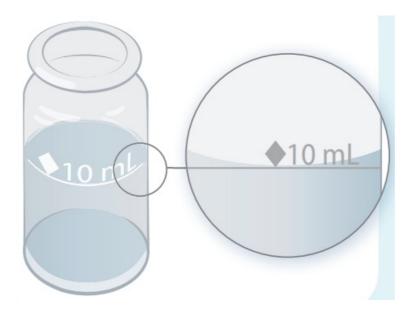


Use the correct Palintest reagents and observe all safety information supplied with them.




Check for bubbles and remove as shown.

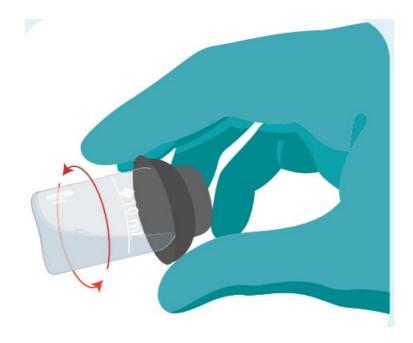

Align tubes as shown. Always use the tube cap as it excludes light from the optics.

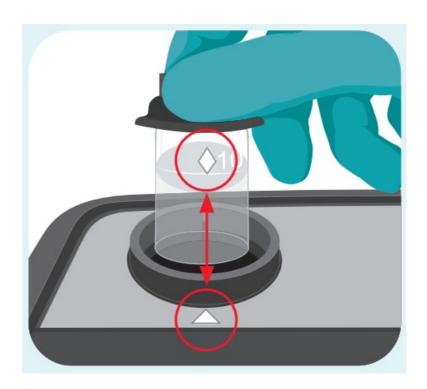

Keep instrument clean. Dry all tubes before placing in the optical cell.

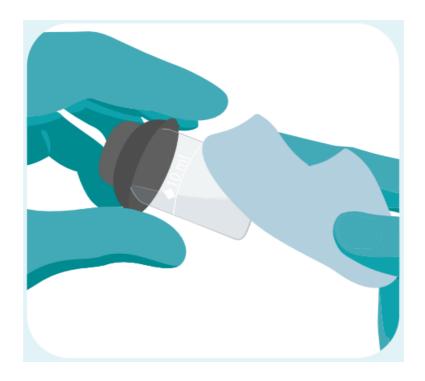

Verify instrument calibration

Check Standards are available for your Lumiso Expert. Use these regularly to confirm that the instrument is within calibration.





A Halma company
©Palintest Limited 2022. All rights reserved.
Registered office: Palintest House, Kingsway, Team Valley,
Tyne & Wear, England NE11 0NS.

Registered in England, Number 1204118 ZIINST7210 V1

Documents / Resources



Palintest Lumiso Expert Benchtop Kit [pdf] User Guide Lumiso Expert Benchtop Kit, Expert Benchtop Kit, Benchtop Kit

References

- P Palintest Water Analysis Technologies
- Palintest Water Analysis Technologies
- P Water Analysis Products, Water Safety Testing Products Palintest
- P Lumiso Expert Palintest

Manuals+,