

spokeo
spokeo
Bluetooth
Streaming
From Various
Headset



spokeo Bluetooth Streaming From Various Headset User Guide

[Home](#) » [spokeo](#) » spokeo Bluetooth Streaming From Various Headset User Guide 

Contents

- 1 [spokeo Bluetooth Streaming From Various Headset](#)
- 2 [Product Usage Instructions](#)
- 3 [ABOUT SPOKEO](#)
- 4 [BOX CONTENTS](#)
- 5 [GETTING STARTED WITH SPOKEO](#)
- 6 [PARTS AND CONTROLS](#)
- 7 [HANDLING, SECURITY, WARNINGS AND WARRANTY](#)
- 8 [Cleaning and maintenance](#)
- 9 [Volume level warning and hearing damages](#)
- 10 [TROUBLESHOOTING](#)
- 11 [TECHNICAL SPECIFICATIONS](#)
- 12 [Documents / Resources](#)
 - 12.1 [References](#)
- 13 [Related Posts](#)

spokeo

spokeo Bluetooth Streaming From Various Headset



Specifications

- Product: Spokeo Assistive Listening Device
- Manufacturer: Audiome S.A.S.
- Features: Supports direct conversations and Bluetooth streaming, high-definition sound quality
- Technology: Patented bi-conduction technology for air and bone sound transmission
- Components: Smart Mic, Sound Boost, Charging Cable, Wall Chargers, Protection Case, Instruction Guides

Product Usage Instructions

Getting Started

To charge the Spokeo Smart Mic and Sound Boost together:

1. Press and hold the orange button on the Smart Mic for 3 seconds until the blue LED light turns on.
2. Slide the headset's orange button to the Bluetooth setting.

Parts and Controls

The Spokeo device includes the following controls:

- Power indicator LED light
- Reset button
- Bluetooth/Microphone mode switch
- ON/OFF switch – Volume OFF/1/2 switch
- Microphone magnetic charging pogo pin
- Type C USB charging port
- 3.5 mm Jack line audio input
- Charging indicator LED light

Handling, Security, Warnings, and Warranty

When using the Spokeo Sound Boost headset

- Identify left (L) and right (R) indicators for correct positioning.
- Adjust for a comfortable fit on the interlocutor's head.

For the Spokeo Smart Mic microphone

- Attach to your shirt collar using the built-in clip.
- Position close to your mouth for optimal audio capture.

Cleaning and maintenance

- Use electronic cleaning wipes for sanitizing.
- Avoid corrosive materials or harsh chemicals.
- Store in a dry, moderate-temperature environment for longevity.

FAQ

Q: How do I know if Spokeo is fully charged?

A: The charging indicator LED light will turn green when fully charged.

ABOUT SPOKEO

- Product overview and key features
- Spokeo is a versatile assistive listening device designed for effortless communication. It supports both direct one-on-one conversations and Bluetooth streaming from various devices such as TVs, phones, tablets, and radios.
- Spokeo is a remote microphone/headset bundle offering high-definition sound quality through a patented bi-conduction technology combining air and bone sound transmission.
- Supplier information and customer service contact
- Spokeo is manufactured by Audiome S.A.S.

For customer service inquiries, please contact Audiome headquarter at:

- Tel.: +33 (0)1 83 64 72 27
- Email: Info@myspokeo.com
- Address: 60 rue de Fontenay, 92350 Le Plessis-Robinson, France For further information, visit: www.myspokeo.com

BOX CONTENTS

- 1 Spokeo Smart Mic (Microphone)
- 1 Spokeo Sound Boost (Headset)
- 1 USB-C / USB-A Charging Cable (1.8 m) 1 US Wall Charger block
- 1 EU Wall Charger block

- 1 EVA Protection Case
- 1 Instruction Guide Manual
- 1 Quick Starter Guide Manual

GETTING STARTED WITH SPOKEO

Charging Instructions:


Spokeo Smart Mic (microphone) and Spokeo Sound Boost (headset) have been engineered to be charged together seamlessly.

- First, attach Smart Mic to Sound Boost's pogo pin socket (see 5 and 11 on page 10).
- Once properly inserted, plug the USB C cable to the headset charging port (see 6 on page 10) and to a power source on the other end (wall charger blocks or any other power source).
- While charging, a red indicator light will be visible, which will turn green once fully charged (see 8 on page 10).


Using Spokeo in microphone mode 3 easy steps

- Powering ON Smart Mic (microphone)

Press and hold the orange button for 3 seconds until the blue LED light turns on. (see 9 on page 10).

- Switching the headset to microphone mode 

Slide the headset's orange button to the microphone setting (see 3 on page 10).


Once connected to Smart Mic, Sound Boost notifies  System connected

- Powering ON Sound Boost (headset)


Slide the black button from OFF to the desired volume preset 1 or 2. (see 4 on page 10).

Once powered ON, Sound Boost notifies  Power ON  and System connected once connected to the Smart Mic?

Switching OFF Spokeo

- Sound Boost: Turn Sound Boost off by sliding the black button to OFF (see 4 on page 10). Sound Boost notifies  Power off and power indicator LED light (see 1 on page 10) switches from blue to red to off.
- Smart Mic: Press and hold the microphone ON/OFF switch button (see 9 on page 10) for 3 seconds and release. The microphone LED light indicator switches from blue to red to off. (see 10 on page 10).

Using Spokeo in Bluetooth mode 3 easy steps

- Switching the headset to Bluetooth mode 

Slide the headset's orange button to the Bluetooth setting  (see 3 on page 10).

- Powering ON Sound Boost (headset)

Slide the black button from OFF to the desired volume preset 1 or 2.


(see 10 on page 10). Once powered ON, Sound Boost notifies  Power ON.

- Connecting to an external device:

Turn on Bluetooth on the external device and search for new device to connect to. Select «Spokeo» from the list of available devices and press “Pair”.


Once connected to the external Bluetooth device, Sound Boost (headset) notifies System connected.

Switching OFF Spokeo Sound Boost:

Turn Sound Boost by sliding the black button to OFF (see 4 on page 10). Sound Boost notifies  Power off, and power indicator LED light (see 1 on page 10) switches from blue to red to off.

PARTS AND CONTROLS

1.

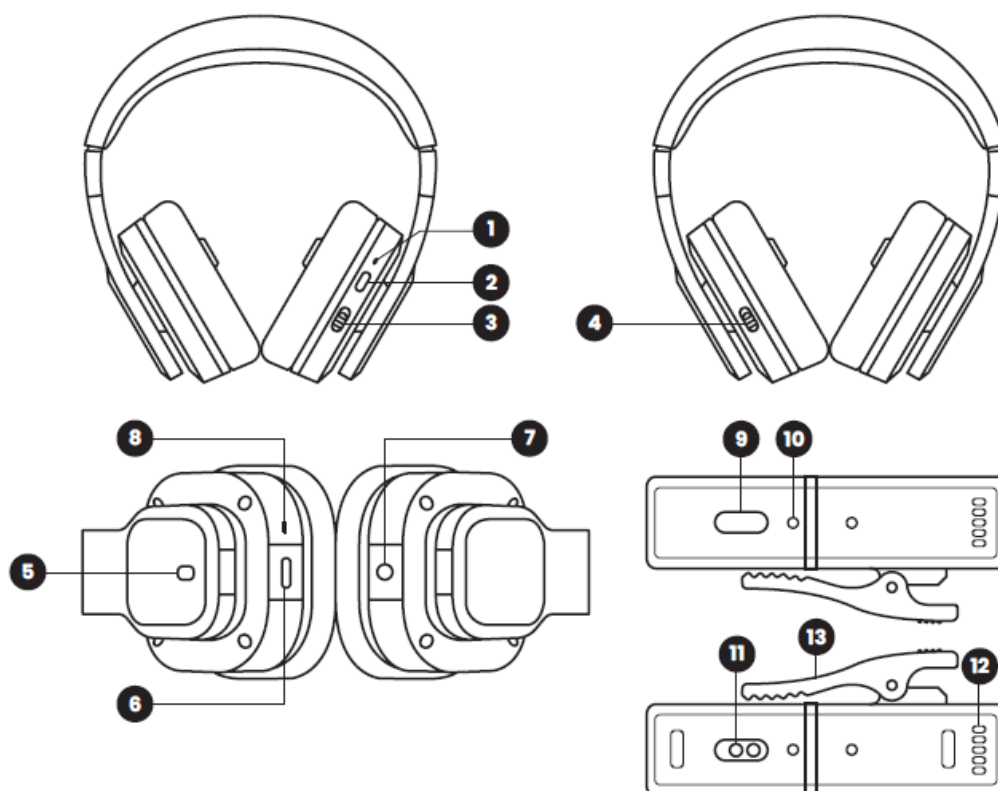

Power indicator LED light

2. Reset button
3. Bluetooth /Microphone mode switch
4. ON/OFF switch – Volume OFF/1/2 switch
5. Microphone magnetic charging pogo pin socket female
6. Type C USB charging port
7. 3.5 mm Jack line audio input
8. Charging indicator LED light
9. ON/OFF switch
10. . Microphone LED light indicator
11. Microphone magnetic charging pogo pin male

12. Microphone sound input

13. Clip

HANDLING, SECURITY, WARNINGS AND WARRANTY

Spokeo Sound Boost (Headset)

Identify the left (L) and right (R) indicators on the headset. Carefully position the headset on your interlocutor's head according to the L and R markers.

Adjust the slide-sliders for a comfortable fit.

Spokeo Smart Mic (Microphone)

Attach the microphone to your shirt collar using the built in clip (see 13 on page 10). Position the microphone close to your mouth for optimal audio capture.

Cleaning and maintenance

Use cleaning wipes designed for electronics to sanitize both Sound Boost (headset) and Smart Mic (microphone). Avoid using any corrosive materials or harsh chemicals. Store Spokeo in a dry, moderate-temperature environment for optimal longevity.

Initial Setup

Perform a full charging cycle before initial use for both Smart Mic (microphone) and Sound Boost (headset).

Battery life and charging cycle

To optimize the battery lifespan, use Spokeo at least once every two weeks and recharge it after use. The charging time for a full charge is a maximum of 2 hours, providing up to 5 hours of continuous use. If needed, you can use another CCC certified compatible charger.

Volume level warning and hearing damages

For safety reasons, always start by setting your Sound Boost headset to volume preset 1. If you find the sound intensity insufficient, switch to volume preset 2.

When using Spokeo in Bluetooth mode, make sure to begin with the external Bluetooth source volume set to 50%, and adjust it as needed.

Warranty and after-sales service

Spokeo comes with a 2-years manufacturer's warranty against any manufacturing defects.

To request a warranty repair, visit myspokeo.com/policies/refund-policy

For any other inquiries, visit: www.myspokeo.com

TROUBLESHOOTING

Issues	Probable Causes	Proposed actions
Unable to boot	The battery level could be too low. When it reaches a low level, the power indicator LED flashes every 30 seconds, and Sound Boost audibly notifies 🎵 Battery low.	Charge the bundle headset/microphone.
Can't connect to a Bluetooth device	The external device is not Bluetooth compatible or not on discoverable mode.	Make sure external device Bluetooth setting is turned on and on discoverable mode.
	The distance between Sound Boost (headset) and the external device is too far.	keep the distance within 1 meter for pairing purpose and within 10 meters for streaming.
	The battery level might be too low.	Charge Sound Boost (headset).
Loud and disruptive noise	The battery level might be too low.	Charge Sound Boost (headset).
Staccato / interruption problems	Sound boost might be too far from Smart Mic or from the external Bluetooth device.	Bring the two devices closer together.
System does not respond /frozen	An electric incident might have occurred.	Turn both Smart Mic and Sound Boost OFF and ON.
Need more help ?	Follow the reset to manufacturer setting protocol. The reset procedure is performed through the headset only.	Ensure that the headset is sufficiently charged. Set the headset volume button to position 1. Press and hold the Reset button on the headset for 10 seconds until you hear a confirmation beep indicating the reset. Once reset, the headset indicator will flash rapidly blue.

TECHNICAL SPECIFICATIONS

Spokeo Sound Boost (Headset)

- Bluetooth version: V 5.0
- Wireless connection range: about 10 m
- Operating frequency range:
- -2480MHZ
- Build-in rechargeable lithium battery capacity: 500 mAh
- Charging socket(in): Type C USB
- Charging socket(out): Pogo-Pin 5.0V-400mA
- Normal charging time: about 2 hours
- Headset working time: ≥5H
- Moving-coilspeaker: Ø40mm
- Moving-coil speaker resistance: 32Ω±10%

- Moving-coil speaker rated power: 30mW
- Bone conduction speaker: Ø10.8mm
- Bone conduction speaker resistance: $8\Omega \pm 15\%$
- Bone conduction speaker rated power: 230mW
- Frequency response: 20Hz-20KHz
- Harmonic distortion: THD3% 1KHz
 - Preset 1 BT maximum sound pressure: 106 dB (at 1kHz)
 - Preset 2 BT maximum sound pressure: 111 dB (at 1kHz)

Spokeo Sound Mic (Microphone)

- IC solution: QCC3034
- Bluetooth version: V 5.0
- Transmission frequency: 2.4G-2.48G
- Transmission power: Class 2
- Bluetooth profile : A2DP, HFP
- Wireless range: $\geq 10\text{m}$
- Charging voltage: 5V
- Charging socket: Pogo-Pin
- Working time: $\geq 5\text{H}$
- Charging time: $\leq 2\text{H}$
- Battery: 3.7V 110mAh (Polymer Lithium Battery)
- Microphone sensitivity: $-45 \pm 2\text{dB}$
- Polar pattern: Unidirectional


Certifications

Spokeo has successfully cleared extensive quality tests and holds CE, FCC Union, United State and UK markets.



www.myspokeo.com

Documents / Resources

 <small>User guide Manuel d'utilisation</small>	spokeo Bluetooth Streaming From Various Headset [pdf] User Guide Bluetooth Streaming From Various Headset, Streaming From Various Headset, Various Headset , Headset
---	--

References

- [User Manual](#)

[Manuals+](#), [Privacy Policy](#)

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.